



Manuals.plus /

› nosiboo /

› Nosiboo Pro Electric Nasal Aspirator Instruction Manual

nosiboo Nosiboo Pro

Nosiboo Pro Electric Nasal Aspirator Instruction Manual

Brand: nosiboo | Model: Nosiboo Pro

INTRODUCTION

The Nosiboo Pro Electric Nasal Aspirator is designed to provide effective and gentle relief for congested babies and toddlers. Developed with special care, this device aims to make nose cleaning sessions easier for both children and parents. Its child-friendly design helps distract the child during use, while adjustable suction power ensures comfort and efficiency for different ages and levels of congestion.

The device is BPA-free, utilizing high-quality, soft materials that have undergone strict testing. Its unique Colibri head is engineered for easy cleaning, ensuring hygiene after each use.

Number 1 nasal aspirator in Korea

– from the heart of EU!



Image: The Nosiboo Pro Electric Nasal Aspirator, featuring its distinctive green and white design with the clear Colibri head and connecting tube.

SETUP

Before first use, ensure all components are clean. The Nosiboo Pro comes with a main unit, a Colibri head, connecting tube, and a power cord. The Colibri head consists of several parts that need to be assembled.

1. **Assemble the Colibri Head:** Insert the small internal piece into one half of the Colibri head. Then, attach the other half of the Colibri head, ensuring it snaps securely together.
2. **Attach the Nose Tip:** Select the appropriate nose tip (newborn or standard) and firmly attach it to the assembled Colibri head.
3. **Connect the Tube:** Attach one end of the flexible connecting tube to the Colibri head and the other end to the designated port on the main Nosiboo Pro unit.
4. **Power Connection:** Connect the power cord to the device, then plug it into a 110V electrical socket.

Air flow technology.

Effective. Reliable. Powerful. Gentle.

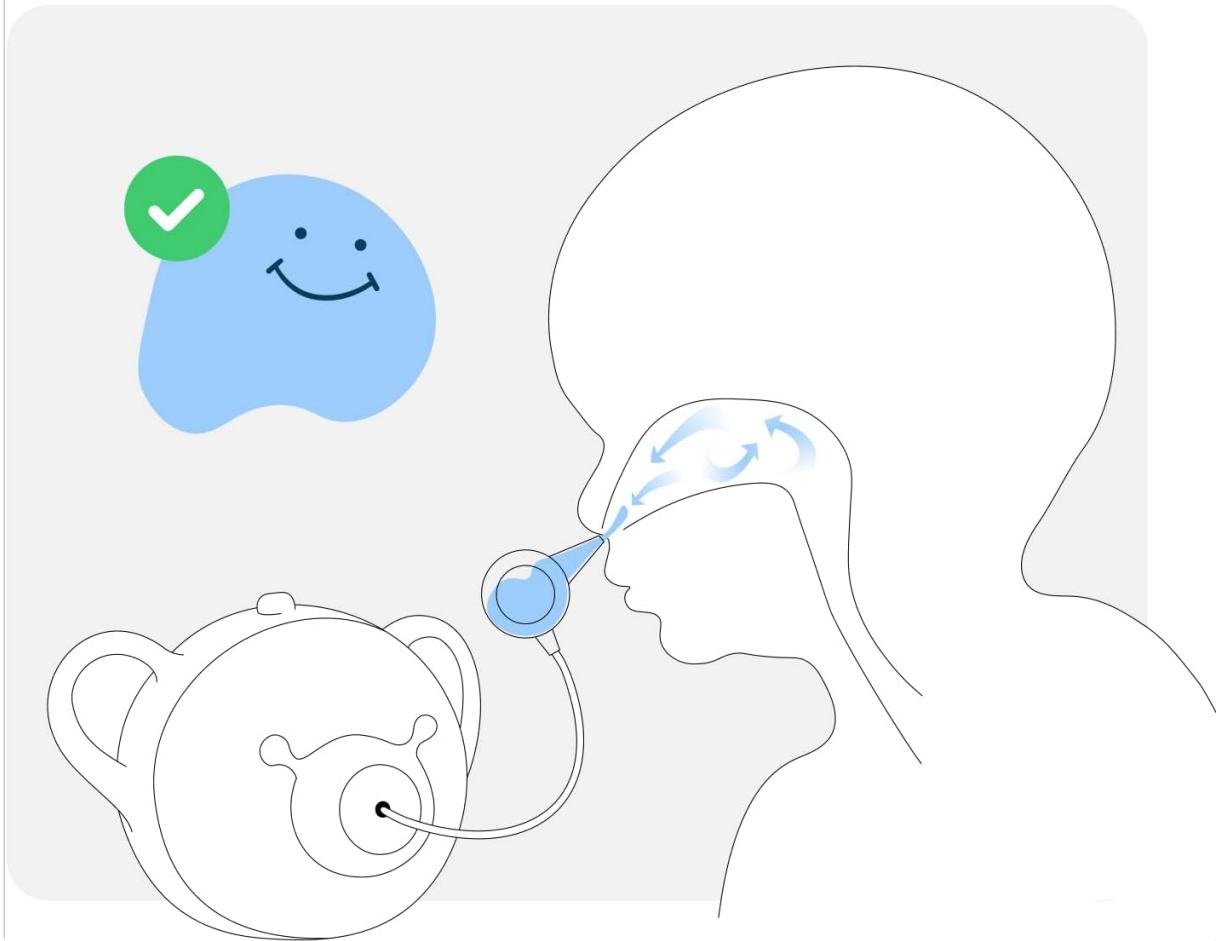


Image: All components of the Nosiboo Pro Electric Nasal Aspirator laid out, including the main unit, Colibri head, connecting tube, power cord, cleaning brush, extra air mesh pack, and standard nose tip, illustrating the contents of the box.

OPERATING INSTRUCTIONS

The Nosiboo Pro is designed for ease of use and effective mucus removal.

1. **Prepare the Child:** For dry, very thick mucus, it is recommended to dissolve it first with a mild saline solution. Wait a few minutes after applying the solution, then gently wipe any excess discharge with a tissue before using the aspirator.
2. **Adjust Suction Power:** Turn the control button on the main unit to adjust the airflow power. The suction can be set from minimum to maximum, depending on the child's age and the thickness of the mucus. The maximum power is designed to be efficient yet gentle.
3. **Position the Nose Tip:** Gently place the nose tip into one of the child's nostrils, ensuring a good seal.
4. **Activate Suction:** Turn on the device. The continuous suction will draw mucus into the Colibri head. Repeat for the other nostril as needed.
5. **Automatic Shut-off:** The device will automatically turn off after 13 minutes of continuous use. It can be restarted if further use is required. In cases of overuse or insufficient ventilation, a thermal switch will temporarily shut off the device for approximately 45 minutes for safety.

in the box

one piece of each and manuals



Nosiboo Pro



power cord



Colibri head with newborn's nose tip



connecting tube



standard nose tip



cleaning brush



extra air mesh pack

Image: A close-up view of the Nosiboo Pro's main unit, highlighting the adjustable suction power dial, which allows users to select between minimum and maximum settings for tailored use.

Video: A demonstration by a seller showing the Nosiboo Pro in action with a baby, illustrating how to use the aspirator for effective nasal clearing.

MAINTENANCE

Regular cleaning of the Nosiboo Pro is essential for hygiene and optimal performance. The Colibri head is the only part that requires cleaning after each use.

1. **Disassemble the Colibri Head:** After use, detach the Colibri head from the connecting tube. Separate the two halves of the Colibri head and remove the internal piece and nose tip.
2. **Clean Components:** Wash all disassembled parts of the Colibri head under running water. A cleaning brush is provided to ensure thorough cleaning of all crevices.
3. **Clean the Connecting Tube:** If mucus enters the connecting tube, detach it from the main unit and wash it thoroughly. Ensure the tube is completely dry before reassembly to prevent moisture buildup.
4. **Filter Maintenance:** The main unit contains a filter to prevent mucus from entering the device's motor. Check this filter periodically and replace it if it appears soiled or damaged. Replacement filters are available as part of an accessory set.
5. **Dry and Reassemble:** Allow all parts to air dry completely before reassembling the Colibri head and

connecting it back to the main unit.

Adjustable suction power

– two nose tips.

Special care for newborns.



Image: The Colibri head of the Nosiboo Pro disassembled into its individual parts, demonstrating how easily it can be taken apart for thorough cleaning after each use.

Video: A user-generated video demonstrating how to replace the fabric filter in the Nosiboo Pro, including a DIY method for creating new filters from suitable material.

TROUBLESHOOTING

If the Nosiboo Pro is not functioning as expected, please check the following:

- **No Suction or Weak Suction:** Ensure all parts of the Colibri head are correctly assembled and securely attached. Check that the connecting tube is firmly inserted into both the Colibri head and the main unit. Verify that the power cord is securely plugged into both the device and the electrical outlet. Inspect the internal filter for any blockages or damage; replace if necessary.
- **Device Not Turning On:** Confirm the power cord is properly connected. If the device recently shut off automatically due to overuse, wait approximately 45 minutes for the thermal switch to reset before attempting to restart.
- **Unusual Noises:** Ensure all components are correctly assembled. If unusual noises persist, discontinue use and contact customer support.

SPECIFICATIONS

Product Dimensions	7.56 x 7.64 x 7.56 inches; 3.62 Pounds
Item Model Number	green
Manufacturer	Attract Kft.
Brand Name	nosiboo
Model Name	Nosiboo Pro
Age Range Description	from birth
Material Type	Plastic
Material Type Free	Bisphenol-A (BPA) Free
Power Source	Corded Electric
Water Resistance Level	Not Water Resistant

WARRANTY AND SUPPORT

Each Nosiboo Pro is packaged in a protective foil inside the box. This foil should be intact upon unboxing. If the foil is damaged or missing, or if the box seal is broken, the device may have been returned previously and should be sent back to the Amazon fulfillment centers.

For any questions, concerns, or support regarding your Nosiboo Pro, please contact the NOSIBOO USA seller on Amazon. They will be glad to assist you with product inquiries, troubleshooting, or replacement parts.