

[Manuals.plus](#) /

› [AURAGLOW](#) /

› Auraglow Moldable Teeth Whitening Trays Instruction Manual

AURAGLOW AGTWK1

Auraglow Moldable Teeth Whitening Trays

Model: AGTWK1

INTRODUCTION

The Auraglow Moldable Teeth Whitening Trays are designed for comfortable and effective at-home teeth whitening. These trays are made from a 20% slimmer EVA material, allowing for a closer fit compared to standard trays. They are reusable, BPA-free, and latex-free, providing a custom impression of your teeth when molded correctly. This manual provides detailed instructions for molding, using, and maintaining your trays.

WHAT'S INCLUDED

Your Auraglow Moldable Teeth Whitening Trays package includes the following components:


- 4 Slim Moldable Trays (2 sets for upper and lower teeth)
- 1 Retainer Storage Case

Note: Whitening gel is not included with this product and must be purchased separately.

What's included



1 Retainer Case

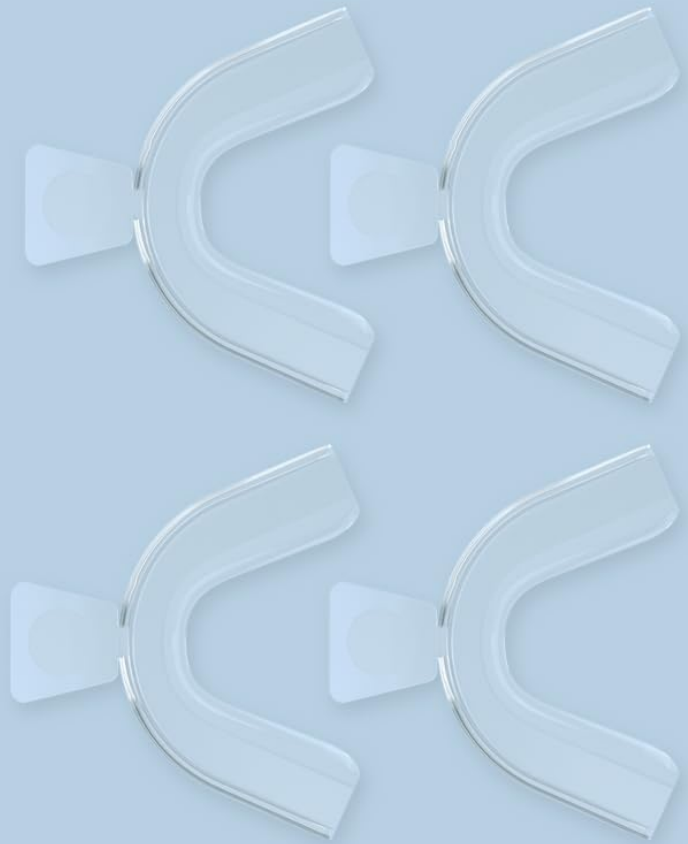

4 Slim Moldable Trays

Contents of the Auraglow Teeth Whitening Trays package.

SETUP: MOLDING YOUR TRAYS

Follow these steps carefully to mold your Auraglow trays for a custom fit:

1. **Step 1: Prepare Hot Water**

Bring a cup of water to a low boil. Ensure there is enough water to fully submerge one tray.

2. **Step 2: Soften the Tray**

Hold the first tray by the tab and dip it into the hot water for approximately 3-4 seconds. The tray will become soft and pliable.

3. **Step 3: Mold to Your Teeth**

Remove the softened tray from the water. Carefully place it in your mouth and bite down gently, pressing the tray against your teeth with your fingers to create a precise impression.

4. **Step 4: Repeat for Second Tray**

Repeat Step 2 and Step 3 for the second tray (for the other arch). If the fit is not satisfactory, you can re-dip the tray

in hot water and mold it again. An extra set of trays is provided in case of errors.

5. **Step 5: Trim for Comfort**

Once both trays are molded and cooled, you can trim off the front tabs using scissors. You may also trim any excess material around the edges for a more comfortable and precise fit.



Visual guide for molding the teeth whitening trays.

Instructions for molding



Step 1



Bring a cup of water to a low boil, deep enough to fully submerge the trays.

Step 2



Hold the first tray by the tab and dip it into the hot water for approximately 3-4 seconds.

Step 3



Remove the tray from the water, place it in your mouth and mold with your teeth. Repeat step 2 for the other tray.

Step 4



Try both trays on together to see how they fit. If necessary, dip in the water to mold again or use the extra set of trays provided.

Step 5



Trim off the front tabs. You can also trim the tray area if there is excess material for a better fit.

Step 6



Fill the tray with Auraglow whitening gel and wear based on the instructions.

Example of molded trays and key features.

OPERATING: USING YOUR TRAYS

Once your trays are molded, they are ready for use with your preferred teeth whitening gel (such as Auraglow teeth whitening gel, sold separately).

1. Apply Whitening Gel

Carefully apply a thin line of whitening gel into the front surface of each molded tray, ensuring even distribution. Avoid overfilling to prevent gel from contacting gums.

2. Insert Trays

Place the trays onto your upper and lower teeth, pressing gently to ensure the gel makes full contact with your tooth surfaces.

3. Wear Time

Wear the trays for the duration recommended by your specific whitening gel product. Typically, this ranges from 15 to 30 minutes per session.

4. Remove and Rinse

After the recommended time, remove the trays and rinse your mouth thoroughly with water to remove any residual gel.

MAINTENANCE

Proper care of your trays will ensure their longevity and hygiene:

- **Cleaning**

After each use, rinse the trays thoroughly with cool water to remove any remaining whitening gel. You may use a soft toothbrush and mild soap if necessary, but avoid abrasive cleaners.

- **Drying**

Allow the trays to air dry completely before storing them.

- **Storage**

Store the clean, dry trays in the provided retainer case to protect them from dust and damage. Keep the case in a cool, dry place.

TROUBLESHOOTING

- **Trays do not fit properly after molding:** If the trays do not provide a snug fit, re-dip them in hot water for 3-4 seconds and repeat the molding process. Ensure you press firmly against your teeth to capture a clear impression. You have an extra set of trays if needed.
- **Excessive gel contact with gums or sensitivity:** This product is for the trays only. If you experience gum irritation or tooth sensitivity, reduce the amount of whitening gel used or consult the instructions provided with your specific whitening gel product. If symptoms persist, discontinue use and consult a dental professional.
- **Trays feel bulky:** The Auraglow trays feature a 20% slimmer design. If they still feel bulky, ensure you have trimmed any excess material as described in Step 5 of the molding instructions.

SPECIFICATIONS

Model Number	AGTWK1
Product Dimensions	5.9 x 4.3 x 0.7 inches
Item Weight	3.2 ounces
Manufacturer	Auraglow
Material	EVA (BPA Free, Latex Free, Silicone Free)
Age Range Description	Adult

WARRANTY AND SUPPORT

Specific warranty information for the Auraglow Moldable Teeth Whitening Trays is not provided in the product details. For any product-related inquiries, support, or warranty claims, please refer to the manufacturer's official website or contact their customer service directly.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question