

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [ELENKER](#) /

› [ELENKER Knee Walker with 10-inch Front Wheels and Dual Braking System \(Model B06XB7GHFH\) Instruction Manual](#)

ELENKER B06XB7GHFH

ELENKER Knee Walker Instruction Manual

Model: B06XB7GHFH

INTRODUCTION

This manual provides essential instructions for the safe and effective use of your ELENKER Knee Walker. Please read this manual thoroughly before assembly and operation to ensure proper function and to prevent injury. The ELENKER Knee Walker is designed to provide mobility support for individuals recovering from foot, ankle, or leg injuries, or those undergoing foot or ankle surgery, serving as an alternative to crutches.

SAFETY INFORMATION

Your safety is paramount. Adhere to the following guidelines:

- Always ensure all locking mechanisms are securely engaged before use.
- Do not exceed the maximum weight capacity of 300 lbs.
- Maintain a clear path and be aware of your surroundings to avoid obstacles.
- Exercise caution when navigating uneven surfaces, inclines, or transitions between different flooring types.
- Regularly inspect the knee walker for any signs of wear, damage, or loose components.
- Consult a healthcare professional for guidance on proper usage and any specific limitations.

PRODUCT OVERVIEW

The ELENKER Knee Walker is a steerable medical scooter featuring 10-inch PVC front wheels and 8-inch PVC rear wheels, equipped with a dual braking system for enhanced control. It includes fully adjustable handlebars and a comfortable knee pad.



Figure 1: ELENKER Knee Walker, showcasing its overall design and components.

SETUP AND ASSEMBLY

Follow these steps for initial setup and assembly:

1. Unpack all components and verify against the packing list.
2. Attach the handlebar assembly to the main frame, ensuring it is securely fastened.
3. Install the knee pad onto its support post.
4. Adjust the handlebar height and knee pad height to a comfortable and safe position (refer to the 'Adjustments' section).
5. Ensure the dual braking system is correctly connected and functional.
6. Attach the detachable cloth bag to the front of the handlebar post, if desired.

Save Space

Compact design for easy transport and storage



Figure 2: Illustration of the folding mechanism for compact transport and storage.

OPERATING INSTRUCTIONS

To operate your ELENKER Knee Walker:

- Place the injured leg with the knee resting comfortably on the knee pad.
- Use your uninjured leg to propel the scooter forward.
- Steer using the handlebars, similar to a bicycle.
- To slow down or stop, gently squeeze the handbrakes. For parking, engage the locking handbrakes by pushing them down.
- The 10-inch front wheels and 8-inch rear wheels are designed for both indoor and outdoor use, providing stability on various terrains.

Suitable for Various Terrains

Suitable for living room, office, outdoors and various roads.



Living room



Office



Outdoors



Figure 3: The knee walker in use on an outdoor surface, highlighting its versatility.



Dual Braking System

Supports immediate stop and slowdown

Figure 4: Detail of the dual braking system, illustrating its function for immediate stopping and slowing.

ADJUSTMENTS

The ELENKER Knee Walker offers several adjustable features for personalized comfort and fit:

- **Handlebar Height:** Adjust the handlebars to a height between 28.15 inches and 37.2 inches. Loosen the adjustment knob, slide the handlebars to the desired height, and securely tighten the knob.
- **Knee Pad Height:** The knee pad can be adjusted from 16.93 inches to 20.86 inches. Loosen the adjustment knob beneath the knee pad, set to the appropriate height, and re-tighten.
- **Knee Pad Position:** The knee pad can also be adjusted forward or backward to optimize comfort and balance.

Adjustable Design

Adjustable knee pad and handlebars to fit you well



Figure 5: Adjustable design with key dimensions for handlebars and knee pad.

Knee Pad

Comfortable and soft

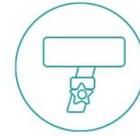


Figure 6: The comfortable knee pad, designed for extended use.

MAINTENANCE

To ensure the longevity and safe operation of your knee walker:

- Keep the frame clean with a damp cloth. Avoid abrasive cleaners.
- Periodically check all nuts, bolts, and screws for tightness. Tighten as necessary.
- Inspect the wheels for wear and tear. Ensure they rotate freely.
- Test the brakes regularly to confirm they are responsive and lock securely.
- Store the knee walker in a dry place when not in use.

TROUBLESHOOTING

If you encounter issues with your ELENKER Knee Walker, consider the following:

- **Difficulty Steering:** Check for debris in the wheels or steering mechanism. Ensure the front axle is clear.
- **Brakes Not Engaging:** Inspect brake cables for damage or looseness. Adjust cable tension if needed. Ensure the locking mechanism is fully depressed.

- **Unstable Feel:** Verify that all frame connections are tight and the knee pad and handlebar heights are correctly adjusted for your body.
- **Squeaking Noises:** Lubricate moving parts with a silicone-based lubricant.



Textured 10" Front Wheels

Featuring textured wheels, it is more non-slip, durable and provides greater stability for those with poor balance

Figure 7: Detail of the textured 10-inch front wheels for improved grip and stability.

SPECIFICATIONS

Feature	Specification
Brand	ELENKER
Model Number	B06XB7GHFH
Color	Green
Material	Polyvinyl Chloride
Front Wheel Size	10 inches
Rear Wheel Size	8 inches

Handlebar Height Range	28.15" - 37.2"
Knee Pad Height Range	16.93" - 20.86"
Weight Capacity	300 lbs
Product Weight	22 lbs
Age Range	Adult

WARRANTY AND SUPPORT

For warranty information, customer support, or to purchase replacement parts, please refer to the official ELENKER website or contact their customer service directly. Keep your purchase receipt for warranty claims.