

Geemarc Amplidect 295 SOS Pro

Geemarc Amplidect 295 SOS Pro User Manual

Model: Amplidect 295 SOS Pro (DECT295SOSPRO_WH_VDE)

Brand: Geemarc

1. INTRODUCTION

The Geemarc Amplidect 295 SOS Pro is a cordless telephone system specifically designed to meet the needs of seniors and individuals with mild to moderate hearing loss. This system includes an amplified cordless handset, a base unit with an integrated answering machine, and a portable SOS emergency call medallion. Its features, such as large buttons, a backlit display, and powerful amplification, ensure ease of use and enhanced safety for daily communication.

2. SAFETY INFORMATION

Please read and understand all safety instructions before using this product. Failure to follow these instructions may result in injury or damage to the device.

- **Hearing Protection:** Due to the amplified volume capabilities, exercise caution when adjusting the reception volume. Prolonged exposure to high volumes may cause hearing damage.
- **Water and Moisture:** Do not expose the phone or medallion to water or moisture. Do not use near bathtubs, sinks, or other wet areas.
- **Heat Sources:** Keep the device away from heat sources such as radiators, heat registers, stoves, or other appliances that produce heat.
- **Batteries:** Use only the specified batteries for the handset and medallion. Dispose of batteries properly according to local regulations. Do not attempt to open or modify batteries.
- **Cleaning:** Unplug the product from the wall outlet before cleaning. Use a soft, damp cloth for cleaning. Do not use liquid cleaners or aerosol cleaners.
- **Servicing:** Do not attempt to service this product yourself. Refer all servicing to qualified service personnel.

3. PACKAGE CONTENTS

Verify that all items are present in your package:

- 1 x Cordless Handset
- 1 x Base Unit with Answering Machine
- 1 x SOS Emergency Medallion
- 1 x Power Adapter for Base Unit
- 1 x Telephone Line Cord
- Rechargeable Batteries for Handset and Medallion
- 1 x User Manual (this document)

4. PRODUCT OVERVIEW

Familiarize yourself with the main components of your Geemarc Amplidect 295 SOS Pro system.



Figure 4.1: The complete Geemarc Amplidect 295 SOS Pro system, showing the cordless handset, base unit, and SOS emergency medallion.



Figure 4.2: Front view of the cordless handset and base unit, highlighting the large buttons and display.



Figure 4.3: Rear view of the base unit, showing connection ports and battery compartment

access.



Figure 4.4: Close-up of the SOS emergency medallion, featuring the prominent red SOS button.



Figure 4.5: Rear view of the SOS medallion, showing the volume controls and M1/M2 direct dial buttons.

4.1. Cordless Handset

The handset features a large, backlit display and oversized buttons for easy dialing and navigation. It includes a dedicated SOS button on the side for emergency calls and a 'boost' button for temporary volume amplification during calls. The handset is hearing aid compatible (T-position).

4.2. Base Unit

The base unit serves as the charging station for the handset and houses the digital answering machine. It provides direct access to answering machine functions such as message playback and deletion. The base also features a powerful ringer (up to 80dB).

4.3. SOS Emergency Medallion

This portable medallion can be worn around the neck and allows for emergency calls and answering incoming calls remotely. It has a built-in microphone and speaker, an SOS button, and two direct memory keys (M1, M2) for pre-programmed contacts. The medallion activates the SOS function, which automatically dials up to four pre-recorded numbers sequentially until a connection is made.

5. SETUP

Follow these steps to set up your Geemarc Amplidect 295 SOS Pro system.

1. Connect the Base Unit:

- Connect the power adapter to the base unit and plug it into a wall outlet.
- Connect the telephone line cord to the base unit and plug the other end into your telephone wall jack.

2. Install Handset Batteries:

- Open the battery compartment on the back of the handset.
- Insert the supplied rechargeable batteries, observing the correct polarity.
- Place the handset on the base unit to charge for at least 16 hours before first use.

3. Charge SOS Medallion:

- Place the SOS medallion in its charging cradle. Ensure it is properly seated and the charging indicator lights up.

4. Initial Language Setting:

- Upon first power-on, the system may default to German. Navigate through the menu to find the language settings (look for "Sprache" if in German) and select English. Refer to the detailed manual for specific menu navigation.

5. Set Date and Time:

- Access the handset menu to set the current date and time. This is important for accurate call logging and answering machine timestamps.

6. Program SOS Numbers:

- The SOS function can dial up to four pre-programmed emergency numbers. Access the SOS settings in the handset menu to enter these numbers. Ensure they are reliable contacts who can respond in an emergency.

7. Personalize Answering Machine Message:

- You can record your own outgoing message for the answering machine. Follow the instructions in the handset menu or base unit controls to record your personalized greeting.

6. OPERATING INSTRUCTIONS

6.1. Making and Receiving Calls

- **Making a Call:** Dial the number using the large keypad and press the green call button.
- **Receiving a Call:** When the phone rings, press the green call button on the handset or the answer button on the SOS medallion.
- **Adjusting Volume:** Use the volume up/down buttons on the side of the handset during a call. Press the 'boost' button for an additional amplification of up to +30dB.
- **Hearing Aid Compatibility:** If you use a hearing aid, switch it to the 'T' position to utilize the inductive coupling feature for clearer sound.

6.2. Answering Machine Functions

The answering machine can record messages up to 15 minutes in total. Messages can be accessed from the base unit or the handset.

- **Playback Messages:** Press the 'Play/Stop' button on the base unit or navigate to the answering machine menu on the handset.

- **Delete Messages:** Use the 'Delete' button on the base unit or the corresponding option in the handset menu.
- **Record Outgoing Message (OGM):** Follow the menu prompts to record or re-record your personalized outgoing message.

6.3. SOS Emergency Function

The SOS function can be activated from both the handset and the SOS medallion.

- **Activating SOS:** Press and hold the red SOS button on the side of the handset or the central red SOS button on the medallion.
- **SOS Call Sequence:** The system will automatically dial the pre-programmed emergency numbers one by one until a contact answers and acknowledges the call. A pre-recorded message will play to the recipient.
- **Communicating via Medallion:** When an SOS call is connected via the medallion, you can speak directly through its built-in microphone and listen through its speaker.

6.4. Medallion Direct Dial (M1/M2)

The SOS medallion features two additional memory buttons (M1 and M2) on its back, allowing for quick dialing of two specific contacts.

- **Programming M1/M2:** Program these numbers via the handset menu.
- **Making Calls:** Press and hold the M1 or M2 button on the medallion to dial the corresponding pre-programmed number. Communication will be hands-free through the medallion.

7. MAINTENANCE

Proper care ensures the longevity of your device.

- **Cleaning:** Wipe the handset, base, and medallion with a soft, slightly damp cloth. Do not use harsh chemicals or abrasive cleaners.
- **Battery Care:** For optimal battery life, allow the handset and medallion to fully discharge occasionally before recharging. Avoid leaving them on the charger indefinitely if not in use for extended periods.
- **Storage:** If storing the device for a long time, remove the batteries from the handset and medallion. Store in a cool, dry place.

8. TROUBLESHOOTING

If you encounter issues with your Geemarc Amplidect 295 SOS Pro, try these solutions:

Problem	Possible Solution
No dial tone / Cannot make calls	Check if the telephone line cord is securely connected to both the base unit and the wall jack. Ensure the power adapter is plugged in and the base unit is powered on.
Handset not charging	Ensure the handset is correctly placed on the base unit. Check that the charging contacts on both the handset and base are clean. Verify the power adapter is connected.

Problem	Possible Solution
SOS function not working	Verify that emergency numbers are correctly programmed in the handset menu. Ensure the SOS medallion is charged and within range of the base unit.
Poor sound quality / Cannot hear clearly	Adjust the reception volume using the handset's volume buttons. Press the 'boost' button for maximum amplification. If using a hearing aid, ensure it is in 'T' position. Check for interference from other electronic devices.
Answering machine not recording messages	Check if the answering machine is turned ON. Ensure there is sufficient recording time available.
Medallion not responding / out of range	Ensure the medallion is fully charged. Move closer to the base unit. Obstacles like thick walls can reduce range.
Display shows incorrect language (e.g., German)	Access the handset menu, navigate to 'Settings' or 'System', and then 'Language' (or 'Sprache' if in German) to select English.

9. SPECIFICATIONS

- **Brand:** Geemarc
- **Model Name:** Amplidect 295 SOS Pro
- **Model Number:** DECT295SOSPRO_WH_VDE
- **Color:** White
- **Screen Size:** 2.4 Inches
- **Operating System:** Geemarc (Proprietary)
- **Connectivity Technology:** DECT (Cordless), Bluetooth (for specific features if applicable, though primarily DECT)
- **Special Features:** Amplified Audio, SOS Emergency Call Function, Large Buttons, Backlit Display, Answering Machine, Hearing Aid Compatible (T-position)
- **Battery Type:** Lithium-ion (Rechargeable)
- **Digital Absorption Rate (DAS):** 0 W/kg (French regulation limit: 2 W/kg)


10. WARRANTY AND SUPPORT


Warranty: This product comes with a manufacturer's warranty of 2 years from the date of purchase. Please retain your proof of purchase for warranty claims.

Support: For technical assistance, troubleshooting beyond this manual, or warranty inquiries, please contact Geemarc customer support or your retailer. Refer to the contact information provided in the product packaging or on the official Geemarc website.

For more information, visit the official Geemarc website:www.geemarc.com

Related Documents - Amplidect 295 SOS Pro

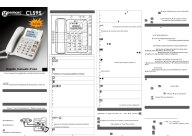	Geemarc AMPLIDECT595 U.L.E Quick User Guide - Setup and Features Comprehensive quick user guide for the Geemarc AMPLIDECT595 U.L.E cordless phone. Learn about installation, setup, features like answering machine, SOS calls, phonebook management, and technical specifications.
	Geemarc CL9000 4G Phone with SOS Emergency Call User Manual Explore the Geemarc CL9000, a 4G amplified phone with SOS emergency features and hearing aid compatibility. This user manual guides you through setup, operation, and safety for clear communication and peace of mind.
	Geemarc CL595 Amplified Corded Phone User Guide User guide for the Geemarc CL595 amplified corded phone, detailing setup, features, and troubleshooting. Learn how to use its amplified sound, large buttons, and various functions for enhanced communication.
	Geemarc CL595 Amplified Corded Phone Quick Start Guide A quick start guide for the Geemarc CL595 amplified corded telephone, detailing its features, setup, operation, SOS function, answering machine, and technical specifications. This hearing aid compatible phone is designed for ease of use with large buttons and clear display.
	Geemarc CL8360 Amplified Mobile Phone User Manual Comprehensive user guide for the Geemarc CL8360 amplified mobile phone. Features include dual SIM, camera, SOS button, large display, and amplified sound for seniors and those with hearing impairments. Covers setup, operation, SMS, phonebook, and troubleshooting.
	Geemarc Liberty 200 TRIO Amplified Digital Cordless Phone User Guide Discover the features of the Geemarc Liberty 200 TRIO amplified digital cordless phone, including call management, call blocking, multi-handset capabilities, and technical specifications. Learn about its hearing aid compatibility and user-friendly design.

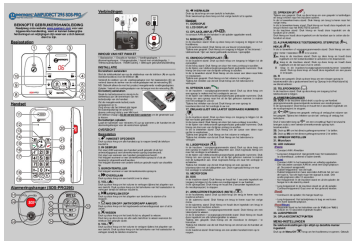

[\[pdf\]](#) Guide

AMPLIDECT295 SOS PRO quick guide Geemarc qug ad295 sos pro du a3 v1 5 hearplanet nl uc f7f5103e2010383ed080095f707020590392a57e91d8200

BEKNOPTTE GEBRUIKERSHANDLEIDING Raadpleeg onze website

www.geemarc.com voor een bijgewerkte handleiding, want er kunnen belangrijke herzieningen en wijzigingen zijn waarvan u zich bewust dient te zijn. Basisstation Handset Alarmeringshanger SOS-PRO295 Verbindingen INHOUD VAN HET PAKKET 1 Basisstat...

lang:nl **score:36** filesize: 1019.85 K page_count: 2 document date: 2020-12-11



[\[pdf\]](#) User Manual Installation Guide Instructions

Sébastien Evrard Geemarc Amplified Telephones 1 IMPORTANT SAFETY INSTRUCTIONS When using your telephone equipment basic safety precautions should always be followed to reduce the risk of fire electric shock and injury including following amplidect 295 SOS pro user manual amplifiedtelephones co uk products large |||

English Geemarc User s manual TABLE OF CONTENTS 1. IMPORTANT SAFETY INSTRUCTIONS 1 2. SET UP THE TELEPHONE..... 2 2.1 Package Contents.....2 2.2 Locat...

lang:en **score:34** filesize: 1.72 M page_count: 56 document date: 2018-09-04

Model	dB Level
AMPLIPOWER 50	60 dB
CL595	50 dB
AMPLIDECT COMBI 595	50 dB
AMPLIDECT 595 U.L.E	50 dB
AMPLIDECT 595 PHOTO	50 dB
AMPLIPOWER 40	40 dB
AMPLIDECT 295	30 dB
AMPLIDECT 295 PHOTO	30 dB
AMPLIDECT 295 SOS-PRO	30 dB
AMPLIDECT COMBI 295	30 dB
AMPLIDECT COMBI-295	30 dB
PHOTOPHONE 100	30 dB
SERENITIES	30 dB
CL100	30 dB
CL1100	30 dB
DALLAS 10	15 dB

[\[pdf\]](#)

Gamme téléphones Geemarc Classement FRV3 téléphones par Volume de Reception FR V1 geemarc fr sites 3 2013 02 |||

CLASSEMENT TLPHONES GEEMARC PAR VOLUME DE RCEPTION

AMPLIPOWER 50 CL595 AMPLIDECT COMBI 595 AMPLIDECT 595 U.L.E

AMPLIDECT 595 PHOTO AMPLIPOWER 40 AMPLIDECT 295 AMPLIDECT 295

PHOTO AMPLIDECT 295 SOS-PRO AMPLIDECT COMBI 295 AMPLIDECT

COMBIPHOTO 295 PHOTOPHONE 100 SERENITIES CL100 CL1100 DALLAS 10

60 ...

lang:en **score:30** filesize: 987.1 K page_count: 1 document date: 2021-07-22



[\[pdf\]](#) User Manual Guide Accessories Catalog

GEEMARC easier TELECOM About Geemarc was founded in 1974 In 1984 we were one of the first to obtain a licence import cordless telephones for British market Today Telecom specialises as Catalogue 2021 21x297 uk interactive revision 1 2 SDH geemarc sites 08 Manual and automatic alarm trigger Up 200m range outdoors Waterproof bracelet IP67 Trigger confirmed by vibration LED Direct vibe handset Placing on your temple will allow sound vibrations travel from pulsator speech center |||

GEEMARC makes your daily life easier R E V I S I O N N 1 2 0 2 1 P R O D U C T C A T A L O G U E GEEMARC TELECOM About Geemarc was founded in 1974. In 1984, we were one of the first to obtain a licence to import cordless telephones for the British market. Today, Geemarc Telecom specialises, as ...

lang:en **score:30** filesize: 7.1 M page_count: 15 document date: 2021-08-02



[\[pdf\]](#) Guide Catalog

Seniors rester connecté vous réussit ! Guide des offres mat C3 A9riels 2019 LB WEB min bazile fr upload

Catalogue matériels interactif Erreur 404 Bazile |||

Seniors, rester connecté vous réussit Roger, 72 ans, fan d arts martiaux adore donner des coups... ... lons compatibles sur la mme base. Les Fixes avec mdailon
CARACTRISTIQUES TECHNIQUES : Geemarc Amplidect 295 SOS Pro et Geemarc Serenities 295 SOS PRO Le Geemarc Amplidect 295 SOS Pro est un tlp...

lang:fr score:29 filesize: 2.68 M page_count: 60 document date: 2020-07-08

geemarc™ Hear the Difference			KLASSIFIKATION DER GEEMARC-TELEFONE NACH EMPFANGSVOLUMEN		
AMPLIPOWER 50		60 dB			
CL595		50 dB			
AMPLIDECT COMBI 595		50 dB			
AMPLIDECT 595 U.L.E		50 dB			
AMPLIDECT 595 PHOTO		50 dB			
AMPLIPOWER 40		40 dB			
AMPLIDECT 295		30 dB			
AMPLIDECT 295 PHOTO		30 dB			
AMPLIDECT 295 SOS-PRO		30 dB			
AMPLIDECT COMBI 295		30 dB			
AMPLIDECT COMBI- PHOTO 295		30 dB			
PHOTOPHONE 100		30 dB			
SERENITIES		30 dB			
CL100		30 dB			
CL1100		30 dB			
DALLAS 10		15 dB			

[\[pdf\]](#)

Gamme téléphones Geemarc Classement GE AMPLIDECT COMBI PHOTO 295 60 dB 50 klassifikation der geemarc telefone nach empfangsvolumen amplidect combi photo 295ampli power db60 db cl595 595 u l e Klassifikation Telefone Empfangsvolumen V1 de sites 5 2018 06 |||

AMPLIPOWER 50 CL595 AMPLIDECT COMBI 595 AMPLIDECT 595 U.L.E
AMPLIDECT 595 PHOTO AMPLIPOWER 40 AMPLIDECT 295 AMPLIDECT 295
PHOTO AMPLIDECT 295 SOS-PRO AMPLIDECT COMBI 295 AMPLIDECT
COMBIPHOTO 295 PHOTOPHONE 100 SERENITIES CL100 CL1100 DALLAS 10
KLASSIFIKATION DER GEEMARC-TELEFONE NACH EMPFANGSVOLUME...

lang:fr score:29 filesize: 987.52 K page_count: 1 document date: 2021-12-06

geemarc™ Hear the Difference			KLASSIFIKATION DER GEEMARC-TELEFONE NACH EMPFANGSVOLUMEN		
AMPLIPOWER 50		60 dB			
CL595		50 dB			
AMPLIDECT COMBI 595		50 dB			
AMPLIDECT 595 U.L.E		50 dB			
AMPLIDECT 595 PHOTO		50 dB			
NEW CL610		40 dB			
AMPLIPOWER 40		40 dB			
NEW PHOTOPHONE 110		30 dB			
NEW DALLAS 20		30 dB			
AMPLIDECT 295		30 dB			
AMPLIDECT 295 PHOTO		30 dB			
AMPLIDECT 295 SOS-PRO		30 dB			
AMPLIDECT COMBI 295		30 dB			
AMPLIDECT COMBI- PHOTO 295		30 dB			
PHOTOPHONE 100		30 dB			
SERENITIES		30 dB			
CL100		30 dB			
CL1100		30 dB			
DALLAS 10		15 dB			

[\[pdf\]](#)

Gamme téléphones Geemarc Classement GE 2023 V2 Tehnični podatki 3305771 CL610 vrvični telefon za starejše foto tipke prostoročno telefoniranje združljivo slušne aparate vključno MegaShop spletna trgovina Megashop si tehnicni geemarc cl610 vrvicni starejše prostoročno združljivo slusne apa asset conrad media10 add 160267 c1 de 003305771DS01 vključno |||

AMPLIPOWER 50 CL595 AMPLIDECT COMBI 595 AMPLIDECT 595 U.L.E
AMPLIDECT 595 PHOTO NEU CL610 AMPLIPOWER 40 NEU PHOTOPHONE 110
NEU DALLAS 20 AMPLIDECT 295 AMPLIDECT 295 PHOTO AMPLIDECT 295 SOS-
PRO AMPLIDECT COMBI 295 AMPLIDECT COMBIPHOTO 295 PHOTOPHONE
100 SERENITIES CL100 CL1100 DALLAS 10 KLASSIFIKATI...

lang:fr score:28 filesize: 1.42 M page_count: 1 document date: 2023-05-16

geemarc™ Hear the Difference			KLASSIFIKATION DER GEEMARC-TELEFONE NACH EMPFANGSVOLUMEN	
AMPLIPOWER 50		60 dB		
CL595		50 dB		
AMPLIDECT COMBI 595		50 dB		
AMPLIDECT 595 U.L.E		50 dB		
AMPLIDECT 595 PHOTO		50 dB		
CL610		40 dB		
AMPLIPOWER 40		40 dB		
PHOTOPHONE 110		30 dB		
DALLAS 20		30 dB		
AMPLIDECT 295		30 dB		
AMPLIDECT 295 PHOTO		30 dB		
AMPLIDECT 295 SOS-PRO		30 dB		
AMPLIDECT COMBI 295		30 dB		
AMPLIDECT COMBI- PHOTO 295		30 dB		
PHOTOPHONE 100		30 dB		
SERENITIES		30 dB		
CL100		30 dB		
CL1100		30 dB		
DALLAS 10		15 dB		

[pdf]

Gamme téléphones Geemarc Classement GE 2023 V2 Klassifikation nach Empfangslautstärke Amplidect 595 Photo online kaufen bei auric auric24 de media documents |||

AMPLIPOWER 50 CL595 AMPLIDECT COMBI 595 AMPLIDECT 595 U.L.E
AMPLIDECT 595 PHOTO NEU CL610 AMPLIPOWER 40 NEU PHOTOPHONE 110
NEU DALLAS 20 AMPLIDECT 295 AMPLIDECT 295 PHOTO AMPLIDECT 295 SOS-
PRO AMPLIDECT COMBI 295 AMPLIDECT COMBIPHOTO 295 PHOTOPHONE
100 SERENITIES CL100 CL1100 DALLAS 10 KLASSIFIKATI...
lang:fr score:28 filesize: 1.35 M page_count: 1 document date: 2023-05-16



[pdf] User Manual

GEEMARC Katalog Geemarc 2022 GE Ausgabe 1 geemarc de sites 5 09 |||
GEEMARC, begleitet Sie sicher in Ihrem Alltag A U S G A B E Nr.1 2 0 2 2 P R O D U
K T K A T A L O G GEEMARC TELECOM ber uns 1974 Geemarc Telecom wurde
gegrndet. 1984 waren wir eine der ersten Firmen, die die Zulassung fr den Vertrieb
von schnurlosen Telefonen auf dem britischen Markt erhie...
lang:de score:28 filesize: 7.23 M page_count: 15 document date: 2022-09-12



[pdf] Catalog

Sébastien Evrard Geemarc AmpliDect 295 SOS Pro ug amplidect295 sos pro fr v1 2 officeeasy media Oe
CatalogPdfAttribute |||
Franais Geemarc Notice TABLE DES MATIRES 1. CONSIGNES DE SCURIT
IMPORTANTES 1 2. MISE EN PLACE DU TLPHONE
..... 2 2.1 Contenu de l emballage
.....2 2.2 Emplacement.....
lang:fr score:28 filesize: 1.65 M page_count: 58 document date: 2020-08-25



[\[pdf\]](#) User Manual

PRODUKTKATALOG Katalog Geemarc 2023 GE Ausgabe 1 geemarc de sites 5 05 |||

GEEMARC, begleitet Sie sicher in Ihrem Alltag A U S G A B E Nr.1 2 0 2 3 P R O D U K T K A T A L O G GEEMARC TELECOM ber uns Folgen Sie uns auf den sozialen Netzwerken #GeemarcTelecomSA Geemarc Telecom wurde 1974 gegrndet. 1984 waren wir eine der ersten Firmen, die die Zulassung fr den Vertri...

lang:de score:27 filesize: 7.02 M page_count: 15 document date: 2023-04-26



[\[pdf\]](#) Guide Accessories Catalog

easier Catalogue Geemarc 2021 21x297 uk interactive revision 1 SDH geemarc sites 2 06

GEEMARC makes your daily life easier R E V I S I O N N1 2 0 2 1 P R O D U C T C A T A L O G U E GEEMARC TELECOM About Geemarc was founded in 1974. In 1984, we were one of the first to obtain a licence to import cordless telephones for the British market. Today, Geemarc Telecom specialises, as ...

lang:en score:27 filesize: 7.12 M page_count: 15 document date: 2021-06-03



[\[pdf\]](#) Guide Accessories Catalog

Catalogue Geemarc UK 2022 21x297 uk interactive revision 1 SDH geemarc sites 2 08 |||

GEEMARC makes your daily life easier R E V I S I O N N1 2 0 2 2 P R O D U C T C A T A L O G U E GEEMARC TELECOM About Geemarc was founded in 1974. In 1984, we were one of the first to obtain a licence to import cordless telephones for the British market. Today, Geemarc Telecom specialises, as ...

lang:en score:27 filesize: 7.37 M page_count: 15 document date: 2022-08-08



[\[pdf\]](#) User Manual

PRODUKTKATALOG Katalog Geemarc 2023 GE Ausgabe 2 geemarc de sites 5 10 |||

GEEMARC, begleitet Sie sicher in Ihrem Alltag A U S G A B E Nr.2 2 0 2 3 P R O D U K T K A T A L O G GEEMARC TELECOM ber uns Folgen Sie uns auf den sozialen Netzwerken #GeemarcTelecomSA Geemarc Telecom wurde 1974 gegrndet. 1984 waren wir eine der ersten Firmen, die die Zulassung fr den Vertri...

lang:de score:26 filesize: 7.22 M page_count: 15 document date: 2023-10-18



[Geemarc Product Catalogue 2023: Amplified Telephones, Accessibility Devices, and More](#)

Explore the comprehensive Geemarc Product Catalogue for 2023, featuring a wide range of amplified telephones, accessibility solutions, alarm clocks, headsets, and mobile devices designed to make daily life easier for everyone.

lang:en score:26 filesize: 7.03 M page_count: 15 document date: 2023-04-19



[\[pdf\]](#) Guide Accessories Catalog

easier Catalogue Geemarc 2024 21x297 uk interactive revision 1 geemarc sites 2 02 |||

GEEMARC makes your daily life easier R E V I S I O N N1 2 0 2 4 P R O D U C T C A T A L O G U E GEEMARC TELECOM Join us on our social media networks #GeemarcTelecomSA About Geemarc was founded in 1974. In 1984, we were one of the first to obtain a licence to import cordless telephones for the B...

lang:en score:26 filesize: 7.44 M page_count: 15 document date: 2024-02-06



[\[pdf\]](#) User Manual

GEEMARC 10 CL1100 äCLA4 Hörerlautstärke von bis zu 30dB Klingeltonlautstärke 80dB 12

Direktwahlkasten Pause Funktion Knochenleitungshandgrerät Katalog Geemarc 2021 GE Ausgabe 1 2

geemarc de sites 5 08 |||

GEEMARC, Ihr tglicher Partner A U S G A B E Nr.1 2 0 2 1 P R O D U K T K A T A L O G GEEMARC TELECOM ber uns 1974 Geemarc Telecom wurde gegrndet. 1984 waren wir eine der ersten Firmen, die die Zulassung fr den Vertrieb von schnurlosen Telefonen auf dem britischen Markt erhielt. Heute konze...

lang:de **score:26** filesize: 7.37 M page_count: 15 document date: 2021-08-02



[\[pdf\]](#) Guide Accessories Catalog

GEEMARC easier TELECOM About Geemarc was founded in 1974 In 1984 we were one of the first to obtain a licence import cordless telephones for British market Today Telecom specialises as Catalogue 2023 21x297 uk interactive revision 2 geemarc sites 10 |||

GEEMARC makes your daily life easier R E V I S I O N N 2 0 2 3 P R O D U C T C A T A L O G U E GEEMARC TELECOM Join us on our social media networks #GeemarcTelecomSA About Geemarc was founded in 1974. In 1984, we were one of the first to obtain a licence to import cordless telephones for the B...

lang:en **score:26** filesize: 7.23 M page_count: 15 document date: 2023-10-18



[\[pdf\]](#) Catalog

CATALOGUE PRODUITS Catalogue Geemarc 2023 FR interactif edition 2 geemarc fr sites 3 10 |||

GEEMARC, pour un quotidien plus serein E D I T I O N N 2 0 2 3 C A T A L O G U E PRODUITS GEEMARC TELECOM Rejoignez-nous sur les rseaux sociaux

#GeemarcTelecomSA propos 1974 La socit Geemarc a t cre en . En 1984, nous tions l une des premires compagnies obtenir l agrment pour l import...

lang:fr **score:25** filesize: 7.2 M page_count: 15 document date: 2023-10-24



[\[pdf\]](#) User Manual Catalog

Ihr täglicher Geemarc 2021 07 01 — DOORBELL 595 U L E RINGER AMPLIDECT595 Läuft wie eine Gegensprechanlage zwischen Doorbell und Telefon Catalogue 21x297 GE interactif version 1 geemarc de sites 5 |||

GEEMARC, Ihr täglicher Partner V E R S I O N Nr.1 2 0 2 1 P R O D U K T K A T A L O G GEEMARC TELECOM ber uns 1974 Geemarc Telecom wurde gegründet. 1984 waren wir eine der ersten Firmen, die die Zulassung für den Vertrieb von schnurlosen Telefonen auf dem britischen Markt erhielt. Heute konze...

lang:de **score:25** filesize: 7.43 M page_count: 15 document date: 2021-07-02



[\[pdf\]](#) Catalog

GEEMARC La société Geemarc a été créée en 1974 En 1984 nous étions l'une des premières compagnies à obtenir agrément pour importation téléphones sans fil sur le marché britannique Catalogue 2021 FR interactif édition 2 geemarc fr sites 3 08 |||

GEEMARC, pour un quotidien plus serein E D I T I O N N 2 0 2 1 C A T A L O G U E PRODUITS GEEMARC TELECOM propos 1974 La société Geemarc a été créée en . En 1984, nous étions l'une des premières compagnies à obtenir l'agrément pour l'importation des téléphones sans fil sur le marché britannique. Geema...

lang:fr **score:25** filesize: 7.3 M page_count: 15 document date: 2021-08-02



[\[pdf\]](#) Guide Accessories Catalog

PRODUCT CATALOGUE 2 sept 2024 — hearing aid headsets clocks alarm and accessories By placing the AMPLICAL140 s microphone next to your doorbell or smoke detector you can have Catalogue Geemarc 21x297 uk interactive revision 3 geemarc sites 09 |||

We make your daily life easier R E V I S I O N N 3 2 0 2 4 P R O D U C T C A T A L O G U E GEEMARC TELECOM About Join us on our social media networks

#GeemarcTelecomSA Geemarc was founded in 1974. In 1984, we were one of the first to obtain a licence to import cordless telephones for the British ...

lang:en **score:25** filesize: 7.65 M page_count: 15 document date: 2024-08-06



[Geemarc Product Catalogue 2026: Amplified Telephones, Accessibility Devices, and Audio Solutions](#)

Discover the comprehensive Geemarc Product Catalogue 2026, featuring a wide range of amplified telephones, hearing aids, accessibility devices, alarm clocks, TV headsets, and audio solutions designed to make daily life easier for seniors, hearing-impaired, and visually-impaired individuals.

lang:en **score:24** filesize: 9.02 M page_count: 17 document date: 2025-09-01



[\[pdf\] Catalog](#)

pour un quotidien Geemarc 2021 03 31 — DOORBELL 595 U L E RINGER AMPLIDECT595 Utilisation du téléphone comme interphone alerteur sonne et flashe lorsque quelqu'un Catalogue FR interactif édition 2 geemarc uk sites |||

GEEMARC, pour un quotidien plus serein E D I T I O N N 2 2 0 2 1 C A T A L O G U E PRODUITS GEEMARC TELECOM propos 1974 La socit Geemarc a t cre en . En 1984, nous tions l'une des premières compagnies obtenir l'agrément pour l'importation des téléphones sans fil sur le marché britannique. Geema...

lang:fr **score:24** filesize: 7.39 M page_count: 15 document date: 2021-03-04




[\[pdf\] Catalog](#)

CATALOGUE PRODUITS Catalogue Geemarc 2025 FR interactif édition 1 geemarc fr sites 3 01 Catalogue 1 |||

Nous facilitons votre quotidien E D I T I O N N 1 2 0 2 5 C A T A L O G U E PRODUITS GEEMARC TELECOM Rejoignez-nous sur les réseaux sociaux #GeemarcTelecomSA propos 1974 La socit Geemarc a t cre en . En 1984, nous tions l'une des premières compagnies obtenir l'agrément pour l'importation d...

lang:fr **score:24** filesize: 7.54 M page_count: 15 document date: 2025-01-16

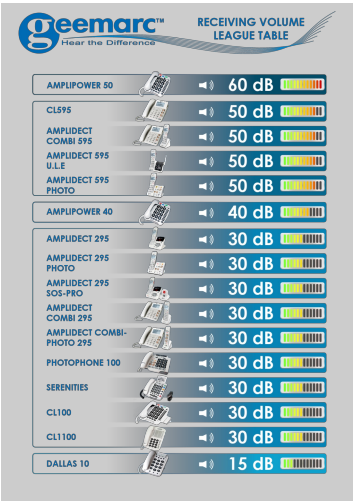

[pdf]

Gamme téléphones Geemarc Classement UK Receiving Volume League Table AMPLIDECT295 PHOTO

Telephones range V1 geemarc uk sites 2 2018 05 |||

AMPLIPOWER 50 CL595 AMPLIDECT COMBI 595 AMPLIDECT 595 U.L.E


AMPLIDECT 595 PHOTO AMPLIPOWER 40 AMPLIDECT 295 AMPLIDECT 295

PHOTO AMPLIDECT 295 SOS-PRO AMPLIDECT COMBI 295 AMPLIDECT

COMBIPHOTO 295 PHOTOPHONE 100 SERENITIES CL100 CL1100 DALLAS 10

RECEIVING VOLUME LEAGUE TABLE 60 dB 50 dB 50 dB 50 dB 5...

lang:fr score:22 filesize: 985.9 K page_count: 1 document date: 2021-12-06

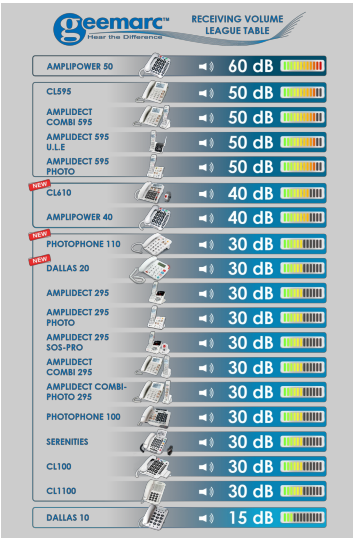

[pdf]

Gamme téléphones Geemarc Classement UK 2023 V2 Receiving Volume League Table CL610 Amplified

Telephone AmplifiedTelephones co uk rv at geemarc cl610 amplified emergency telephone with sos

remote control amplifiedtelephones user products large |||

AMPLIPOWER 50 CL595 AMPLIDECT COMBI 595 AMPLIDECT 595 U.L.E


AMPLIDECT 595 PHOTO NEW CL610 AMPLIPOWER 40 NEW PHOTOPHONE 110

NEW DALLAS 20 AMPLIDECT 295 AMPLIDECT 295 PHOTO AMPLIDECT 295

SOS-PRO AMPLIDECT COMBI 295 AMPLIDECT COMBIPHOTO 295

PHOTOPHONE 100 SERENITIES CL100 CL1100 DALLAS 10 RECEIVING VO...

lang:fr score:21 filesize: 1.35 M page_count: 1 document date: 2023-05-16

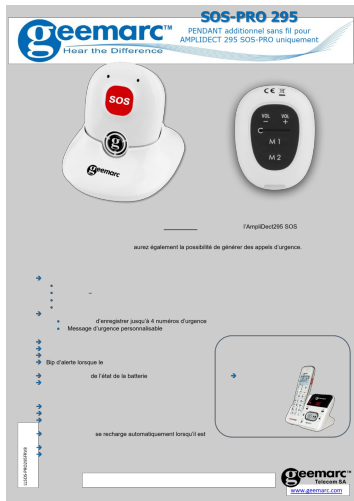

[pdf]

Wam SOS PRO 295 FRV8 geemarc fr sites 3 2018 06 |||

SOS-PRO 295 PENDANT additionnel sans fil pour AMPLIDECT 295 SOS-PRO uniquement DESCRIPTIF Le SOS-PRO 295 est un Dect additionnel uniquement compatible avec l'AmpliDect295 SOS-PRO. Il permet de recevoir vos appels entrants et de converser en mode mains libres. Deux mémoires directes M1 M2 permet...

lang:fr score:21 filesize: 733.19 K page_count: 1 document date: 2021-03-04



[pdf] User Manual Installation Guide Instructions

Sébastien Evrard Geemarc files ecommercecdns uk User's manual 2 Location For maximum coverage and reduce the interference here are some guidelines you should consider when place base unit it at a convenient high central location with no obstructions between handset in an indoor environment 349d71b3341c37b4b7e8c5b2cebb644d 212198 When using your telephone equipment basic safety precautions always be followed to risk of fire electric shock injury including Phone MasterUser environmentWhen the When theGEEMARC AMPLIDECT 295 SOS PRO MANUALfiles

349d71b3341c37b4b7e8c5b2cebb644dphone master nl Specials GEEMARC

349d71b3341c37b4b7e8c5b2cebb644dfiles |||

English Geemarc User's manual TABLE OF CONTENTS 1. IMPORTANT SAFETY INSTRUCTIONS..... 1 2. SET UP THE TELEPHONE..... 2 2.1 Package Contents

..... 2 2.2 Loca...

lang:en score:20 filesize: 1.74 M page_count: 56 document date: 2017-07-25



[pdf]

Wam Download Prospektblatt Seniorentelefon mit Notruf Mobilteil Geemarc AmpliDECT 295 SOS Pro günstig kaufen Hörwerkstatt ProspektblattRtJcx5A4Kj7B9 hoerwerkstatt.de media 6b bb 41 ||| AMPLIDECT295 SOS PRO Verstärktes schnurloses Telefon mit integriertem Anrufbeantworter und mitgeliefertem Anrufer ANRUFEN/SPRECHEN Hrgertekompatibel Einstellbare Hrerlautstärke bis zu 30dB Freisprechfunktion Lautsprecher Einstellung Telefonbuch 50 : Einträge Direktwahl Taste zum Telefonbuch R-...

lang:de score:20 filesize: 827.21 K page_count: 1 document date: 2018-08-16



[pdf] Accessories Borchure

Layout 1 Untitled Earis Wireless Listener - The premium TV listener Using dynamic sound technology to deliver the highest quality currently available in listening equipment Brilliant for people with or Sarabec A4 Brochure Spring 2020 sarabec site assets files 03 03 Wireless TV/Audio Systems 09 Loop Systems 12 Audio/Video Accessories 14 14 Personal Listeners 15 Listening Options 20 Personal Audio Accessories 23 25 Signolux Home Alerting 31 Smoke/CO Alarms 37 Telephone Ring Indicators 39 Baby Monitors 42 42 Doorbells 48 48 Clocks 55 71 55 Mo...

lang:en score:19 filesize: 14 M page_count: 80 document date: 2020-01-06

	Retail Price List	Sarabec Ltd 15 High Force Road MIDDLESBROUGH TS2 1RH e-mail sales@sarabec.co.uk		Telephone 01642 247789 Fax 01642 230827 Telephone 01642 251310 e-mail sales@sarabec.co.uk	
Effective Date: 01 November 2022					
DESCRIPTION	Order Code	Price	VAT	Total	
TV & Audio					
Digital Wireless Personal Listening Systems					
EARS Premium Digital Headset System	238-188	249.00	45.80	294.80	
EARS Premium Digital PR Headset System	238-189	249.00	52.80	301.80	
PR system with EARS Headloop included	238-182	264.00	52.80	316.80	
EARS Premium Digital PR Headset System	238-183	264.00	52.80	316.80	
PR system with Headloop included	238-184	129.00	25.80	154.80	
EARS Premium Digital PR Headset Receiver	238-185	129.00	25.80	154.80	
EARS Headloop	238-186	19.99	4.00	23.99	
Somnux 2.4 V2 Headset System	238-491	116.99	23.38	140.37	
Somnux 2.4 V2 Headset Receiver	238-492	57.00	11.40	68.40	
Somnux 2.4 V2 PR Headset System	238-440	116.99	23.38	140.37	
Somnux 2.4 V2 PR Headloop System	238-441	116.99	23.38	140.37	
Somnux 2.4 V2 PR Headloop System	238-442	123.00	24.78	147.78	
Somnux 2.4 V2 PR Headset System	238-443	123.00	24.78	147.78	
Somnux 2.4 V2 PR Headset Receiver	238-439	57.00	11.40	68.40	
Somnux 2.4 V2 PR Headset	238-440	116.99	23.38	140.37	
Geemar CL7500 OPT1 Wireless Headphones	252-228	124.99	25.00	149.99	
Geemar CL7500 OPT1 Wireless Headphones	252-223	143.66	26.33	169.99	
Bluetooth Personal Listening Systems					
Antona BT2 Bluetooth Headset	233-885	99.99	16.00	115.99	
Antona 3 Max Bluetooth Headset	233-885	99.99	20.00	119.99	
Antona HSC Wireless Microphone	233-828	79.00	15.80	94.80	
Antona TVB TV Transmitter (Optical)	233-822	49.00	9.80	58.80	
LA240 Loop Systems					
LA240 Room Loop System UK	195-258	124.99	25.00	149.99	
LA240 Loop Pad System UK	195-255	124.99	25.00	149.99	
LA240 Loop Amplifier only - UK PSU	195-178	99.99	19.99	119.98	
LA240 Loop Amplifier and 24VDC 1.5A PSU	195-179	89.99	18.00	107.99	
LA240 Loop Amplifier only - DC Load	195-173	89.99	18.00	107.99	
LA240 Loop Amplifier and 24VDC 1.5A PSU	195-199	3.50	0.70	4.20	
LA251/2320/244 Mounting brackets (2)	195-168	3.50	0.70	4.20	
LA90 Loop Systems					
LA90 Portable Loop System	197-182	179.99	35.20	215.19	
LA90 Loop Amplifier, Headset & Display Inc	197-182	179.99	35.20	215.19	

Page 1 of 1

Page 1 of 17

[pdf] Price List

RetailPriceList Report colin Sarabec Price List Retail 01 11 22 sarabec |||

a Description TV Audio Retail Price List Effective Date: 01 November 2022 Sarabec Ltd Telephone 01642 247789 15 High Force Road Fax 01642 230827

MIDDLESBROUGH Textphone 01642 251310 TS2 1RH e-mail sales@sarabec.co.uk

United Kingdom web www.sarabec.com Order Code Price VAT Total EARL...

lang:en score:17 filesize: 882.65 K page_count: 17 document date: 2022-11-03

AMPLIDECT295 SOS PRO

Verstärktes schnurloses Telefon mit integriertem Anrufbeantworter und integrierter Anhang

Hear the Difference

Mitgeliefertes Zubehör:

- 2 Direktwahlknoten
- SOS-Taste/Anrufer für eingehende Anrufe

Geemarc

www.geemarc.com

[pdf] Datasheet

Wam Geemarc AMPLIDECT 295 SOS Pro Cordless Big Button Answerphone incl emergency call

transmitter Backlit White Conrad datasheet 1677412 geemarc amplidect sos pro cordless big button

answerphone backlit white pd asset conrad media10 add 160267 c1 de 001677412DS01 |||

AMPLIDECT295 SOS PRO Verstärktes schnurloses Telefon mit integriertem

Anrufbeantworter und mitgeliefertem Anhang ANRUFEN/SPRECHEN

Hörgerterkompatibel Einstellbare Hörerlautstärke bis zu 30dB Freisprechfunktion

Lautsprecher Einstellung Telefonbuch 50 : Eingänge Direktwahl Taste zum Telefonbuch

R...

lang:de score:17 filesize: 825.88 K page_count: 1 document date: 2018-01-17

		Retail Price List		Sarabec Ltd 15 High Force Road TS2 1RH e-mail: sales@sarabec.co.uk		Telephone 01642 247789 Fax 01642 230827 MIDDLESBROUGH Telephone 01642 251310	
Effective Date: 01 July 2021							
Description	Order Code	Price	VAT	Total			
TV & Audio							
Digital Wireless Personal Listening Systems							
EARS Premium Digital PR Headset System							
EARS Premium Digital PR Headloop System	238-188	209.00	41.80	250.80			
EARS Premium Digital PR Headloop System	238-181	209.00	41.80	250.80			
PR system with EARS Headloop included	238-182	209.00	41.80	250.80			
EARS Premium Digital PR Headset System	238-183	209.00	41.80	250.80			
PR system with Headloop included	238-184	99.00	19.80	118.80			
EARS Premium Digital PR Headset Receiver	238-185	99.00	19.80	118.80			
EARS Headloop	238-186	19.99	3.20	23.19			
EARS HD Digital Headset System	238-183	199.99	38.80	238.79			
Somnux 2.4 Headset System	238-458	105.99	22.00	127.99			
Somnux 2.4 Headset Receiver	238-451	89.00	11.00	100.00			
Somnux 2.4 PR Headset System	238-460	110.99	22.20	133.19			
Somnux 2.4 PR Headset Receiver	238-461	110.99	22.20	133.19			
Somnux 2.4 PR Headloop System	238-462	116.99	23.40	140.39			
Somnux 2.4 PR Headphone System	238-463	110.99	22.20	133.19			
Somnux 2.4 PR Headset	238-458	105.99	22.00	127.99			
Geemar CL7500 OPT1 Wireless Headphones	252-228	124.99	25.00	149.99			
Bluetooth Personal Listening Systems							
Antona 3 Max Bluetooth Headset							
Antona 3 Max Bluetooth Headset	233-885	99.99	20.00	119.99			
Antona HSC Wireless Microphone	233-828	79.00	15.80	94.80			
Antona TVB TV Transmitter (Optical)	233-822	49.00	9.80	58.80			
LA240 Loop Systems							
LA240 Room Loop System UK							
LA240 Room Loop System UK	195-258	124.99	24.00	148.99			
LA240 Loop Pad System UK	195-255	119.99	24.00	143.99			
LA240 Portable Loop System UK (70m2)	195-260	394.99	79.00	473.99			
LA240 Loop Amplifier only - UK PSU	195-178	99.99	18.00	117.99			
LA240 Loop Amplifier and 24VDC 1.5A PSU	195-179	87.29	17.46	104.75			
LA240 Loop Amplifier only - DC Load	195-173	87.29	17.46	104.75			
LA240 Loop Amplifier and 24VDC 1.5A PSU	195-199	3.50	0.70	4.20			
LA90 Loop Systems							
LA90 Portable Loop System							
LA90 Loop Amplifier, Headset & Display Inc	197-182	149.99	30.00	179.99			
LA90 Portable Loop	197-180	104.99	21.00	125.99			
LA90 Portable Loop Spare Battery	197-183	29.99	6.00	35.99			
LA90 Headset Inc	197-186	31.99	6.40	38.39			
TV & Audio Accessories							

Page 1 of 17

Page 1 of 17

[pdf] Price List

RetailPriceList Report Colin Current Price List Retail Sarabec Hearing Products 01 07 21 sarabec |||

a Retail Price List Effective Date: 01 July 2021 Sarabec Ltd Telephone 01642 247789

15 High Force Road Fax 01642 230827 MIDDLESBROUGH Textphone 01642

251310 TS2 1RH e-mail sales@sarabec.co.uk United Kingdom web www.sarabec.com

Description TV Audio Order Code Price VAT Total Digital...

lang:en score:17 filesize: 505.1 K page_count: 17 document date: 2021-07-05



Retail Price List

Effective Date: 05 April 2022

Sarabec Ltd
15 High Force Road
MIDDLESBROUGH
TS2 1RH
United Kingdom
web: www.sarabec.com

Telephone: 01642 247789
Fax: 01642 230827
Textphone: 01642 251310
e-mail: sales@sarabec.co.uk

Descriptions

Order Code Price VAT Total

TV & Audio

Digital Wireless Personal Listening Systems

EARIS Premium Digital Headset System	238-188	249.00	45.80	294.80
EARIS Premium Digital RR Headloop System	238-191	264.00	52.80	316.80
RR system with EARIS Headloop included				
EARIS Premium Digital RR Headset System	238-192	264.00	52.80	316.80
RR system with Earloopcase Headset included				
EARIS Premium Digital Headset Receiver	238-184	129.00	25.80	154.80
EARIS Premium Digital RR Receiver	238-195	129.00	25.80	154.80
EARIS Headloop	238-190	19.99	4.00	23.99
Sonumaxx 2.4 Headset System	238-456	114.99	23.38	148.37
Sonumaxx 2.4 Headset Receiver	238-451	57.00	11.40	68.40
Sonumaxx 2.4 RR Headset System	238-460	116.99	23.38	148.37
Sonumaxx 2.4 RR Headloop System	238-461	116.99	23.38	148.37
Sonumaxx 2.4 RR Headset Headloop Sys	238-462	123.99	24.78	148.68
Sonumaxx 2.4 RR Headphone System	238-463	123.99	24.78	148.68
Sonumaxx 2.4 RR Receiver	238-459	59.00	11.80	70.80
Seemee CL7400 OPT Wireless Headphones	252-228	124.99	25.00	149.99
Bluetooth Personal Listening Systems				
Antena 3 Plus Bluetooth Loudset	233-665	95.99	20.00	115.99
Antena 3MC Wireless Microphone	233-626	79.00	15.80	94.80
Antena TVB TV Transmitter (Optical)	233-622	49.00	9.80	58.80
LA240 Loop Systems				
LA240 Room Loop System UK	135-258	324.99	25.00	349.99
LA240 Loop Pad System UK	135-255	124.99	25.00	149.99
LA240 Portable Loop System UK (70m2)	135-260	299.99	60.00	359.99
LA240 Loop amplifier only - UK PSU	135-178	99.99	19.99	119.98
Included mains lead and 2000C 1.5A PSU				
LA240 Loop amplifier only - DC Lead	135-173	89.99	18.00	107.99
Included Power lead, mounting brackets				
LA251/220-140 Mounting Brackets (2)	135-188	3.50	0.70	4.20

LA90 Loop Systems

LA90 Portable Loop System	137-182	179.99	35.20	215.19
LA90 Loop amplifier Headset & Headset Rev				
LA90 Portable Loop	137-180	106.99	21.40	128.39
LA90 Portable Loop Spare Battery	137-185	29.99	6.00	35.99
LA90 Headset Kit	137-186	32.49	6.50	38.99

TV & Audio Accessories

Digital Connection Converters

Page 1 of 17

Page 1 of 17

[pdf] Price List

RetailPriceList Report Colin Current Price List Retail Sarabec Hearing Products 05 04 22 sarabec ||| a Retail Price List Effective Date: 05 April 2022 Sarabec Ltd Telephone 01642 247789 15 High Force Road Fax 01642 230827 MIDDLESBROUGH Textphone 01642 251310 TS2 1RH e-mail sales sarabec.co.uk United Kingdom web www.sarabec.com Description TV Audio Order Code Price VAT Total Digita... lang:en score:17 filesize: 505.17 K page_count: 17 document date: 2022-04-05

		Retail Price List		Sarabec Ltd 15 High Force Road MIDDLESBROUGH TS2 1RH United Kingdom web: www.sarabec.com		Telephone: 01642 247789 Fax: 01642 230827 Textphone: 01642 251310 e-mail: sales@sarabec.co.uk	
Description		Order Code		Price	VAT	Total	
TV & Audio							
Digital Wireless Personal Listening Systems							
EARIS Digital Headset System		238-188	209.00	41.80	250.80		
EARIS Digital Headloop System		238-191	209.00	41.80	250.80		
RR system with EARIS Headset included							
EARIS Digital Podset Receiver System		238-192	209.00	41.80	250.80		
RR system with Earloopcase Headset included							
EARIS Digital Headset Receiver		238-184	98.00	19.60	117.60		
EARIS Digital RR Receiver		238-195	98.00	19.60	117.60		
Sonumaxx 2.4 Headset System		238-456	109.99	22.00	131.99		
Sonumaxx 2.4 Headset Receiver		238-451	55.00	11.00	66.00		
Sonumaxx 2.4 RR Headset System		238-460	109.99	22.00	131.99		
Sonumaxx 2.4 RR Headloop System		238-461	109.99	22.00	131.99		
Sonumaxx 2.4 RR Headset Headloop System		238-462	114.99	23.00	137.99		
Sonumaxx 2.4 RR Receiver		238-459	55.00	10.40	65.40		
Seemee CL7400 Wireless Headphones		252-233	116.66	23.33	139.99		
LA240 Loop Systems							
LA240 Room Loop System UK		135-258	119.99	24.00	143.99		
LA240 Loop Pad System UK		135-255	119.99	24.00	143.99		
LA240 Portable Loop System UK (70m2)		135-260	394.99	79.00	473.99		
LA240 Loop amplifier only - UK PSU		135-178	92.99	18.60	111.59		
Incl. mains lead and 2000C 1.5A PSU							
LA240 Loop amplifier only - DC Lead		135-173	87.29	17.46	104.75		
Incl. Power lead & mounting brackets							
LA251/220-140 Mounting Brackets (2)		135-188	3.50	0.70	4.20		
LA90 Loop Systems							
LA90 Portable Loop System		137-188	104.99	21.00	125.99		
LA90 Portable Loop Spare Battery		137-185	29.99	6.00	35.99		
LA90 Headset Kit		137-186	28.99	6.00	35.99		
TV & Audio Accessories							
Digital Connection Converters							
SPDF Digital Audio Converter & Leads		238-189	34.99	7.00	41.99		
Converts Coaxial (RCA) to Optical (Digi) to Stereo (RCA)							
Includes Power Supply, Coaxial and Optical leads							
Digital Audio Optical Splitter		954-125	8.99	1.80	10.79	not for sale	
Digital System Accessories							
SW151 SR Digital 200mAh Li-Pol Battery		238-146	15.99	3.40	19.39		
EARIS 400mAh Li-Pol Battery		238-187	17.99	3.60	21.59		
Sonumaxx 2.4 300mAh Li-Pol Batt		238-469	12.99	2.60	15.59		

Page 1 of 4

Page 4 of 18

[pdf] Price List

a Retail Price List 01642 247789 230827 251310 Description Order Code VAT Total TelephoneFax Textphone e mail web sales@sarabec co uk Effective Date 01 August 2019 sarabec Sarabec Ltd 15 High Force Road 08 19 site assets files ||| a Retail Price List Effective Date: 01 August 2019 Sarabec Ltd Telephone 01642 247789 15 High Force Road Fax 01642 230827 MIDDLESBROUGH Textphone 01642 251310 TS2 1RH e-mail sales sarabec.co.uk United Kingdom web www.sarabec.com Description TV Audio Order Code Price VAT Total Digit... lang:en score:17 filesize: 247.54 K page_count: 19 document date: 2019-07-31