

Banglijian Ziv-201A

Banglijian Ziv-201A Rechargeable Digital Hearing Aid User Manual

Model: Ziv-201A


1. INTRODUCTION AND OVERVIEW

The Banglijian Ziv-201A is a compact, rechargeable digital hearing aid designed to provide clear sound with noise-cancelling features. It is suitable for individuals with mild to moderate hearing loss.

Key Features:

- **Digital Noise Cancelling:** Helps reduce background noise for clearer hearing.
- **Rechargeable Battery:** Offers approximately 20-24 hours of continuous use on a full charge.
- **Adjustable Volume Control:** Five volume levels to suit different listening environments.
- **Two Program Modes:** Select between quiet and noisy environment settings.
- **Lightweight and Comfortable:** Designed for ease of use and discreet wear.

Main Functions

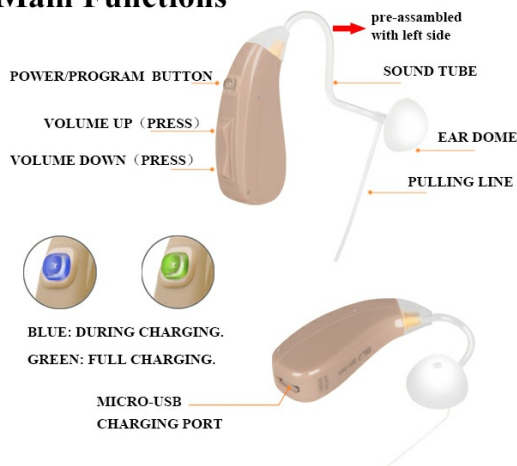

Features

Lightweight--only 3.7g with portable small shockproof case

Digital Chip--smart noise reduction

Prompt tone--clearly know what setting you have it on

Two Programs and Five Volume Levels--adjust to your listening environment

One Finger Operation--easily change the volume and program, adjust to your listening environment.

Automatic Memory--if the device is turned off, the program and volume selection always remains the same until changed by the user.

challenges, indicating that the Ziv-201A is suitable for mild and moderate hearing loss ranges.

Your browser does not support the video tag.

Video: An overview of the Banglijian Ziv-201A hearing aid, highlighting its portable design and clear sound for users experiencing hearing difficulties.

2. PRODUCT COMPONENTS

Your Banglijian Ziv-201A package includes the following items:

- 1 x Hearing aid (pre-assembled with left slim tube)
- 1 x Right slim tube with red mark
- 10 x Ear domes for slim tube (various sizes)
- Thick tube accessories (1 ear hook, 3 thick tubes, 5 ear plugs)
- 1 x Cleaning brush
- 1 x Power adapter
- 1 x Micro-USB cable
- 1 x Shockproof case
- 1 x User manual

COMPONENT



Power Adapter




USB Cable




Shockproof Case



Cleaning Brush

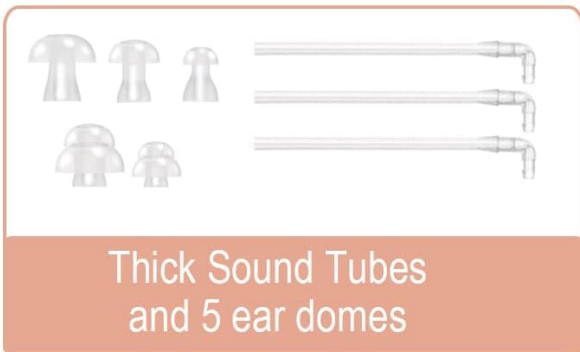

Thick Sound Tubes
and 5 ear domes

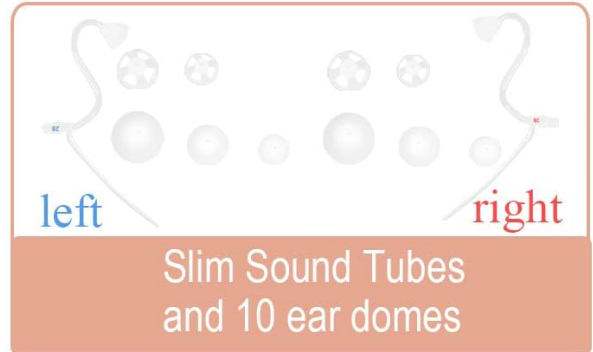

Slim Sound Tubes
and 10 ear domes

Image: A visual representation of all items included in the Banglijian Ziv-201A hearing aid package, such as the hearing aid unit, various ear domes, sound tubes, charging accessories, and cleaning tools.

3. SETUP INSTRUCTIONS

3.1 Initial Charging

Before first use, fully charge the hearing aid. Connect the micro-USB cable to the charging port on the device and the power adapter. The indicator light will be blue during charging and turn green when fully charged.



Image: Close-up of the Banglijian Ziv-201A hearing aid showing the micro-USB charging port and the indicator light, which glows blue during charging and green when fully charged.

3.2 Attaching Sound Tubes and Ear Domes

The hearing aid comes pre-assembled with a left slim tube. A right slim tube is also included. You can choose between slim and thick tubes, and various ear dome sizes, to ensure a comfortable and effective fit.

- **Selecting the Right Tube:** The slim tube is generally suitable for glasses wearers, those needing lower gain, or smaller ear domes. The thick tube is for adjusting length or requiring higher gain.
- **Selecting the Right Ear Dome:** Try different sizes of ear domes to find one that fits snugly inside your ear canal without causing discomfort or feedback.
- **Assembly:** Ensure tight connections between the sound tube, the ear dome, and the device to prevent sound leakage and whistling.

Clockwise rotation can be tightened

Counterclockwise rotation can be loosened

We provide two kinds of Sound Tube (left and right) and ten Ear Domes. By selecting the Sound Tubes and Ear Domes, it can fit to either ear.


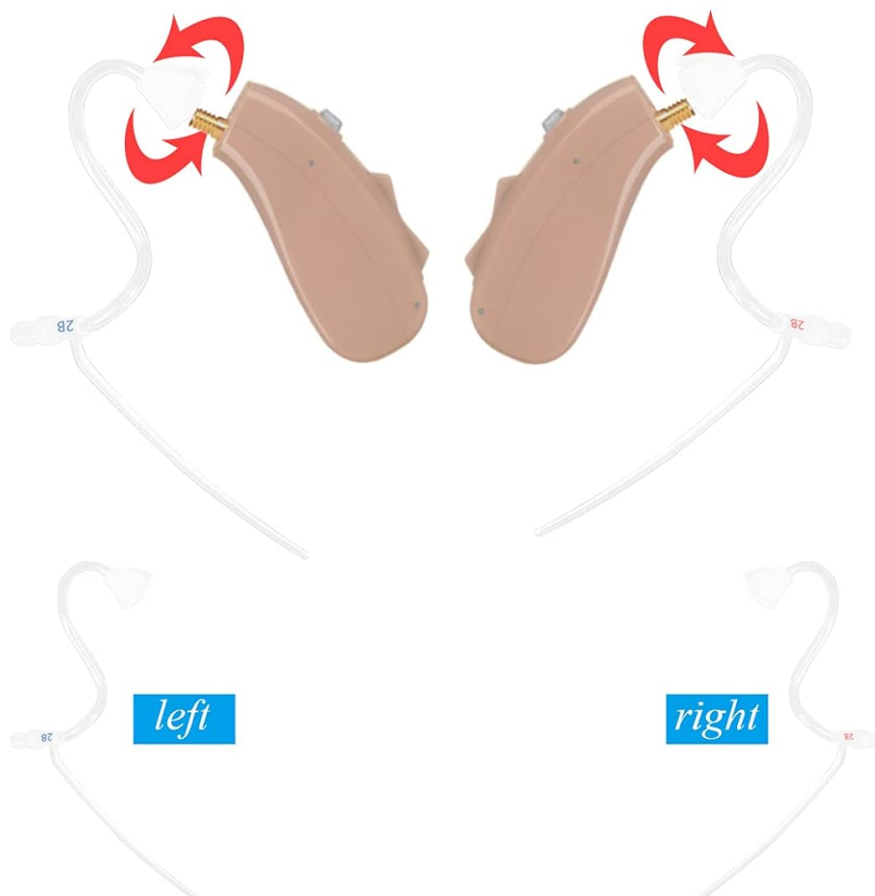


Image: Diagram illustrating how to attach and adjust the sound tubes (left and right) and ear domes to the Banglijian Ziv-201A hearing aid, emphasizing clockwise rotation to tighten and counter-clockwise to loosen.

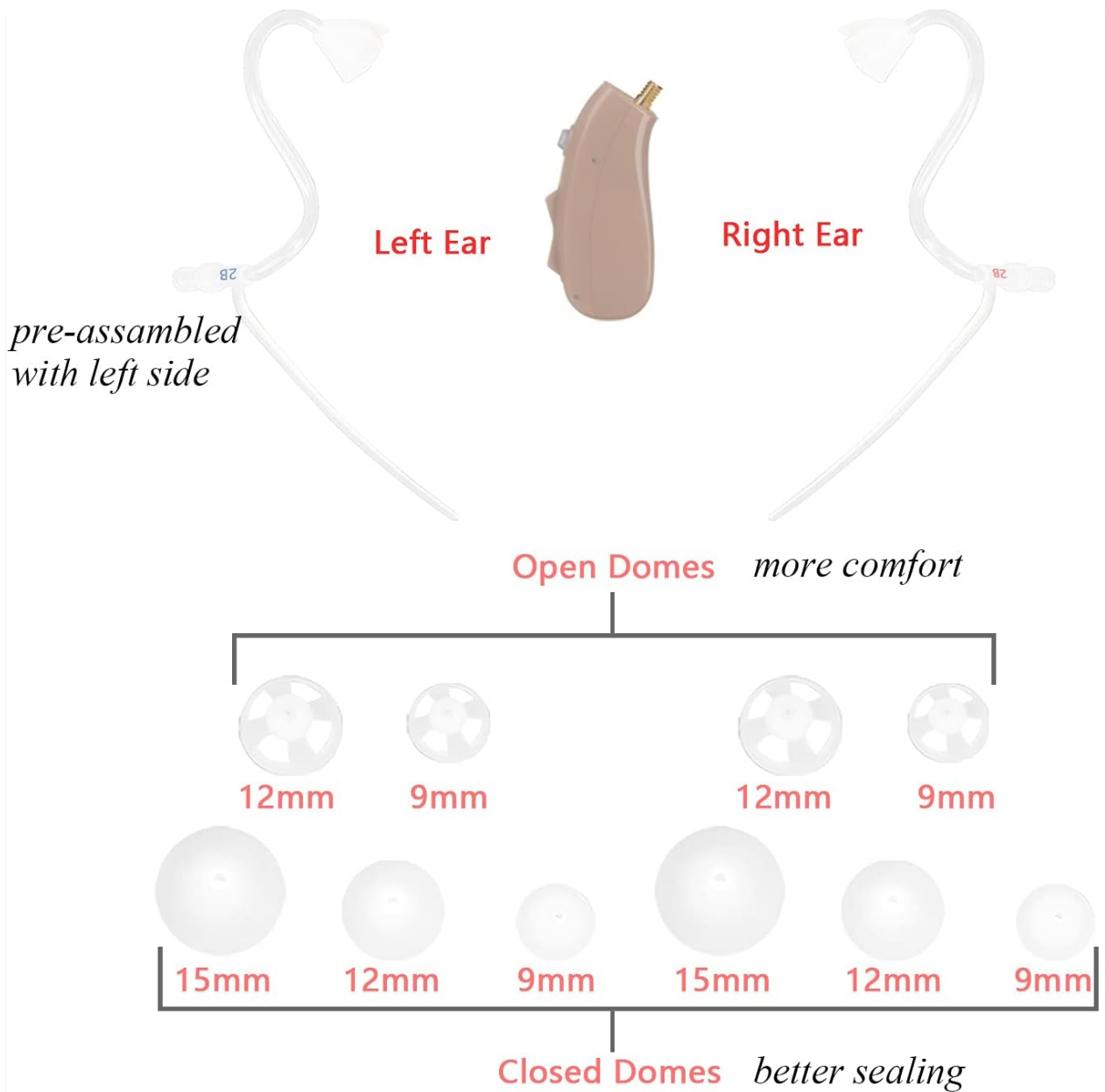

Warm Prompt: Ear domes in thick tube pack can fit slim tube.

Image: A selection of open and closed ear domes in various sizes (9mm, 12mm, 15mm) for the Banglijian Ziv-201A, demonstrating options for comfort and better sealing.

Your browser does not support the video tag.

Video: Demonstrates the pre-assembled left slim sound tube on the 201A unit and how to prepare it for use.

Your browser does not support the video tag.

Video: A comprehensive guide on how to use and wear the Banglijian 201A hearing aid, including power on/off, program switching, and proper placement.

Your browser does not support the video tag.

Video: Instructions on how to replace the slim tube with a thick tube on the 201A hearing aid, demonstrating the process of detaching and reattaching components.

3.3 Wearing the Hearing Aid

Carefully place the hearing aid behind your ear, ensuring the sound tube and ear dome are comfortably and securely positioned in your ear canal. If you wear glasses, guide the sound tube below the glasses temple, keeping it parallel and not pinched.



Image: A person demonstrating the correct way to wear the Banglijian Ziv-201A hearing aid behind the ear, with the sound tube and ear dome properly inserted.

Your browser does not support the video tag.

Video: Detailed instructions on how to properly wear the 201A hearing aid with glasses, ensuring comfort and optimal sound delivery.

4. OPERATING INSTRUCTIONS

4.1 Power On/Off

Press and hold the 'Power/Program' button for 3+ seconds until the indicator light turns green (for power on) or blue (for power off).

4.2 Volume Control

Use the 'Volume Up' and 'Volume Down' buttons to adjust the sound level. There are five volume levels. You will hear a single beep with each press. Three beeps indicate maximum or minimum volume.

4.3 Program Selection

Short press the 'Power/Program' button to switch between the two sound environment programs:

- **Quiet Mode (One Beep):** Ideal for low-frequency, relatively quiet everyday environments.
- **Noisy Mode (Two Beeps):** Suitable for very noisy environments with fewer low-frequency sounds.

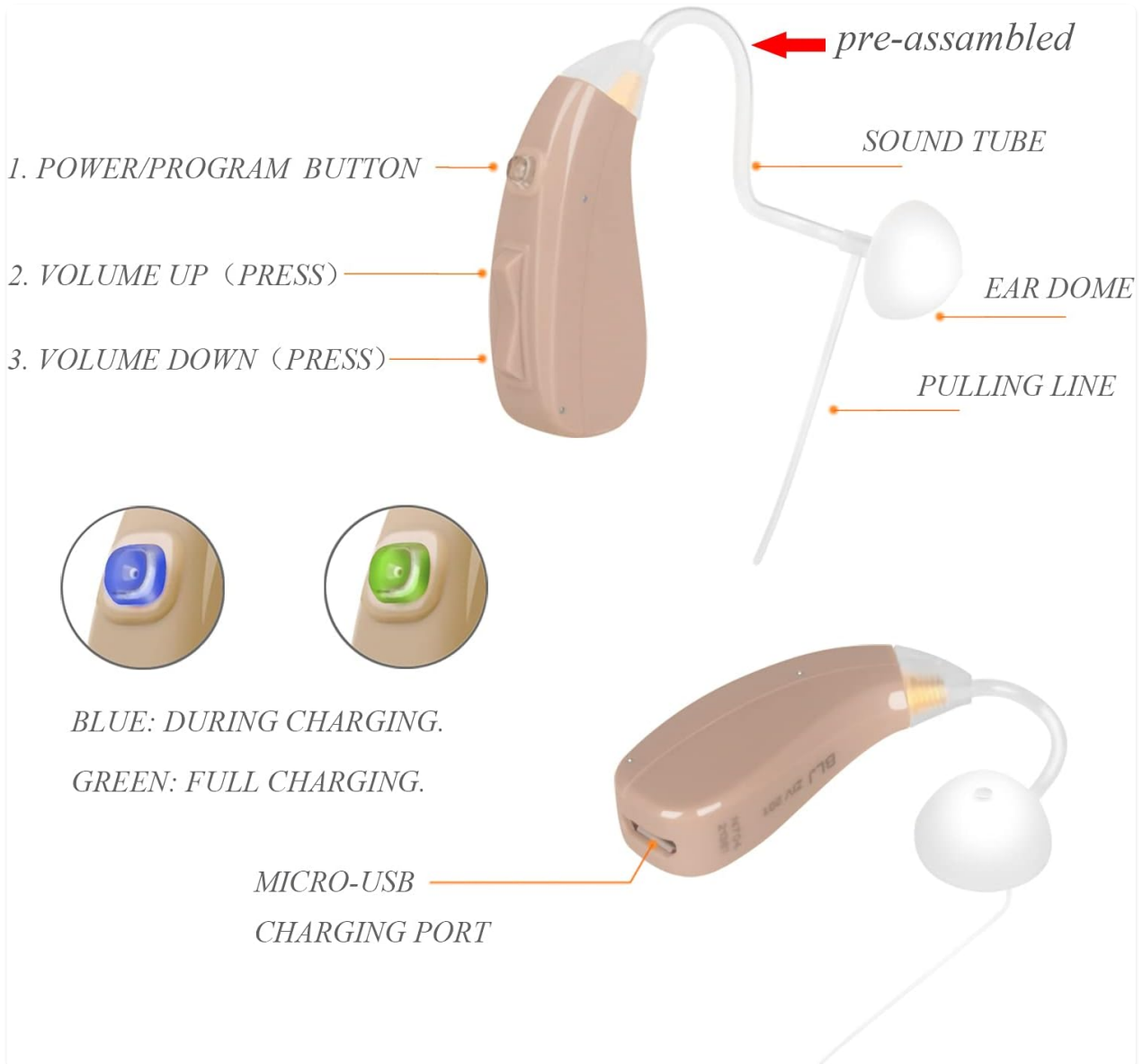


Image: Detailed view of the Banglijian Ziv-201A showing the Power/Program button, Volume Up/Down buttons, and the Micro-USB charging port. Indicator light colors for charging status are also shown.

5. MAINTENANCE

Proper maintenance ensures the longevity and optimal performance of your hearing aid.

5.1 Cleaning

- Regularly clean the surface of the device, sound tube, and ear domes using the provided cleaning brush.
- Ensure no wax or debris accumulates in the ear dome or sound tube.

5.2 Storage

- When not in use, store the device in the provided shockproof case or a hearing aid dry box to prevent

moisture damage.

- Keep the device away from high temperatures, high humidity, and direct sunlight.
- Never expose the hearing aid to fire.

5.3 Battery Care

- Use only the original provided charger and micro-USB cable.
- Charge the device in a timely manner. Prolonged low power states can damage the battery.

6. TROUBLESHOOTING

6.1 Feedback or Whistling

- Check if the ear dome fits snugly inside the ear canal. Replace it with a different size if necessary.
- Ensure tight connections among the sound tube, the ear dome, and the device.

6.2 Hearing Noise

- Allow time for adaptation. It may take 2-4 weeks to get used to the sound.
- Try switching between the program modes (Quiet/Noisy) to find the optimal setting for your environment.
- Start with the lowest possible volume and gradually increase it.

6.3 General Adaptation

It may take 20-30 days to fully adapt to using the hearing aid. Be patient and consistent with its use.

7. SPECIFICATIONS

Feature	Detail
Model Name	Ziv-201A
Power Source	Rechargeable (Lithium Polymer battery)
Battery Life	Approx. 20-24 hours continuous use
Control Type	Button Control
Noise Control	Environmental Noise Cancellation
Fit Type	Behind-the-Ear
Item Weight	3.7 Grams
Product Dimensions	1.45"L x 0.4"W
Material	Acrylonitrile Butadiene Styrene (ABS), Silicone

8. WARRANTY AND SUPPORT

8.1 Warranty Information

The Banglijian Ziv-201A comes with an 18-month warranty from the date of purchase. During this period, maintenance or replacement will be provided for all quality problems under normal conditions of use. The warranty does not cover losses or damages caused by dropping, exposure to water, or damage by animals or humans.

8.2 Customer Support

For any problems or inquiries, please contact the Banglijian support team. You can find contact information on the first page of your user manual or by following these steps:

1. Go to the product listing page on Amazon.
2. Look for the 'Banglijian' icon (usually near the 'Sold by' information) and click it.
3. On the seller's page, look for an 'Ask a question' button and click it.
4. Select 'An item for sale' and then write your message, including your 17-digit Amazon order ID and a clear description of your problem.

Your browser does not support the video tag.

Video: A short video explaining how to contact the seller for support if you have issues with your product.

© 2023 Banglijian. All rights reserved.