

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [Ansell](#) /

› [Ansell Microflex ISO 93-243 Nitrile Disposable Gloves User Manual](#)

## Ansell 93243

# Ansell Microflex ISO 93-243 Nitrile Disposable Gloves User Manual

Model: 93243

## 1. INTRODUCTION

This manual provides essential information for the proper use, care, and disposal of Ansell Microflex ISO 93-243 Nitrile Disposable Gloves. Please read these instructions carefully before use to ensure optimal performance and safety.

## 2. PRODUCT OVERVIEW

The Ansell Microflex ISO 93-243 gloves are high-quality, powder-free disposable gloves designed for various applications requiring hand protection. Made from durable nitrile material, these gloves offer a comfortable fit and reliable barrier protection.

- **Material:** Nitrile
- **Color:** Blue
- **Feature:** Powder-Free, Disposable
- **Length:** 12 inches
- **Thickness:** 4.7 mil



Image 1: A person demonstrating the donning of a long blue nitrile glove, extending protection up the arm.

### **3. SAFETY INFORMATION**

These gloves are intended for single use only. Do not reuse. While nitrile gloves offer good chemical resistance, they are not universally resistant to all chemicals. Always consult a chemical resistance guide if working with hazardous substances. Discontinue use immediately if any signs of allergic reaction or skin irritation occur. This product is accelerator-free to minimize the risk of Type IV allergic reactions.

# 加硫促進剤フリーで アレルギーに配慮



アレルギーの原因となる加硫促進剤を含まない設計で、**敏感肌の方でも安心して着用可能**です。

## 4. PREPARATION FOR USE

---

1. **Inspect Packaging:** Before opening, check the glove packaging for any signs of damage or compromise. Do not use if the packaging is damaged.
2. **Select Correct Size:** Ensure you select the appropriate glove size for a comfortable and secure fit. An ill-fitting glove can compromise dexterity and protection. The Microflex ISO 93-243 is available in various sizes, including Small.
3. **Hand Hygiene:** Wash and dry hands thoroughly before donning gloves.

## 5. USAGE INSTRUCTIONS

---

### 5.1 Donning Gloves

1. Carefully remove one glove from the dispenser box.
2. Slide your hand into the glove, ensuring fingers are fully seated.
3. Pull the cuff up your arm to ensure the long cuff provides maximum coverage, extending up to the elbow for enhanced safety.
4. Repeat for the other hand.



*point*

安全性アップ

肘まで覆うロングカフ設計で

特別に長いカフ（39.5cm）が手首から肘までを保護し、薬品や液体の侵入リスクを低減します。

Image 3: Visual representation of the long cuff design, providing extended protection up to the elbow.

1. Grasp the outside of one glove near the cuff with your other gloved hand.
2. Peel the glove off, turning it inside out, and hold it in your gloved hand.
3. Slide an ungloved finger under the cuff of the remaining glove.
4. Peel the second glove off, turning it inside out over the first glove.
5. Dispose of the gloves immediately.

### **5.3 Recommended Applications**

These gloves are versatile and suitable for a wide range of tasks, including:

- Chemical handling
- Pharmaceutical manufacturing
- Cleaning tasks
- Laboratory work and research and development
- Quality control

汎用性



多様な作業現場で使える



化学薬品取り扱い、製薬、清掃  
作業など、幅広い作業環境で活躍  
する多用途手袋です。

# 化学薬品飛沫から手を守る 高い保護性能



耐久性に優れたニトリル素材が  
化学薬品の飛沫からしっかり手を守り  
安心して作業が行えます。

**Ansell**

# 柔らかく快適な フィット感



高品質なニトリル素材が  
柔らかく手にフィットし、  
突っ張り感が少なく  
長時間作業でも快適です。

## 6. CARE AND STORAGE

Store gloves in their original packaging in a cool, dry place, away from direct sunlight, heat, and fluorescent lighting. Avoid exposure to ozone. Improper storage can degrade the material and reduce the gloves' protective properties.

## 7. DISPOSAL INFORMATION

Dispose of used gloves according to local regulations for waste disposal. If gloves have been exposed to hazardous materials, follow specific protocols for contaminated waste.

## 8. SPECIFICATIONS

Attribute	Detail
Brand	Ansell (Microflex)
Model Number	93243
Series	ISO 93-243
Material	Nitrile
Color	Blue
Size	Small
Length	12 inches
Thickness	4.7 mil
Special Feature	Powder-Free, Disposable

# こんな場面でおすすめ

*Recommendation*



医薬品  
バイオ製造  
研究開発



品質管理

展開  
サイズ

サイズ	XS	S	M	L	XL

Image 7: Recommended applications and available glove sizes (XS, S, M, L, XL).

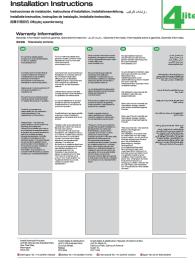
## 9. TROUBLESHOOTING

- **Glove Tears During Donning:** Ensure correct size is selected. Avoid sharp objects. Discard and use a new glove.
- **Skin Irritation:** Although accelerator-free, if irritation occurs, discontinue use and consult a medical professional.
- **Loss of Dexterity:** Verify glove size. An overly large or small glove can impair tactile sensitivity.

© 2023 Ansell. All rights reserved.

### Related Documents - 93243

	<p><a href="#"><u>Ansell Chemical Resistance Guide: Permeation &amp; Degradation Data (8th Edition)</u></a></p> <p>Comprehensive 8th edition guide from Ansell detailing permeation and degradation resistance data for chemical-resistant gloves against a wide range of hazardous chemicals, based on ASTM F 739 standards.</p>
	<p><a href="#"><u>Ansell AMOLED/LI/3M/ST LED Bulkhead Installation Instructions</u></a></p> <p>Comprehensive installation guide for the Ansell AMOLED/LI/3M/ST LED bulkhead, covering warranty details, safety precautions, and step-by-step fitting instructions.</p>
	<p><a href="#"><u>Ansell Installation Instructions and Warranty Information</u></a></p> <p>Official installation instructions and warranty details for Ansell luminaires, covering installation recommendations, warranty terms, and contact information.</p>

	<p><a href="#"><u>Ansell 4lite Outdoor Wall Lantern Installation and Warranty Guide</u></a></p> <p>Detailed installation instructions and warranty information for Ansell 4lite outdoor wall lanterns (Product Codes: 4L2/3412 &amp; 42L/3413). Covers safety, setup steps, product certifications, and contact details for Ansell.</p>
	<p><a href="#"><u>Ansell 4lite Outdoor Lantern Installation Instructions and Warranty Information</u></a></p> <p>Comprehensive installation guide and warranty details for Ansell 4lite outdoor lanterns (models 4L2/3410 &amp; 4L2/3411), including safety precautions and multi-language support.</p>
	<p><a href="#"><u>Ansell 4L2/2306 &amp; 2308 Lighting Installation and Warranty Guide</u></a></p> <p>Comprehensive guide for the installation, maintenance, and warranty information of Ansell 4L2/2306 &amp; 2308 luminaires. Includes multilingual instructions and product details.</p>