

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [ALCO](#) /

› [ALCO Heavy Duty Wheelchair Arm Pads \(Model AL-66055BK-A3\) - Instruction Manual](#)

ALCO AL-66055BK-A3

ALCO Heavy Duty Wheelchair Arm Pads (Model AL-66055BK-A3) - Instruction Manual

Model: AL-66055BK-A3 | Brand: ALCO

INTRODUCTION

This manual provides detailed instructions for the installation, proper use, and maintenance of your ALCO Heavy Duty Desk Length Wheelchair Arm Pads. These arm pads are designed for universal replacement on most manual wheelchairs, offering enhanced comfort and support. Please read these instructions carefully before installation and use.



Image 1: ALCO Heavy Duty Desk Length Wheelchair Arm Pads (Pair) with mounting hardware.

SETUP AND INSTALLATION

The ALCO Heavy Duty Wheelchair Arm Pads are designed for straightforward installation. Ensure you have the correct tools, typically a screwdriver, before beginning.

1. **Preparation:** Remove the existing arm pads from your wheelchair. Keep the old screws if they are in good condition and compatible, or use the new screws provided.
2. **Hole Alignment:** Position the new ALCO arm pad on the armrest frame. The arm pads feature two sets of pre-drilled holes to accommodate different wheelchair models. These hole spacings are approximately 5-1/2 inches (13.97 cm) or 3-5/16 inches (8.41 cm) center-to-center. Select the set of holes that aligns with your wheelchair's armrest frame.
3. **Secure Attachment:** Insert the screws through the aligned holes in the arm pad and into the armrest frame. Tighten the screws securely but do not overtighten, as this may damage the arm pad or the wheelchair frame.
4. **Repeat for Second Pad:** Follow the same steps to install the second arm pad on the opposite side of the wheelchair.
5. **Final Check:** Gently test the stability of both arm pads to ensure they are firmly attached before using the wheelchair.

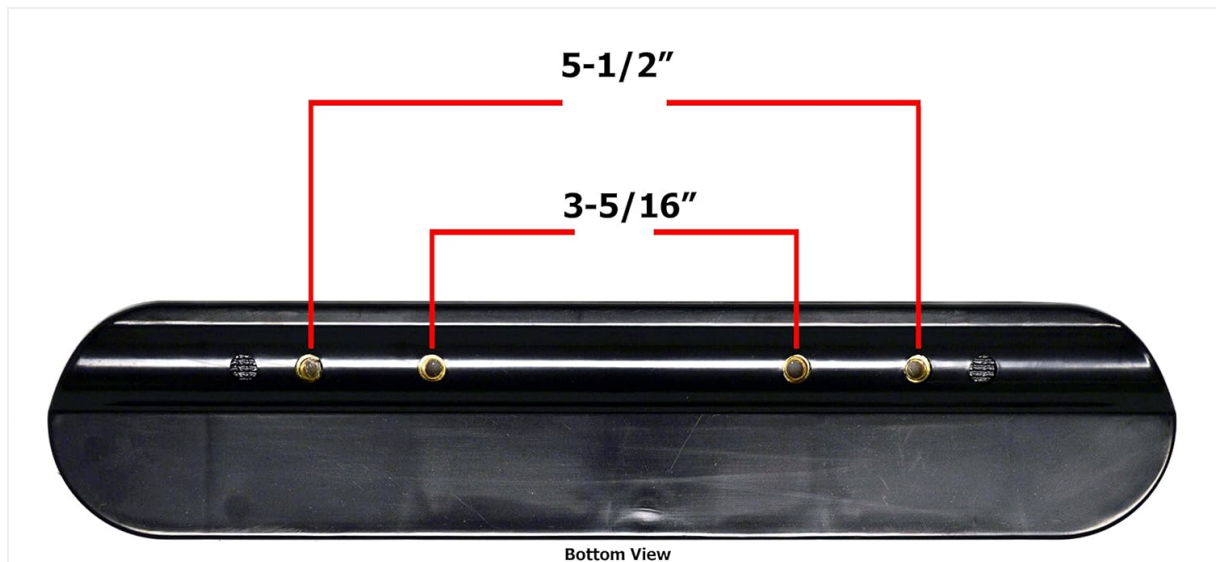


Image 2: Bottom view of arm pad illustrating dual hole spacing options for universal fit.

PROPER USE

Once installed, the ALCO arm pads provide a comfortable and supportive surface for the user's arms. Ensure that the user's arms are resting fully on the pads for optimal support and pressure distribution. These arm pads are designed for manual wheelchairs and are not suitable for power chairs.

MAINTENANCE AND CARE

- **Cleaning:** To clean the arm pads, wipe them with a damp cloth using mild soap and water. Avoid harsh chemical cleaners or abrasive materials, as these can damage the surface.
- **Inspection:** Regularly inspect the arm pads for any signs of wear, tears, or damage. Check the mounting screws periodically to ensure they remain tight.
- **Storage:** If the arm pads are removed for storage, keep them in a clean, dry place away from direct

sunlight and extreme temperatures.

TROUBLESHOOTING

Issue	Possible Cause	Solution
Arm pads do not fit or align with holes.	Incorrect hole spacing selected or incompatible wheelchair model.	Ensure you are using the correct set of holes (5-1/2" or 3-5/16"). Verify your wheelchair model is compatible with universal arm pads.
Arm pads feel loose after installation.	Screws are not tightened sufficiently.	Gently tighten all mounting screws. Do not overtighten.
Damage to arm pad surface.	Exposure to harsh chemicals or abrasive cleaning.	Clean only with mild soap and water. Replace if damage affects functionality or comfort.

SPECIFICATIONS

- **Model Number:** AL-66055BK-A3
- **Brand:** ALCO
- **Type:** Heavy Duty Desk Length Wheelchair Arm Pads
- **Color:** Black
- **Length:** 10-1/4 inches (26.04 cm)
- **Hole Spacing Options (Center-to-Center):** 5-1/2 inches (13.97 cm) or 3-5/16 inches (8.41 cm)
- **Product Dimensions:** 25.73 x 5.41 x 3.81 cm (approximately 10.13 x 2.13 x 1.5 inches)
- **Weight:** 260 g (approximately 0.57 lbs) per pair
- **Compatibility:** Universal fit for most manual wheelchairs (e.g., Medline, Drive, Invacare, E&J, Guardian, Lumex, Tuffcare). Not for power chairs.

WARRANTY INFORMATION

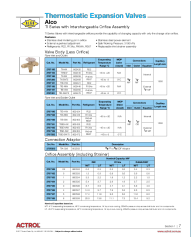
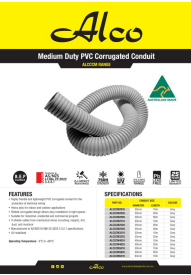

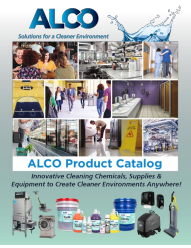


ALCO products are manufactured to high-quality standards. For specific warranty details regarding your ALCO Heavy Duty Wheelchair Arm Pads, please refer to the warranty documentation included with your purchase or contact the retailer or manufacturer directly. Keep your proof of purchase for any warranty claims.

CUSTOMER SUPPORT

If you have any questions about the installation, use, or maintenance of your ALCO Heavy Duty Wheelchair Arm Pads, or if you encounter any issues not covered in this manual, please contact your retailer or the ALCO customer support team for assistance. Provide your product model number (AL-66055BK-A3) when seeking support.

For further information, you may visit the official ALCO website or contact the seller from whom you purchased the product.

Related Documents - AL-66055BK-A3

	<p>Alco & Danfoss Thermostatic Expansion Valves: Product Catalog</p> <p>Comprehensive catalog of Alco and Danfoss Thermostatic Expansion Valves (TX Valves). Details TI, BA, C, HF, TRAE+, TCLE, TJLE, TJRE, TER, TIR, THR, T2, TE2, TUA, TUAE series with features, specifications, refrigerant compatibility, and capacity data for HVAC and refrigeration.</p>
	<p>Alco Medium Duty PVC Corrugated Conduit (ALCCCM Range) - Datasheet</p> <p>Datasheet for Alco's Medium Duty PVC Corrugated Conduit (ALCCCM Range). Features include flexibility, durability, UV stabilization, and compliance with AS/NZS 61386.22:2022. Includes detailed specifications for various conduit sizes.</p>
	<p>Alco Hex Head Screw Bolt Anchors - ALCHSBA Range Specifications</p> <p>Detailed specifications for Alco's ALCHSBA range of hex head screw bolt anchors, featuring dual thread, friction anchoring, and easy installation.</p>
	<p>ALCO Product Catalog: Cleaning Chemicals, Supplies & Equipment</p> <p>ALCO offers a comprehensive range of cleaning chemicals, supplies, and equipment, including Mandala Brand products. Discover disinfectants, floor care, laundry solutions, kitchen hygiene, and more, backed by expert service and support for various industries.</p>
	<p>AlcoLimit RBT Mini User Manual: Setup, Operation, and Specifications</p> <p>Comprehensive user manual for the AlcoLimit RBT Mini breathalyser. Learn about initial setup, performing tests, specifications, troubleshooting, recalibration, warranty, and important safety information.</p>
	<p>Alpina Watches: Warranty and Operating Instructions Manual</p> <p>Comprehensive warranty and operating instructions for Alpina watches, covering model details, care, maintenance, and international warranty information. Includes guides for various models like AL-235, AL-240, and more.</p>