



Manuals.plus /

- › Joerns Healthcare /
- › EasyCare Electric Adjustable Bed Frame User Manual

Joerns Healthcare B0787JMBW3

EasyCare Electric Adjustable Bed Frame User Manual

Brand: Joerns Healthcare | Model: B0787JMBW3

INTRODUCTION

The EasyCare Electric Adjustable Bed Frame is a versatile healthcare bed designed to meet the diverse needs of patients and caregivers in various care facilities. This manual provides essential information for the safe and effective setup, operation, and maintenance of your EasyCare bed frame.



Figure 1: EasyCare Electric Adjustable Bed Frame overview. This image shows the complete bed frame with its adjustable sections and base.

SAFETY INFORMATION

Adherence to the following safety guidelines is crucial for preventing injury and ensuring the longevity of the product:

- Always use the pedal locking system to secure the bed in place before making any adjustments to its position or when transferring a patient.
- Exercise caution when expanding the bed's width to maintain stability and prevent potential tip-overs.
- Assemble the bed frame meticulously according to the provided instructions to avoid improper setup and potential hazards.
- Handle all surfaces and components of the bed with care to prevent damage and ensure proper functionality.

SETUP AND FEATURES

The EasyCare bed frame is designed for flexible configuration and ease of use. Key features that facilitate setup and adaptability include:

- **Integrated Width Expansion:** The bed frame allows for easy width adjustment to accommodate various patient needs and room layouts. This feature is designed for quick and secure modification.
- **Pedal Locking System:** A robust pedal-activated locking system ensures the bed remains stationary when

required, providing stability during patient transfers or adjustments.

- **Redesigned Panel Retentions:** Improved panel retention mechanisms ensure secure attachment of head and foot boards, enhancing overall bed stability.
- **New Mattress Corner Stops:** Adjustable corner stops are integrated to secure the mattress in place, preventing shifting as the bed's width is adjusted.

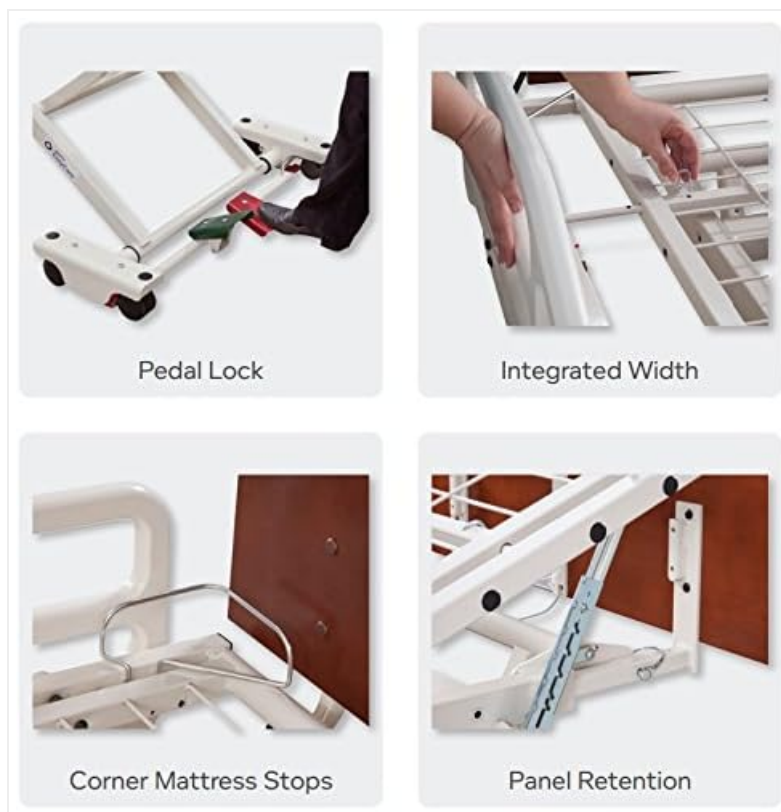


Figure 2: Detailed view of key EasyCare bed frame features, including the pedal lock mechanism, integrated width adjustment, corner mattress stops, and panel retention system.

Standard Features Overview

The EasyCare bed frame incorporates several standard features designed for patient comfort, caregiver convenience, and durability:

- **A. 7-inch to 30-inch Travel Range:** Offers a wide range of height adjustments, from an ultra-low height for safety to a high position for caregiver access.
- **B. Revolutionary Arched Legs:** Deliver a strong and stable wide stance, providing ample equipment clearance beneath the bed.
- **C. Four-section Mattress Support Platform:** Provides full range of comfort and therapeutic positions, designed to reduce patient and caregiver injuries.
- **D. Therapeutic Grid Surface:** Improves airflow and offers optimal support for the sleep surface.
- **E. Hand Pendant and Optional Staff Control:** Allows both patient and caregiver to adjust the bed to desired positions for comfort or care needs.
- **F. Toroidal Architecture:** The bed's steel frame engineering is significantly stronger than competitor's angle iron bed frames.
- **G. EasyLift Leg Lift:** Provides quick positioning of the knee section, ideal for patients with blood flow disorders.
- **H. Four DC Motors:** Supply quiet and smooth operation during frame articulation.



Standard Features

- A. 7-inch to 30-inch travel range** offers safety, whether resting at an ultra low height or providing care at the highest height.
- B. Revolutionary arched legs** deliver a super-strong wide stance for equipment clearance.
- C. Four-section mattress support platform** offers a full range of comfort and therapeutic positions designed to reduce patient and caregiver injuries and potentially reduce shear.
- D. Therapeutic grid surface** offers improved airflow and ultimate support for the newly reinforced sleep surface.
- E. Hand pendant and optional staff control** allow both patient and caregiver the flexibility to adjust the bed into the desired position for comfort or care needs.
- F. Toroidal architecture** engineering of the bed's steel frame is up to 30 times stronger than competitor's angle iron bed frames.
- G. EasyLift Leg Lift** provides quick positioning of the knee section of the bed for improved circulation – ideal for patients with blood flow disorders.
- H. Four DC motors** supply amazingly quiet and smooth operation during frame articulation. *The EasyCare now features a 4-year warranty on the frame electronics.*

Figure 3: Diagram illustrating the standard features of the EasyCare bed frame, labeled A through H, corresponding to the descriptions above.

OPERATING INSTRUCTIONS

The EasyCare bed frame is operated using an ergonomically designed hand pendant, which provides intuitive control over the bed's various functions.

Using the Hand Pendant

The hand pendant features soft-touch edges, finger grips, and raised buttons for ease of use, especially for geriatric needs. It may also include optional enhanced backlighting and illuminated lock-out status indicators.



Simplified Controls

Ergonomically designed pendant meets geriatric needs featuring soft touch edges, finger grips and optional enhanced back lighting, raised buttons and illuminated lock-out status.

Figure 4: Close-up of the simplified control hand pendant, showing its ergonomic design and button layout for adjusting bed positions.

To adjust the bed, press and hold the corresponding buttons on the hand pendant. The four DC motors ensure quiet and smooth articulation of the bed sections. Common adjustments include:

- **Head Section Adjustment:** Raises or lowers the head portion of the bed.
- **Foot Section Adjustment:** Raises or lowers the foot portion of the bed, including the EasyLift Leg Lift for knee section positioning.
- **Overall Bed Height Adjustment:** Changes the bed's height from its lowest (7 inches) to highest (30 inches) position.
- **Trendelenburg/Reverse Trendelenburg:** Allows tilting the entire bed frame.

Always ensure the patient is safely positioned before making any adjustments. For detailed button functions, refer to the labels on your specific hand pendant.

MAINTENANCE

Proper maintenance ensures the longevity and hygienic operation of your EasyCare bed frame.

- **Cleaning:** Wipe all surfaces of the bed frame with a dry cloth regularly to remove dust and debris. For deeper cleaning, use a mild, non-abrasive cleaner and a damp cloth, then dry thoroughly. Avoid harsh chemicals that may damage the finish.
- **Inspection:** Periodically inspect all moving parts, electrical cords, and connections for any signs of wear, damage, or loose fittings. Ensure the pedal locking system engages securely.
- **Lubrication:** Consult a qualified technician for any necessary lubrication of moving parts, as specified by Joerns Healthcare.

TROUBLESHOOTING

If you encounter issues with your EasyCare bed frame, consider the following basic troubleshooting steps before contacting support:

- **Bed Not Responding:**
 - Check if the power cord is securely plugged into a functional electrical outlet.
 - Ensure all connections between the hand pendant, motors, and control box are firm.

- Verify that any lock-out features on the hand pendant or control box are disengaged.
- **Unusual Noises During Operation:**
 - Ensure no objects are obstructing the moving parts of the bed frame.
 - Check for any loose bolts or components that may be causing rattling.
- **Bed Not Holding Position:**
 - Confirm that the pedal locking system is fully engaged.
 - Ensure the bed is on a level surface.

If the issue persists after performing these checks, please contact Joerns Healthcare customer support for further assistance.

SPECIFICATIONS

Attribute	Detail
Brand	Joerns Healthcare
Model Number (ASIN)	B0787JMBW3
Manufacturer	Joerns Healthcare
Color	Brown
Size	Full
Package Dimensions	80 x 35 x 30 inches
Item Weight	220 Pounds
Product Care Instructions	Wipe with Dry Cloth
Date First Available	December 12, 2017

WARRANTY AND SUPPORT

Warranty Information: The EasyCare Electric Adjustable Bed Frame comes with a **4-year warranty on the frame electronics**, ensuring reliability and peace of mind regarding its operational components.

Customer Support: For any inquiries, technical assistance, or warranty claims, please contact Joerns Healthcare directly. You can find their contact information on their official website or through your purchase documentation.

For more information, visit the Joerns Healthcare Store on Amazon.