

## Skil-Care 751680

# Skil-Care Contour Plus Cushion with Gel Pad User Manual

Model: 751680 | Brand: Skil-Care

## INTRODUCTION

This manual provides essential information for the proper use, care, and maintenance of your Skil-Care Contour Plus Cushion with Gel Pad. Designed for enhanced comfort and postural support, this cushion features a wedged design and an integrated gel pad to assist in pressure distribution and stability for individuals using wheelchairs or requiring specialized seating.

## PRODUCT OVERVIEW

The Skil-Care Contour Plus Cushion is engineered to provide superior seating support and pressure relief. Key features include:

- **Advanced Contoured Support:** Features a wedged base to prevent forward sliding, promoting proper pelvic alignment and reducing sacral sitting for improved postural stability.
- **Pressure-Relieving Gel Pod:** An integrated gel insert in the coccyx area distributes weight evenly, alleviating pressure points and minimizing the risk of pressure ulcers during prolonged sitting.
- **Durable High-Density Foam Core:** Provides firm support that retains its shape over time, ensuring consistent comfort and stability for long-term use.
- **Low-Shear LSII Cover:** A moisture-proof, wipe-clean cover reduces friction and shear, protecting sensitive skin and facilitating easy maintenance.

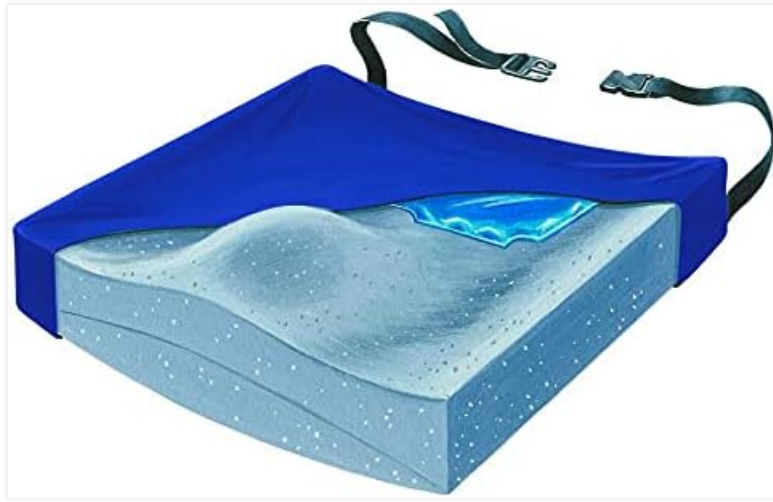


Figure 1: Cutaway view of the Skil-Care Contour Plus Cushion, illustrating the high-density foam core and the integrated gel pad.

**SkiL-Care™**  
CORPORATION

**Skil-Care products  
make life better  
and safer for  
patients**



Figure 2: The Skil-Care Contour Plus Cushion, showing its full profile and securing straps.

1. **Unpack the Cushion:** Carefully remove the cushion from its packaging. Inspect for any signs of damage.
2. **Orient the Cushion:** Identify the front and back of the cushion. The wedged design means the front (higher side) should face forward, towards the user's knees, and the back (lower side) towards the backrest.
3. **Place in Wheelchair/Chair:** Position the cushion firmly on the seat of the wheelchair or desired chair. Ensure it is centered and sits flat.
4. **Secure the Cushion (if applicable):** If your cushion includes securing straps, attach them around the backrest or frame of the chair to prevent the cushion from shifting during use. Adjust straps for a snug fit without compressing the cushion.



**SkiL-Care is a one-stop shop for all of your healthcare needs.**



**From bed and wheelchair positioning to mattresses to skin protection, we have everything you need to live a comfortable and independent life.**



Figure 3: Proper placement of the SkiL-Care Contour Plus Cushion within a wheelchair.

## USAGE GUIDELINES

For optimal comfort and therapeutic benefit, observe the following guidelines:

- **Proper Seating Posture:** Ensure the user sits with their back fully against the backrest of the chair, allowing the cushion's wedged design to promote correct pelvic tilt and prevent forward sliding.



- **Regular Repositioning:** While the cushion aids in pressure relief, regular repositioning and pressure redistribution are still recommended for individuals at high risk of pressure injuries.
- **Monitor Skin Integrity:** Periodically check the user's skin for any signs of redness, irritation, or pressure points, especially in areas contacting the gel pad.
- **Avoid Punctures:** Do not use sharp objects near the cushion, as this may damage the gel pad or foam core.

## Skil-Care™ CORPORATION



**Skil-Care products are for people of all walks of life.**

Figure 4: A user demonstrating proper seating on the Skil-Care Contour Plus Cushion.

## CARE AND MAINTENANCE

Proper care extends the life of your cushion and maintains its hygienic condition.

### Cleaning Instructions:

- **Cover Cleaning:** The Low-Shear LSII cover is moisture-proof and wipe-clean. Use a mild soap and water solution or a hospital-grade disinfectant to wipe down the cover. Rinse thoroughly with clean water and allow to air dry completely before re-use. Do not machine wash or tumble dry the cover.
- **Cushion Core:** The foam core and gel pad should not be submerged in water. If the core becomes soiled,

spot clean with a damp cloth and mild detergent. Ensure the core is completely dry before replacing the cover.

- **Product Care Instructions:** Hand Wash Only.

Storage:

Store the cushion in a clean, dry place away from direct sunlight and extreme temperatures. Avoid placing heavy objects on the cushion during storage to prevent deformation.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Cushion feels uncomfortable or too firm.	Incorrect orientation; initial break-in period; user's specific needs.	Ensure the cushion is oriented correctly (wedged side forward). Allow a few days for the foam to adapt. Consult a healthcare professional if discomfort persists.
Cushion slides or shifts on the chair.	Straps not secured or adjusted properly; smooth chair surface.	Ensure securing straps are tightly fastened. If no straps, consider a non-slip mat underneath the cushion.
Cover appears dirty or stained.	Regular use and spills.	Wipe clean with mild soap and water or disinfectant. Do not machine wash. Refer to "Care and Maintenance" section.
Signs of skin irritation or redness.	Prolonged sitting without repositioning; improper cushion fit; underlying medical condition.	Increase frequency of repositioning. Re-evaluate cushion fit. Consult a healthcare professional immediately.

SPECIFICATIONS

Attribute	Detail
Brand	Skil-Care
Model Number	751680
Material	Gel Pad, High-Density Foam
Cover Material	Low-Shear LSII (Moisture-proof, Wipe-clean)
Color	Blue
Dimensions (W x D x H)	18" W x 16" D x 5" H (Front) x 2" H (Back)
Item Weight	4.59 pounds
Firmness Description	Firm
Recommended Use	Wheelchair, General Seating
Care Instructions	Hand Wash Only (Cover: Wipe Clean)

Attribute	Detail
UPC	671509207100


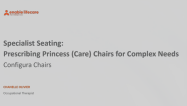


## WARRANTY AND SUPPORT



Specific warranty details for the Skil-Care Contour Plus Cushion are not provided within this manual. For warranty information, technical support, or to inquire about replacement parts, please contact Skil-Care directly through their official website or customer service channels. Please refer to the product packaging or the manufacturer's website for the most current contact information.

You can visit the Skil-Care Store on Amazon for more information:[Skil-Care Store](#)

© 2024 Skil-Care Corporation. All rights reserved. Information subject to change without notice.

### Related Documents - 751680

	<p><a href="#">5 Stars United Seat Cushion Product Use and Care Guide</a></p> <p>Instructions for using and cleaning the 5 Stars United Seat Cushion, including how to restore its shape, proper placement, and washing instructions for the cover.</p>
	<p><a href="#">Specialist Seating: Prescribing Princess (Care) Chairs for Complex Needs</a></p> <p>A guide to understanding and prescribing specialist seating, focusing on the Configura Advance chair for individuals with complex postural, neurological, or mobility needs. Covers learning objectives, indications for use, chair features, clinical justification, and case studies.</p>
	<p><a href="#">2025 OTC Benefits Catalog   Health and Wellness Products</a></p> <p>Explore the 2025 Over-The-Counter (OTC) Benefits Catalog from Security Health Plan, featuring a wide range of health, wellness, and personal care products from leading brands like CURAD, Medline, Tylenol, and more.</p>
	<p><a href="#">Contour Legacy Leg Pillow: Instructions for Use and Care</a></p> <p>Detailed instructions for the Contour Legacy Leg Pillow by Tristar Products, Inc. Learn how to use it, manage the break-in period, address new foam smell, understand foam discoloration, and follow detailed cleaning guidelines for your Tristar Products leg pillow.</p>

 <p>The image shows the Stryker IsoTour Gel Support Surface, a blue and orange medical device. Below it is the cover of the operations manual, which features the Stryker logo and the product name.</p>	<p><a href="#">Stryker IsoTour Gel Support Surface Operations Manual</a></p> <p>This operations manual provides detailed instructions for the Stryker IsoTour Gel Support Surface, a medical device for pressure injury prevention and treatment. It covers setup, operation, maintenance, and replacement procedures, highlighting features like TruTurn® technology and low air loss for acute care settings.</p>
 <p>The image shows the cover of the WHILL Model RK manual. It features a blue and white design with the WHILL logo and the product name. The cover also includes a small image of the device and some text in Japanese.</p>	<p><a href="#">WHILL Model RK: Enhancing Lifestyle and Independence with Personal Mobility</a></p> <p>Explore how the WHILL Model RK personal mobility device empowers users to regain independence, enjoy daily activities, and pursue hobbies, featuring user testimonials and a detailed care plan.</p>