



[Manuals.plus](#) /

> [Vive Mobility](#) /

> Vive Mobility 2 in 1 Walker Wheelchair Combo User Manual

Vive Mobility MOB1018RED

Vive Mobility 2 in 1 Walker Wheelchair Combo User Manual

Model: MOB1018RED

INTRODUCTION

The Vive Mobility 2 in 1 Walker Wheelchair Combo is a versatile mobility aid designed to provide both walking support and comfortable transport. This innovative device seamlessly transitions between a rollator walker and a transport wheelchair, offering maximum mobility and flexibility for users. Its foldable aluminum frame ensures portability and ease of storage, making it an ideal solution for various situations.



Image: The Vive Mobility 2 in 1 Walker Wheelchair Combo in Ruby Red, showcasing its dual functionality and 300lb capacity.

KEY FEATURES

- **2-in-1 Design:** Easily converts from a rollator walker to a transport wheelchair, providing adaptable mobility support.
- **Adjustable Handle Height:** Ergonomically designed handles can be adjusted to suit individual user height for optimal comfort and control.
- **Dual Braking System:** Features an integrated dual braking system for enhanced safety and control during use.
- **Comfortable Seating:** Equipped with a spacious 18-inch mesh seat and a soft padded backrest for comfortable resting or transport.
- **Removable Footrests:** Includes adjustable and removable footrests for comfortable transportation in wheelchair mode.
- **Compact & Foldable:** Lightweight aluminum frame (approx. 21 lbs) folds easily for convenient storage and

transport in car trunks or small spaces.

- **High Weight Capacity:** Safely supports users up to 300 pounds.
- **Integrated Storage:** Comes with a spacious detachable storage bag for personal items and an integrated cane holder.
- **Durable Construction:** Made with a corrosion-resistant aluminum frame and 8-inch all-terrain wheels for stability and longevity.



Walker & Wheelchair in One

2 mobility aids in 1

- Easily Foldable
- Adjustable Footrests
- Customizable Backrest

Image: A person using the Vive Mobility 2 in 1 Walker Wheelchair Combo in walker mode, highlighting its versatility.

Made to Fit your Needs

Customizable Experience



Image: A detailed diagram illustrating the various features of the combo unit, including adjustable handles, removable footrests, spacious bag, and cane holder.

SETUP AND ASSEMBLY

The Vive Mobility 2 in 1 Walker Wheelchair Combo comes largely pre-assembled for quick and easy setup. Assembly can typically be completed in around 5 minutes.

Initial Unpacking and Unfolding:

1. Carefully remove the unit from its packaging.
2. Place the folded unit on a flat, stable surface.
3. Gently pull the sides of the frame apart until the unit is fully extended and the seat locks into place. Ensure all locking mechanisms are secure before use.

Attaching Footrests:

The footrests are removable and can be attached when using the unit as a transport chair.

1. Locate the attachment points on the front lower frame of the unit.
2. Align the footrest pegs with the corresponding holes and insert them.
3. Push down firmly until the footrests click securely into place. Ensure they are stable before use.

Adjusting Handle Height:

The handles can be adjusted to the user's height for comfortable and ergonomic use.

1. Locate the height adjustment pins on the handle tubes.
2. Press the pins and slide the handles up or down to the desired height.
3. Release the pins, ensuring they fully engage in the adjustment holes to lock the handles in position.

Your browser does not support the video tag.

Video: An official product video from Vive Health Direct demonstrating the features and setup of the Vive Rollator Wheelchair (Red).

Your browser does not support the video tag.

Video: An official product video from Vive Health Direct introducing the Vive 2 in 1 Rollator Transport Wheelchair and its functionalities.

OPERATING INSTRUCTIONS

Using as a Rollator Walker:

1. Ensure the footrests are folded up or removed if not needed.
2. Adjust the handle height so your arms are slightly bent at the elbow when holding the grips.
3. To move, push the rollator forward a comfortable distance, then step forward. Maintain a steady pace.
4. To engage the brakes, squeeze the brake levers on the handles. To lock the brakes for sitting, push the brake levers down until they click into the locked position. Pull them up to release.
5. To sit, ensure the brakes are locked. Turn around so your back is to the seat, and slowly lower yourself onto the seat.


Converting to a Transport Chair:

1. Ensure the unit is on a flat, stable surface and the brakes are locked.
2. If not already attached, securely attach the footrests as described in the Setup section.
3. The backrest can be adjusted for comfort in transport mode.
4. The user should sit facing forward with their feet on the footrests.
5. The caregiver should stand behind the unit, holding the handles, and release the brakes to push the user.

Folding for Storage/Transport:

1. Ensure the unit is on a flat surface and the brakes are unlocked.
2. Lift the strap or handle located on the center of the seat. The unit will begin to fold inward.
3. Continue lifting until the unit is fully folded. You may need to gently push the sides together.

4. The footrests can be folded up or removed for a more compact size.



Easy to transport
Designed for convenience

- Foldable • Weighs 19 lbs • Travel Friendly

Image: The combo unit in its folded state, demonstrating its compact design for easy transport and storage.



Premium Comfort

Adjustable to your needs

- Comfortable & Spacious Seat
- Padded Backrest • Travel Friendly

Image: A person seated in the unit in transport chair mode, highlighting the comfortable and spacious seating.

MAINTENANCE AND CARE

Regular maintenance ensures the longevity and safe operation of your Vive Mobility 2 in 1 Walker Wheelchair Combo.

- **Cleaning:** Wipe down the frame and seat regularly with a damp cloth and mild detergent. Avoid abrasive cleaners.
- **Wheel Inspection:** Periodically check the wheels for debris, wear, and proper rotation. Remove any hair or lint that may impede movement.
- **Brake Check:** Test the brakes regularly to ensure they engage and lock properly. If brakes feel loose or ineffective, they may need adjustment.
- **Frame Integrity:** Inspect the aluminum frame for any signs of damage, cracks, or loose fasteners. Do not use

the unit if any structural damage is observed.

- **Storage:** Store the unit in a dry, clean place away from extreme temperatures and direct sunlight when not in use.

TROUBLESHOOTING

If you encounter any issues with your Vive Mobility 2 in 1 Walker Wheelchair Combo, refer to the following common solutions:

- **Unit does not fold or unfold smoothly:** Check for any obstructions in the folding mechanism or seat area. Ensure the seat is lifted completely when folding.
- **Brakes are not engaging or locking:** Inspect the brake cables for kinks or damage. The brake tension may need adjustment. Consult a qualified technician if you are unsure how to adjust them.
- **Wheels are not rolling freely:** Check for debris (hair, lint, dirt) wrapped around the wheel axles. Clean the wheels thoroughly. Ensure the wheels are properly attached and not loose.
- **Unit feels unstable or wobbly:** Verify that all frame connections are secure and that the unit is fully unfolded and locked into place. Check for any loose screws or bolts and tighten them if necessary.

If the issue persists after attempting these troubleshooting steps, please contact Vive Mobility customer support for assistance.

SPECIFICATIONS

Brand	Vive Mobility
Model Number	MOB1018RED
Color	Ruby Red
Material	Aluminum
Item Weight	21 Pounds
Product Dimensions (Unfolded)	28.5"D x 25.5"W x 38"H
Weight Capacity	300 lbs
Seat Width	18 inches
Wheel Size	8 inches (All-Terrain)

Easy to Fold & Take With You

Weighs only 19 pounds

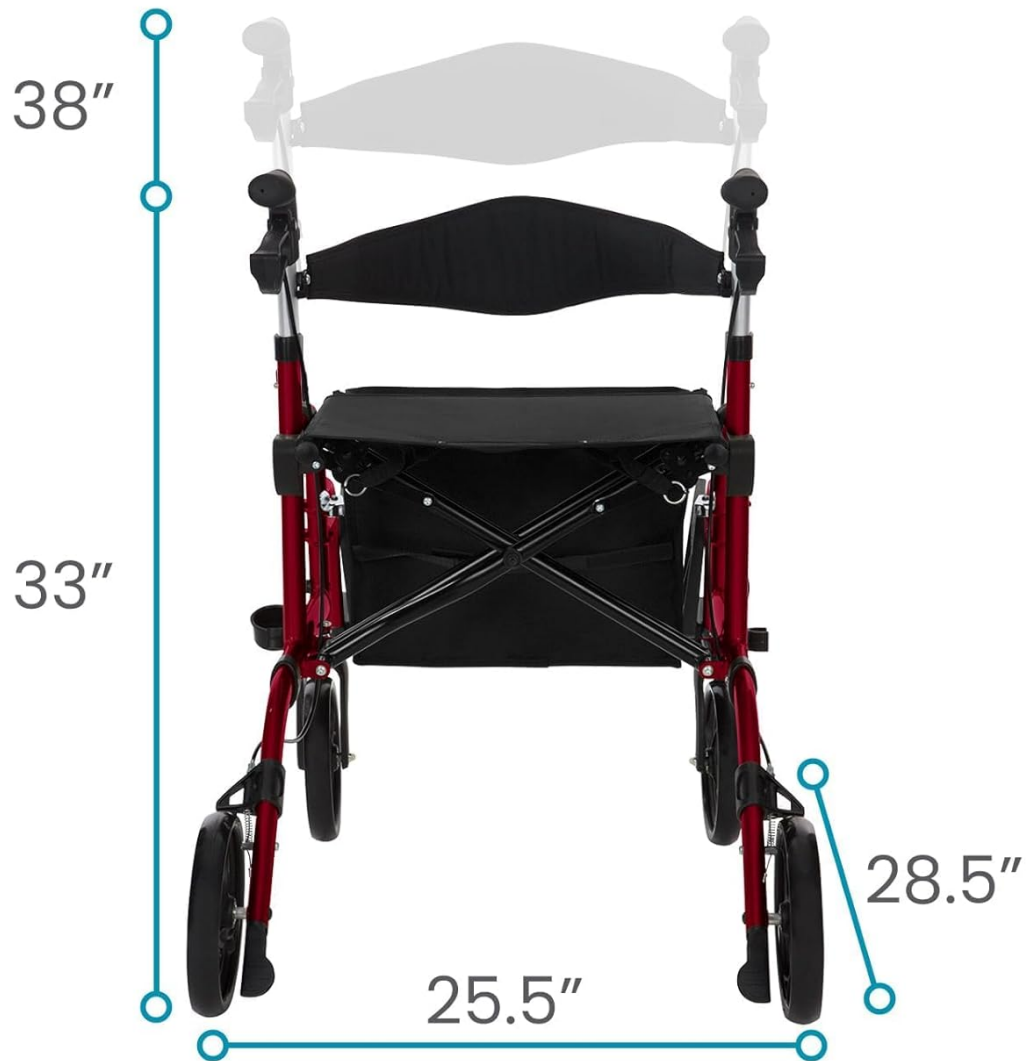


Image: A visual representation of the key dimensions of the Vive Mobility 2 in 1 Walker Wheelchair Combo.

WARRANTY AND SUPPORT

The Vive Mobility 2 in 1 Walker Wheelchair Combo comes with a **Lifetime Guarantee**, reflecting Vive Health's commitment to quality and durability. We offer 24/7 support and high availability of replacement parts to ensure your continued satisfaction.

For any questions, assistance, or to inquire about replacement parts, please contact Vive Health Direct customer service.

Lasts For Years

Stable. Sturdy. Lifetime Guarantee



Image: Vive Health's commitment to quality, highlighted by their Lifetime Guarantee.

© 2024 Vive Mobility. All rights reserved.

For more information, visit [Vive Mobility Product Page](#)