

axion STIM-PRO 2S

axion Vaginal Probe STIM-PRO 2S Instruction Manual

For use with Sanitas EMS Devices SEM 42, SEM 43, SEM 44

1. INTRODUCTION

This manual provides essential instructions for the safe and effective use of the axion Vaginal Probe, model STIM-PRO 2S. This probe is designed for electrostimulation of the pelvic floor muscles, primarily for incontinence treatment and post-childbirth regression, when used with compatible EMS (Electrical Muscle Stimulation) devices.

The axion Vaginal Probe STIM-PRO 2S is specifically compatible with Sanitas EMS devices SEM 42, SEM 43, and SEM 44. Always refer to the instruction manual of your EMS device for detailed information on program selection and general device operation.

2. SAFETY INFORMATION

Please read all safety warnings and instructions carefully before using the axion Vaginal Probe. Failure to do so may result in injury or improper use.

- **Intended Use:** This probe is intended for vaginal electrostimulation of the pelvic floor muscles only. Do not use for any other purpose.
- **Compatibility:** Use this probe exclusively with Sanitas EMS devices SEM 42, SEM 43, and SEM 44. Using it with incompatible devices may cause malfunction or injury.
- **Hygiene:** This is a personal hygiene product. For health protection and hygiene reasons, the probe cannot be returned once its original sealed packaging has been opened.
- **Material Sensitivity:** The probe's electrodes are made of stainless steel, which may contain nickel. If you have a known nickel allergy, consult your doctor before use.
- **Medical Advice:** Always consult a healthcare professional before starting any electrostimulation therapy, especially if you have pre-existing medical conditions, are pregnant, or have an IUD.
- **Discomfort or Pain:** If you experience any discomfort, pain, or unusual sensations during use, immediately stop the stimulation and consult your doctor.
- **Storage:** Keep out of reach of children. Store in a clean, dry place.

3. SETUP

Follow these steps to prepare your axion Vaginal Probe for use:

1. **Inspect the Probe:** Before each use, visually inspect the probe and its connecting cables for any signs of damage, wear, or corrosion. Do not use if damaged.
2. **Clean the Probe:** Thoroughly clean the probe before its first use and after each subsequent use as described in the "Maintenance and Care" section.
3. **Apply Contact Gel:** For optimal conductivity and comfort, apply a sufficient amount of a water-based contact gel (such as axion contact gel) to the electrode surfaces of the probe.
4. **Connect to EMS Device:** Connect the probe's cable to the appropriate output channel of your compatible Sanitas EMS device (SEM 42, 43, or 44). Ensure the connection is secure.

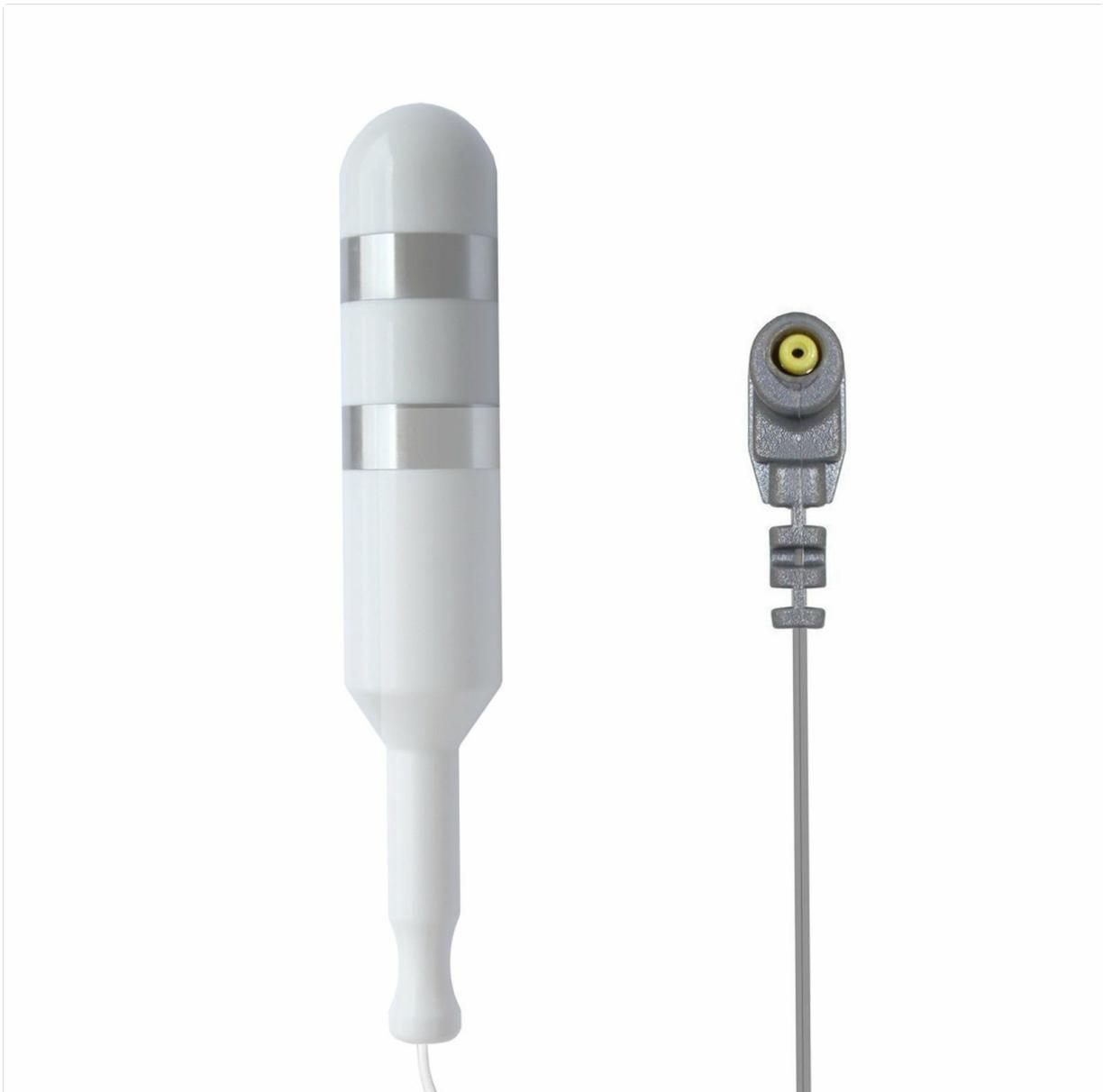


Image 1: The axion Vaginal Probe STIM-PRO 2S with its connecting cables. Note that the EMS device is sold separately.



Image 2: Detailed view of the axion Vaginal Probe, highlighting the electrode rings and cable attachment point.

4. OPERATING INSTRUCTIONS

Once the probe is connected to your EMS device and prepared, follow these general operating guidelines:

1. **Insertion:** Carefully insert the probe into the vagina, similar to a tampon, ensuring the electrode surfaces are in contact with the vaginal walls.
2. **Select Program:** Turn on your Sanitas EMS device and select the appropriate program for pelvic floor stimulation as recommended by your healthcare professional or as detailed in your EMS device's manual.
3. **Adjust Intensity:** Gradually increase the intensity of the stimulation on your EMS device until you feel a comfortable, yet effective, muscle contraction. Avoid setting the intensity too high, which could cause discomfort or pain.
4. **Duration:** Follow the recommended treatment duration as specified by your healthcare professional or the EMS device program.
5. **Removal:** After the session, turn off the EMS device and gently remove the probe.

5. MAINTENANCE AND CARE

Proper cleaning and storage will ensure the longevity and hygienic use of your probe:

- **Cleaning:** After each use, disconnect the probe from the EMS device. Clean the probe thoroughly with warm water and a mild, non-abrasive soap. Rinse well to remove all soap residue.
- **Disinfection:** For additional hygiene, the probe can be disinfected using a suitable medical device disinfectant spray or wipes. Follow the disinfectant manufacturer's instructions.
- **Drying:** Allow the probe to air dry completely before storing. Do not use heat to dry.
- **Storage:** Store the clean, dry probe in its original packaging or a clean, dry, and dust-free container. Keep it away from direct sunlight and extreme temperatures.

- **Cable Care:** Avoid kinking or sharply bending the connecting cables. Store them loosely coiled.

6. TROUBLESHOOTING

If you encounter issues during use, consider the following:

- **No Sensation or Weak Stimulation:**
 - Ensure the probe is correctly inserted and making good contact with the vaginal walls.
 - Verify that sufficient contact gel has been applied to the electrodes.
 - Check that the probe's cable is securely connected to the EMS device.
 - Confirm that the EMS device is turned on, has sufficient battery power, and the correct program is selected.
 - Gradually increase the intensity on the EMS device.
- **Discomfort or Irritation:**
 - Reduce the intensity of the stimulation.
 - Ensure adequate contact gel is used.
 - Check for any visible damage to the probe. Discontinue use if damaged.
 - If irritation persists, stop use and consult a healthcare professional.
- **Device Malfunction:** If the EMS device itself is not functioning, refer to its specific instruction manual for troubleshooting.

If problems persist after following these steps, contact the manufacturer or your supplier for further assistance.

7. SPECIFICATIONS

Feature	Specification
Model Number	STIM-PRO 2S
Brand	axion
Length	14.7 cm (147 mm)
Maximum Diameter	2.5 cm (25 mm)
Cable Length	30 cm
Electrode Material	Stainless Steel (may contain nickel)
Compatibility	Sanitas EMS devices SEM 42, SEM 43, SEM 44

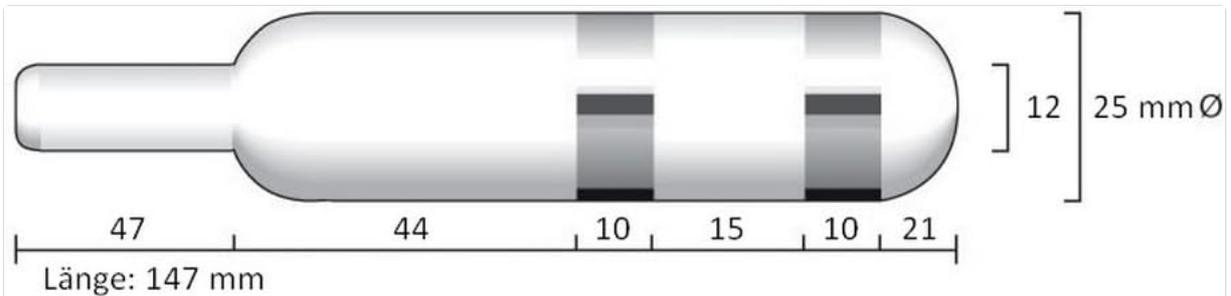


Image 3: Dimensional diagram of the axion Vaginal Probe STIM-PRO 2S.

8. WARRANTY AND SUPPORT

The axion Vaginal Probe STIM-PRO 2S is a product from a certified German manufacturer of medical devices, axion GmbH, ensuring high quality and reliability.

For specific warranty terms and conditions, please refer to the documentation included with your purchase or visit the official axion website. Warranty coverage typically applies to manufacturing defects.

If you require technical support, have questions about the product, or need to report an issue, please contact your retailer or axion GmbH directly through their official customer service channels. Ensure you have your product model number (STIM-PRO 2S) and purchase details available when contacting support.