



[Manuals.plus](#) /

› [Oral-B](#) /

› [Oral-B Glide Pro-Health Comfort Plus Floss User Manual](#)

Oral-B 037000724841

Oral-B Glide Pro-Health Comfort Plus Floss User Manual

Brand: Oral-B | Model: 037000724841

[Overview](#)

[Setup](#)

[Operating](#)

[Maintenance](#)

[Troubleshooting](#)

[Specifications](#)

[Warranty &](#)

[Support](#)

1. PRODUCT OVERVIEW

The Oral-B Glide Pro-Health Comfort Plus Floss is engineered to provide a superior cleaning experience. Its advanced design allows it to slide up to 50% easier in tight spaces compared to regular floss, effectively removing plaque from between teeth and just below the gumline. This floss is strong, smooth, and shred-resistant, leaving your mouth with a fresh minty clean feeling.

Key features include:

- **#1 Dentist Recommended:** Trusted by dental professionals.
- **Easy Sliding:** Designed to glide effortlessly in narrow gaps between teeth.
- **Powerful Cleaning:** Combines comfort with effective plaque removal.
- **Prevents Gingivitis:** Aids in preventing gum inflammation as part of a comprehensive oral hygiene routine.
- **Natural Wax Coating:** Features a light coating for improved grip and handling.
- **Travel-Friendly:** Compact size for convenient use on the go.

Oral-B

Glide FLOSS

VALUE
3
PACK

TOUGH ON
PLAQUE
GENTLE ON
GUMS

MICRO-TEXTURED



Oral-B

Glide

Oral-B

Glide

Oral-B

Glide

COMFORT
PLUS

MINT FLOSS

TOTAL 131.2 YD (120 m)

Image: A three-pack of Oral-B Glide Pro-Health Comfort Plus Mint Floss, highlighting its 'Tough on Plaque, Gentle on Gums' and 'Micro-Textured' features.

2. INITIAL SETUP

The Oral-B Glide Pro-Health Comfort Plus Floss comes ready for immediate use. No assembly or special preparation is required. Simply open the packaging and extract a single floss dispenser.





Image: Three individual Oral-B Glide Comfort Plus Mint Floss cases, as they appear when removed from the multi-pack packaging.

3. HOW TO USE ORAL-B GLIDE FLOSS

For optimal oral hygiene, follow these steps for effective flossing:

1. **Dispense Floss:** Pull out approximately 18-24 inches (45-60 cm) of floss from the dispenser.
2. **Wrap Floss:** Wind most of the floss around your middle fingers, leaving about 1-2 inches (2.5-5 cm) to work with.
3. **Grip Floss:** Hold the floss tightly between your thumbs and forefingers.
4. **Guide Floss:** Gently guide the floss between your teeth using a rubbing motion. Do not snap the floss into your gums.
5. **Clean Tooth Surface:** Once the floss reaches the gumline, curve it into a C-shape against one tooth. Gently slide it into the space between the gum and the tooth.
6. **Move Up and Down:** Move the floss up and down against the tooth surface and just below the gumline. Repeat this for each tooth, using a clean section of floss for every tooth.
7. **Remove Floss:** To remove the floss, use the same back-and-forth motion to bring it up and out from between the teeth.

Glides smoothly with **Micro-Textured Fiber**



Image: A person demonstrating the proper technique of flossing with Oral-B Glide floss, showing how it glides smoothly with micro-textured fiber.

Removes Plaque where brushes don't reach

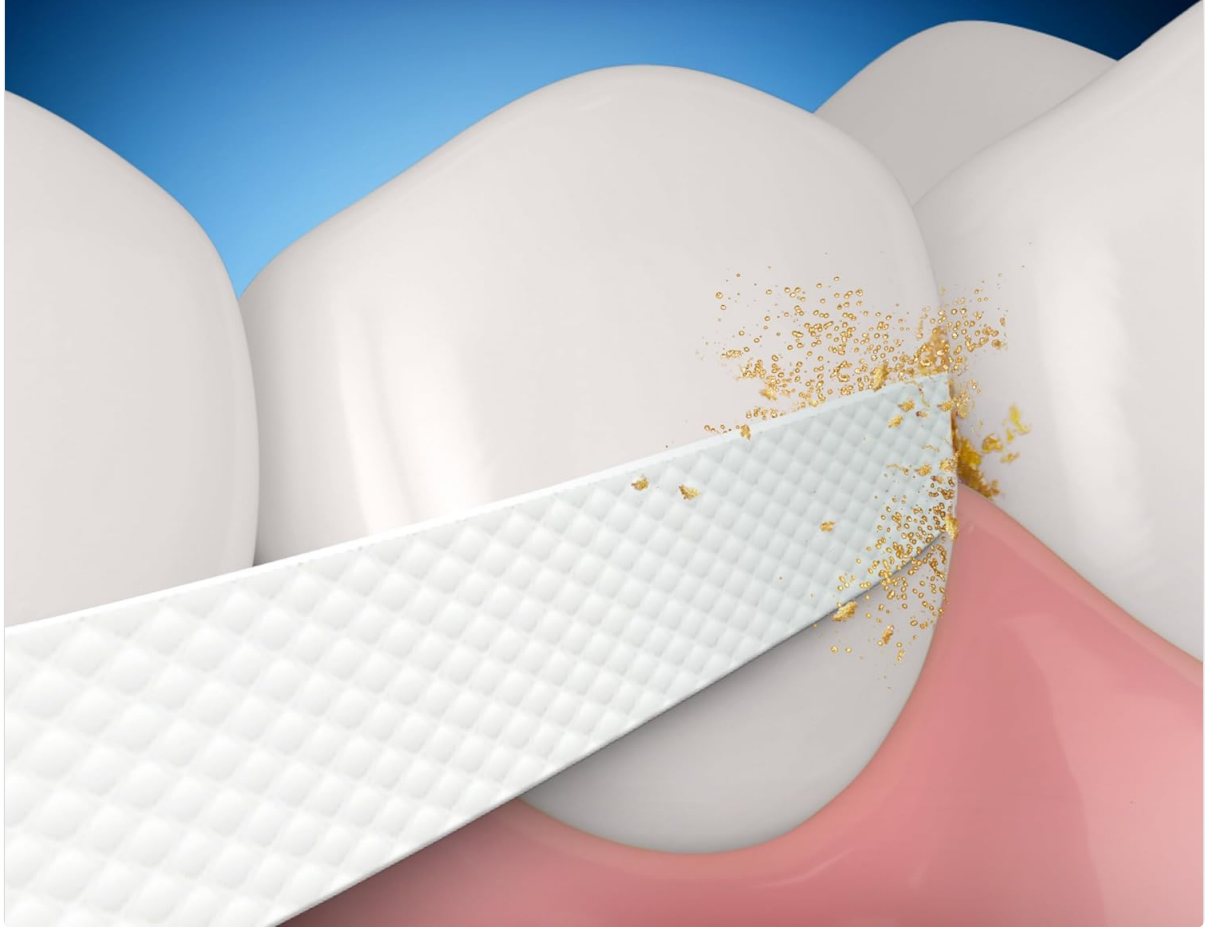


Image: An animated illustration showing dental floss effectively removing plaque from between teeth, highlighting its ability to reach areas brushes cannot.

Slides Easily into tight spaces

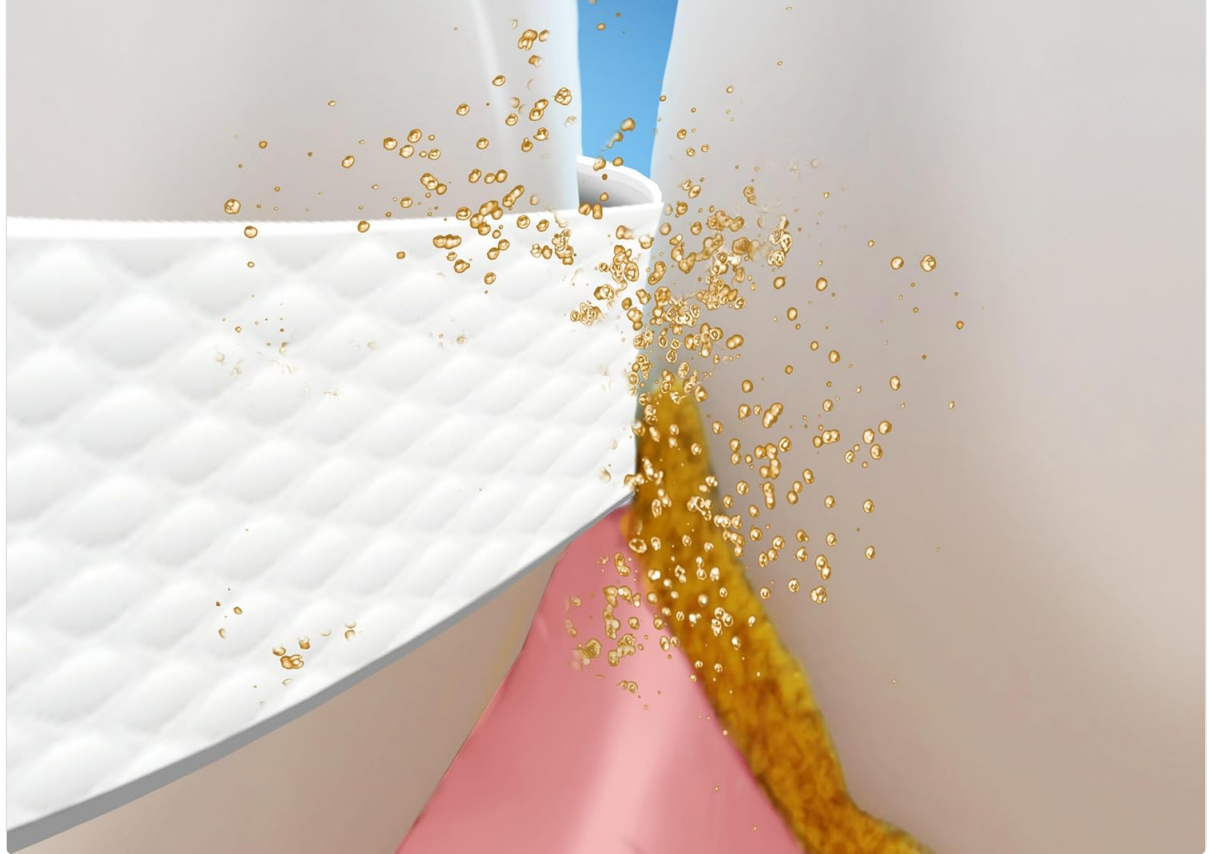


Image: An animated illustration demonstrating how Oral-B Glide floss easily slides into tight spaces between teeth, dislodging debris.

4. MAINTENANCE AND CARE

Oral-B Glide Pro-Health Comfort Plus Floss is designed for single-use per section. After flossing, dispose of the used floss responsibly. The floss dispenser itself requires minimal maintenance.

- **Storage:** Store the floss dispenser in a cool, dry place away from direct sunlight and extreme temperatures.
- **Cleaning:** The exterior of the floss case can be wiped clean with a damp cloth if necessary. Do not immerse the case in water.
- **Recycling:** The plastic floss case is recyclable through the Oral-B 'Recycle on Us' program. Visit recycling.oralb.com for more information. The carton packaging is also recyclable.



Image: The back of the Oral-B Glide floss packaging, showing the UPC barcode and a QR code for the 'Recycle on Us' program, indicating the recyclability of the plastic case and carton.

5. COMMON QUESTIONS

While Oral-B Glide Pro-Health Comfort Plus Floss is designed for ease of use, here are answers to some common concerns:

- **Floss Shredding:** This floss is designed to be shred-resistant. If you experience shredding, ensure you are not applying excessive force or snapping the floss. Gently guide it between teeth. The micro-textured fiber is engineered for durability.
- **Difficulty Sliding:** While designed for easy sliding, very tight contacts between teeth can still present a challenge. Ensure you are using a gentle back-and-forth motion to ease the floss through.
- **Gum Sensitivity:** If your gums are sensitive or bleed, this may indicate gum inflammation (gingivitis). Continue flossing gently daily, and consult your dentist if symptoms persist. This floss is gentle on gums.

Durable, Gentle Fibers slide without shredding



Image: A close-up view of the Oral-B Glide floss, emphasizing its 'Micro-Textured Fiber' which contributes to its durability and gentle sliding properties.

6. PRODUCT SPECIFICATIONS

Feature	Detail
Brand	Oral-B
Model Number	037000724841
Product Dimensions	1.39 x 3.14 x 5.99 inches (packaging)
Item Weight	0.64 ounces
Floss Type	Dental Floss (Tape)
Flavor	Mint
Special Feature	Flavored, Micro-Textured Fiber

Manufacturer	Procter & Gamble - HABA Hub
ASIN	B07FLBBWJR

7. WARRANTY AND CUSTOMER SUPPORT

Oral-B stands behind the quality of its products. Your satisfaction with Oral-B Glide Pro-Health Comfort Plus Floss is guaranteed.

- **Satisfaction Guarantee:** If you are not satisfied with your purchase, you may be eligible for a money-back guarantee.
- **Claiming Guarantee:** For guarantee claims, please call **1-877-769-8791** within 60 days of purchase. You will need to provide the UPC from the packaging and your purchase receipt.
- **General Inquiries:** For any other questions or support regarding Oral-B products, please visit the official Oral-B website or refer to the contact information provided on the product packaging.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is a dental hygiene tool.