

[Manuals.plus /](#)

› [MOMMED /](#)

› MomMed Ovulation and Pregnancy Test Kit User Manual

## MOMMED 1993

# MomMed Ovulation and Pregnancy Test Kit User Manual

Model: 1993

## PRODUCT OVERVIEW

The MomMed Ovulation and Pregnancy Test Kit (Model 1993) provides a comprehensive at-home solution for tracking ovulation and detecting early pregnancy. This kit is designed for ease of use and high accuracy, empowering individuals to understand their fertility cycles better.



## WHAT'S INCLUDED

Your MomMed Ovulation and Pregnancy Test Kit includes the following components:

- 50 Ovulation (LH) Test Strips
- 20 Home Pregnancy (HCG) Test Strips
- 70 Urine Collection Cups
- Time Tracking Chart
- Instruction Manual



Image: Detailed view of the kit contents, showing individually packaged test strips, urine cups, and accompanying documentation.

## SETUP AND PREPARATION

Before performing any test, ensure you have all necessary components ready and understand the instructions. For best results, collect your sample using the provided urine cups.

1. Unpack one individually sealed test strip (LH or HCG) from its pouch.
2. Prepare a clean, dry urine collection cup.
3. Collect a urine sample in the cup. For LH tests, it is recommended to collect samples at the same time each day, avoiding first morning urine for optimal LH detection. For HCG tests, first morning urine is recommended for early detection.

**Be With You.**  
**MomMed**



*Image: Visual guide showing the preparation steps: having test strips ready, collecting urine sample in a cup.*

## **OPERATING INSTRUCTIONS**

---

### **Ovulation (LH) Test Instructions**

The LH test detects the luteinizing hormone surge, which typically occurs 24-36 hours before ovulation. This indicates your most fertile window.


1. Immerse the test strip into the urine sample in the cup for 5 seconds. Do not immerse past the MAX line indicated on the strip.
2. Remove the strip and lay it flat on a clean, dry, non-absorbent surface.
3. Read the results in 10-15 minutes. Do not interpret results after 15 minutes.




# LH

## Instructions

- 1 Immerse the strip into the urine
- 2 Take the strip out after 5 seconds

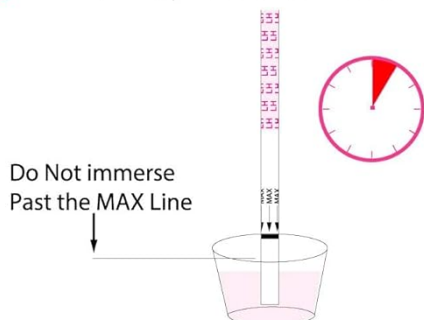

- 3 Lay the strip flat

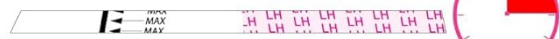

- 4 Read results in 5 minutes

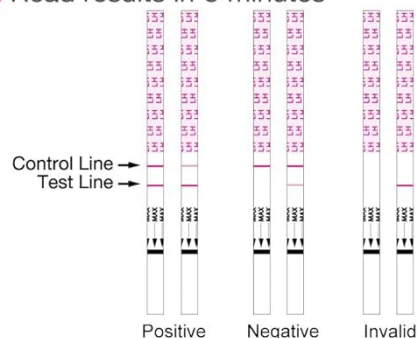


Image: Step-by-step visual instructions for using the LH ovulation test strip, including dipping, waiting, and result interpretation.

## Pregnancy (HCG) Test Instructions

The HCG test detects the presence of Human Chorionic Gonadotropin, a hormone produced during pregnancy. This test can detect pregnancy as early as 11 days past ovulation.


1. Immerse the test strip into the urine sample in the cup for 5 seconds. Do not immerse past the MAX line indicated on the strip.
2. Remove the strip and lay it flat on a clean, dry, non-absorbent surface.
3. Read the results in 5 minutes. Do not interpret results after 5 minutes.




# HCG

## Instructions

- 1 Immerse the strip into the urine
- 2 Take the strip out after 5 seconds

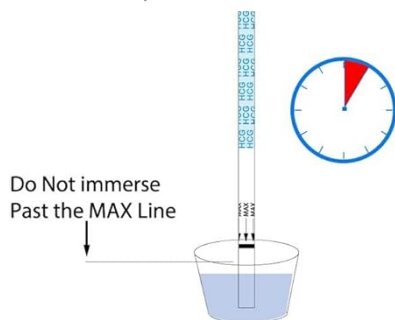

- 3 Lay the strip flat

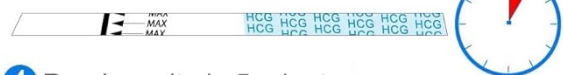

- 4 Read results in 5 minutes

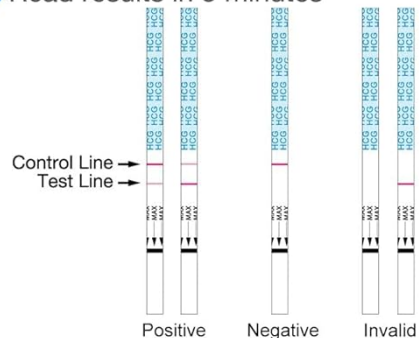


Image: Step-by-step visual instructions for using the HCG pregnancy test strip, including dipping, waiting, and result interpretation.

## Interpreting Results

Both LH and HCG test strips will display a Control Line (C) and a Test Line (T).

- **Positive Result:** Two lines appear, and the Test Line (T) is as dark as or darker than the Control Line (C). This indicates an LH surge (for ovulation tests) or pregnancy (for HCG tests).
- **Negative Result:** Only the Control Line (C) appears, or the Test Line (T) is lighter than the Control Line (C). This indicates no LH surge or no pregnancy.
- **Invalid Result:** No Control Line (C) appears, or the Control Line (C) is very faint. The test is invalid and should be repeated with a new strip.



25  
mIU  
25 mIU



## MomMed Ovulation and Pregnancy Test Kit

Image: Visual guide to interpreting test results as Positive, Negative, or Invalid.

### Time Tracking Chart

Use the provided Time Tracking Chart to log your daily test results. This helps visualize your cycle progression and identify patterns for future planning.

# LH

Using the Time Tracking Chart we provided, you can track your ovulation efficiently and easily.

### Ovulation Test Strip Instruction

**1** Immerse the strip into the urine for 5 seconds

**2** Lay the test on a flat and read result in 5 minutes

**3** Read the results

Control Line  
Test Line

Positive Negative Invalid

**TIPS:**  
If two color bands are visible, they do not indicate the same result. Only when the test line is as **DARK** as the control line or **DARKER** than the control line is a positive result, indicating that you will ovulate within 24-48 hours.

### Time Tracking Chart

CD	DATE-TIME	TEST STRIP	RESULT
9	09/10 7:50 am		N
10	10/10 7:50 am		N
	11/10 7:50 am		N
	12/10 7:50 am		N
	13/10 7:50 am		N
	14/10 7:50 am		✓
	14/10 9:50 am		✓
	14/10 4:00 pm		✓
	14/10 6:00 pm		✓
15	15/10 7:50 am		N
16	16/10 7:50 am		N

CD: Day of Cycle to Begin Testing. Please refer to the Instruction Insert for more details.



Image: A hand marking results on the MomMed Time Tracking Chart, showing how to log date, time, test strip type, and result.

## Easy to Read Results

We are ready to solve your confusion If you have any problems with the test results

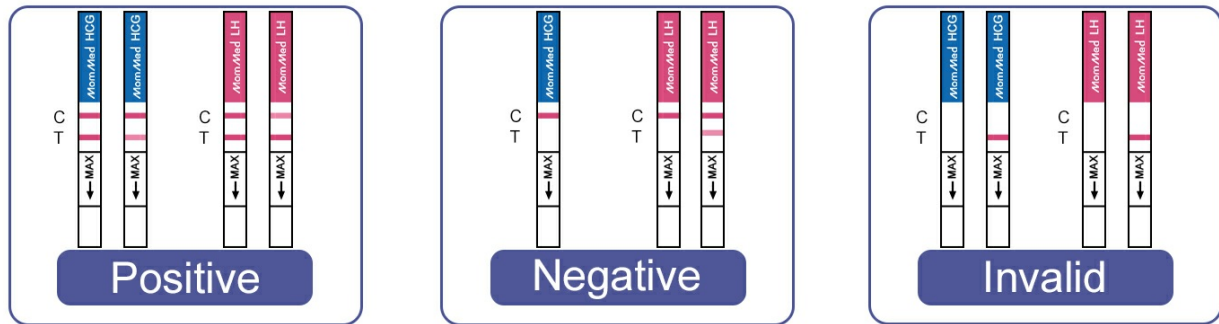


Image: The Time Tracking Chart overlaid on an image of a mother holding a baby, emphasizing its use for family planning.

### Official Product Video: How to Use

Your browser does not support the video tag.

Video: An official MomMed video demonstrating the usage of both ovulation and pregnancy tests, including sample collection, dipping, and result interpretation. It also highlights the time tracking chart.

### MAINTENANCE AND STORAGE

To ensure the longevity and accuracy of your MomMed test strips:

- Store the test strips in their original sealed pouches at room temperature (36-86°F or 2-30°C).
- Do not store in direct sunlight, moisture, or heat.
- Do not freeze the test strips.
- Use the test strip immediately after opening the sealed pouch.
- Dispose of used test strips and urine cups according to local regulations.

# LH & HCG

Wider 3TT (0.12in) Strips

Impressive 99% Accuracy Rate

Easy to Carry and Keep



Image: Close-up of MomMed LH and HCG test strips, highlighting their wider 3mm size and the convenience of a small pouch for storage.

## TROUBLESHOOTING

If you encounter issues with your test results, consider the following:

- **Invalid Result (No Control Line):** This indicates the test did not function correctly. Ensure the strip was immersed for the correct duration and not past the MAX line. Repeat the test with a new strip.
- **Faint Lines:** For LH tests, a faint test line indicates no LH surge. For HCG tests, a faint test line may indicate very early pregnancy or a low HCG level. Re-test in 48-72 hours with first morning urine for HCG.
- **Irregular Cycles:** If you have irregular cycles, consistent daily testing with LH strips and using the time tracking chart or a fertility tracking app (like the MomMed App) can help identify your unique ovulation pattern.

# Tracking by MomMed Free App



MomMed



Image: A smartphone displaying the MomMed app interface, showing features for cycle tracking and community support. For further assistance, please refer to the detailed instruction manual included in your kit or contact MomMed customer support.

## SPECIFICATIONS

Feature	Detail
Brand	MOMMED
Model Number	1993
Test Type	Ovulation (LH) & Pregnancy (HCG)
Measurement Accuracy	99%
Number of Pieces	70 (50 LH, 20 HCG, 70 Urine Cups)
Special Feature	Accurate

Package Dimensions	6.46 x 4.72 x 2.52 inches; 8.16 ounces
--------------------	--

## **IMPORTANT INFORMATION & LEGAL DISCLAIMER**

---

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

For any product-related inquiries or support, please refer to the contact information provided in the product packaging or visit the official MomMed website.