

Pic Solution 02004990000000

Pic Solution AirCube Kit - Piston Nebulizer Replacement Kit User Manual

This manual provides essential instructions for the safe and effective use of your Pic Solution AirCube Kit. Please read carefully before use and keep for future reference.

1. PRODUCT OVERVIEW

The Pic Solution AirCube Kit is a comprehensive replacement accessory set designed for use with compatible Pic piston nebulizers. It ensures continued effective aerosol therapy at home, providing all necessary components for a complete and hygienic treatment experience.



Image 1.1: Overview of the Pic Solution AirCube Kit components.

- **Complete Kit:** Contains all necessary accessories for aerosol therapy, compatible with Pic piston nebulizer systems.
- **Compact Design:** Allows for replacement of components or provides a dedicated kit for each family member.
- **Simple and Practical:** Accessories kit for convenient aerosol therapy, stored in an easy-to-organize transparent pouch.
- **Ergonomic Shape:** Soft and ergonomic masks designed to fit both adult and child faces comfortably.
- **Includes:** Soft adult mask, soft child mask, mouthpiece, nasal fork, air tube, nebulizer kit, and a practical transparent pouch.



Image 1.2: Front view of the complete kit contents.



Image 1.3: Nebulizer cup and air tube.



Image 1.4: Adult nebulizer mask.





Image 1.5: Child nebulizer mask.



Image 1.6: Mouthpiece and nasal fork.



Image 1.7: Kit stored in its transparent pouch.

2. SETUP AND ASSEMBLY

Follow these steps to correctly assemble your Pic Solution AirCube Kit with a compatible Pic piston nebulizer device. Refer to the diagram below for visual guidance.

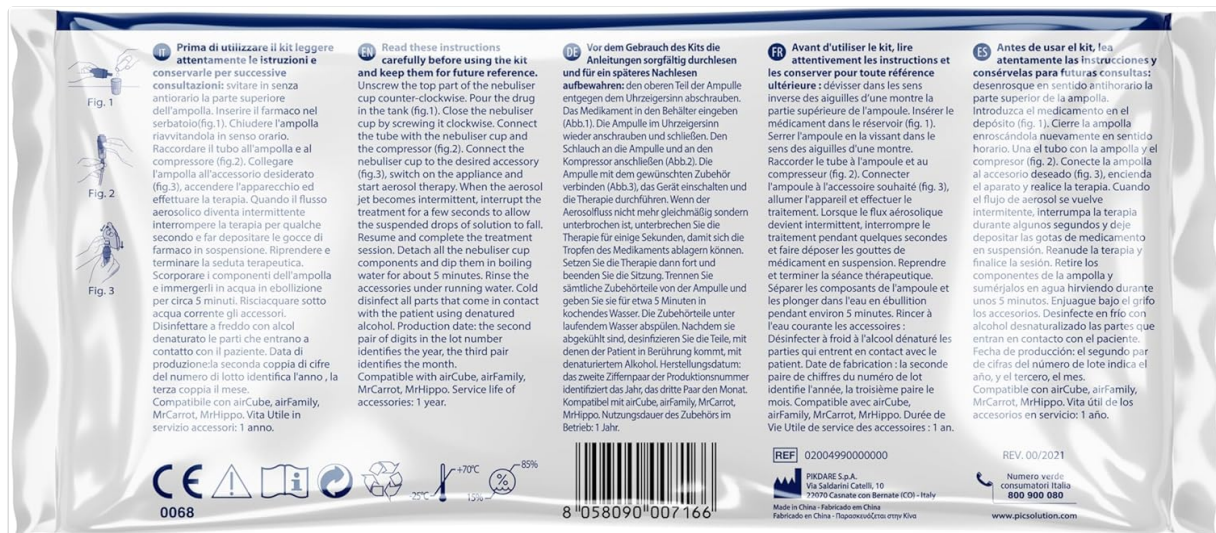


Image 2.1: Detailed instruction diagram for assembly and use.

- 1. Prepare the Nebulizer Cup:** Unscrew the top part of the nebulizer cup counter-clockwise. Pour the prescribed medicine into the tank. Close the nebulizer cup by screwing it clockwise until secure.
- 2. Connect the Air Tube:** Connect one end of the air tube to the nebulizer cup and the other end to the compressor outlet of your Pic piston nebulizer device.
- 3. Attach the Accessory:** Connect the desired accessory (adult mask, child mask, mouthpiece, or nasal fork) to the nebulizer cup. Ensure all connections are firm and secure.

3. OPERATING INSTRUCTIONS

Once the AirCube Kit is assembled and connected to your nebulizer device, proceed with the therapy as follows:

- Start Therapy:** Switch on your Pic piston nebulizer device. The device will begin to produce an aerosol mist.

- **Inhale:** Breathe normally through the chosen accessory (mask, mouthpiece, or nasal fork).
- **Intermittent Aerosol:** If the aerosol jet becomes intermittent during treatment, briefly interrupt the therapy for a few seconds. This allows any suspended drops of solution to fall back into the nebulizer cup. Resume breathing until the medicine is fully nebulized.
- **Complete Treatment:** Continue the therapy session as prescribed by your doctor.



Image 3.1: User demonstrating nebulizer use with a mask.

Important: Always follow your doctor's instructions regarding medication dosage and treatment duration. Do not exceed the recommended fill level for the nebulizer cup.

4. MAINTENANCE AND CLEANING

Proper cleaning and maintenance are crucial for the hygiene, effectiveness, and longevity of your AirCube Kit components. Clean the accessories after each use.

1. **Disassemble:** After each therapy session, detach all components of the nebulizer cup and the accessories (masks, mouthpiece, nasal fork).
2. **Sterilize:** Dip the disassembled components in boiling water for approximately 5 minutes.
3. **Rinse:** Rinse all parts thoroughly under running water to remove any residue.
4. **Disinfect:** Disinfect all parts that come into contact with the patient using denatured alcohol.
5. **Dry and Store:** Ensure all components are completely dry before reassembling or storing them in the transparent pouch.

Replacement: The nebulizer kit and its accessories have a service life of approximately 1 year. It is recommended to replace the kit annually or sooner if any signs of wear, damage, or reduced performance are observed.

5. TROUBLESHOOTING

If you encounter any issues while using your Pic Solution AirCube Kit, consider the following common troubleshooting steps:

- **No Aerosol Production:**
 - Ensure the nebulizer cup is correctly assembled and sealed tightly.

- Verify that the air tube is securely connected to both the nebulizer cup and the compressor outlet.
- Check if the compressor device is powered on and functioning correctly.

- **Weak Aerosol Flow:**

- Inspect the air tube for any kinks, blockages, or damage.
- Ensure the nebulizer cup and accessories are clean and free from medication residue.
- If the kit is old, it might require replacement as per the recommended service life (1 year).

- **Leaking from Connections:**

- Confirm that all connections (nebulizer cup, air tube, accessories) are tight and properly seated.
- Ensure the nebulizer cup is closed securely after adding medication.

If problems persist after attempting these steps, please refer to the main instruction manual of your Pic piston nebulizer device or contact customer support for further assistance.

6. SPECIFICATIONS

Feature	Detail
Product Name	Air Cube Aerosol Kit
Units	1 Total
Brand	Pic Solution
Country of Origin	Italy
ASIN	B07G478WJ9
Product Dimensions (L x W x H)	22 x 14 x 6 cm
Product Weight	100 grams
Manufacturer	Pikdare S.p.A.
Item Model Number	02004990000000
Date First Available	July 3, 2022

7. SUPPORT AND WARRANTY

The Pic Solution AirCube Kit is a medical device (CE 0068). Always read the warnings and instructions for use carefully before operating.

This product is manufactured by Pikdare S.p.A.

For customer support or inquiries, please contact:

- **Customer Service Phone:** [800 900 080](tel:800900080)

- **Official Website:** www.picsolution.com

The AirCube Kit is compatible with Pic AirCube, AirFamily, MrCarrot, and MrHippo nebulizer devices. The recommended service life for the accessories in this kit is 1 year from the date of first use.

© 2026 Pic Solution. All rights reserved.