

Zoll 710-0731-01

Zoll Power Supply and Power Cord for Zoll Impact 731 Ventilators - User Manual

INTRODUCTION

This manual provides essential information for the safe and effective use of the Zoll Power Supply and Power Cord, designed specifically for Zoll Impact 731 Ventilators.

The package includes:

- Zoll Power Supply (Part Number: 708-0042-00)
- Zoll Power Cord (Part Number: 024-0012-00)



Figure 1: The Zoll Power Supply and Power Cord. This image shows the black rectangular power supply unit connected to a power cord with a standard 3-prong plug on one end and a specialized 2-pin connector on the other, designed for the Zoll Impact 731 Ventilator.

SETUP AND INSTALLATION

Follow these steps to properly set up your Zoll Power Supply and Power Cord:

1. **Unpack Components:** Carefully remove the power supply unit and power cord from their packaging. Inspect both items for any visible damage. Do not use if damaged.
2. **Connect Power Cord to Supply:** Insert the standard 3-prong end of the power cord into a compatible wall outlet. Connect the other end of the power cord to the input port on the Zoll Power Supply unit. Ensure a secure connection.
3. **Connect Power Supply to Ventilator:** Locate the power input port on your Zoll Impact 731 Ventilator. Connect the 2-pin connector from the Zoll Power Supply to this port. Ensure the connection is firm and fully seated.
4. **Verify Connection:** Once all connections are made, ensure that the power supply unit is placed on a stable, flat

surface, away from excessive heat or moisture.

Important: Always ensure the ventilator is turned off before connecting or disconnecting the power supply.

OPERATING INSTRUCTIONS

The Zoll Power Supply and Power Cord are designed for continuous operation with the Zoll Impact 731 Ventilator. Once properly connected, the power supply will provide the necessary electrical current to operate and charge the ventilator.

- **Power On:** After connecting the power supply, turn on your Zoll Impact 731 Ventilator according to its specific operating instructions.
- **Power Indicator:** The power supply unit may have an indicator light to show that it is receiving power and functioning correctly. Refer to the ventilator's manual for specific power status indicators.
- **Disconnection:** To disconnect, first power off the ventilator, then unplug the power cord from the wall outlet, and finally disconnect the 2-pin connector from the ventilator.

MAINTENANCE AND CARE

Proper maintenance ensures the longevity and reliable performance of your power supply:

- **Cleaning:** Clean the exterior of the power supply and cords with a soft, dry cloth. Do not use liquid cleaners or abrasive materials. Ensure the unit is unplugged before cleaning.
- **Storage:** When not in use, store the power supply and cord in a cool, dry place, away from direct sunlight and extreme temperatures. Coil the cords loosely to prevent kinking or damage.
- **Inspection:** Periodically inspect the power cord and connectors for any signs of wear, fraying, or damage. If any damage is observed, discontinue use immediately and replace the unit.
- **Ventilation:** Ensure the power supply unit has adequate ventilation during operation to prevent overheating. Do not cover the unit.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No power to ventilator	Power cord not fully plugged in; Wall outlet not active; Faulty power supply/cord.	Ensure all connections are secure. Test the wall outlet with another device. If issues persist, the power supply or cord may need replacement.
Power supply unit is hot	Normal operation (slight warmth); Poor ventilation; Overload.	Ensure adequate airflow around the unit. If excessively hot, disconnect and contact support.
Ventilator not charging	Power supply not providing sufficient power; Ventilator battery issue.	Verify the power supply is correctly connected. Refer to the Zoll Impact 731 Ventilator manual for battery troubleshooting.

SPECIFICATIONS

Model Number: 710-0731-01

Included Components: Zoll Power Supply (708-0042-00), Zoll Power Cord (024-0012-00)

Compatible Devices: Zoll Impact 731 Ventilators

Output Voltage: 24 Volts (DC)

Connector Type: 2-Pin

Color: Black

Manufacturer: Zoll

Specification Met: UL

Date First Available: August 17, 2018

WARRANTY INFORMATION

Specific warranty details for the Zoll Power Supply and Power Cord (Part Number: 710-0731-01) are typically provided at the point of purchase or within the product packaging. Please retain your proof of purchase for warranty claims. For detailed warranty terms and conditions, please refer to the official Zoll website or contact Zoll customer support directly.

CUSTOMER SUPPORT

For technical assistance, product inquiries, or to report issues, please contact Zoll customer support:


- **Brand:** Zoll
- **Manufacturer:** Zoll
- **Website:** www.zoll.com
- **Model Number:** 710-0731-01

When contacting support, please have your product model number (710-0731-01) and any relevant purchase information ready.



© 2024 Zoll. All rights reserved.

Documents - Zoll – 710-0731-01

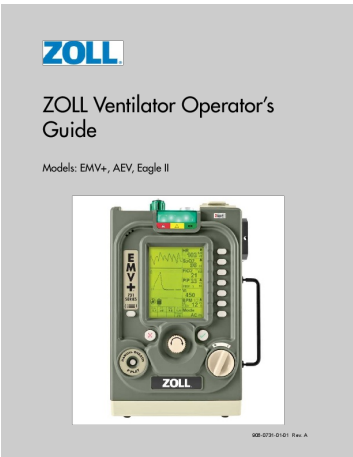

[\[pdf\]](#) User Manual Guide

Ventilator Operator s Guide book abizzarro ZOLL Medical Information about this manual 2020 Low O2 Supply Fault Alarm triggers when the supply pressure is 35 psig 241 Zoll User frankshospitalworkshop equipment documents ventilators user manuals

ZOLL Ventilator Operator s Guide Models: EMV , AEV, Eagle II 906-0731-01-01 Rev. A The issue date ... ese items, contact ZOLL or your local distributor. Part Number 024-0012-00 708-0042-00 703-0731-01 **710-0731-01** 704-0EMV-XX 708-0041-XX 402-0032-00 465-0024-00 465-0025-00 465-0027-01 465-0028-01 820-...

lang:en score:31 filesize: 3.67 M page_count: 164 document date: 2017-12-07



ZOLL Ventilator Operator's Guide

Models: EMV+, AEV, Eagle II



906-0731-01-01 Rev. A

[\[pdf\]](#) Guide

Ventilator Operator s Guide book abizzarro 906 0731 01 sf a zoll media public site products ventilators ashx |||

ZOLL Ventilator Operator s Guide Models: EMV , AEV, Eagle II 906-0731-01-01 Rev.


A The issue date ... ese items, contact ZOLL or your local distributor. Part Number 024-0012-00 708-0042-00 703-0731-01 **710-0731-01** 704-0EMV-XX 708-0041-XX 402-0032-00 465-0024-00 465-0025-00 465-0027-01 465-0028-01 820-...

lang:en **score:24** filesize: 9.33 M page_count: 164 document date: 2017-12-07



ZOLL Ventilator Operator's Guide

Models: EMV+, AEV, Eagle II

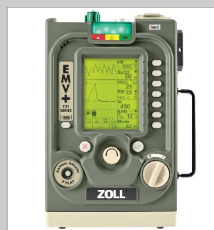

906-0731-01-05 Rev. B

[\[pdf\]](#) Guide

906 0731 01 05 B book DLELACHEUR ZOLL Ventilator Operator s Guide Zoll Gui sequence 1 isAllowed y bmet ewh org bitstream handle 20 500 12091 2105 |||

ZOLL Ventilator Operator s Guide Models: EMV , AEV, Eagle II 906-0731-01-05 Rev.


B The issue date ... ese items, contact ZOLL or your local distributor. Part Number 024-0012-00 708-0042-00 703-0731-01 **710-0731-01** 704-0EMV-XX 708-0041-XX 402-0032-00 465-0024-00 465-0025-00 465-0027-01 465-0028-01 820-...

lang:en **score:22** filesize: 5.6 M page_count: 156 document date: 2017-05-17



ZOLL Ventilator Operator's Guide

Models: EMV+, AEV, Eagle II

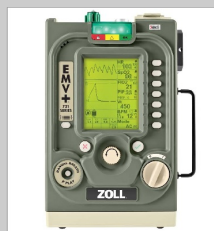

906-0731-01-05 Rev. B

[\[pdf\]](#) Guide Accessories

906 0731 01 05 B book DLELACHEUR ZOLL Ventilator Operator s Guide iFixitDo not connect to an electrical outlet controlled by a wall switch or dimmer The protection against defibrillator depends on the use of accessories Zoll Models EMV AEV Eagle IIdocuments cdn ifixit 4N3HmLnUpsppyuSH |||

ZOLL Ventilator Operator s Guide Models: EMV , AEV, Eagle II 906-0731-01-05 Rev.

B The issue date ... ese items, contact ZOLL or your local distributor. Part Number 024-0012-00 708-0042-00 703-0731-01 **710-0731-01** 704-0EMV-XX 708-0041-XX 402-0032-00 465-0024-00 465-0025-00 465-0027-01 465-0028-01 820-...

lang:en **score:21** filesize: 5.6 M page_count: 156 document date: 2017-05-17



Guide de l'utilisateur du respirateur ZOLL

Modèles : EMV+, AEV et Eagle II



906-0731-01-02 Rev. A

[\[pdf\]](#) Guide

906 0731 01 02 A book DLELACHEUR Guide de l utilisateur du respirateur ZOLL 15 mag 2017 — Chapitre 1 Informations générales DESCRIPTION DU PRODUIT rev a portable critical care ventilator french zoll media public site products ventilators ashx |||

Guide de l utilisateur du respirateur ZOLL Modles : EMV , AEV et Eagle II 906-0731-01-02 Rev. A Ce ... 1,83 m version pour les tats-Unis 703-0731-01 Bloc-batterie, 6,6 Ah, 14,8 V, lithium-ion, 4S3P **710-0731-01** Bloc d alimentation CA/CC et cordon d alimentation avec embout NEMA 5-15P 704-0EMV-XX ...

lang:fr **score:20** filesize: 5.71 M page_count: 156 document date: 2017-05-16

ZOLL

Guida all'uso del ventilatore ZOLL



906-0731-01 Rev F

[\[pdf\]](#)

906 0731 01 11 SF F book DLELACHEUR zoll media Product Materials product manuals ventilators ashx
sc lang en |||

Guida all uso del ventilatore ZOLL 906-0731-01-11 Rev. F La data di emissione della
Guida all uso ... 1,8 m Versione per Stati Uniti 703-0731-01 Batterie: 6,6 Ah, 14,8 V,
agli ioni di litio, 4S3P **710-0731-01** Alimentatore CA/CC e cavo con terminazione
NEMA 5-15P 704-0EMV-XX Prolunga da 230 cm ...

lang:it score:18 filesize: 8.49 M page_count: 162 document date: 2024-12-27

Ventilators

ZOLL



[\[pdf\]](#) Accessories Borchure

Vent Accessories Brochure indd tgibson vent accessories brochure pn 335 zoll media public site products
ventilators ashx la en gb hash 68C237C5867C60A2F09BA82FF8B2EAF9AB04AD86

Ventilators ACCESSORIES AND DISPOSABLES ACCESSORIES FOR ALL
MARKETS Ensures safe and effective tre ... or ZOLL ventilators. Power Supply 024-
0012-00 Power Cord 708-0042-00 Power Supply and Power Cord **710-0731-01**
OTHER Test lung for demonstrations and product education. One-Liter Dark Blue
Test Lung...

lang:en score:18 filesize: 1.11 M page_count: 7 document date: 2017-09-05

ZOLL

ZOLL Ventilator Operator's Guide

Models: EMV+, AEV, Eagle II



906-0731-01 Rev C

[\[pdf\]](#) Guide

906 0731 01 05 SF C book abizzarro zoll media public site products ventilators sf c ashx

ZOLL Ventilator Operator s Guide Models: EMV , AEV, Eagle II 906-0731-01-05 Rev.
C The issue date ... ese items, contact ZOLL or your local distributor. Part Number
024-0012-00 708-0042-00 703-0731-01 **710-0731-01** 704-0EMV-XX 708-0041-XX
402-0032-00 465-0024-00 465-0025-00 465-0027-01 465-0028-01 820-...

lang:en score:10 filesize: 10.22 M page_count: 156 document date: 2017-12-28



Guida all'uso del ventilatore ZOLL

Modelli: EMV+, AEV, Eagle II



906-0731-01-11 Rev. B

[\[pdf\]](#)

906 0731 01 11 B book DLELACHEUR Guida all uso del ventilatore ZOLL Dec 15 2017 — Capitolo 1
Informazioni generali DESCRIZIONE DEL PRODOTTO SF zoll media public site products ventilators sf b
ashx |||

Guida all uso del ventilatore ZOLL Modelli: EMV , AEV, Eagle II 906-0731-01-11 Rev.
B La data di e ... 1,8 m Versione per Stati Uniti 703-0731-01 Batterie: 6,6 Ah, 14,8 V,
agli ioni di litio, 4S3P **710-0731-01** Alimentatore CA/CC e cavo con terminazione
NEMA 5-15P 704-0EMV-XX Prolunga da 230 cm ...

lang:it **score:10** filesize: 5.58 M page_count: 154 document date: 2017-12-19



Portable Critical Care Ventilator Operator's Guide

Models: EMV+, AEV, Eagle II



731 Series

906-0731-01 Rev. J

[\[pdf\]](#) Guide

906 0731 01 Rev J book MHirshson SF api zoll media public site products ventilators sf j ashx |||
Portable Critical Care Ventilator Operator s Guide Models: EMV , AEV, Eagle II 731
Series 906-0731 ... ese items, contact ZOLL or your local distributor. Part Number
024-0012-00 708-0042-00 703-0731-01 **710-0731-01** 704-0EMV-XX 708-0041-XX
402-0032-00 465-0024-00 465-0025-00 465-0027-00 465-0028-00 820-...

lang:en **score:10** filesize: 10.71 M page_count: 134 document date: 2017-06-02



Guía del operador del ventilador ZOLL

Modelos: EMV+, AEV, Eagle II



906-0731-01-10 Rev. D

[\[pdf\]](#)

906 0731 01 10 SF D book DLELACHEUR zoll media public site products ventilators es sf d ashx
Gua del operador del ventilador ZOLL Modelos: EMV , AEV, Eagle II 906-0731-01-10
Rev. D La fecha d ... ulos, comunquese con su distribuidor local de ZOLL. Nmero ref.
024-0012-00 708-0042-00 703-0731-01 **710-0731-01** 704-0EMV-XX 708-0041-XX
402-0032-00 465-0024-00 465-0025-00 465-0027-00 465-0028-00 820-...

lang:es **score:10** filesize: 5.81 M page_count: 154 document date: 2020-05-21



Bedienungsanleitung für das ZOLL-Beatmungsgerät

Modelle: EMV+, AEV, Eagle II



906-0731-01-08 Rev. B

[\[pdf\]](#)

906 0731 01 08 SF B book DLELACHEUR rev b portable critical care ventilator german zoll media public
site products ventilators ashx |||

Bedienungsanleitung für das ZOLL-Beatmungsgerät Modelle: EMV , AEV, Eagle II 906-0731-01-08 Rev. B D ... an ZOLL oder den nächstgelegenen Vertriebspartner.

Teilenummer 024-0012-00 708-0042-00 703-0731-01 **710-0731-01** 704-0EMV-XX
708-0041-XX 402-0032-00 465-0024-00 465-0025-00 465-0027-01 465-0028-01 820-
...

lang:de **score:10** filesize: 5.49 M page_count: 156 document date: 2017-05-17



Guía del operador del ventilador ZOLL

Modelos: EMV+, AEV, Eagle II



906-0731-01-10 Rev. C

[\[pdf\]](#)

906 0731 01 10 SF C book DLELACHEUR rev c portable critical care ventilator spanish zoll media public
site products ventilators ashx

Guía del operador del ventilador ZOLL Modelos: EMV , AEV, Eagle II 906-0731-01-10
Rev. C La fecha de ... ublicación, comuníquese con su distribuidor local de ZOLL. Número ref.

024-0012-00 708-0042-00 703-0731-01 **710-0731-01** 704-0EMV-XX 708-0041-XX
402-0032-00 465-0024-00 465-0025-00 465-0027-00 465-0028-00 820-...

lang:es **score:9** filesize: 5.75 M page_count: 154 document date: 2017-09-22