



Manuals.plus /

› Rayovac /

› Rayovac Hearing Aid Batteries Size 312 Instruction Manual

## Rayovac BD15516

# Rayovac Hearing Aid Batteries Size 312 Instruction Manual

Model: BD15516

## INTRODUCTION

---

This manual provides essential instructions for the safe and effective use of Rayovac Size 312 Hearing Aid Batteries. These zinc-air batteries are designed for advanced hearing aid devices, offering reliable power and a long shelf life. Please read this manual thoroughly before use.

## SAFETY INFORMATION

---

- **Keep out of reach of children:** Batteries can be swallowed. Ingesting batteries can cause serious injury or death. Seek immediate medical attention if a battery is swallowed.
- **Single Use:** These are single-use batteries and are not rechargeable. Attempting to recharge them may lead to leakage or rupture.
- **Disposal:** Dispose of used batteries according to local regulations. Do not dispose of in fire.
- **Storage:** Store batteries in their original packaging in a cool, dry place away from direct sunlight and extreme temperatures.
- **Mercury-Free:** Rayovac Size 312 batteries feature a mercury-free formula.

## SETUP AND INSTALLATION

---

Proper handling and installation ensure optimal performance and safety.

### 1. Opening the Battery Package

The Rayovac battery package is designed for ease of use and secure storage. It features a lock dial system on the back to prevent accidental spinning and a front-facing door for easy access to individual batteries.

## Front Facing Door



*Image: The front-facing door allows for easy removal of a single battery.*

### 2. Removing a Battery

1. Locate the front-facing door on the package.
2. Open the door and carefully slide out one battery. The long, wide tabs on the battery are designed for easier handling.
3. Close the door to keep the remaining batteries clean and secure.

## Long, Wide Battery Tabs



*Image: The long, wide tab on the battery facilitates easier handling.*

### 3. Air Activation

Rayovac hearing aid batteries are zinc-air batteries, which means they are activated by air. Once you remove the protective tab, air enters the battery and activates it. For optimal performance, allow the battery to air activate for approximately 1 minute before inserting it into your hearing aid.

### 4. Inserting into Hearing Aid

1. Open the battery compartment of your hearing aid.
2. Insert the activated battery with the flat, positive (+) side facing up, or as indicated by your hearing aid manufacturer's instructions.
3. Close the battery compartment securely. Your hearing aid should now power on.

### Setup and Battery Change Demonstration Videos

**Meet Our Easy to Use, Award-Winning Package**

Your browser does not support the video tag.

*Video: This video demonstrates the features of the Rayovac battery packaging, including the lock dial system and front-facing door for easy battery access.*


### How to Change Your Hearing Aid Battery

Your browser does not support the video tag.

*Video: This video provides a step-by-step guide on how to replace a used hearing aid battery with a new Rayovac battery.*


## OPERATING INSTRUCTIONS

Once installed, your Rayovac Size 312 batteries will power your hearing aid. These batteries are engineered for advanced hearing aid devices, supporting features such as audio streaming and phone calls.

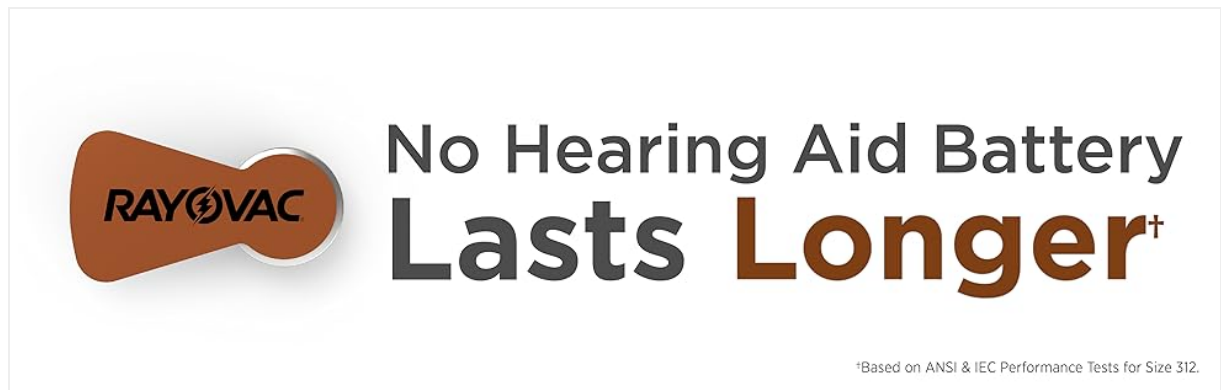

*Image: Rayovac batteries are designed for long-lasting performance.*

Rayovac Size 312 batteries are designed to deliver consistent power throughout their lifespan. No other Size 312 hearing aid battery lasts longer on ANSI & IEC performance tests.

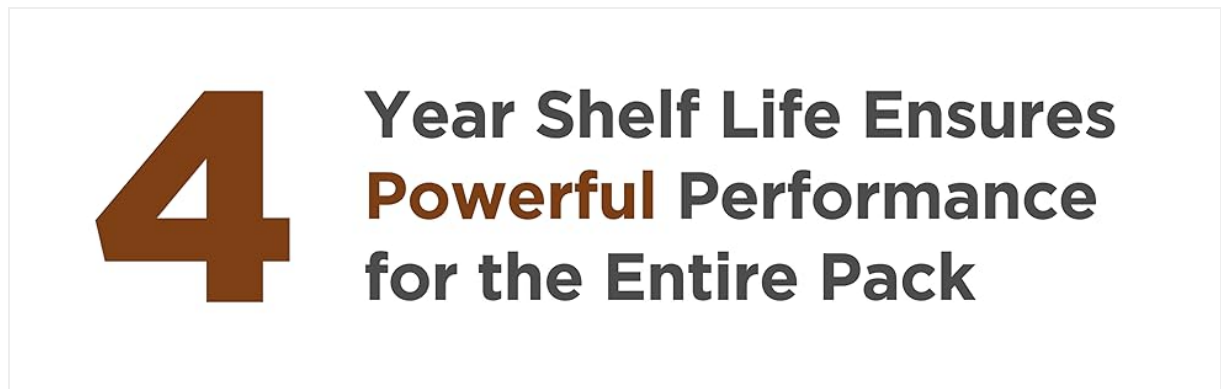

*Image: These batteries boast a 4-year shelf life, ensuring readiness when needed.*

## MAINTENANCE

Proper maintenance of your batteries and their storage can extend their effectiveness.

- **Secure Storage:** Always reclose the front-facing door of the battery package after removing a battery. This keeps the remaining batteries clean, secure, and protected from environmental factors.
- **Handling:** The long, wide tabs on the batteries are designed to make handling easier, especially for individuals with dexterity challenges.
- **Avoid Contamination:** Do not touch the battery with dirty or greasy hands, as this can affect performance.

## TROUBLESHOOTING

---

If you experience issues with your Rayovac hearing aid batteries, consider the following common solutions:

- **Hearing Aid Not Powering On:**

- Ensure the protective tab has been completely removed from the battery.
- Verify the battery is inserted correctly into the hearing aid compartment, with the correct polarity.
- Allow 1 minute for air activation after removing the tab before insertion.

- **Short Battery Life:**

- Ensure the battery tab was removed just prior to use, not significantly in advance.
- Check the expiration date on the battery package.
- Environmental factors like humidity and temperature can affect battery life.
- High-drain usage (e.g., continuous streaming) will naturally reduce battery life.

- **Difficulty Removing Batteries from Package:**

- Ensure the lock dial on the back of the package is disengaged to allow the batteries to rotate freely to the front-facing door.
- The package features a perforated rounded edge, allowing the dial to be completely removed for added convenience if preferred.

If problems persist, consult your hearing aid professional or contact Rayovac customer support.

## SPECIFICATIONS

---

Feature	Detail
Model Number	BD15516
Battery Size	312
Battery Cell Type	Zinc-Air
Voltage	1.45 Volts
Reusability	Single Use
Unit Count (per pack)	8 Count (Pack of 7)
Shelf Life	4 Years
Product Dimensions	0.31 x 0.31 x 0.14 inches
Item Weight	0.06 Pounds (per pack)
Manufacturer	Spectrum Brands Inc.

## WARRANTY AND SUPPORT

---

Rayovac products are manufactured by Spectrum Brands Inc. For any questions regarding product performance, warranty claims, or technical support, please contact Spectrum Brands Inc. directly or refer to the contact information provided on the product packaging. Keep your purchase receipt for warranty purposes.

