Widex TV-Dex

Widex TV-Dex Audio Streamer User Manual

1. INTRODUCTION

The Widex TV-Dex is a dedicated audio streaming system designed to enhance the television and HiFi listening experience for users of compatible Widex hearing aids. This system streams high-quality audio directly to your hearing aids with minimal delay, providing a clear and immersive sound experience. The TV-Dex system consists of two primary components: the TV-Controller II and the TV-Base.

2. PACKAGE CONTENTS

Please verify that all the following items are included in your package:

- TV-Controller II
- TV-Base
- Mini USB power supply
- Lanyard
- 3.5mm jack to 3.5mm jack audio cable
- 3.5mm jack to phono (RCA) audio cable
- SCART adaptor
- User Guide (this document)

Image 2.1: Overview of the Widex TV-Dex package contents.

3. PRODUCT OVERVIEW

3.1. TV-Controller II

The TV-Controller II is the portable unit that receives audio from the TV-Base and streams it to your hearing aids. It features controls for power and volume adjustment. It is designed to be worn around the neck using the provided lanyard during use.

3.2. TV-Base

The TV-Base serves as the charging station for the TV-Controller II and connects to your television or HiFi system to transmit audio. It also connects to a power source.

Image 3.1: The TV-Controller II docked in the TV-Base.

Image 3.2: The TV-Base unit.

4. SETUP INSTRUCTIONS

Follow these steps to set up your Widex TV-Dex system:

1. **Connect the TV-Base to Power:** Plug the mini USB power supply into the TV-Base and then into a wall outlet.
2. **Connect the TV-Base to your Television/HiFi System:**
 - Identify an available audio output port on your television or HiFi system. Common options include a 3.5mm headphone jack, RCA (red and white) audio output, or a SCART port.

- Use the appropriate cable provided (3.5mm jack to 3.5mm jack, 3.5mm jack to phono, or SCART adaptor) to connect the TV-Base to your audio source.
 - Ensure the connection is secure.
3. **Charge the TV-Controller II:** Place the TV-Controller II into the TV-Base. The controller will begin charging automatically. A full charge provides approximately 10 hours of uninterrupted sound.
 4. **Pairing (if necessary):** The TV-Dex system is typically pre-paired. If pairing is required, refer to your hearing aid's user manual or consult your hearing care professional for specific instructions.

5. OPERATING INSTRUCTIONS

Once the TV-Base is connected and the TV-Controller II is charged, you can begin streaming audio:

1. **Wear the TV-Controller II:** Lift the TV-Controller II from the TV-Base and place the lanyard around your neck. Position the controller comfortably.
2. **Power On:** Press and hold the power button on the TV-Controller II until the indicator light illuminates.
3. **Activate Streaming:** Ensure your television or HiFi system is playing audio. The TV-Dex should automatically detect the audio signal and begin streaming to your compatible Widex hearing aids.
4. **Adjust Volume:** Use the volume up (+) and volume down (-) buttons on the TV-Controller II to adjust the audio level to your preference.
5. **Power Off:** When you are finished, press and hold the power button on the TV-Controller II until the indicator light turns off. Return the controller to the TV-Base for charging.

6. CHARGING

The TV-Controller II recharges when placed in the TV-Base. Ensure the TV-Base is connected to a power source. A fully charged TV-Controller II provides approximately 10 hours of continuous audio streaming. It is recommended to place the controller in the base after each use to ensure it is ready for your next listening session.

7. MAINTENANCE AND CARE

- **Cleaning:** Wipe the TV-Dex components with a soft, dry cloth. Do not use abrasive cleaners, solvents, or immerse the devices in water.
- **Storage:** Store the TV-Dex in a cool, dry place when not in use. Avoid extreme temperatures and direct sunlight.
- **Handling:** Handle the devices with care. Avoid dropping them or subjecting them to strong impacts.

8. TROUBLESHOOTING

If you encounter issues with your Widex TV-Dex, try the following:

- **No Sound:**
 - Ensure the TV-Base is securely connected to both power and the television's audio output.
 - Verify that the television or HiFi system's audio output is active and not muted.
 - Check that the TV-Controller II is powered on and charged.
 - Confirm that your Widex hearing aids are properly paired and functioning.
- **Poor Sound Quality:**
 - Ensure all cables are securely connected.

- Check for interference from other electronic devices.
 - Adjust the volume on both the TV-Controller II and your television/HiFi system.
- **TV-Controller II Not Charging:**
 - Ensure the TV-Base is connected to a working power outlet.
 - Verify that the TV-Controller II is correctly seated in the TV-Base.

If problems persist, contact your hearing care professional or Widex customer support for assistance.

9. SPECIFICATIONS

Feature	Specification
Model Name	TV-Dex
Item Weight	15.5 ounces
Package Dimensions	6.54 x 4.65 x 3.94 inches
Connectivity Technology	USB (for power)
Mounting Type	Table Mount (for TV-Base)
Battery Life (Controller)	Up to 10 hours of streaming

10. SAFETY INFORMATION


Please read and follow these safety guidelines:


- Do not attempt to open or modify the TV-Dex components. This may cause damage and void any warranty.
- Keep the device away from water and moisture.
- Use only the provided power supply and cables.
- Keep out of reach of children and pets.
- Dispose of electronic waste responsibly according to local regulations.

11. WARRANTY AND SUPPORT

For information regarding warranty coverage, technical support, or service, please refer to the warranty card included with your product or contact your authorized Widex dealer or hearing care professional. You may also visit the official Widex website for further assistance.

© 2023 Widex. All rights reserved.

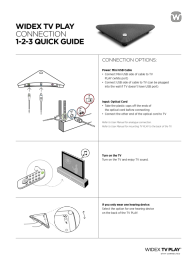	Widex MOMENT™ Hearing Aid User Instructions & Guide Comprehensive user guide for the Widex MOMENT™ MBR3D hearing aid. Learn about features, operation, customization, connectivity, and safety. Includes instructions for the MOMENT app and Bluetooth pairing.
	Widex TV PLAY Connection: 1-2-3 Quick Guide A concise guide to connecting the Widex TV PLAY device for seamless audio streaming to hearing aids. Learn about power and optical connections, and how to stream audio using your hearing aids, Widex RC-DEX remote, or the Widex EVOKE app.
	Widex COM-DEX User Guide: Wireless Audio Streaming for Hearing Aids Comprehensive user guide for the Widex COM-DEX, a wireless accessory that streams audio from mobile phones and other devices directly to Widex hearing aids, enhancing communication and listening experiences.
	Widex RC-DEX Remote Control User Instructions User instructions for the Widex RC-DEX remote control, detailing its features, operation, care, troubleshooting, and regulatory information for use with Widex tinnitus devices.
	WIDEX Sound Assist User Instructions User instructions for the WIDEX Sound Assist, a device that streams audio, acts as a remote control for hearing aids, and enhances speech understanding in noisy environments. Learn about pairing, features, maintenance, and troubleshooting.
	Widex CLEAR™ Hearing Aid User Instructions Comprehensive user guide for Widex CLEAR™ hearing aids (models 440, 330, 220). Covers setup, operation, care, troubleshooting, safety, and features like the Zen program.

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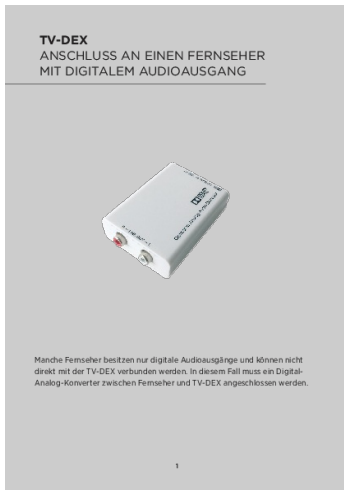
lang:tl **score:69** filesize: 2.71 M page_count: 52 document date: 2020-11-18

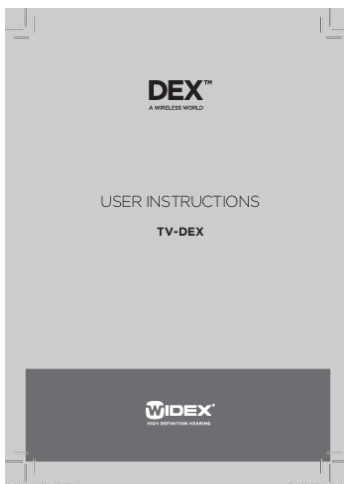[\[pdf\]](#) User Manual Quick Start Guide Instructions

9 514 0130 041 01 indd Widex TV Dex User manual 836945 archive org manualzz id

User instructions TV-Dex 9 514 0130 041 #01.indd 1 20-04-2011 15:27:30 Package Contents User manual TV-Controller Box 1: Mini USB power supply Box 2: Lanyard Box 3: a. 3.5 mm jack to 3.5 mm jack cable b. 3.5 mm jack to phono cable Box 4: a. TV-Base b. SCART adaptor For European markets only ...

lang:en **score:66** filesize: 2.7 M page_count: 48 document date: 2011-04-20

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lang:no score:60 filesize: 1.99 M page_count: 40 document date: 2014-11-11

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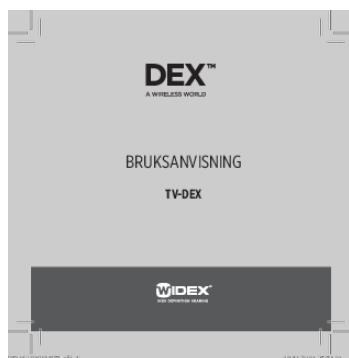
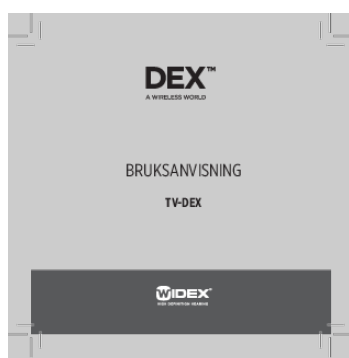
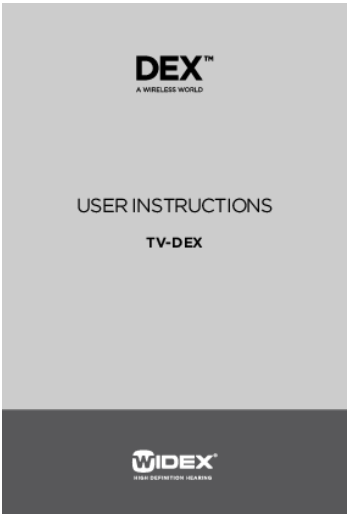
[pdf] User Manual Quick Start Guide Instructions Guide
d2zd0gbjm6m8nn cloudfront net 2018 06 WidexTV DEX user guide Widex Ihr Partner für perfekten
Hörsystem Klang USER INSTRUCTIONS User Instructions manual TV Controller Box 1 Mini USB power
supply 2 Lanyard 3 a 5 mm jack to cable b phono 4 PACKAGE CONTENTS HEARING AID AIDS BASE On
the back of base you will find following in puts Power Audio underside there is switch with 39 WIDEX A S
Nymoellevej 6 DK 3540 Lyngby Denmark widex Manual no CoselgiTV C inetpub wwwroot BIZ Axapta
documents 9 514 0130 001 03 webfiles WebFiles
USER INSTRUCTIONS **TV-DEX** 2 CONTENTS PACKAGE CONTENTS
..... 4 DESCRIPTION OF DEVICE
5 TV-CONTROLLER
lang:i-klingon score:55 filesize: 2 M page_count: 40 document date: 2014-10-27

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	iPhone	Android	Windows Phone
Smartphone	✓	✓	✓
Tablet	✓	✓	✓
Smartwatch	✓	✓	✓
Smart TV	✓	✓	✓
Smart Home	✓	✓	✓
Smart Car	✓	✓	✓
Smart City	✓	✓	✓
Smart Industry	✓	✓	✓
Smart Retail	✓	✓	✓
Smart Logistics	✓	✓	✓
Smart Energy	✓	✓	✓
Smart Agriculture	✓	✓	✓
Smart Healthcare	✓	✓	✓
Smart Education	✓	✓	✓
Smart Government	✓	✓	✓
Smart Defense	✓	✓	✓
Smart Aerospace	✓	✓	✓
Smart Maritime	✓	✓	✓
Smart Aviation	✓	✓	✓
Smart Space	✓	✓	✓

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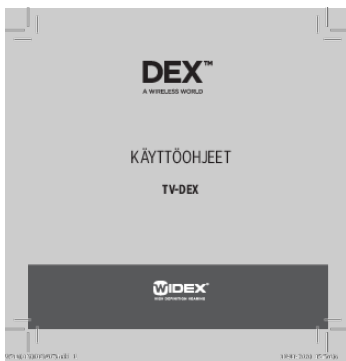
lang:no score:52 filesize: 2.5 M page_count: 52 document date: 2021-04-20

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lang:de score:52 filesize: 1.82 M page_count: 2 document date: 2010-06-03

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lang:no score:52 filesize: 1.67 M page_count: 2 document date: 2010-06-03

[\[pdf\]](#)

TV DEX WIDEX COM Printed by KLS 2016 05 9 502 1839 001 01 j9 001c# j 01v# FOR ENJOYING 1 2 3 4 5
6 Controller Base C inetpub wwwroot BIZ Axapta documents webfiles widex WebFiles |||

TV-DEX CONVENIENT STEREO SOUND 1. 2. 1. TV-C ontroller 2. TV-Base 3.
Transmission on/off 4. Volume up 5. Volume down 6. Room Off 3. 4. 5. 6.
WWW.WIDEX.COM 9 502 1839 001c #01v Printed by KLS/2016-05 9 502 1839 001 #01
TV-DEX TV-DEX FOR ENJOYING TV A WHOLE NEW WORLD OF WIRELESS
SOUND The ...

lang:en score:52 filesize: 361.13 K page_count: 2 document date: 2016-06-08



[\[pdf\]](#) Flyer

TV DEX The Hear ClinicThe new DEX™ assistive listening devices from Widex let you take advantage of a
whole world sounds They make communication between your hearing aidsWidex Flyerthehearclinic co uk
2017 02 Flyer |||

TV-DEX CONVENIENT STEREO SOUND 1. 2. 1. TV-C ontroller 2. TV-Base 3.
Transmission on/off 4. Volume up 5. Volume down 6. Room Off www.scola.info 3. 4. 5.
6. www.widex.com WWW.WIDEX.COM 9 502 1839 001c #01v Printed by FB/04-10 9
502 1839 001 #01 TV-DEX TV-DEX FOR ENJOYING TV A WHOLE...

lang:en score:51 filesize: 1.56 M page_count: 2 document date: 2010-06-03

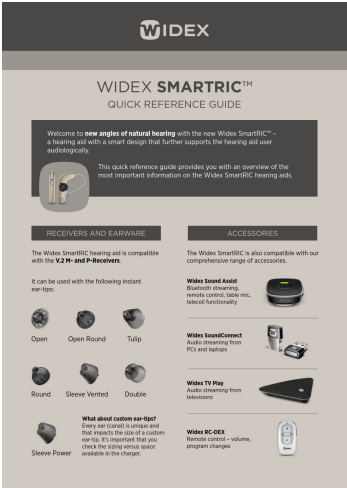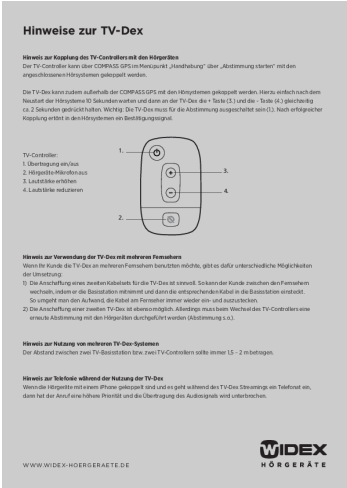[pdf] User Manual Quick Start Guide Instructions

USER INSTRUCTIONS Hearing aids from Widex User manual TV Controller II Box 1 Mini USB power supply

Manual no 9 514 0130 001 07 Issue 2016 08 WIDEX A S Nymoellevvej 6 DK 3540 Lyngby PACKAGE

CONTENTS 2 5 C inetpub wwwroot BIZ Axapta documents webfiles widex WebFiles

USER INSTRUCTIONS **TV-DEX** 2 CONTENTS PACKAGE CONTENTS
. 5 DESCRIPTION OF DEVICE

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quick reference guide digital rev 66254ff6a4a9450fa038f0795cb71624 |||

WIDEX SMARTRICTM QUICK REFERENCE GUIDE Welcome to new angles of natural
hearing with the new Widex S ... Zen/SoundRelax Audibility Extender *Incl. Widex Sound
Assist, TV-PLAY, RC-DEX, **TV-DEX**, COM-DEX, FM DEX. ...

lang:en score:49 filesize: 388.61 K page_count: 2 document date: 2024-02-27

[pdf]

9 514 0130 016 08 indd TV DEX Widex Pro widex biz axapta documents |||

TV-DEX 9 514 0130 016 #08.indd 1 18-11-2020 16:00:22 2 9 514 0130 016 #08.indd 2 18-11-2020 16:00:22 5 7 TV- II
..... 9 TV- 10 LED
lang:i-klington **score:47** filesize: 2.96 M page_count: 52 document date: 2020-11-18

[pdf] Accessories

Tomohiko Ito RC TV DEX QuickMatching widexpro media widex ja jp document library accessory dex quickmatching |||

RC/TV-DEX RC-DEX/**TV-DEX** 5cmDEX 5 5cm RC-DEX TV- DEX5 RC-DEX **TV-DEX** TVTV-5 5TV- TV-TV- TV- TV- 5 **TV-DEX** TV- TV-TV- TV- TV-TV- ...
lang:i-klington **score:47** filesize: 363.96 K page_count: 1 document date: 2017-05-09

[pdf] User Manual Quick Start Guide Instructions

USER INSTRUCTIONS Hearing aids from Widex User manual TV Controller II Box 1 Mini USB power supply widex 5qr0q3 fa f Manual no 9 514 0130 041 07 Issue 2016 08 CIB201 1214 2 5 Lanyard 3 a mm jack to cable b phono 4 Base SCART adaptor For European markets only 6 DESCRIPTION OF DEVICE The DEX comprises and PACKAGE CONTENTS C inetpub wwwroot BIZ Axapta documents webfiles WebFiles USER INSTRUCTIONS **TV-DEX** 2 CONTENTS PACKAGE CONTENTS
..... 5 DESCRIPTION OF DEVICE
.....

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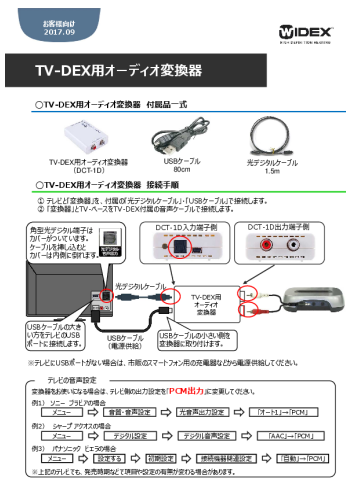
TV DEX

Widex Audio Converter for DEX201709 widex media ja jp document library

accessory dex audio converter tv dex201709 |||

2017.09 **TV-DEX TV-DEX TV-DEX** DCT-1D USB 80cm 1.5m **TV-DEX** USB TV-**TV-DEX** DCT-1D DCT-1D USB USB **TV-DEX** USB USB USB PCM 1 1PCM 2 AACPCM 3 PCM ... lang:fr score:46 filesize: 142.66 K page_count: 1 document date: 2017-09-14

[pdf] Datasheet Guide

FEATURE AND PRODUCT QUICK GUIDE Widex Evoke Data Sheet 750f 1 c7e6e2e10e985b59dd85

14c71175e65161d63fa2ccab74fc782b ssl cf2 rackcdn

Call Absolute Hearing Solutions 614-452-4280 WIDEX EVOKETM FEATURE AND PRODUCT QUICK GUIDE ACHIEVE! ... devices* Max gain MPO Bandwidth IP rating * COM-DEX, COM-DEX REMOTE MIC, CALL-DEX, UNI-DEX, RC-DEX. **TV-DEX**, FM DEX, PHONE-DEX 2 Available in some countries only. S-receiver Yes 10 100 hrs No Via DE... lang:en score:46 filesize: 3.74 M page_count: 12 document date: 2020-05-13

[pdf] Guide

C inetpub wwwroot BIZ Axapta documents 9 502 4231 001 02 Widex Unique Feature and Product Guide

webfiles widex 80 WebFiles thehearingclinic co uk 2017

WIDEX UNIQUE QUICK GUIDE AND COLOURS STYLE Receiver Open fit option Battery size Battery life Preference control Program shift Telecoil Widex CROS compatible DEX devices* Max gain MPO Bandwidth IP-rating * COM-DEX incl. COM-DEX Remote Mic, CALL-DEX, UNI-DEX, RC-DEX. **TV-DEX**, FM DEX, PHONE-DEX Avail... lang:it score:45 filesize: 1.49 M page_count: 2 document date: 2016-10-12

[pdf] Instructions Accessories

WIDEX MOMENT™ RIC BTE ORDER FORM Ear Texture q Soft Average Firm Impression Open Jaw

Closed Previous Model if applicable User Instruction language EN FR Other Widex MOMENT 0321 fillable

widexpro media widex en ca document library order form moment bte |||

SOUND LIKE NO OTHER Service Priority: q Next Day 50 Bill To: q Standard WIDEX MOMENT™ RIC BT ... CHARGER7 Anthracite Grey q q REMOTE LINK q WIDEX CHARGE n CLEAN CHARGER7 Champagne White q q **TV-DEX** Emerald Green q q **TV-DEX** D to A CONVERTER 3. Requires COM-DEX for use. 4. Requires FM rec... lang:en score:41 filesize: 553.46 K page_count: 2 document date: 2021-06-24