

## Stander 4350-CB

# Stander Wonder Rollator Plus Instruction Manual

Model: 4350-CB

Brand: Stander

## ABOUT YOUR WONDER ROLLATOR PLUS

The Stander Wonder Rollator Plus is a lightweight, compact folding mobility aid designed for seniors and adults seeking reliable support and independence. Its patented EZ Fold-N-Go technology allows for easy portability and storage, making it an ideal companion for daily activities and travel.

## IMPORTANT SAFETY GUIDELINES

- Always engage the locking brakes when sitting on the seat or when the rollator is stationary.
- Ensure all parts are securely locked in place before use.
- Do not exceed the maximum weight capacity of 250 pounds.
- Use caution when navigating uneven surfaces or inclines.
- Regularly inspect the rollator for any signs of wear or damage.
- Consult a healthcare professional to ensure this device is appropriate for your specific mobility needs.
- Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

## SETUP AND ASSEMBLY

The Stander Wonder Rollator Plus comes pre-assembled, requiring no complex setup. Follow these steps for initial preparation:

1. Unbox the rollator.
2. Adjust the handle height to fit your wrist.

3. Unlock the brakes to move.
4. Lock the brakes when ready to rest.
5. Fold the rollator by pulling the handle upwards and back.
6. Unfold the rollator by pressing down on the seat.

## NO ASSEMBLY REQUIRED



This image shows the Stander Wonder Rollator Plus in its unfolded state, highlighting its dimensions (23" x 26.5") and adjustable handle height (33"-36"). It also displays the compact folded size (10" width).

# HANDLE HEIGHT RANGE



A close-up view of the rollator's handle, demonstrating how to adjust the height from 33 to 36 inches for optimal user comfort and posture.

## OPERATING INSTRUCTIONS

### Walking with the Rollator

To walk, ensure the brakes are disengaged. Push the rollator forward, maintaining a comfortable distance. Keep your body upright and centered between the handles for stability.

### Engaging and Disengaging Brakes

The rollator features rear locking brakes for secure stopping. To engage the brakes, push the brake levers down until they click into the locked position. To disengage, pull the levers upwards. Always engage brakes when sitting or when the rollator is stationary.

### Using the Perch Seat

The 18.5 by 6-inch non-slip travel seat provides a convenient resting spot. Before sitting, ensure the rollator is on a flat, stable surface and both brakes are fully engaged. Sit facing the front of the rollator. The compact seat and back strap offer a place to sit for a quick rest while traveling or at the store.

### Folding and Unfolding

To fold the rollator, pull the strap located on the seat upwards and back. This will collapse the frame into a compact size. To unfold, simply push down on the seat until the frame locks into its open position.



**EZ FOLD-N-GO**

27" - 36"
14.5 lbs
NO
10"
POUCH & BASKET
22"
SINGLE-STEP

HEIGHT ADJUSTMENT
PRODUCT WEIGHT
ASSEMBLY REQUIRED
FOLDED WIDTH
STORAGE
SEAT HEIGHT
FOLDING STYLE



**BASIC ROLLATOR**

31" - 35"
19 lbs
YES
17.5"
BAG
20.5"
MULTI-STEP

A woman is shown easily lifting and placing the folded Stander Wonder Rollator Plus into the back seat of a car, demonstrating its portability and compact design for travel.





A close-up image illustrating the simple one-step folding mechanism of the rollator, where a user pulls a strap to collapse the device.

## FEATURES AND COMPONENTS


---




A full view of the Stander Wonder Rollator Plus in Cobalt Blue, showcasing its four wheels, handles, seat, and under-seat storage basket.

- **Compact Folding Design:** Patented EZ Fold-N-Go technology allows it to fold 3 times smaller than the average rollator to a 10 x 11 inch footprint, saving space while on the go.
- **Perch Seat:** The 18.5 by 6-inch non-slip travel seat safely supports 250 pounds.
- **Adjustable Handles:** Padded handles can be adjusted from 33 to 36 inches tall for optimal user height.
- **6-inch Wheels:** Provide stability and security for both indoor and outdoor use.
- **Rear Locking Brakes:** Allow users to control their own movement and secure the rollator when stationary.
- **Organizer Pouch:** Conveniently located for personal items.
- **Storage Basket:** Located under the seat for carrying larger items.




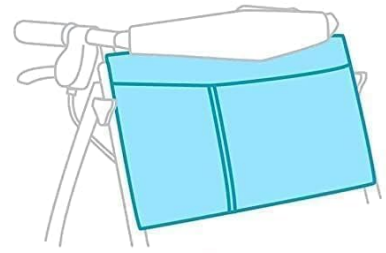

**ORGANIZER POUCH**

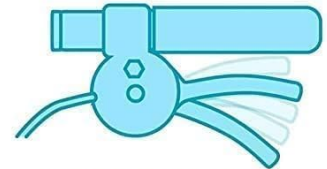

**LOCKING BREAKS**

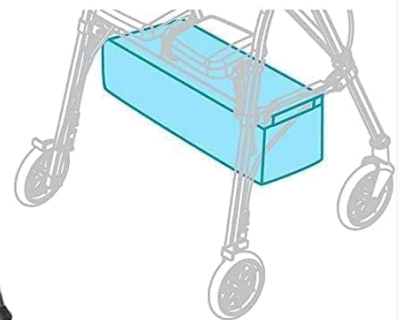

**BASKET**



**6" WHEELS**

This image displays a man using the rollator, with callouts pointing to its key features: an organizer pouch, locking brakes, a storage basket, and 6-inch wheels.



This image shows a woman comfortably seated on the rollator's perch seat, demonstrating its 250 lbs weight capacity. Another view shows her easily holding the folded rollator, emphasizing its lightweight design at 14.5 lbs.

## MAINTENANCE AND CARE

Proper maintenance ensures the longevity and safe operation of your rollator.

- **Cleaning:** Wipe down the frame and seat with a damp cloth and mild detergent. Avoid abrasive cleaners.
- **Wheel Inspection:** Periodically check the wheels for debris or wear. Clean as needed to ensure smooth rolling.
- **Brake Check:** Test the brakes regularly to ensure they engage and disengage properly. If brakes feel loose or ineffective, consult the manufacturer or a qualified technician.
- **Fasteners:** Ensure all screws, bolts, and fasteners are tight. Do not overtighten.
- **Storage:** Store the rollator in a dry place away from extreme temperatures.



# TROUBLESHOOTING

Refer to the table below for common issues and their solutions.

Problem	Possible Cause	Solution
Rollator does not roll smoothly	Debris in wheels or loose wheels	Clean wheels, check for obstructions, ensure wheels are securely attached.
Brakes are not engaging/disengaging	Brake cable too loose/tight or worn brake pads	Adjust brake cable tension (consult user manual for specific instructions) or contact support for replacement parts.
Rollator feels unstable	Handles not properly adjusted or uneven surface	Ensure handles are at the correct height and locked. Use on flat, stable surfaces.
Difficulty folding/unfolding	Mechanism obstructed or not fully engaged	Check for obstructions around the folding mechanism. Ensure the seat is fully pushed down when unfolding, or the strap is fully pulled when folding.

For issues not listed here, please refer to the official user manual or contact Stander customer support.

# PRODUCT SPECIFICATIONS

- **Product Dimensions:** 23 x 25 x 33 inches (unfolded)
- **Folded Footprint:** 10 x 11 inches
- **Item Weight:** 14.5 Pounds
- **Weight Capacity:** Up to 250 Pounds
- **Handle Height Adjustment:** 33 to 36 inches
- **Seat Dimensions:** 18.5 x 6 inches (non-slip)
- **Wheel Size:** 6-inch
- **Material:** Aluminum
- **Model Number:** 4350-CB
- **Color:** Cobalt Blue
- **Manufacturer:** Standers, Inc
- **Country of Origin:** China






# WARRANTY INFORMATION AND CUSTOMER SUPPORT

For detailed warranty information, please refer to the official user manual provided with your product or visit the [Stander Store](#).

A downloadable PDF user manual is available [here](#).

For further assistance, technical support, or to inquire about replacement parts, please contact Stander customer service.

## Related Documents - 4350-CB

	<p><a href="#">Wonder Rollator User Manual: Features, Safety, and Operation</a></p> <p>Comprehensive guide for the Stander Wonder Rollator, covering features, safety precautions, warranty, and operating instructions for this compact folding mobility aid.</p>
	<p><a href="#">Stander Bed Cane Assembly and Safety Instructions</a></p> <p>Comprehensive guide to the safe assembly and installation of the Stander Bed Cane, a mobility aid designed to assist individuals in standing from bed. Includes important safety warnings, warranty information, and step-by-step instructions.</p>
	<p><a href="#">Stander 7630 Wonder Tray Premium Owners Manual and Assembly Instructions</a></p> <p>Comprehensive owners manual and assembly instructions for the Stander 7630 Wonder Tray Premium, a mobility aid designed to assist users when standing from chairs or couches. Includes safety warnings, specifications, and step-by-step assembly guidance.</p>
	<p><a href="#">Stander Security Pole 1150: Safety Warnings, Specifications, and Assembly Instructions</a></p> <p>Comprehensive guide to the Stander Security Pole (Model 1150), including essential safety warnings, product specifications, limited warranty details, and step-by-step assembly instructions for safe and effective installation.</p>
	<p><a href="#">Stander 8140 Wonder Pole Lite Owner's Manual and Assembly Instructions</a></p> <p>Detailed owner's manual and assembly guide for the Stander 8140 Wonder Pole Lite. Includes comprehensive safety warnings, product specifications, and step-by-step instructions for installation. Designed as a mobility and support aid for safe standing and sitting.</p>

[illegible]

## Stander 8900 Series Bed Rail: Assembly, Installation, and Safety Guide

Comprehensive guide for Stander 8900 Series bed rails (models 8930, 8940), covering installation, safety warnings, warranty, and accessories. Ensure safe and proper use of your Stander bed rail.