

Activeshop Basic 202

Activeshop Basic 202 Cosmetic Treatment Chair User Manual

1. PRODUCT OVERVIEW

The Activeshop Basic 202 Cosmetic Treatment Chair is a versatile and adjustable chair designed for various professional applications, including cosmetic procedures, makeup application, and massage therapy. This chair features a comfortable design with a face cut-out in the headrest, practical storage drawers, and durable, easy-to-clean ecological leather upholstery.

Key components include:

- Seat
- Armrests (2 pieces)
- Storage Drawers (2 pieces)
- Bracket for disposable covers
- Guide rails for storage drawers (2 pieces)
- Steel frame components (2 pieces)
- Stop buffers for the steel frame



Figure 1: The Activeshop Basic 202 Cosmetic Treatment Chair in a reclined position, highlighting its adjustable features and integrated storage.

2. SAFETY INFORMATION

Please read all safety instructions carefully before assembly and use.

- **Weight Limit:** Do not exceed a maximum weight of 150 kg (approximately 330 lbs) on the chair. Exceeding this limit may cause damage to the chair or injury to the user.
- **Sharp Corners:** Exercise caution during assembly and use, as some components may have sharp corners.
- **Assembly Required:** This product requires assembly. Ensure all parts are correctly installed and tightened before use.
- **Stable Surface:** Always place the chair on a flat, stable, and level surface to prevent tipping.
- **Children and Pets:** Keep children and pets away from the chair during assembly and ensure they do not play on or around the chair unsupervised.

3. SETUP AND ASSEMBLY

The Activeshop Basic 202 Cosmetic Treatment Chair requires assembly. Please follow these general steps. A detailed assembly manual is typically included with the product.

1. **Unpack Components:** Carefully remove all parts from the packaging. Verify that all listed components (seat, armrests, drawers, frame parts, etc.) are present.
2. **Assemble the Steel Frame:** Connect the steel frame components according to the provided diagrams. Ensure all bolts and fasteners are securely tightened.
3. **Attach Guide Rails and Drawers:** Install the guide rails for the storage drawers onto the frame. Slide the two storage drawers into their respective compartments.
4. **Mount the Seat and Backrest:** Secure the main seat and backrest sections to the assembled frame.

Pay attention to the orientation of the headrest with the face cut-out.

5. **Install Armrests:** Attach the two armrests to the sides of the chair.
6. **Attach Bracket:** Secure the bracket for disposable covers to the backrest.
7. **Final Check:** Before first use, ensure all connections are tight and the chair is stable. Test all adjustable mechanisms.

If you encounter any difficulties during assembly, refer to the specific assembly instructions provided in your product packaging or contact customer support.

4. OPERATING INSTRUCTIONS

The Activeshop Basic 202 chair is designed for ease of use and adjustability to suit various treatment needs.

4.1 Adjusting the Backrest and Legrest

The backrest and legrest sections of the chair are manually adjustable. Locate the adjustment levers or knobs, typically found on the sides or underneath the chair sections. Lift or lower the sections to the desired angle and secure them in place using the locking mechanism.



Figure 2: The chair adjusted to an upright position, ideal for various cosmetic procedures.



Figure 3: The chair fully flattened, suitable for massage or full-body treatments.

4.2 Using the Face Cut-out

For treatments requiring the client to lie face down, the headrest features a face cut-out. When not in use, this opening can be covered with the provided cushion or plug for client comfort.

4.3 Utilizing Storage Drawers

The two integrated storage drawers provide convenient access to tools, supplies, or personal items during treatments. Pull the drawers out gently to open and push them in to close securely.

4.4 Using the Disposable Cover Bracket

The bracket located on the backrest is designed to hold rolls of disposable covers, napkins, or towels, ensuring hygiene and easy access during procedures.



Figure 4: Rear view of the chair, illustrating the disposable cover bracket and backrest adjustment points.

5. MAINTENANCE AND CARE

Proper maintenance ensures the longevity and hygienic condition of your Activeshop Basic 202 Cosmetic

Treatment Chair.

- **Cleaning Upholstery:** The ecological leather upholstery is easy to clean. Wipe surfaces regularly with a soft, damp cloth and mild, non-abrasive soap. Avoid harsh chemicals, solvents, or abrasive cleaners, as these can damage the material.
- **Frame Cleaning:** Clean the steel frame with a damp cloth. Dry thoroughly to prevent rust.
- **Moving Parts:** Periodically check all moving parts and adjustment mechanisms. If any squeaking or stiffness occurs, apply a small amount of silicone-based lubricant to the pivot points.
- **Fastener Check:** Regularly inspect all bolts, screws, and fasteners to ensure they remain tight. Retighten as necessary to maintain stability and safety.
- **Storage:** When not in use for extended periods, store the chair in a clean, dry environment away from direct sunlight and extreme temperatures.

6. TROUBLESHOOTING

This section addresses common issues you might encounter with your cosmetic treatment chair.

- **Chair is Unstable/Wobbly:**
 - Ensure the chair is placed on a flat, level surface.
 - Check all assembly bolts and fasteners. Tighten any loose connections.
 - Verify that the floor glides/feet are properly adjusted and making full contact with the floor.
- **Adjustment Mechanisms are Stiff:**
 - Apply a small amount of silicone-based lubricant to the pivot points and moving parts of the backrest and legrest mechanisms.
 - Ensure no obstructions are preventing movement.
- **Drawers are Sticking:**
 - Check the drawer guide rails for any debris or obstructions.
 - Ensure the drawers are correctly aligned on their rails.
- **Upholstery Damage:**
 - For minor scuffs, clean gently with a mild cleaner.
 - For tears or significant damage, contact Activeshop customer support for repair options or replacement parts.

If you experience issues not covered here, please contact Activeshop customer support for assistance.

7. SPECIFICATIONS

Feature	Detail
Model	Basic 202
Brand	Activeshop
Dimensions (L x W x H)	180 x 60 x 75 cm (70.9 x 23.6 x 29.5 inches)
Weight	39 kg (86 lbs)

Maximum Load Capacity	150 kg (330 lbs)
Color	White
Material	Ecological Leather (PU Leather), Steel Frame
Assembly Required	Yes
Features	Adjustable backrest and legrest, face cut-out in headrest, 2 storage drawers, disposable cover bracket

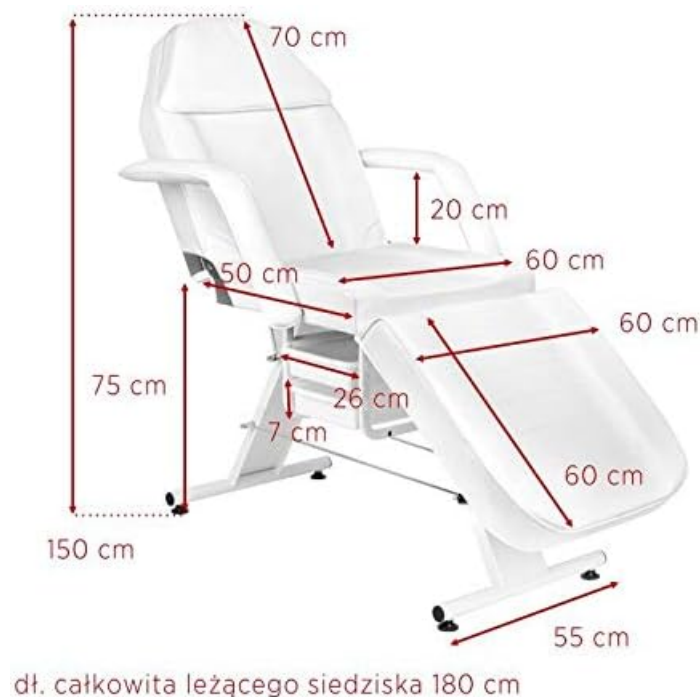


Figure 5: Dimensional diagram of the Activeshop Basic 202 Cosmetic Treatment Chair.

8. WARRANTY AND SUPPORT

For warranty information, please refer to the documentation included with your purchase or visit the official Activeshop website. If you require technical support, replacement parts, or have questions regarding your Activeshop Basic 202 Cosmetic Treatment Chair, please contact Activeshop customer service.

Contact information can typically be found on the manufacturer's website or in your product packaging.