

Signia Streamline TV 03

Signia Streamline TV 03 Instruction Manual

Model: Streamline TV 03 | Brand: Signia

INTRODUCTION

The Signia Streamline TV 03 is a dedicated accessory designed to stream audio directly from your television to compatible Signia NX and BT 13 hearing aids. This device enables a personalized listening experience, allowing you to enjoy TV programs with enhanced clarity and individual volume control, without affecting the listening experience of others in the room.



Image 1: The Signia Streamline TV 03 device. This compact, rectangular unit features a black base and a silver top panel, with two small light indicators on its front edge, signifying power and connection status.

SETUP

Follow these steps to set up your Signia Streamline TV 03 for optimal performance:

- 1. Unpack and Identify Components:** Ensure all necessary cables are present. The Streamline TV 03 typically includes various audio cables to accommodate different television setups.
- 2. Connect to Television Audio Output:**
 - Locate the audio output ports on your television. Most modern TVs feature several options.
 - The **Optical audio out** port is generally recommended for the best sound quality. If available, use the optical cable provided to connect the Streamline TV 03 to your TV's optical output.
 - Alternatively, if your TV has RCA audio out jacks (red and white), use the appropriate RCA cable to connect the Streamline TV 03.
 - Ensure the cable is securely connected to both the TV and the Streamline TV 03 unit.
- 3. Power Connection:** Plug the power adapter into the Streamline TV 03 unit, then connect it to a power outlet. The unit will power on automatically.

4. Pairing with Hearing Aids:

- Ensure your Signia NX or BT 13 hearing aids are turned on and within range (up to 30 feet) of the Streamline TV 03.
- The pairing process is typically automatic or initiated through the Signia app. Refer to your hearing aid's specific instructions for pairing with accessories.
- Your hearing care professional may need to add a specific "TV program" to your hearing aids for optimal functionality, especially if you intend to use the device without the Signia app.

OPERATING INSTRUCTIONS

Once set up, the Streamline TV 03 provides a seamless audio experience:

- **Automatic Connection:** After initial pairing, your hearing aids should automatically connect to the Streamline TV 03 when it is powered on and within range.
- **Individual Volume Control:** You can adjust the volume of the streamed audio directly through your hearing aids or via the Signia app on your smartphone. This allows for personalized volume settings without altering the TV's main volume.
- **Enhanced Clarity:** The Streamline TV 03 bypasses ambient room noise, delivering clear, direct audio to your hearing aids, which can significantly improve speech understanding.
- **Signia App Integration:** If you use the Signia app, controls for the Streamline TV 03 may appear within the app interface, offering convenient management of streaming functions and volume.
- **Range:** The device can stream audio to your hearing aids from up to 30 feet away.

MAINTENANCE

The Signia Streamline TV 03 requires minimal maintenance to ensure continued performance:

- **Cleaning:** Gently wipe the exterior of the unit with a soft, dry cloth to remove dust. Avoid using liquid cleaners or abrasive materials.
- **Storage:** Store the device in a cool, dry place when not in use.
- **Cable Management:** Ensure all connected cables are not kinked or strained to prevent damage.

TROUBLESHOOTING

If you encounter issues with your Streamline TV 03, consider the following:

- **No Sound:**
 - Verify that the Streamline TV 03 is powered on and its indicator lights are active.
 - Check all cable connections between the TV and the Streamline TV 03. Ensure they are secure and correctly inserted into the audio output ports of the TV.
 - Confirm that your TV's audio output settings are correctly configured to send audio through the selected port (e.g., Optical, RCA).
 - Ensure your hearing aids are paired and connected to the Streamline TV 03.
 - Adjust the volume on your hearing aids or through the Signia app.
- **Intermittent Connection/Poor Sound Quality:**
 - Ensure your hearing aids are within the 30-foot operating range of the Streamline TV 03.

- Minimize obstructions between the Streamline TV 03 and your hearing aids.
- Check for potential interference from other electronic devices.
- Restart both the Streamline TV 03 and your hearing aids.
- **Difficulty Pairing:**
 - Ensure your hearing aids are compatible with the Streamline TV 03 (Signia NX & BT 13 models).
 - Follow the specific pairing instructions for your hearing aid model.
 - If issues persist, contact your hearing care professional for assistance, as they may need to adjust your hearing aid's programming.
- **Unable to Pause Streaming:** Some users have noted difficulty pausing streaming when conversing. This functionality may depend on your specific hearing aid model and app settings. Consult your hearing care professional for personalized advice.

SPECIFICATIONS



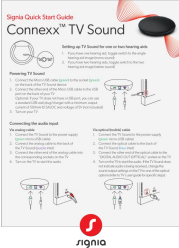

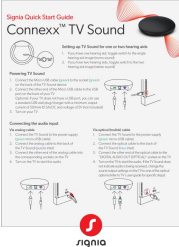

Brand Name	Signia
Model Name	Streamline TV 03
Item Model Number	03
ASIN	B07N35HMCD
Color	Silver & Black
Item Weight	1.15 pounds
Package Dimensions	6.4 x 4.6 x 3.7 inches
Connectivity Technology	Bluetooth (for hearing aids), HDMI, USB (for TV connections, though optical/RCA are primary audio inputs)
Recommended Uses	Home Viewing
Voltage	1.1E+2 Volts (AC)
Included Components	Streamline TV 03 unit, various audio cables, power adapter. (Note: "Stand, TV" from original data refers to the TV itself, not included with Streamline TV 03)

WARRANTY AND SUPPORT

For specific warranty information, please refer to the documentation included with your purchase or contact Signia customer support directly. Support resources, including FAQs and contact details, are typically available on the official Signia website.

It is recommended to consult with your hearing care professional for any programming adjustments or advanced troubleshooting related to your hearing aids and their interaction with the Streamline TV 03.

Related Documents - Streamline TV 03

	<p>Signia Pocket Guide: Hearing Aid Solutions for Government Services & Veterans</p> <p>Comprehensive guide to Signia hearing aids for Government Services, featuring the Signia Xperience platform, rechargeable models, and product details for veterans. Includes fitting and support information.</p>
	<p>Signia Pocket Guide for Government Services</p> <p>Signia's Pocket Guide details advanced hearing aid solutions, including the Augmented Xperience and Xperience platforms, for veterans and military personnel. It covers product features, fitting information, and support resources.</p>
	<p>Signia Connexx TV Sound: Quick Start Guide for Seamless Audio Streaming</p> <p>Comprehensive guide to setting up, powering, and pairing your Signia Connexx TV Sound device with hearing aids. Learn about analog/optical connections, Bluetooth pairing, and app control for an enhanced TV audio experience.</p>
	<p>Signia Motion C&G SP IX Hearing Aid Quick Start Guide: Setup, Charging, and Care</p> <p>A comprehensive quick start guide for Signia Motion C&G SP IX hearing aids, covering component identification, turning on/off, insertion, removal, charging, maintenance, user controls, and support resources.</p>
	<p>Signia Connexx TV Sound Quick Start Guide</p> <p>A quick start guide for the Signia Connexx TV Sound, providing instructions on setup, powering, audio connection via analog or optical cables, pairing hearing aids, starting audio streaming using hearing aid controls or the Signia app, connecting sound bars, switching devices, and removing pairings.</p>
	<p>Signia Portable Charger RIC User Guide - Charging and Maintenance</p> <p>Official user guide for the Signia Portable Charger RIC. Learn how to charge your hearing instruments, understand components, perform maintenance, and troubleshoot issues. Includes safety information.</p>