

## Vive Mobility 3 Wheel Walker

# Vive Mobility 3 Wheel Walker Instruction Manual

Model: 3 Wheel Walker

Brand: Vive Mobility

## 1. INTRODUCTION

This manual provides essential information for the safe and effective use of your Vive Mobility 3 Wheel Walker. Please read all instructions carefully before assembly and operation. Retain this manual for future reference.

The Vive Mobility 3 Wheel Walker is a lightweight, foldable, three-wheeled rollator designed to provide support and enhance mobility for individuals requiring assistance with walking. It features a stable design, adjustable handles, and a convenient storage bag.

## 2. SAFETY INFORMATION

Always prioritize safety when using the 3 Wheel Walker. Failure to follow safety guidelines may result in injury or damage to the product.

- Ensure all parts are securely assembled and locked before use.
- Regularly check the brakes for proper function.
- Do not exceed the maximum weight capacity of 250 pounds.
- Avoid leaning excessively forward or sideways on the walker, as this may compromise stability.
- Use caution on uneven surfaces. The walker is designed for flat surfaces such as cement or brick. It may not be suitable for rough terrain like gravel or grass.
- Do not use the walker as a seat. This model does not include a seating function.

## 3. SETUP AND ASSEMBLY

The Vive Mobility 3 Wheel Walker comes largely pre-assembled, requiring minimal setup.

1. **Unfolding the Walker:** Carefully unfold the walker until it locks into its open position. Ensure all frame components are fully extended and secure.
2. **Attaching Handles:** Insert the handle tubes into the main frame. Ensure the push pins engage securely into the adjustment holes.

3. **Adjusting Handle Height:** Press the push pins on each handle and slide the handles up or down to your desired height. Release the pins to lock the handles in place. Ensure both handles are at the same height. Your hands should be at a natural angle for comfort and easy access to the brakes.
4. **Attaching the Accessory Bag:** Secure the included stylish grey accessory bag to the designated attachment points on the walker frame.

## Sturdy & Long-Lasting

Supports up to 250 pounds



**Figure 3.1:** The Vive Mobility 3 Wheel Walker in use, highlighting its pre-assembled frame and adjustable handles for user comfort.

# Lightweight & Foldable

Stays out of the way until you need it

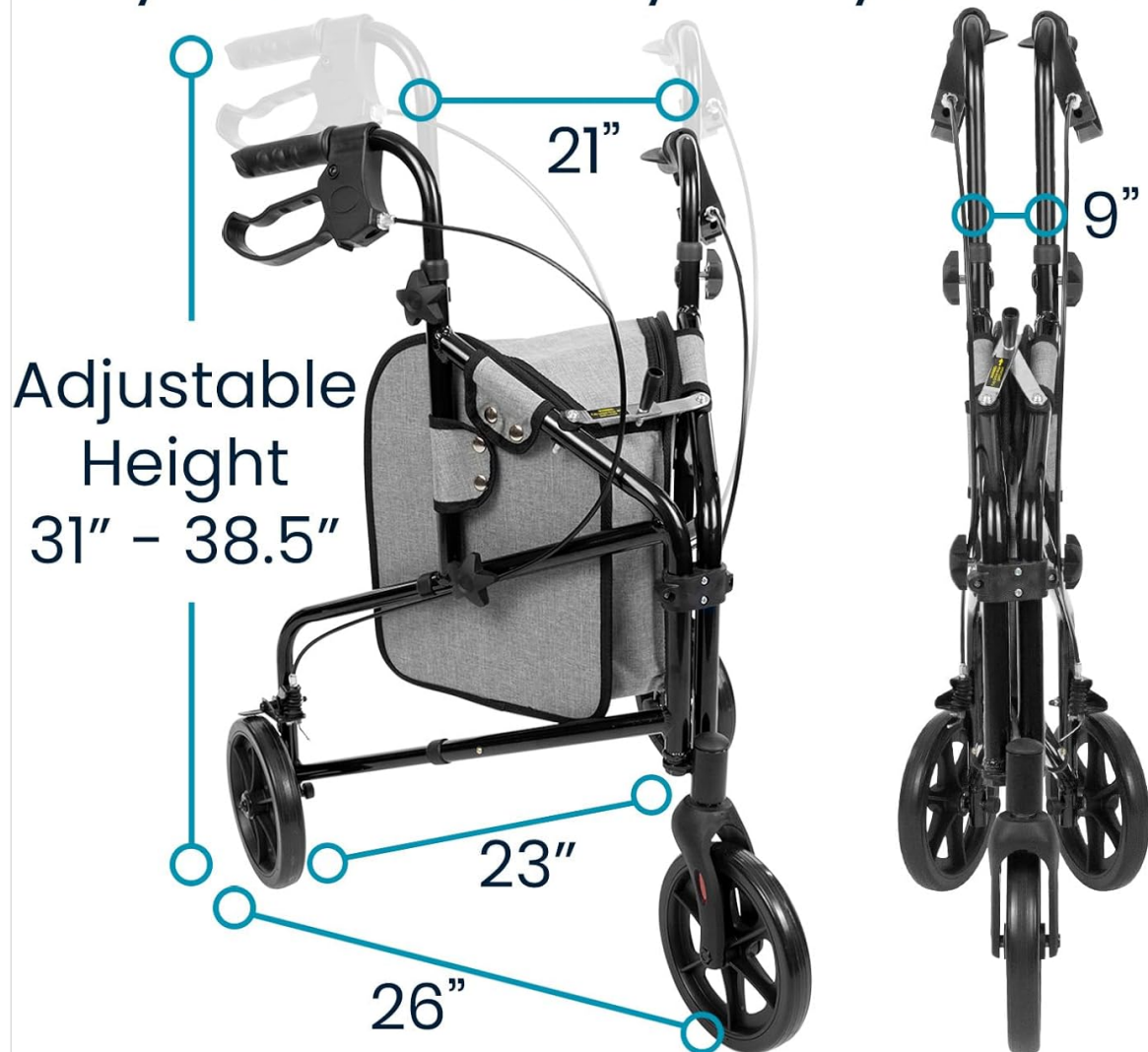


Figure 3.2: The walker's dimensions and its compact folded state, illustrating its lightweight and foldable design.

## 4. OPERATING INSTRUCTIONS

Follow these steps for safe and effective operation of your 3 Wheel Walker.

1. **Walking with the Walker:** Position yourself behind the walker, holding both handles firmly. Push the walker forward a short distance, then step forward, maintaining your balance. Repeat this motion, keeping the walker slightly ahead of you.
2. **Using the Brakes:** To slow down or stop, gently squeeze the brake levers located on the handles. To engage the parking brake, press the lever on the back of the walker downwards until it locks. Ensure the parking brake is engaged when standing still or resting.
3. **Maneuvering:** The three-wheel design allows for enhanced maneuverability in small spaces and crowded areas. Use the handles to guide the front wheel in the desired direction.
4. **Using the Accessory Bag:** The removable grey accessory bag provides convenient storage for personal items. Place items securely inside the bag before moving.
5. **Folding for Storage/Transport:** To fold the walker, disengage any parking brakes. Follow the



folding mechanism instructions (typically involves pulling a strap or lever) to collapse the frame. Once folded, the walker is compact for storage or transportation.



**Figure 4.1:** The Vive Mobility 3 Wheel Walker, illustrating its design and 250lb weight capacity.

# Easy to Transport & Maneuver



- Weighs only 11 lbs
- 360 degree swivel wheels

**Figure 4.2:** The walker being easily transported, demonstrating its lightweight and foldable design for convenient storage in a vehicle.

# Travel Bag Included

So you never have to leave your essentials behind



**Figure 4.3:** The included travel bag, showing its capacity and how it attaches to the walker for carrying essentials.

## 5. MAINTENANCE

Regular maintenance ensures the longevity and safe operation of your walker.

- **Cleaning:** Wipe the frame with a damp cloth and mild soap as needed. The accessory bag is removable and easy to clean. Avoid abrasive cleaners.
- **Brake Inspection:** Periodically check the brake cables and levers for wear or damage. Ensure brakes engage smoothly and effectively.
- **Wheel Inspection:** Check wheels for debris and ensure they rotate freely. Inspect tires for wear.
- **Fastener Check:** Regularly inspect all nuts, bolts, and fasteners to ensure they are tight and secure.
- **Storage:** Store the walker in a dry, clean place when not in use.

## 6. TROUBLESHOOTING

This section addresses common issues you might encounter with your 3 Wheel Walker.

- **Walker feels unstable:** Ensure all frame components are fully extended and locked. Verify that the handle height is adjusted correctly and evenly on both sides. Do not exceed the weight capacity.
- **Brakes are not engaging properly:** Check the brake levers and cables for any obstructions or damage. If the issue persists, consult a qualified technician.
- **Difficulty folding/unfolding:** Ensure no items are obstructing the folding mechanism. Check for any loose or missing fasteners that might affect the frame's movement.
- **Wheels are not rolling smoothly:** Inspect the wheels for debris (e.g., hair, lint) that might be caught in the axles. Clean as necessary.

## 7. SPECIFICATIONS

Feature	Specification
Model	3 Wheel Walker
Brand	Vive Mobility
Material	Aluminum
Color	Black
Product Dimensions (Unfolded)	26"D x 21"W x 31"H (Adjustable height up to 38.5")
Item Weight	11 Pounds
Weight Capacity	250 Pounds
Wheel Type	360-degree swivel wheels (front)

## 8. WARRANTY INFORMATION

The Vive Mobility 3 Wheel Walker is backed by a **LIFETIME GUARANTEE**. This guarantee reflects our commitment to the durability and longevity of our product. For specific details regarding the terms and conditions of the lifetime guarantee, please refer to the official Vive Mobility website or contact customer support.

## 9. SUPPORT

For further assistance, questions, or to report any issues with your Vive Mobility 3 Wheel Walker, please contact Vive Mobility customer support. Contact information can typically be found on the product packaging or the official Vive Mobility website.

