



Manuals.plus /

- › Bemis /
- › Bemis Assurance 3" Raised Elongated Toilet Seat Instruction Manual

## Bemis 7YE85320TSS

# Bemis Assurance 3" Raised Elongated Toilet Seat Instruction Manual

Model: 7YE85320TSS

## INTRODUCTION

---

The Bemis Assurance 3-inch Raised Elongated Toilet Seat is designed to provide enhanced comfort and safety for individuals with limited mobility. This manual provides comprehensive instructions for the proper installation, use, and maintenance of your new toilet seat, ensuring a secure and hygienic experience.

Key features include a 3-inch elevation for easier sitting and standing, a secure Snap 2 Secure™ hinge system to prevent wobbling, and an innovative Clean Shield™ guard for improved hygiene. The seat is constructed from durable plastic and is tested to support up to 1,000 lbs.



**BEST OVERALL**

RAISED TOILET SEAT

BY

**bob vila**™

### **3" Raised Seat**

For easier sitting & standing

### **Securely Bolts to Toilet**

Hinges will never loosen

### **Clean Shield Guard**

Keeps everything in the bowl

### **Stable Grip Bumpers**

Ensure seat won't shift or move

**MOST RECOMMENDED**

By Occupational Therapists

Image: Key features of the Bemis Assurance Raised Toilet Seat, including its 3-inch elevation, secure bolting mechanism, Clean Shield guard, stable grip bumpers, and occupational therapist recommendation.


## **SETUP & INSTALLATION**

---

Proper installation is crucial for the safety and stability of your Bemis Assurance toilet seat. Please follow these steps carefully. An installation tool is included with your product.

### **1. Verify Toilet Type**

Ensure your toilet is an elongated model. This seat is designed to fit all elongated toilets, including popular brands like Kohler, American Standard, and TOTO. The product dimensions are 19.13" L x 16.25" W x 4.25" H.

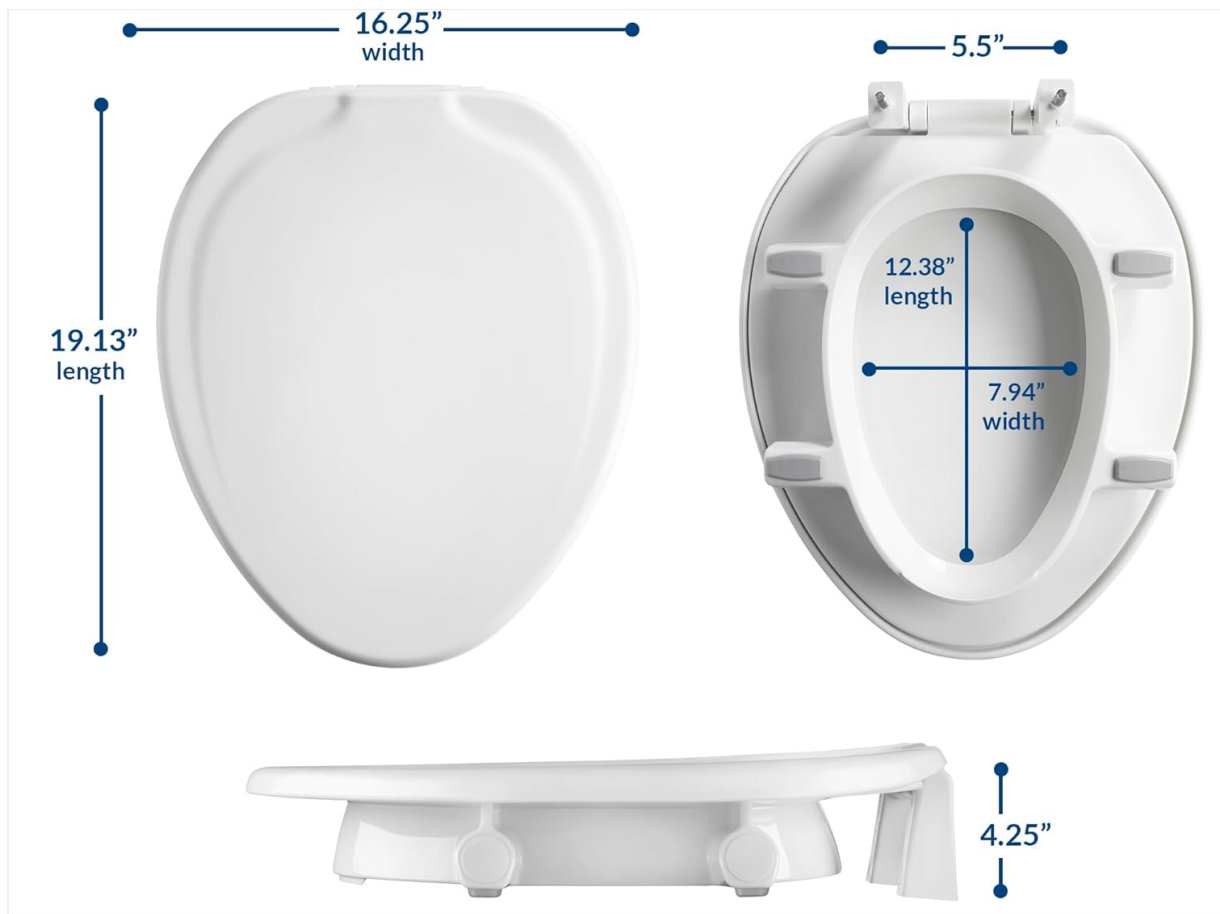

# Elongated Product Dimensions

Image: Elongated product dimensions for the Bemis Assurance Raised Toilet Seat, confirming compatibility with various elongated toilet models.

## 2. Remove Existing Toilet Seat

Carefully remove your current toilet seat by unscrewing the bolts that secure it to the toilet bowl.

## 3. Place New Seat

Position the Bemis Assurance raised toilet seat onto the toilet bowl, aligning the mounting bolts with the holes in the bowl.

## 4. Secure with Snap 2 Secure™ System

From underneath the toilet bowl, lightly finger-tighten both Snap 2 Secure™ nuts onto the mounting bolts. While holding the seat in place with one hand, use the provided installation tool to tighten the nuts until both installation heads snap off. This indicates a secure and proper installation.

## 5. Optional Support Arms/Caps

If using optional support arms (sold separately), slide them into the designated holes until the button lock makes contact and clicks. If not using support arms, insert the provided caps into the holes to cover them.

## Installation Video

Your browser does not support the video tag.

Video: Detailed step-by-step instructions for installing the Bemis Independence Clean Shield Toilet Seat, demonstrating the Snap 2 Secure system and optional support arm attachment.

## OPERATING INSTRUCTIONS

---

The Bemis Assurance Raised Toilet Seat is designed for straightforward use. Once securely installed, it provides a stable and elevated seating surface.

### 1. Sitting and Standing

The 3-inch elevated height reduces the distance needed to sit down and stand up, minimizing strain on joints and muscles. Use the optional support arms (if installed) for additional leverage and stability.

### 2. Secure Fit

The Snap 2 Secure™ system ensures the seat remains firmly attached to the toilet bowl, preventing wobbling or shifting during use. This secure fit is critical for user safety.

Your browser does not support the video tag.

Video: Demonstration of the secure attachment of the Bemis Independence Clean Shield Toilet Seat, highlighting its stability compared to other elevated seats.

## MAINTENANCE

---

Maintaining your Bemis Assurance toilet seat is simple and helps ensure its longevity and hygiene.

### 1. Cleaning

The smooth plastic surface and Clean Shield™ design make cleaning easy. To preserve the seat's glossy finish, clean with mild soap and water, anti-bacterial wipes, or anti-bacterial hand soap. Avoid using abrasive cleaners, bleach, or harsh chemicals, as these can damage the finish.



## CLEANING YOUR TOILET SEAT

To preserve your seat's glossy finish, we recommend cleaning with either:

- Anti-bacterial wipes
- Dish soap
- Anti-bacterial hand soap

**Do Not Use  
Bleach or Harsh Chemicals**

Image: Recommended cleaning methods for the Bemis Assurance Toilet Seat, emphasizing the use of gentle cleaners and avoiding harsh chemicals.

## 2. Inspection

Periodically check the seat for any signs of wear or damage. Ensure the mounting bolts remain tight and the seat does not wobble. If any issues are found, refer to the troubleshooting section or contact customer support.

## TROUBLESHOOTING

---

If you encounter any issues with your Bemis Assurance toilet seat, refer to the common problems and solutions below.

### Seat Wobbles or Shifts

- **Solution:** Re-tighten the Snap 2 Secure™ nuts using the provided installation tool until the heads snap off. Ensure the seat is properly aligned before tightening.

### Gap Between Seat and Toilet Bowl

- **Solution:** Ensure the seat is correctly positioned and the mounting bolts are fully engaged. Some minor gaps may occur depending on the toilet bowl's specific curvature, but the secure fastening should prevent movement.

## Lid Slams Shut

- **Note:** This model does not feature a slow-close mechanism. The lid will close at a normal speed.

Your browser does not support the video tag.

Video: Testing the strength and durability of the Bemis Independence Clean Shield Toilet Seat, demonstrating its ability to withstand significant loads and movement.

## SPECIFICATIONS

<b>Brand</b>	Bemis
<b>Model Number</b>	7YE85320TSS
<b>Color</b>	White
<b>Material</b>	Plastic
<b>Style</b>	Raised Seat
<b>Shape</b>	Elongated (Oval)
<b>Product Dimensions (L x W x H)</b>	19.13" x 16.25" x 4.25"
<b>Item Weight</b>	8.56 Pounds
<b>Weight Capacity</b>	Up to 1,000 lbs
<b>Included Components</b>	(1) Raised Toilet Seat & Hardware
<b>UPC</b>	073088158246



Image: Guide on how to measure your toilet seat to ensure correct sizing for Bemis products.

## WARRANTY & SUPPORT

Bemis is committed to providing high-quality products and excellent customer service.

## Manufacturer's Warranty

Your Bemis Assurance 3" Raised Toilet Seat is covered by a manufacturer's warranty for 1 year from the date of purchase. This warranty applies to customers in the US and Canada and covers defects in materials and workmanship under normal use.

## Customer Support

For any questions, assistance, or warranty claims, please contact Bemis Customer Service:

- **Phone:** 888.722.6488
- **Hours:** Monday - Friday, 8 AM - 5 PM CST



Image: Bemis Customer Service and Support information, including contact number, hours, and warranty details.