

## Nuby 50550

# Nuby Comfort All Silicone Breast Pump Instruction Manual

Model: 50550

## INTRODUCTION

The Nuby Comfort All Silicone Breast Pump is designed for convenient and hands-free milk collection. Its natural suction capability allows for efficient milk expression while breastfeeding or using another pump. This lightweight and portable pump is made from 100% soft silicone, ensuring comfort against the skin and quick, easy expression. It features a suction cup base to prevent spills and a sealing plug for secure milk storage.

Capacity: 4 oz / 120 ml.

## IMPORTANT SAFETY INFORMATION

- Always inspect the product before each use. Discard at the first sign of damage or weakness.
- Keep all components out of reach of children when not in use.
- This product is intended for personal use only. Do not share with others.
- Do not use abrasive cleaners or solvents.
- Ensure all parts are clean and dry before assembly and storage.
- If you experience pain or discomfort, discontinue use and consult a healthcare professional.

## COMPONENTS

The Nuby Comfort All Silicone Breast Pump Nursing Set includes:

- 1 x Nuby Comfort All Silicone Breast Pump (4 oz / 120 ml)
- 1 x Sealing Plug



Image: The Nuby Comfort All Silicone Breast Pump shown with its included sealing plug inserted at the top, ready for milk storage.

## SETUP AND FIRST USE

1. **Disassembly:** The Nuby Comfort All Silicone Breast Pump is a single-piece unit and does not require disassembly for cleaning.
2. **Cleaning:** Before first use, and after every subsequent use, wash the pump and sealing plug thoroughly in warm, soapy water. Rinse well with clean water.
3. **Sterilization:** For initial sterilization, you may boil the silicone pump and plug in water for 5 minutes. Ensure there is enough water to prevent contact with the bottom or sides of the pot.
4. **Drying:** Allow all components to air dry completely or use a clean cloth.

## OPERATING INSTRUCTIONS

The Nuby Comfort All Silicone Breast Pump uses natural suction to collect milk. It is ideal for catching let-down milk while breastfeeding on the opposite breast or for gentle expression.

1. **Preparation:** Ensure the pump is clean and dry. Position yourself comfortably.
2. **Application:** Center the pump opening over your nipple. Gently squeeze the base of the pump to create a vacuum, then release to allow the pump to attach to your breast. Adjust the position until you feel a comfortable, firm suction.
3. **Milk Collection:** Once attached, the pump will draw out milk using continuous suction. Milk will flow naturally into the collection cup. No manual pumping is necessary.
4. **Hands-Free Use:** The strong suction allows for hands-free operation, making it convenient to use while nursing your baby on the other breast.
5. **Releasing Suction:** To remove the pump, gently squeeze the base again or break the seal by inserting a finger between your breast and the pump flange.
6. **Storing Milk:** Once milk is collected, remove the pump from your breast. Insert the sealing plug firmly into the opening of the pump to prevent spills. You can then transfer the milk to a suitable storage container or bottle.

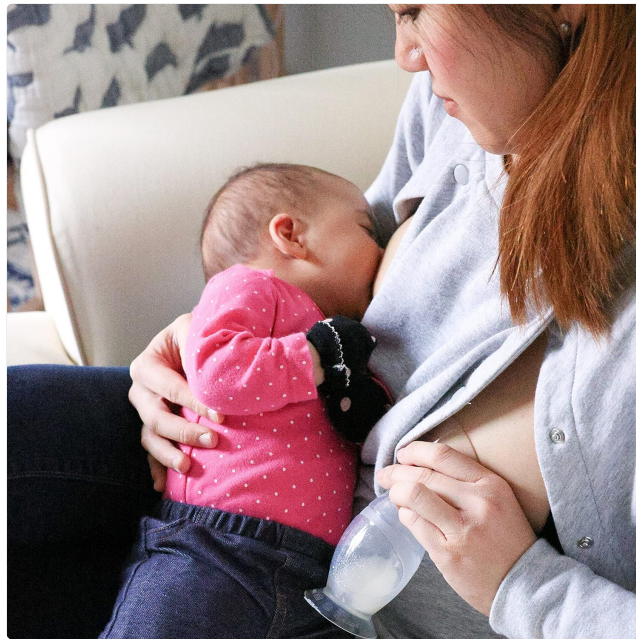


Image: A mother is shown breastfeeding her baby while simultaneously using the Nuby Comfort All Silicone Breast Pump on her other breast to collect milk hands-free.

## MAINTENANCE AND CLEANING

Proper cleaning and maintenance are crucial for hygiene and extending the life of your breast pump.

- **After Each Use:** Wash the pump and sealing plug immediately after each use.
- **Washing:** Use warm, soapy water and a soft brush or sponge to clean all surfaces. Rinse thoroughly under running water.
- **Dishwasher Safety:** This product is **not** dishwasher safe. Hand washing is recommended.
- **Drying:** Allow all components to air dry completely on a clean drying rack. Ensure no moisture remains before storage.
- **Storage:** Store the clean, dry pump and plug in a clean, dry place away from direct sunlight and extreme temperatures.

## TROUBLESHOOTING

- **Weak Suction:**
  - Ensure the pump is centered correctly on your nipple and forms a complete seal with your breast.
  - Try repositioning the pump or adjusting the squeeze pressure to create a stronger vacuum.
  - Check for any damage or tears in the silicone that might affect suction.
- **Discomfort or Pain:**
  - If the suction feels too strong or causes pain, reduce the initial squeeze pressure when applying the pump.
  - Ensure your nipple is comfortably centered within the flange and not rubbing against the sides.
  - If discomfort persists, discontinue use and consult a healthcare professional.
- **Pump Falling Off / Spillage:**
  - Ensure the suction cup base is placed on a flat, stable surface if you set the pump down.
  - Always use the sealing plug when storing milk in the pump or transporting it.
  - Re-apply the pump if the suction weakens or breaks during use.

SPECIFICATIONS


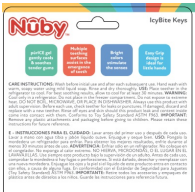
Feature	Detail
Product Dimensions	3 x 3 x 6.5 inches
Item Model Number	50550
Material Type	Silicone
Material Free	BPA Free
Additional Features	Lightweight, Portable
Number Of Items	1
Style	Nursing set
Batteries Required	No
Dishwasher Safe	No
Item Weight	5.9 ounces
Capacity	4 oz / 120 ml





WARRANTY AND SUPPORT

Specific warranty information for the Nuby Comfort All Silicone Breast Pump is not provided in this manual. For warranty details, product support, or to contact customer service, please refer to the product packaging or visit the official Nuby website.

© 2025 Nuby. All rights reserved.

Related Documents - 50550

	<p><a href="#">Nuby RapidCool Flask: User Instructions and Guide</a></p> <p>Comprehensive instructions for using the Nuby RapidCool Flask to prepare baby formula, including preparation steps, digital lid features, cooling times, cleaning, and warranty information.</p>
	<p><a href="#">Nuby IcyBite Keys Teether: Features and Care Instructions</a></p> <p>Discover the features and essential care instructions for the Nuby IcyBite Keys teether, designed to provide soothing relief for babies during teething. Learn how to properly clean and use this safe and effective teething toy.</p>

	<p><a href="#">Nuby RapidCool Bottle Maker Instructions and Guide</a></p> <p>Comprehensive guide to using the Nuby RapidCool bottle maker, including preparation, making a feed, digital LED indicators, cooling times, cleaning, storage, disposal, and NHS guidelines.</p>
	<p><a href="#">Nuby UV Dummy Steriliser: User Guide and Safety Instructions</a></p> <p>A comprehensive guide to the Nuby UV Dummy Steriliser, detailing its features, operation, safety warnings, technical specifications, and warranty information. Learn how to effectively sterilise soothers using UV light.</p>
	<p><a href="#">Nuby My Bathtime Clock &amp; Thermometer - User Manual</a></p> <p>User manual for the Nuby My Bathtime Clock &amp; Thermometer, detailing its functions as a clock, thermometer, and timer for safe and enjoyable bath times. Includes instructions for use, warnings, and battery replacement.</p>
	<p><a href="#">Nuby RapidCool Flask: User Instructions and Guide</a></p> <p>Comprehensive instructions for using the Nuby RapidCool Flask to prepare baby formula, including preparation steps, digital lid features, cooling times, cleaning, and warranty information.</p>

Documents - Nuby – 50550

no relevant documents