

DonJoy 01A-L

DJO 01A-L AIRCAST SP Walker Short Pneumatic Safety Shoe User Manual

Model: 01A-L

1. INTRODUCTION

This manual provides essential information for the proper use, setup, and maintenance of your DJO 01A-L AIRCAST SP Walker Short Pneumatic Safety Shoe. Please read these instructions thoroughly before use to ensure safe and effective application.

The SP Walker features a low rocker sole designed to promote natural ambulation, enhancing mobility. Its generous foot base accommodates dressings comfortably. The lightweight, semi-rigid shell includes a pre-inflated anterior air cell and two adjustable air cells at the malleoli, allowing for individual inflation to achieve a secure and total contact fit.

2. SAFETY INFORMATION

Important: Consult your healthcare professional before using this device. Follow all instructions provided by your doctor or therapist.

- Do not overtighten straps or overinflate air cells, as this may restrict circulation.
- Regularly inspect the device for wear and tear. Discontinue use if damage is observed.
- This device is intended for single patient use.
- Keep out of reach of children.
- If you experience increased pain, numbness, or discomfort, consult your healthcare professional immediately.

3. WHAT'S IN THE BOX

Your package should contain the following item:

- 1 x SP Walker Short Pneumatic Safety Shoe

4. SETUP AND APPLICATION



Figure 1: The DJO 01A-L AIRCAST SP Walker Short Pneumatic Safety Shoe. This image displays the complete grey short pneumatic walker boot with its green adjustable straps and the 'AIRCRAFT' logo visible on the side.

1. **Prepare the Boot:** Unfasten all straps and open the liner.
2. **Position Foot:** Place your foot and leg into the boot, ensuring your heel is seated firmly at the back.
3. **Secure Liner:** Wrap the soft liner around your foot and leg, securing it snugly but not too tightly.
4. **Fasten Straps:** Begin fastening the straps from the bottom (toes) upwards (calf). Ensure each strap is secure but allows for comfortable circulation.
5. **Inflate Air Cells:** The boot features a pre-inflated anterior air cell and two adjustable air cells at the malleoli (ankle bones). Locate the inflation bulbs/valves. Inflate the air cells until you feel comfortable compression and support. Do not overinflate. The air cells should provide a "total contact" fit.
6. **Check Fit:** Stand up and walk a few steps to assess the fit and comfort. Adjust straps and air cell inflation as needed. The boot should feel secure and supportive without causing pressure points or discomfort.

5. OPERATING INSTRUCTIONS

Once the SP Walker is properly applied, you can begin ambulating as directed by your healthcare professional.

- **Walking:** The low rocker sole is designed to promote a more natural gait. Walk slowly and deliberately, focusing on a smooth heel-to-toe motion.
- **Adjustments:** Periodically check the fit of the boot. If swelling changes, you may need to adjust the inflation of the air cells or the tightness of the straps.
- **Support:** Use crutches or other assistive devices as recommended by your healthcare provider, especially during the initial stages of recovery.

6. MAINTENANCE

- **Cleaning:** The soft goods (liner) can be hand washed in cold water with mild soap. Air dry completely before reattaching to the boot. Do not machine wash or tumble dry.
- **Shell Cleaning:** Wipe the rigid shell with a damp cloth and mild disinfectant. Do not immerse the shell in water.
- **Storage:** Store the boot in a cool, dry place away from direct sunlight when not in use.
- **Inspection:** Regularly inspect the boot for any signs of wear, damage, or air leakage from the cells. If any issues are found, contact your supplier or healthcare professional.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Discomfort or pressure points	Straps too tight, air cells overinflated, or improper positioning.	Loosen straps, slightly deflate air cells, or re-apply the boot ensuring proper heel seating. Consult healthcare professional if discomfort persists.
Boot feels loose or unstable	Straps too loose, air cells underinflated, or swelling has decreased.	Tighten straps, inflate air cells further (without overinflating), or consult healthcare professional for re-evaluation of fit.
Air cells not holding air	Valve not properly sealed or air cell damage.	Ensure valve is fully closed. If damage is suspected, discontinue use and contact your supplier or manufacturer.

8. SPECIFICATIONS

Model	01A-L
Product Name	AIRCAST SP Walker Short Pneumatic Safety Shoe
Brand	DonJoy
Size	Men Size 10-13, Women Size 11-15 (Large)
Package Dimensions	1 x 1 x 1 inches
Item Weight	1 Pound
Manufacturer	DJO Global, Inc.
First Available Date	September 19, 2019



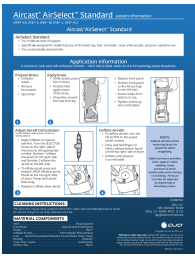
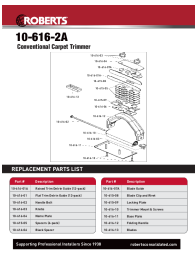
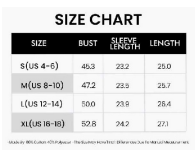
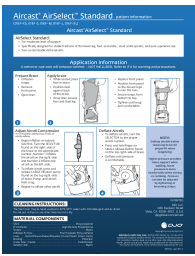
9. WARRANTY AND SUPPORT

For warranty information and product support, please contact the manufacturer directly. Keep your purchase receipt as proof of purchase.

Manufacturer Contact Information:

- Phone: 760-727-1280
- Manufacturer: DJO Global, Inc.

Related Documents - 01A-L

	<p>Comment Traiter une Entorse de Cheville</p> <p>Ce guide patient de DJO fournit des informations complètes sur le traitement des entorses de la cheville, couvrant les blessures courantes, l'instabilité chronique, le protocole GREC/RICE, les plans de traitement en plusieurs phases, la cryothérapie, les solutions orthopédiques, les options chirurgicales et la rééducation post-opératoire. Il détaille divers produits DJO et Aircast conçus pour la gestion et la prévention des entorses.</p>
	<p>Serendipity Latex Sizing Chart for Women</p> <p>A comprehensive sizing chart for Serendipity Latex products, detailing measurements for neck, bust, underbust, bicep, wrist, waist, hip, thigh, knee, calf, and ankle in both inches and centimeters across various sizes from XS to 2XL.</p>
	<p>Aircast AirSelect Standard Walking Brace Instructions</p> <p>Instructions and information for the Aircast AirSelect Standard walking brace, including preparation, application, compression adjustment, deflation, and cleaning.</p>
	<p>Roberts 10-616-2A Conventional Carpet Trimmer Replacement Parts List</p> <p>Official replacement parts list for the Roberts 10-616-2A Conventional Carpet Trimmer. Find part numbers and descriptions for all components.</p>
	<p>KIRUNDO Women's Fleece Hoodie Size Chart Cotton Polyester Sweatshirt Measurements</p> <p>Find the perfect fit with the KIRUNDO Women's Casual Cotton Long Sleeve Fleece Hoodie. This size chart provides detailed measurements for bust, sleeve length, and garment length for sizes S through XL.</p>
	<p>Aircast AirSelect Standard Walking Brace: Patient Information and Application Guide</p> <p>Detailed patient information and application guide for the Aircast AirSelect Standard walking brace. Learn how to prepare, apply, adjust compression, and deflate the brace. Includes cleaning instructions and material components.</p>

