

[Manuals.plus](#) /

› [August](#) /

› August Audar Smart Band KERi Instruction Manual

## August KERi

# August Audar Smart Band KERi Instruction Manual

Model: KERi

## INTRODUCTION

---

The August Audar Smart Band KERi is a versatile health and fitness tracker designed to monitor various aspects of your well-being. This manual provides detailed instructions for setting up, operating, maintaining, and troubleshooting your device to ensure optimal performance and a seamless user experience. Please read this manual thoroughly before using your Smart Band.

## SETUP

---

### 1. Charging the Smart Band

Before initial use, fully charge your Audar Smart Band KERi. The device features a convenient integrated USB charger. Simply remove the strap from one end of the tracker to reveal the USB connector.

# Up to 30 days of battery on a single charge



Image: The Audar Smart Band KERi being plugged into a USB port for charging. The image highlights the integrated USB connector on the device.



Image: A detailed view of the Audar Smart Band KERi's integrated USB charging port, showing how it connects to a power source.

- Connect the USB connector to any standard USB charging port (e.g., computer, USB wall adapter).
- The battery indicator on the screen will show charging progress.

- A full charge can provide up to 30 days of battery life.

## 2. Installing the Audar App

To unlock the full potential of your Smart Band and view your health data, download the official Audar app on your smartphone.

- Scan the QR code provided in the product packaging or search for "Audar" in your device's app store.
- The app is compatible with both Android and iPhone devices.

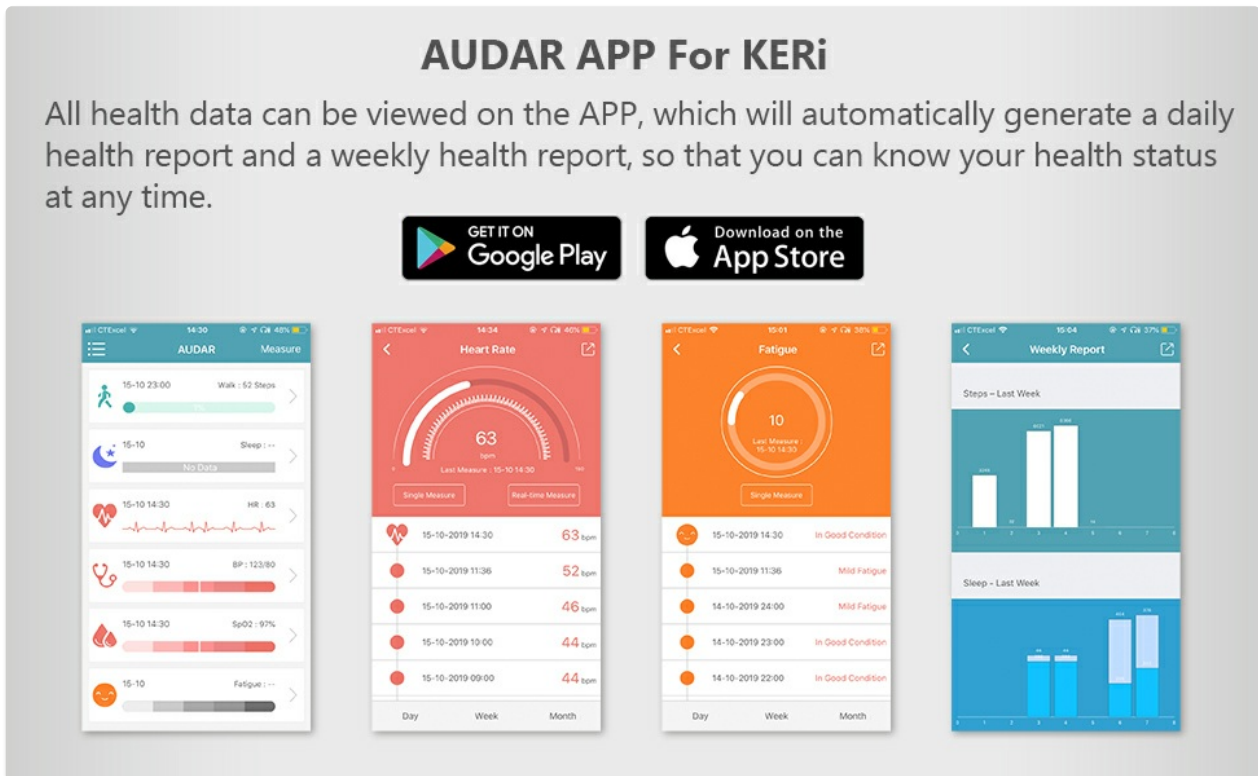


Image: Screenshot of the Audar app interface, displaying health data and showing icons for Google Play and Apple App Store downloads.

## 3. Pairing with Your Smartphone

- Ensure Bluetooth is enabled on your smartphone.
- Open the Audar app and follow the on-screen instructions to search for and connect to your KERi Smart Band.
- Confirm the pairing request on both your phone and the Smart Band if prompted.

## 4. Initial Configuration

- Once paired, the app will guide you through setting up your personal profile, including age, height, weight, and gender, for accurate data tracking.
- Synchronize the time and date with your phone.

## OPERATING INSTRUCTIONS

## 1. Basic Navigation

The Audar Smart Band KERi features a color touchscreen display for easy navigation.

- **Tap:** To select an option or wake the screen.
- **Swipe:** To scroll through different menus and data screens.

## 2. Activity Tracking

The Smart Band automatically tracks your daily activity.

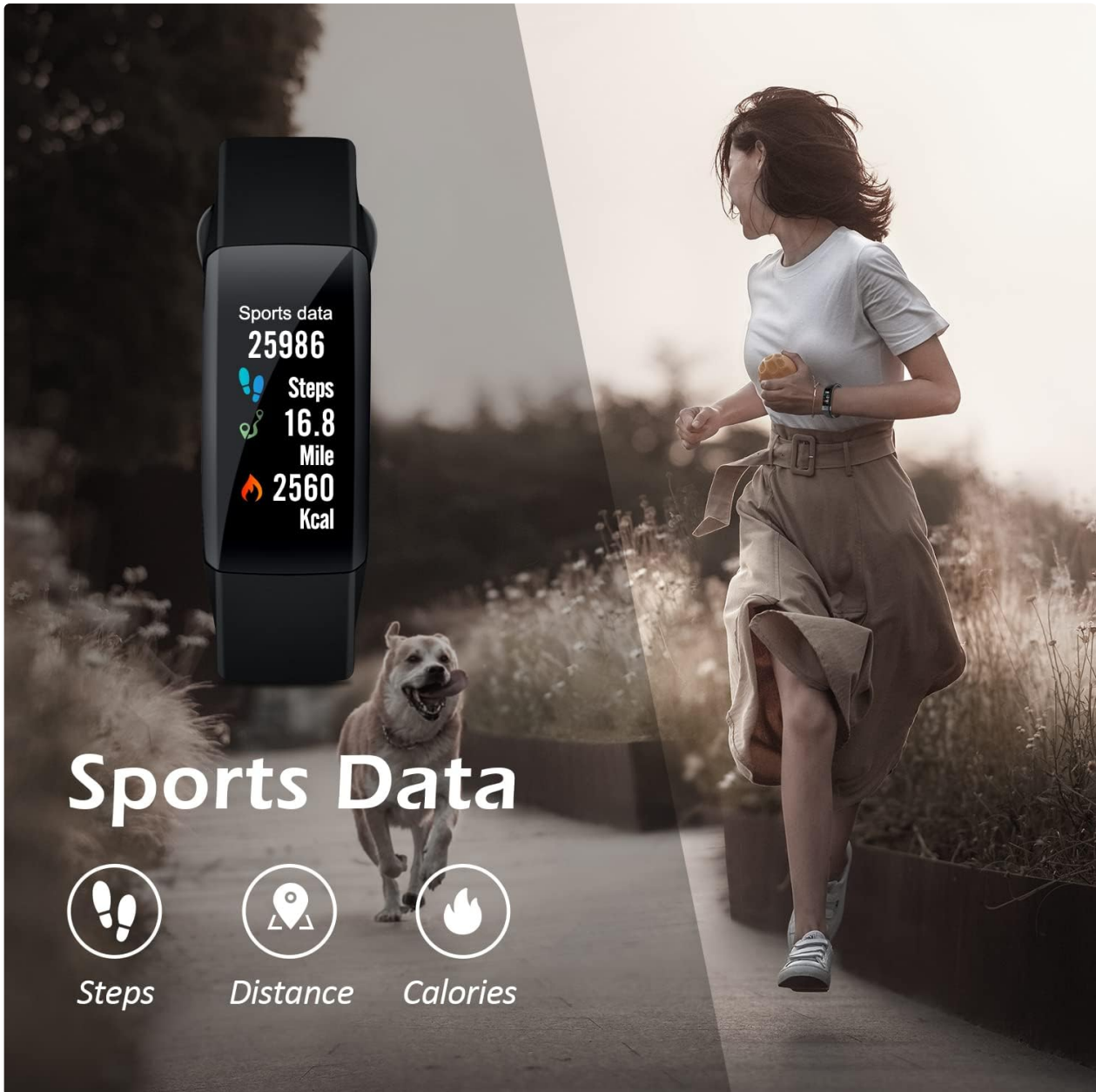


Image: A person running outdoors with the Audar Smart Band KERi on their wrist, showing the display with steps, distance, and calorie data.

- **Steps:** Counts the number of steps taken throughout the day.
- **Distance:** Estimates the distance covered based on your steps and profile data.
- **Calories:** Calculates calories burned during activity.
- **GPS Tracking:** When connected to your smartphone, the app can use your phone's GPS to map your route during outdoor activities.

## 3. Health Monitoring

The KERi Smart Band provides continuous and on-demand health monitoring.

- **Heart Rate Monitor:** Measures your heart rate in real-time. Access this feature from the band's menu or through the app.
- **Blood Pressure Monitor:** Provides an estimate of your blood pressure. For accurate medical readings, consult a healthcare professional.
- **Blood Oxygen Monitor:** Measures your blood oxygen saturation levels.
- **Sleep Monitor:** Automatically tracks your sleep patterns, including deep sleep, light sleep, and wake times. View detailed analysis in the Audar app.

# SLEEP

*Sleep quality and fatigue state management*



Image: A person sleeping with the Audar Smart Band KERi on their wrist, showing the band's display indicating sleep duration and an accompanying smartphone app displaying sleep analysis.

## 4. Notifications

Receive alerts directly on your wrist when connected to your smartphone.



Image: A close-up of a wrist wearing the Audar Smart Band KERi, with the screen showing an incoming call notification.

- **Call Alerts:** Displays incoming call notifications.
- **Message Alerts:** Shows notifications for SMS, WhatsApp, Messenger, and other selected apps.
- Configure which app notifications you wish to receive through the Audar app settings.

## 5. Additional Features

- **Remote Camera Control (Shake Shooting):** Use the band to remotely trigger your smartphone's camera by shaking your wrist.
- **Find My Phone:** Activate this feature from the band to make your paired phone ring, helping you locate it.
- **Alarm Clock:** Set silent vibrating alarms through the app.
- **CareMate Function:** Share your health data with family members through the app, allowing them to monitor your well-being remotely.

**Caremate Function**  
*You can receive data (Heart Rate, Blood pressure, Steps...) from member of your family on your phone.*

Image: Screenshots of the Audar app demonstrating the CareMate function, showing how users can subscribe to and view family members' health data.

## MAINTENANCE

---

### 1. Cleaning

- Regularly clean your Smart Band and strap to prevent skin irritation and maintain device hygiene.
- Wipe the device with a soft, damp cloth. Avoid using harsh chemicals or abrasive cleaners.
- Ensure the charging contacts are clean and dry before charging.

### 2. Water Resistance

The Audar Smart Band KERi is IP67 waterproof.



Image: The Audar Smart Band KERi being splashed with water, illustrating its IP67 waterproof capability. Icons for heart rate, calories, steps, distance, SMS alerts, and blood oxygen are also visible.

- It is resistant to rain, sweat, and can be submerged up to 1 meter underwater for a limited time.
- It is suitable for showering and swimming in shallow water.
- Avoid prolonged exposure to hot water, steam, or high-velocity water.

### 3. Storage

- Store the device in a cool, dry place when not in use.
- Avoid extreme temperatures.

## TROUBLESHOOTING

If you encounter issues with your Audar Smart Band KERi, refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on or respond	Low battery	Charge the device for at least 30 minutes.
Cannot pair with smartphone	Bluetooth off, app issue, device out of range	Ensure Bluetooth is on. Restart phone and band. Reinstall app. Move band closer to phone.
Inaccurate health readings	Improper fit, movement during measurement, sensor obstruction	Ensure band is snug but comfortable. Remain still during measurements. Clean sensor on back of band.
No data synchronization to app	Bluetooth connection lost, app not running in background	Check Bluetooth connection. Ensure app has necessary permissions and is running. Manually sync from app.
Short battery life	Frequent notifications, continuous monitoring, old battery	Reduce notification frequency. Disable continuous monitoring for certain features. Ensure full charge.

If the problem persists, please refer to the support section for further assistance.

## SPECIFICATIONS



Image: A diagram illustrating the dimensions of the Audar Smart Band KERi, including screen size and band length.

Feature	Detail
Model Name	KERi
Brand	August
Screen Size	0.96 Inches
Display Type	LCD
Resolution	240 x 240
Connectivity Technology	Bluetooth, USB
Wireless Communication Standard	Bluetooth
Operating System	Android (Compatible with Android & iPhone)
Metrics Measured	Heart Rate, Blood Pressure, Sleep Duration, Blood Oxygen
Special Features	Alarm Clock, Calorie Tracker, Notifications, GPS (via phone), Remote Camera Control
Water Resistance Level	IP67 (Waterproof up to 1 meter)
Battery Cell Type	Lithium Polymer
Battery Average Life	Up to 30 days
Band Material Type	Rubber or Silicone
Case Material Type	Ceramic
Item Weight	29 Grams
Item Dimensions	25.8 x 2.2 x 1.3 centimetres
Band Length	130 - 190 Millimetres
Color	Black

## WARRANTY INFORMATION

For detailed warranty information regarding your August Audar Smart Band KERi, please refer to the warranty card included in your product packaging or visit the official August website. Warranty terms and conditions may vary by region and retailer.



Image: The contents of the Audar Smart Band KERi product box, showing the smart band, charging cable, and instruction manuals.

## CUSTOMER SUPPORT

If you require further assistance or have questions not covered in this manual, please contact August customer support through one of the following channels:

- **Online Support:** Visit the official August website for FAQs, support articles, and contact forms.

- **Email Support:** Refer to your product packaging or the August website for the customer service email address.
- **Brand Store:** You can also visit the [August Brand Store](#) for more product information and support resources.

© 2024 August. All rights reserved.