

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [EAR BASICS](#) /

› [EAR BASICS Basic Ear Nude Personal Sound Amplifier Instruction Manual](#)

EAR BASICS Basic Ear Nude

EAR BASICS Basic Ear Nude Personal Sound Amplifier Instruction Manual

Model: Basic Ear Nude

[Included](#) [Setup](#) [Operation](#) [Product Overview](#) [Maintenance](#) [What's Troubleshooting](#) [Specifications](#)

1. PRODUCT OVERVIEW

The EAR BASICS Basic Ear Nude Personal Sound Amplifier is designed to enhance sound for clearer hearing in various daily situations. This device is lightweight, discreet, and rechargeable, offering up to 20 hours of continuous use on a single charge. It is not intended as a substitute for a hearing aid for individuals with diagnosed hearing impairment.

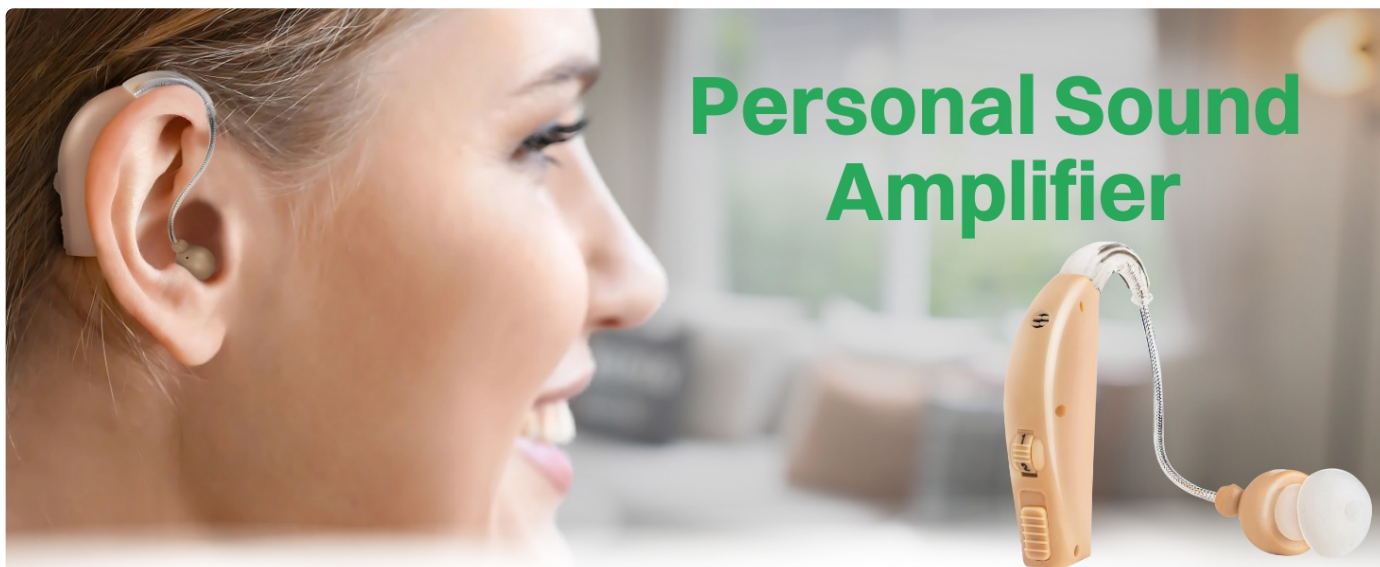


Image: The EAR BASICS Personal Sound Amplifier is designed for discreet wear, fitting comfortably behind the ear.

Key Features:

- **Sound Amplification:** Amplifies sound up to 30 times.
- **Compact Design:** Small, lightweight, and contoured for a comfortable fit, making it virtually unnoticeable.
- **Rechargeable:** Features a rapid recharging station, providing up to 20 hours of continuous use per charge.

- **Customizable Fit:** Includes 6 silicone ear tips of various sizes to ensure a secure and comfortable fit.
- **Easy Maintenance:** Comes with a cleaning brush for simple upkeep.

2. WHAT'S INCLUDED

Your EAR BASICS Personal Sound Amplifier package contains the following items:

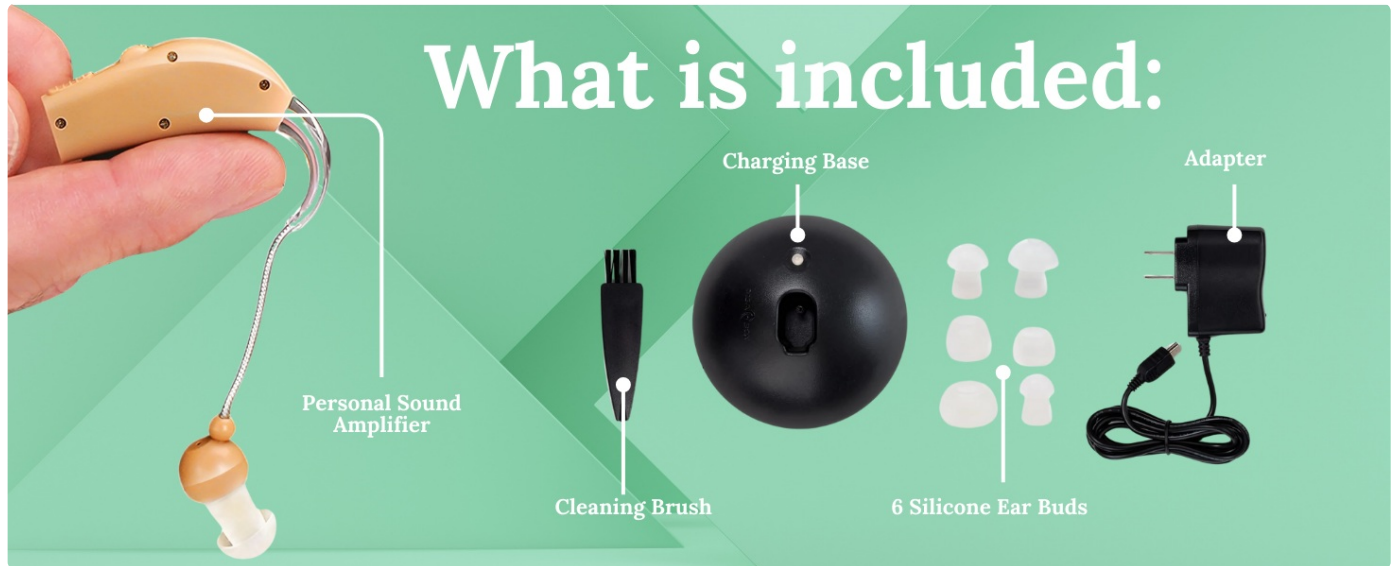


Image: All components included in the package: Personal Sound Amplifier, Charging Base, Adapter, 6 Silicone Ear Buds, and Cleaning Brush.

- 1 x EAR BASICS Personal Sound Amplifier
- 1 x Recharging Station
- 1 x AC Adapter
- 6 x Silicone Ear Tips (various sizes)
- 1 x Cleaning Brush

3. SETUP GUIDE

3.1. Initial Charging

1. Connect the AC adapter to the recharging station.
2. Plug the AC adapter into a standard wall outlet.
3. Place the Personal Sound Amplifier onto the recharging station. Ensure the amplifier is correctly seated for charging. A red indicator light on the charging station will illuminate to confirm charging is in progress.
4. Allow the device to fully charge before first use. The indicator light will change color (e.g., green) or turn off when charging is complete.



Image: The Personal Sound Amplifier placed on its recharging station. The charging indicator light is visible.

3.2. Selecting and Attaching Ear Tips

1. Examine the provided 6 silicone ear tips. They come in various sizes to accommodate different ear canals.
2. Select an ear tip that feels comfortable and provides a snug fit in your ear canal. A proper fit is crucial for optimal sound amplification and to minimize feedback.
3. Gently attach the chosen silicone ear tip to the sound amplifier's ear canal piece. Ensure it is securely fastened.

4. OPERATING INSTRUCTIONS

4.1. Wearing the Amplifier

1. With the ear tip attached, carefully place the amplifier over the top of your ear, so the main unit rests behind your ear.
2. Guide the ear tip into your ear canal until it feels comfortable and secure. The clear tubing should follow the contour of your ear.



Small & Perfectly Contoured

Image: The amplifier is designed to fit snugly and discreetly behind the ear, with the ear tip inserted into the ear canal.

4.2. Powering On/Off and Volume Adjustment

1. Locate the small power switch on the amplifier unit. Slide the switch to the "On" position to activate the device.
2. Once powered on, adjust the volume using the volume control dial or button (refer to your specific model for exact control type). Start with the lowest volume setting and gradually increase it until you reach a comfortable listening level.
3. To turn off the amplifier, slide the power switch to the "Off" position.



Lightweight

Image: The compact design of the amplifier, highlighting the location of the power and volume controls.

5. MAINTENANCE AND CARE

Regular cleaning and proper care will ensure the longevity and optimal performance of your Personal Sound Amplifier.

5.1. Daily Cleaning

- After each use, gently wipe the amplifier and ear tip with a soft, dry cloth to remove any earwax or debris.
- Use the provided cleaning brush to carefully clear any wax buildup from the ear tip and sound outlet.
- Do not use water or cleaning solutions directly on the amplifier unit, as this can damage the electronics.

5.2. Storage

- When not in use, store the amplifier in its recharging station or a clean, dry case.
- Avoid exposing the device to extreme temperatures, humidity, or direct sunlight.

6. TROUBLESHOOTING

Problem	Possible Cause	Solution
No sound or weak sound	Device is off. Low battery. Ear tip is clogged with earwax. Volume is set too low.	Turn the device on. Charge the device. Clean the ear tip and sound outlet with the cleaning brush. Increase the volume gradually.
Whistling or feedback noise	Improperly fitted ear tip. Ear tip not fully inserted. Volume too high.	Try a different size ear tip for a better seal. Ensure the ear tip is fully and snugly inserted into the ear canal. Reduce the volume.
Device not charging	Charging cable or adapter not connected properly. Faulty power outlet.	Ensure all connections are secure. Try a different power outlet.

7. PRODUCT SPECIFICATIONS

Model Name:	Basic Ear Nude
Brand:	EAR BASICS
ASIN:	B083ZJWBH4
Material:	Plastic
Item Weight:	5.6 ounces
Package Dimensions:	7.3 x 4.2 x 2.4 inches
Batteries:	1 Product Specific rechargeable battery required (included)
Manufacturer:	Global TV Concepts
Date First Available:	February 25, 2020

8. WARRANTY AND SUPPORT

For information regarding warranty coverage, technical support, or customer service, please refer to the contact details provided with your product packaging or visit the official EAR BASICS website. Keep your purchase receipt for warranty claims.

