



Manuals.plus /

› HOMCOM /

› HOMCOM 6-Level Adjustable Curved Bath Stool Spa Shower Chair Instruction Manual

## HOMCOM 713-052

# HOMCOM 6-Level Adjustable Curved Bath Stool Instruction Manual

Model: 713-052 | Brand: HOMCOM

## INTRODUCTION

---

This instruction manual provides essential information for the safe and effective use of your HOMCOM 6-Level Adjustable Curved Bath Stool. Please read this manual thoroughly before assembly and use, and retain it for future reference. This bath stool is designed to provide support and comfort during bathing or showering, particularly for individuals who may have difficulty standing for extended periods, including the elderly, injured, or pregnant women.



Image: The HOMCOM bath stool positioned inside a shower area.

## SAFETY INFORMATION

---

- Always ensure all legs are securely locked into place before use.
- Verify the stability of the stool on a flat, non-slip surface.
- Do not exceed the maximum weight capacity of 300 lbs (136 kg).
- Regularly inspect the stool for any signs of wear, damage, or loose components. Discontinue use if any issues are found.
- Clean the stool regularly to prevent mold or mildew buildup, especially the blue EVA cushion.
- Use caution when transferring on and off the stool.

## PACKAGE CONTENTS

---

Please check the package to ensure all components are present and undamaged:

- 1 x Curved Seat Frame (White, high-density plastic)
- 1 x EVA Cushion (Blue, non-slip)
- 2 x U-shaped Leg Supports (Aluminum alloy)

- 4 x Adjustable Leg Extensions with Non-Slip Feet (Aluminum alloy with rubber tips)
- Fasteners and small parts for assembly



Image: Detailed view of the seat features including the non-slip cushion and drainage holes.

## ASSEMBLY INSTRUCTIONS

---

The HOMCOM Bath Stool is designed for tool-free assembly. Follow these steps:

1. Unpack all components and lay them out on a clean, flat surface.
2. Attach the two U-shaped leg supports to the underside of the curved seat frame. Ensure they click securely into the designated slots.
3. Insert the four adjustable leg extensions into the bottom of the U-shaped leg supports.
4. Adjust the height of each leg by pressing the push-pins and sliding the extension to the desired level. Ensure all four legs are set to the same height for stability. You will hear a 'click' sound when the pin locks into a hole.
5. Place the blue EVA cushion onto the seat frame. It is designed to fit snugly and provide additional comfort and slip resistance.

Your browser does not support the video tag.

Video: Watch this video for a visual guide on assembling and using the HOMCOM Bath Stool.

## OPERATING INSTRUCTIONS

### Height Adjustment

The bath stool features 6 levels of adjustable height to accommodate various user needs and bathtub/shower heights. To adjust, simply press the spring-loaded push-pin on each leg and slide the inner leg extension up or down until the desired height is reached. Ensure the pin fully engages in a hole for secure locking. Always adjust all four legs to the same height to maintain stability.

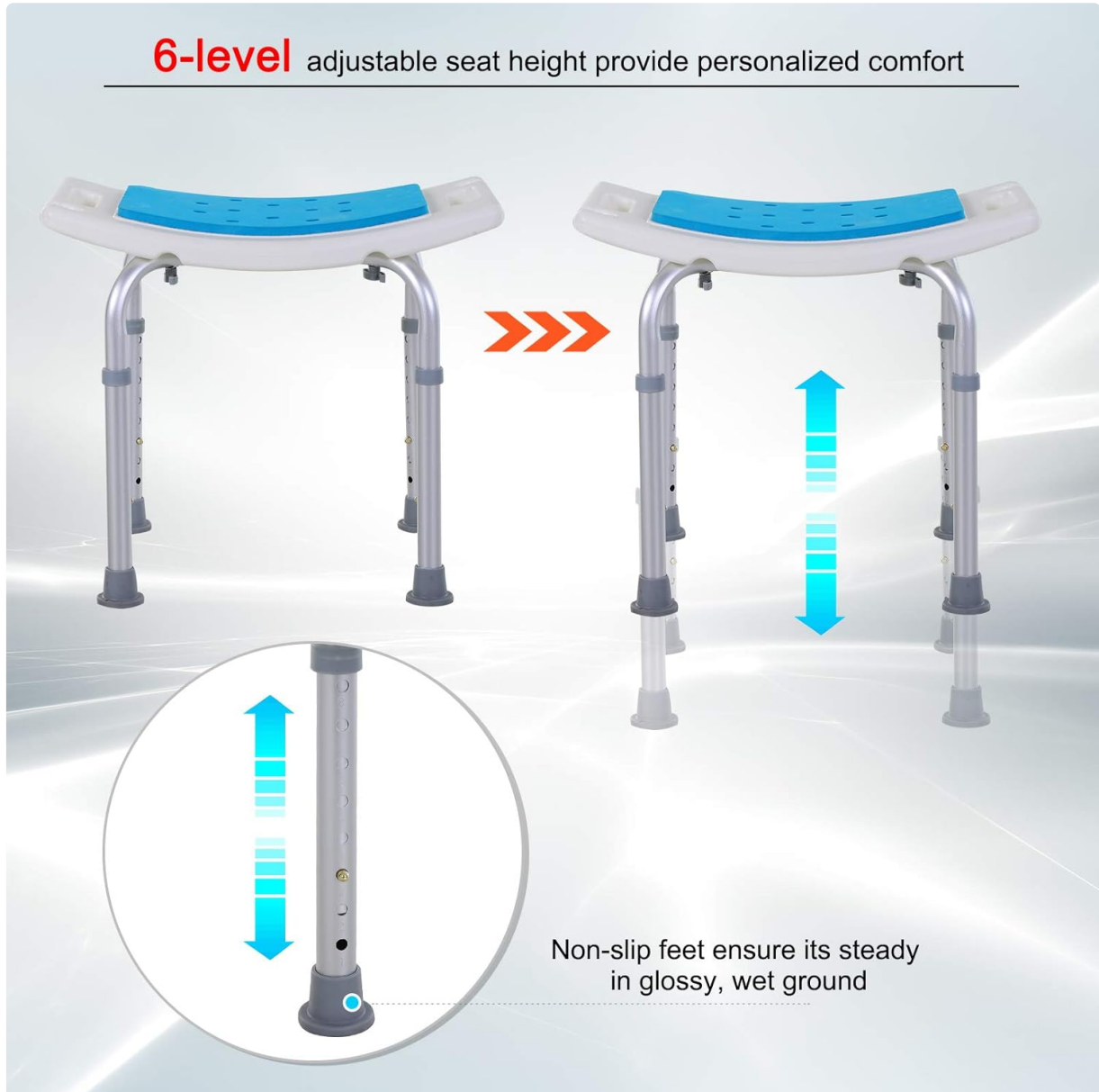


Image: Visual guide for adjusting the stool's height and highlighting non-slip feet.

### Placement and Use

Place the bath stool firmly inside your shower or bathtub. The non-slip feet are designed to provide stability on wet surfaces. The curved seat offers ergonomic comfort, and the smooth edge notch can hold a hand shower. The stool is suitable for use by seniors, individuals with disabilities, those recovering from injuries, or pregnant women seeking a safer and more comfortable bathing experience.



Image: Examples of users benefiting from the bath stool.



Image: Stool in a bathtub.



Image: Stool in a shower stall.

## MAINTENANCE

---

- **Cleaning:** After each use, rinse the stool with clean water. For thorough cleaning, use a mild soap or non-abrasive bathroom cleaner and a soft cloth. Avoid harsh chemicals or abrasive scrubbers that could damage the materials.
- **Drying:** To prevent mustiness and prolong the life of the EVA cushion, dry the blue seat thoroughly after each use.
- **Inspection:** Periodically check all connections and leg adjustments to ensure they remain secure. Inspect the non-slip feet for wear and tear; replace if necessary to maintain stability.
- **Storage:** Store the stool in a dry, well-ventilated area when not in use.

## TROUBLESHOOTING

---

- **Stool feels unstable:** Ensure all four legs are adjusted to the same height and the push-pins are fully engaged. Verify that the non-slip feet are clean and making full contact with the floor.
- **Difficulty adjusting height:** Check for any debris or obstruction in the leg extension holes. Ensure the push-pin is fully depressed before attempting to slide the leg.
- **Cushion slipping:** Ensure the EVA cushion is properly seated on the frame. Clean both the cushion

and the seat surface to remove any soap residue or buildup that might reduce friction.

## SPECIFICATIONS

<b>Brand</b>	HOMCOM
<b>Model Number</b>	713-052
<b>Material</b>	Aluminum Alloy, PE, EVA
<b>Color</b>	White (Frame), Blue (Cushion)
<b>Product Dimensions (Overall)</b>	20.5" W x 12.25" D x 15.25" - 20.5" H (52cm W x 31cm D x 39-52cm H)
<b>Adjustable Height Levels</b>	6 levels
<b>Maximum Weight Capacity</b>	300 lbs (136 kg)
<b>Item Weight</b>	4 lbs (1.8 kg)



Image: Key dimensions of the bath stool.

## WARRANTY AND SUPPORT

---

For warranty information or customer support, please refer to the documentation provided with your purchase or contact HOMCOM customer service directly. Keep your purchase receipt for warranty claims.