

[manuals.plus](#) /

- › [BEBIRD](#) /
- › [BEBIRD R1 Ultra Ear Wax Removal Tool Camera Instruction Manual](#)

## BEBIRD R1 Ultra

# BEBIRD R1 Ultra Ear Wax Removal Tool Camera Instruction Manual

Brand: BEBIRD

Model: R1 Ultra

## INTRODUCTION

The BEBIRD R1 Ultra is an advanced visual ear cleaning device designed for comprehensive ear care. It features a high-resolution camera and various attachments to assist in the safe and effective removal of earwax, as well as visual inspection of the ear canal, nose, and throat. This manual provides essential information for the proper setup, operation, and maintenance of your device.



*Image: The BEBIRD R1 Ultra device, showcasing its sleek design and integrated camera tip.*

## IMPORTANT SAFETY INFORMATION


- **Consult a Professional:** If you experience ear pain, discharge, or significant hearing loss, consult a medical professional before using this device.
- **Gentle Use:** Always use the device gently and avoid inserting it too deeply into the ear canal to prevent injury to the eardrum or ear canal.
- **Children:** Use with children should always be supervised by an adult. The device includes a safety stopper to prevent over-insertion.
- **Cleanliness:** Ensure all tips are thoroughly cleaned and disinfected before and after each use to prevent infections.
- **Power Source:** Use only the specified charging method and cable provided with the device.
- **Water Resistance:** The device features an IP67 waterproof lens, but the main body is not fully submersible. Avoid exposing the main unit to excessive moisture.

## PACKAGE CONTENTS

---

Your BEBIRD R1 Ultra package includes the following items:

- BEBIRD R1 Ultra Device with 10-megapixel lens
- Various ear scoop tips (for wet and dry earwax)
- Painless ear tweezers attachment
- Visual ear dropper attachment
- Canal inspection covers (Ear, Nasal, Throat/Oral)
- Charging cable
- User Manual

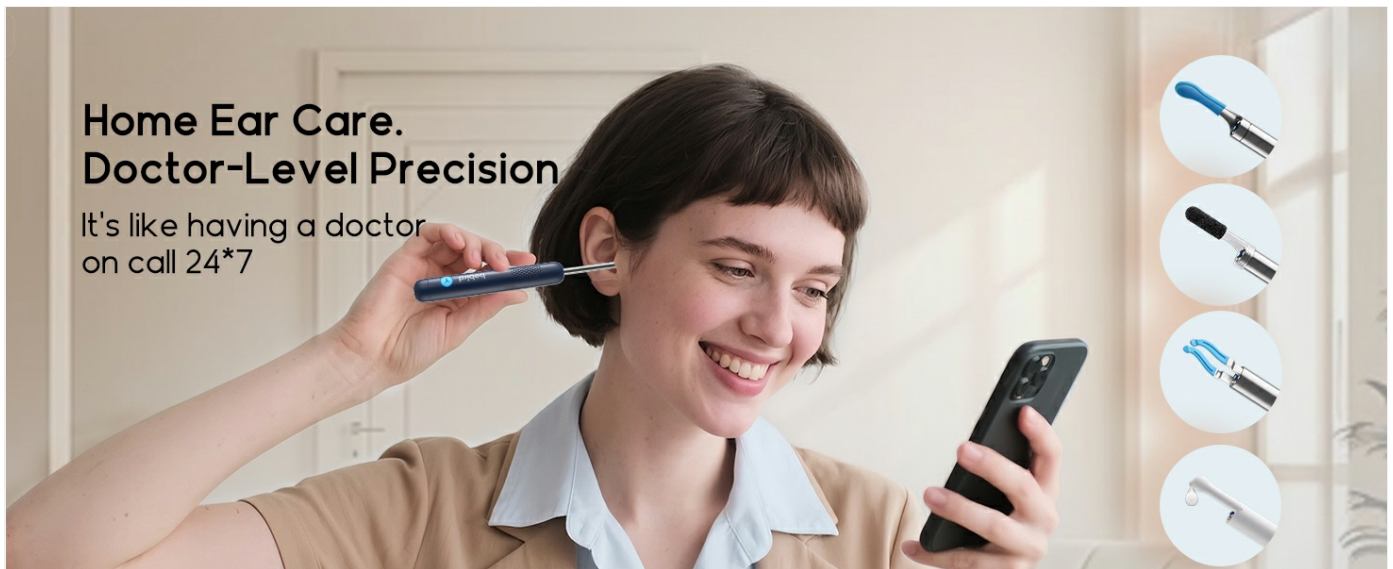

*Image: The BEBIRD R1 Ultra device and its comprehensive set of interchangeable tips and accessories.*

## SETUP GUIDE

---

1. **Charge the Device:** Before first use, fully charge the BEBIRD R1 Ultra using the provided charging cable. The indicator light will show charging status.
2. **Download the App:** Search for the "Bebird" app on your smartphone's app store (iOS or Android) and install it.
3. **Power On:** Press and hold the power button on the device to turn it on. The LED light at the tip will illuminate.
4. **Connect via Wi-Fi:**
  - Open your phone's Wi-Fi settings.
  - Locate and connect to the Wi-Fi network named "bebird-XXXX" (where XXXX is a unique identifier). No password is usually required.
  - Once connected, open the Bebird app. The device should automatically connect to the app.
5. **Attach a Tip:** Select the appropriate tip for your intended use (e.g., ear scoop, tweezers, inspection cover) and gently screw it onto the device's head.



## Visual

Seeing is believing



## Blind

Uncertain cleaning, no clear results



Image: Illustration of the one-tap Wi-Fi pairing process between the BEBIRD device and a smartphone.


## OPERATING INSTRUCTIONS

### General Use & Visual Inspection


1. Ensure the device is connected to the app and a suitable inspection cover is attached.
2. Gently insert the device into the area you wish to inspect (ear canal, nostril, or mouth). Observe the live feed on your smartphone screen.
3. The 10-megapixel lens provides a clear, magnified view for detailed examination.
4. Use the app's features to capture photos or record videos of your examination for personal records or to share with a healthcare professional.

## Upgrade to 10-Megapixel Lens

Get a clear visual on what's going on your ears

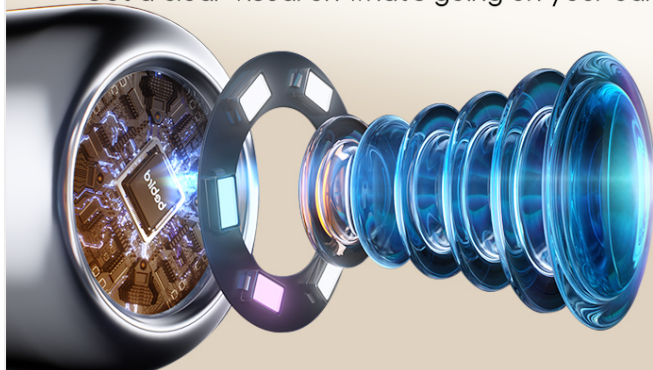


Bebird

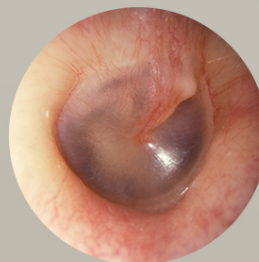

1080P HD Camera

Other

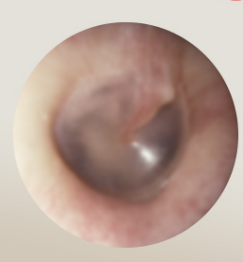

480P HD Camera

Image: Demonstrates the multi-functional use of the BEBIRD R1 Ultra for oral, nasal, and ear examinations.

### Earwax Removal

1. Select the appropriate ear scoop tip (for wet or dry wax) or the ear tweezers. Ensure the safety stopper is correctly positioned to prevent over-insertion.
2. With the device powered on and connected to the app, carefully guide the tip into the ear canal while observing the live view on your screen.



3. Gently scoop, grip, or scrape the earwax. Avoid aggressive movements or pushing wax deeper into the canal.
4. For stubborn earwax, consider using the visual ear dropper with a wax-softening solution (consult a doctor for recommendations).



Image: Various ear cleaning methods using the BEBIRD R1 Ultra, including scooping, tweezing, and using a sponge or dropper.

## Using the Visual Dropper

1. Attach the visual ear dropper tip to the device.
2. Press the sleeve button to draw medication liquid into the dropper.
3. Carefully insert the dropper into the ear canal, guided by the live camera view.
4. Press the button to administer the medication precisely to the target area.



Image: Step-by-step guide on how to use the precision dropper for ear cleaning and medication application.

## Video Demonstrations

Your browser does not support the video tag.

Video: Demonstration of the BEBIRD C3 Ear Camera Otoscope with Light, showing its functionality for ear inspection.

Your browser does not support the video tag.

Video: Overview of the Bebird R1 Ultra's capabilities for ear, nose, and throat examination.

Your browser does not support the video tag.

Video: Illustrates the process of visually guided earwax cleaning using the device.

Your browser does not support the video tag.

Video: A general demonstration of the Bebird Ear Wax Removal Tool with Camera.

## MAINTENANCE

### Cleaning the Tips

- After each use, carefully remove the attached tip by unscrewing it.
- Clean the tip with rubbing alcohol or an appropriate disinfectant. Ensure no debris remains.
- Allow the tip to air dry completely before reattaching or storing.
- The IP67 waterproof lens allows for easy cleaning of the camera tip itself with a damp cloth, but avoid submerging the main device body.

### Storage

- Store the BEBIRD R1 Ultra and its accessories in a clean, dry place, away from direct sunlight and extreme temperatures.
- Keep out of reach of children.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Device not turning on	Low battery	Charge the device fully.
Cannot connect to Wi-Fi	Incorrect Wi-Fi network selected or device not powered on.	Ensure the device is on and select the correct "bebird-XXXX" network. Restart both the device and your phone's Wi-Fi.
Blurry image on app	Camera lens is dirty or focus issue.	Clean the camera lens gently with a soft cloth. The device has a fixed focus, ensure proper distance.
Tips breaking easily	Excessive force or improper storage.	Use gentle pressure. Store tips properly to prevent damage. Replace worn or damaged tips.

## SPECIFICATIONS

- **Model:** BEBIRD R1 Ultra
- **Camera Resolution:** 10 Megapixels (1080P HD)
- **Magnification:** 3X
- **Lens:** 3-layer anti-distortion lenses
- **Waterproof Rating:** IP67 (lens only)

- **Connectivity:** Wi-Fi (20fps WIFI Chip for seamless visuals)
- **Gyroscope:** 9-Axis Gyroscope for steady image orientation
- **Dimensions:** 2 x 2 x 2 inches; 2.2 ounces (device only, approximate)
- **Manufacturer:** Heifeng Zhizao (Shenzhen) Technology Co., Ltd
- **First Available:** May 16, 2020
- **Batteries:** 1 9V batteries required (as per product specifications, typically refers to internal rechargeable battery for this type of device).



*Image: Details on the BEBIRD R1 Ultra's 10-megapixel camera, 1080P HD resolution, and IP67 waterproof lens.*

## WARRANTY AND SUPPORT

---

The BEBIRD R1 Ultra typically comes with a **1-year warranty** from the date of purchase. Please refer to your product packaging or the official BEBIRD website for specific warranty terms and conditions.

### Customer Support

For any questions, technical assistance, or warranty claims, please contact BEBIRD customer support through the following methods:

- **Via Amazon:**
  - Go to "Your Orders" on Amazon.
  - Locate your BEBIRD R1 Ultra order.
  - Select "Get product support" or "Chat" to contact the seller.
- **Official Website:** Visit the official BEBIRD website for further support options and contact information.

# Your Safety Matters: Bebird's Priority

Threaded design for stability, skin-safe silicone for comfort.

Safe Spiral Structure



Soft & Gentle Material

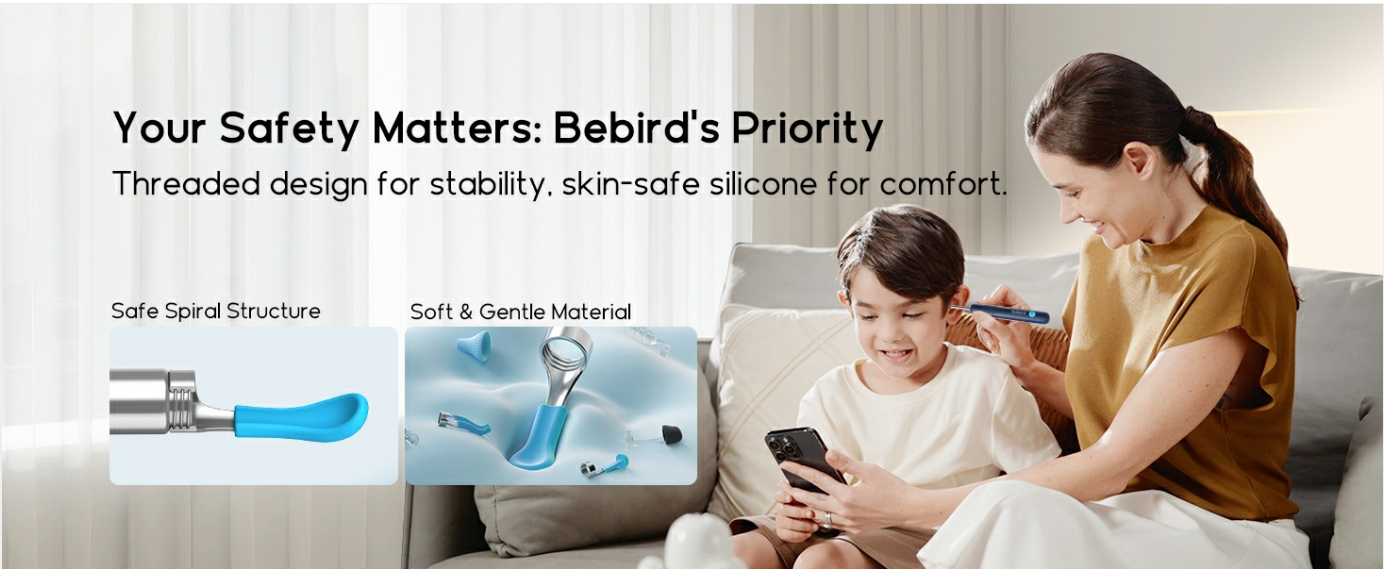


Image: A family utilizing the BEBIRD R1 Ultra for home ear care, highlighting its user-friendly design.



© 2024 BEBIRD. All rights reserved.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.