

## Stander Let's Go Indoor Rollator (Walnut), B086VT84JX

# Stander Let's Go Indoor Rollator User Manual

Model: Let's Go Indoor Rollator (Walnut) | Brand: Stander

## WELCOME TO YOUR STANDER LET'S GO INDOOR ROLLATOR

Thank you for choosing the Stander Let's Go Indoor Rollator. This lightweight, four-wheel Euro style walker is designed to provide stable and convenient mobility within your home. Its narrow frame, integrated tray, and storage basket offer practical assistance for daily activities. This manual provides essential information for the safe and effective use, maintenance, and care of your new rollator.

## IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using the rollator. Failure to follow these instructions could result in injury or damage.

- Ensure all parts are securely fastened and the rollator is stable before each use.
- Do not exceed the maximum weight capacity of the rollator (refer to specifications).
- Always engage the brakes when standing still or transferring to and from the rollator.
- Be aware of your surroundings and avoid obstacles.
- Do not use the rollator on uneven surfaces, stairs, or escalators.
- Keep hands and fingers clear of moving parts, especially during folding and unfolding.
- Regularly inspect the wheels, brakes, and frame for any signs of wear or damage.

**Legal Disclaimer:** Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is a mobility aid and not a medical device for treatment of conditions.

## SETUP AND INITIAL ASSEMBLY

The Stander Let's Go Indoor Rollator is designed for minimal assembly. It typically arrives ready for use with simple unfolding and handle adjustment.

### Unfolding the Rollator

1. Place the rollator on a flat, stable surface.
2. Gently pull the sides of the frame apart until it fully extends and locks into place. Ensure all locking mechanisms are engaged.

### Adjusting Handle Height

The single padded handle is height adjustable to suit individual user needs, ranging from 32.5 to 37.5 inches tall.

1. Locate the easy adjustment knob on the handle stem.
2. Loosen the knob and slide the handle to the desired height. The handle should be at wrist level when standing upright with arms relaxed.
3. Tighten the knob securely to lock the handle in position.



**Figure 1:** Rollator dimensions and adjustable handle height. The diagram illustrates the overall size and the range of handle height adjustment.

## OPERATING YOUR ROLLATOR

### Maneuvering and Braking

The Let's Go Rollator features four 6-inch wheels, with front swivel wheels for enhanced maneuverability in tight spaces. The braking system provides stability and control.

- To move, push the rollator forward while maintaining a comfortable distance behind it.
- To slow down or stop, gently squeeze the brake levers located on the handle.
- To engage the parking brake (lock the wheels), push the brake levers down until they click into place. This is essential when standing still or using the tray.
- To release the parking brake, pull the brake levers up.



**ADJUSTABLE  
FOAM HANDLE**



**UNILEVER BRAKE  
HANDLE**



**6-INCH  
SWIVEL WHEELS**



**Let's Go Indoor Rollator™**

**Figure 2:** Key operational features including the adjustable handle, brake lever, and swivel wheels.

### **Utilizing the Tray and Storage Basket**

The rollator comes equipped with a 16 by 14.5 inch removable tray and a fabric basket, ideal for carrying personal items, food, or beverages.

- The tray can be placed on top of the rollator frame to transport items or serve as a small table, for example, when eating meals from your favorite chair.
- The fabric basket is located beneath the tray and provides additional storage for larger items.
- Always ensure items are securely placed on the tray or in the basket to prevent spills or falls.

# REMOVABLE TRAY



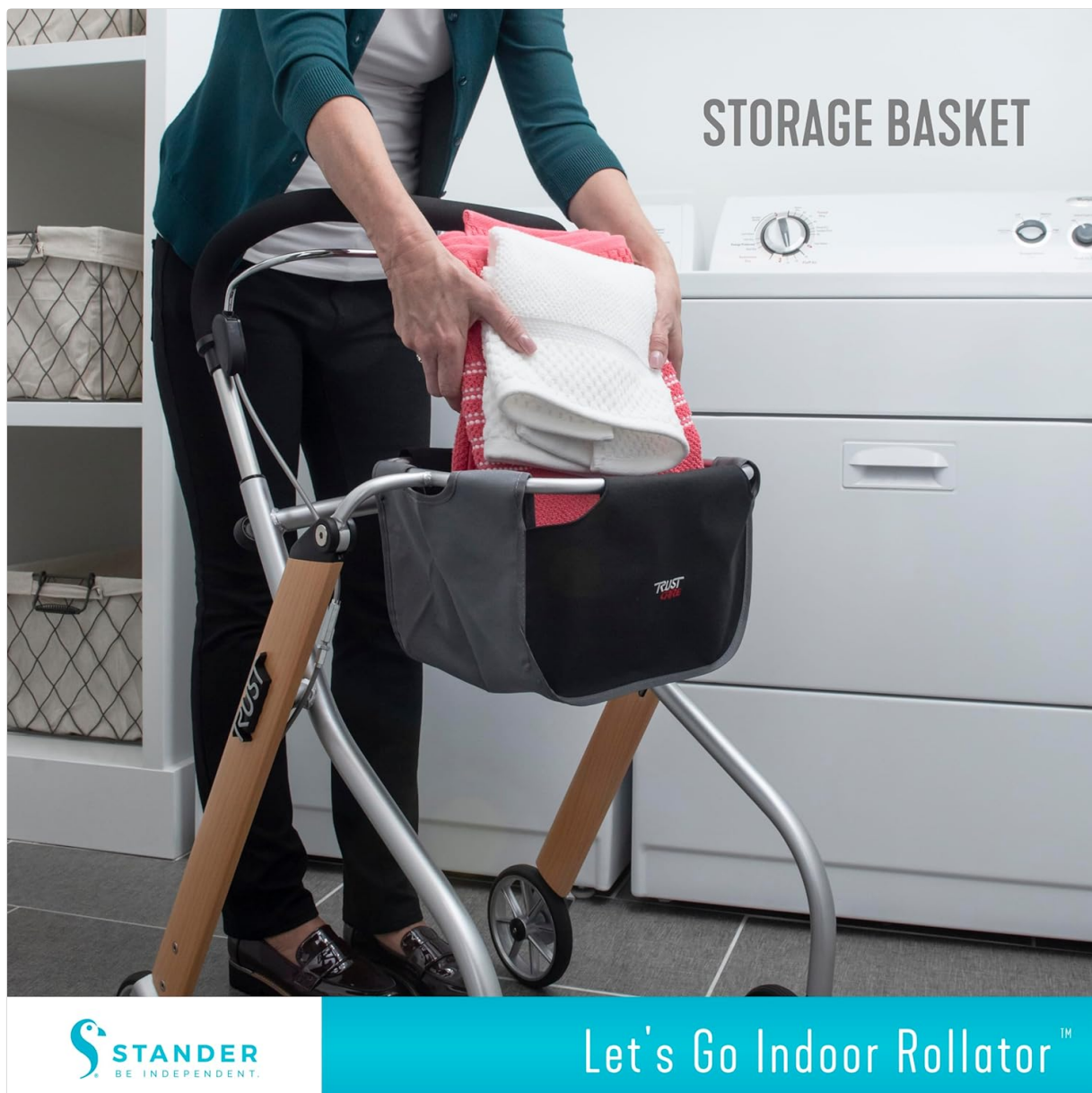
16" x 14.5" TRAY  
FOR EATING

 **STANDER**  
BE INDEPENDENT.

Let's Go Indoor Rollator™

**Figure 3:** The removable tray in use, demonstrating its utility for carrying items like food.





**Figure 4:** The storage basket, showing its capacity for carrying various household items.

## FOLDING FOR STORAGE

The lightweight 14.8-pound frame of the rollator is designed to fold easily for convenient storage and transport.

1. Remove any items from the tray and storage basket.
2. Lift the central folding mechanism (if applicable, refer to specific model instructions for exact location).
3. Gently push the sides of the rollator together until it collapses into its folded position.
4. The rollator folds to an 11-inch width, allowing for compact storage in narrow spaces.



**Figure 5:** The rollator in its folded state, demonstrating its space-saving design for storage.

## MAINTENANCE AND CARE

Proper maintenance will ensure the longevity and safe operation of your rollator.

- **Cleaning:** Wipe down the frame and handle with a damp cloth and mild detergent. Avoid abrasive cleaners. The tray can be washed with soap and water.
- **Wheels:** Periodically check the wheels for debris (hair, lint) and clean as necessary to ensure smooth rolling.
- **Brakes:** Test the brakes regularly to ensure they engage properly and hold the rollator securely. If the brakes feel loose or do not engage effectively, they may need adjustment. Consult a qualified technician or the manufacturer's support for brake adjustments.
- **Frame:** Inspect the aluminum frame for any cracks, bends, or signs of damage. Do not use the rollator if the frame is compromised.
- **Storage:** Store the rollator in a dry place, away from extreme temperatures and direct sunlight.

## TROUBLESHOOTING COMMON ISSUES

This section addresses common issues you might encounter with your rollator.

Problem	Possible Cause	Solution
Rollator does not roll smoothly.	Debris in wheels; loose wheel axle.	Clean wheels of any obstructions. Check if wheel axles are tight.
Brakes are not engaging or holding.	Brake cable too loose; worn brake pads.	Tighten brake cable (consult manual or technician). Inspect and replace worn brake pads if necessary.
Handle feels wobbly after adjustment.	Adjustment knob not fully tightened.	Ensure the handle adjustment knob is securely tightened.
Rollator does not fold completely.	Obstruction; folding mechanism not fully disengaged.	Check for anything blocking the folding mechanism. Ensure the folding release is fully activated.

## PRODUCT SPECIFICATIONS

Detailed specifications for the Stander Let's Go Indoor Rollator (Walnut).

Feature	Specification
Product Dimensions (D x W x H)	23.6"D x 22.8"W x 37.5"H
Item Weight	14.8 Pounds
Material	Aluminum
Color	Walnut
Handle Height Adjustment	32.5 to 37.5 inches
Tray Dimensions	16" x 14.5"
Wheel Size	6-inch
Manufacturer	Stander
ASIN	B086VT84JX



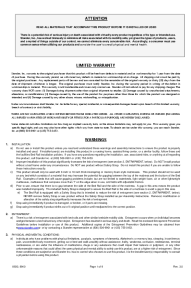

## SUPPORT AND WARRANTY INFORMATION

For specific warranty details, product support, or to purchase replacement parts, please refer to the official Stander website or the documentation included with your product. You may also contact Stander customer service directly for assistance.

**Note on Videos:** No official product videos with 'Seller' as creator type were available in the provided data for embedding in this manual.

**Online Resources:** For additional information and support, please visit the [Stander Store on Amazon](#).

Related Documents

	<p><a href="#">The Stander Step: Bedside Safety and Mobility Aid</a></p> <p>Introducing the Stander Step, a versatile and safe mobility aid designed to assist with getting in and out of bed. Features adjustable height, ergonomic grips, and an organizer pouch for enhanced independence.</p>
	<p><a href="#">Stander Security Pole 1150: Safety Warnings, Specifications, and Assembly Instructions</a></p> <p>Comprehensive guide to the Stander Security Pole (Model 1150), including essential safety warnings, product specifications, limited warranty details, and step-by-step assembly instructions for safe and effective installation.</p>
	<p><a href="#">Stander 8900 Series Bed Rail: Assembly, Installation, and Safety Guide</a></p> <p>Comprehensive guide for Stander 8900 Series bed rails (models 8930, 8940), covering installation, safety warnings, warranty, and accessories. Ensure safe and proper use of your Stander bed rail.</p>
	<p><a href="#">Stander Bed Cane Assembly and Safety Instructions</a></p> <p>Comprehensive guide to the safe assembly and installation of the Stander Bed Cane, a mobility aid designed to assist individuals in standing from bed. Includes important safety warnings, warranty information, and step-by-step instructions.</p>