



[Manuals.plus](#) /

› [GROWNSY](#) /

› GROWNSY Electric Baby Nasal Aspirator User Manual

GROWNSY BC023

GROWNSY Electric Baby Nasal Aspirator User Manual

MODEL: BC023

Brand: GROWNSY

PRODUCT OVERVIEW

The GROWNSY Electric Baby Nasal Aspirator (Model BC023) is designed to provide gentle yet effective relief from nasal congestion for infants and toddlers. This device features adjustable suction levels, soothing music, and a soft light function to make the cleaning process comfortable for your child. It is equipped with a USB-C rechargeable battery for convenience and comes with multiple food-grade silicone tips suitable for different ages.



Image: The GROWNSY Electric Baby Nasal Aspirator shown alongside its product packaging, highlighting its features and design.

Key Components:

- **Main Unit:** Contains the motor, control buttons, battery, and light/music functions.
- **Suction Chamber:** Transparent, detachable chamber for collecting mucus.
- **Silicone Tips:** Multiple sizes (3 food-grade silicone tips) for different age groups and comfort.
- **USB-C Charging Port:** For convenient recharging.

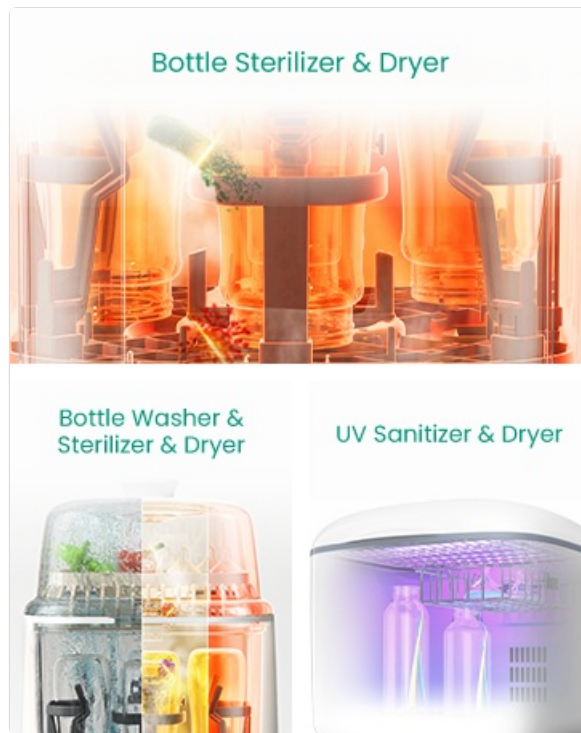


Image: A detailed view of the nasal aspirator's components, including the main unit, suction chamber, and silicone tips, with an illustration of a baby using the device.

SAFETY INFORMATION


It is crucial to adhere to the following safety guidelines to ensure the well-being of your child and the proper functioning of the device:

- Do not insert the nozzle too deep into the baby's nostril to prevent irritation or damage to nasal tissues.
- Keep the device out of reach of children when not in use.
- Ensure all components are clean and dry before and after each use.
- Do not immerse the main unit in water or any other liquid.
- Use only the provided food-grade silicone tips.

SETUP AND FIRST USE

1. **Unpack the Device:** Carefully remove the nasal aspirator and all its accessories from the packaging.
2. **Charge the Device:** Before first use, fully charge the aspirator using the provided USB-C cable. Connect the cable to the USB-C port on the device and a suitable power adapter. The indicator lights will show charging status. A full charge provides up to 30 days of use.
3. **Assemble the Suction Chamber:** Twist the transparent suction chamber clockwise onto the main unit until it clicks securely into place. Ensure a tight fit to maintain suction power.
4. **Install the Sealant Ring:** Fit the sealant rubber ring around the suction chamber base tightly before screwing the suction set into the suction port. This prevents leaks and ensures optimal performance.
5. **Select a Silicone Tip:** Choose the appropriate food-grade silicone tip based on your baby's age and nasal passage size. Gently attach it to the top of the suction chamber.

USB-C Rechargeable Convenience

 Enjoy 120Min of continuous use



No filters required, battery-included:
hassle-free nasal aspirator for life

Image: The nasal aspirator being charged via USB-C, demonstrating its rechargeable convenience, with additional silicone tips visible.

OPERATING INSTRUCTIONS

Follow these steps for effective and comfortable nasal aspiration:

1. **Prepare Nasal Passages:** If the mucus is dry or sticky, it is recommended to soften it with saline water for 30 seconds to 1 minute before aspiration.

2. **Power On:** Press the power button to turn on the device.
3. **Adjust Suction Level:** The aspirator features 3 adjustable suction levels. Press the suction level button to cycle through the options (Soft, Moderate, Strong) to choose the appropriate gear for different congestion situations or according to the baby's age.
4. **Activate Soothing Functions (Optional):** Press the music button to play calming lullabies and the light button to activate the gentle light function. These features are designed to distract and soothe your baby during the process.
5. **Perform Aspiration:** Gently insert the silicone tip into your baby's nostril, ensuring it forms a seal without going too deep. The device will automatically begin suctioning. Move the tip gently to clear mucus.
6. **Turn Off:** Once aspiration is complete, turn off the device immediately to prevent mucus backflow.

Soothe Your Fussy Baby

Cozy Light & Sweet Lullaby

Enjoy fun and hassle-free
nasal care at no extra cost



FSA or HSA eligible



Image: A baby looking at the nasal aspirator, which emits a gentle light and plays music, illustrating the soothing features during use.

Essential Item for A Baby Registry



- ✓ Gets snot out quickly
- ✓ Gentle on baby
- ✓ multifunctional
- ✓ Easy to clean

BEST ELECTRIC NASAL ASPIRATOR

what to expect.

 babycentre.

 Babygaga

 babylist

 Best Products Reviews

Parents the **BUMP**

Image: A mother gently using the GROWNSY Electric Baby Nasal Aspirator on her baby, demonstrating the ease of use and gentle application.

MAINTENANCE AND CLEANING

Proper cleaning and maintenance are essential for hygiene and longevity of your aspirator:

1. **Empty Suction Chamber:** Empty the transparent mucus collection chamber before mucus exceeds 1/2 of its total storage capacity to maintain consistent suction power and prevent overflow.

2. **Disassemble for Cleaning:** After each use, detach the silicone tip and the suction chamber from the main unit.
3. **Clean Components:** Rinse the silicone tip and the suction chamber thoroughly with warm soapy water. Ensure all mucus residue is removed. The silicone nozzle can also be sterilized in boiling water.
4. **Dry Components:** Allow all detachable parts to air dry completely before reassembling. Ensure the storage chamber and device are clean and DRY after use.
5. **Clean Main Unit:** Wipe the main unit with a soft, damp cloth. Do not immerse the main unit in water.



Image: An exploded view diagram illustrating how the nasal aspirator can be disassembled for thorough cleaning of its components.

TROUBLESHOOTING

If you encounter any issues with your GROWNSY Nasal Aspirator, please refer to the following common solutions:

Problem	Possible Cause	Solution
No suction or weak suction	Suction chamber not properly assembled or sealed. Mucus chamber is full. Silicone tip is clogged. Battery is low.	Ensure the suction chamber is twisted clockwise until it clicks and the sealant ring is properly fitted. Empty the mucus chamber. Clean the silicone tip and suction chamber thoroughly. Recharge the device.
Device does not turn on	Battery is depleted. Device malfunction.	Charge the device fully. If charging does not resolve the issue, contact customer support.
Liquid entered the main unit	Improper handling or overfilling of mucus chamber.	Immediately remove the storage chamber. Run the machine for 5-10 minutes to drain the liquid. Ensure the device is completely dry before next use.

For further assistance, please refer to the official [Troubleshooting Guide \(PDF\)](#) or contact GROWNSY customer support.

PRODUCT SPECIFICATIONS

Specification	Detail
Model Number	BC023
Brand	GROWNSY

Specification	Detail
Age Range (Description)	Newborn Baby Infant Toddler
Material Type	Acrylonitrile Butadiene Styrene (ABS)
Material Type Free	Latex Free
Power Source	Battery Powered (USB-C Rechargeable)
Product Dimensions	6.1"W x 1.38"H
Item Weight	11.3 ounces
Water Resistance Level	Not waterproof (main unit)
Style	Automatic
UPC	738632348336, 860011589501

WARRANTY AND SUPPORT

GROWNSY offers a **12-month warranty** on this product, covering manufacturing defects from the date of purchase. Please retain your proof of purchase for warranty claims.

For any questions, technical support, or warranty inquiries, please visit the official GROWNSY store or contact their customer service directly.

Official GROWNSY Store: [GROWNSY Store on Amazon](#)

Additional User Manual (PDF): [User Manual \(PDF\)](#)

© 2024 GROWNSY. All rights reserved.

This manual is for informational purposes only. Product specifications are subject to change without notice.