

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [Yaheetech](#) /

› [Yaheetech Aluminium 3 Folding Massage Table with Rolling Stool User Manual](#)

### Yaheetech Aluminium 3 Folding Massage Table with Rolling Stool

## Yaheetech Aluminium 3 Folding Massage Table with Rolling Stool User Manual

Model: Aluminium 3 Folding Massage Table with Rolling Stool

### INTRODUCTION

This manual provides comprehensive instructions for the safe and effective setup, operation, and maintenance of your Yaheetech Aluminium 3 Folding Massage Table and Rolling Stool. Please read all instructions carefully before use to ensure proper functionality and longevity of your product.



Image: The Yaheetech Aluminium 3 Folding Massage Table and Rolling Stool set.

## SAFETY INFORMATION

---

**WARNING:** This product can expose you to chemical Lead which is known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to [www.P65Warnings.ca.gov/furniture](http://www.P65Warnings.ca.gov/furniture).

- Ensure all components are securely assembled before use.
- Do not exceed the maximum weight capacity of the massage table or stool.
- Keep hands and fingers clear of moving parts during adjustment.
- Use on a flat, stable surface to prevent tipping.
- Regularly inspect the product for any signs of wear or damage. Discontinue use if damage is found.

## PACKAGE CONTENTS

---

Verify that all items are present in your package:

- 1 x Aluminium 3 Folding Massage Table
- 1 x Adjustable Rolling Stool
- 1 x Adjustable Headrest
- 2 x Armrests
- 1 x Storage Bag (integrated with table when folded)

## SETUP INSTRUCTIONS

---

### Massage Table Setup

1. **Unfold the Table:** Place the folded massage table on its side on a flat surface. Release the chrome clasps and gently open the table.
2. **Extend the Legs:** Allow the legs to unfold completely. Ensure all leg support cables are taut and the legs are locked into position.
3. **Adjust Height:** To adjust the table height, locate the small round knobs on each leg. Press the knobs and slide the inner leg section up or down to the desired height. Ensure all four legs are set to the same height and the knobs are securely engaged in the adjustment holes. The table offers 8 different height positions.
4. **Attach Headrest:** Insert the headrest support bars into the designated slots at one end of the table. Adjust the angle of the headrest as needed.
5. **Attach Armrests:** Slide the armrests into the slots located on the sides of the table.



Soft and Resilient High-Density Foam  
Skin-Friendly And Breathable Luxurious PVC Leather



Image: Visual guide for unfolding the massage table and attaching the headrest and armrests.

## Rolling Stool Assembly

The rolling stool typically requires minimal assembly. Follow these general steps:

1. **Attach Casters:** Insert the five casters (wheels) into the holes at the end of each leg of the stool base. Push firmly until they click into place.
2. **Connect Base to Seat:** Place the stool seat upside down. Insert the pneumatic cylinder (if separate) into the center hole of the stool base. Then, place the seat onto the top of the pneumatic cylinder.

## OPERATING INSTRUCTIONS

### Massage Table Operation

- **Height Adjustment:** As described in the setup section, press the round knobs on the legs to adjust the table height. Ensure all legs are at the same level for stability.
- **Headrest Adjustment:** The headrest can be tilted to various angles for client comfort. Loosen the adjustment levers, set the desired angle, and tighten the levers securely.
- **Armrest Usage:** The armrests provide additional comfort and support. Ensure they are fully inserted into their slots before use.



Image: Diagram showing the adjustable dimensions of the massage table, including height range.

## Rolling Stool Operation

- **Height Adjustment:** To adjust the stool's height, pull the lever located under the seat while seated or standing. Release the lever at the desired height.
- **Swivel Function:** The stool seat offers 360-degree rotation, allowing for easy movement and access without repositioning the base.
- **Mobility:** The five smooth-rolling casters allow for quiet and effortless movement across various floor surfaces.



Image: Diagram illustrating the 360-degree rotation and adjustable height features of the rolling stool.

## MAINTENANCE

- **Cleaning the Table:** Wipe the massage table surfaces with a dry cloth. For deeper cleaning, a damp cloth with a few drops of gentle cleanser can be used, then wipe dry immediately.
- **Cleaning the Stool:** The PU leather surface of the stool can be cleaned with a damp cloth and a mild cleanser. The chrome iron base can be wiped with a dry cloth.
- **Regular Inspection:** Periodically check all screws, bolts, and adjustment mechanisms for tightness and proper function. Tighten any loose fasteners.
- **Storage:** When not in use, fold the massage table and store it in a dry, cool place to prolong its lifespan.

## TROUBLESHOOTING

- **Table Instability:** Ensure all four legs are extended to the same height and the adjustment knobs are fully engaged. Check that the leg support cables are taut.
- **Stool Not Adjusting Height:** Verify that the height adjustment lever is fully engaged and the pneumatic cylinder is correctly seated.

- **Difficulty Folding/Unfolding:** Ensure no obstructions are present and the clasps are fully released. Do not force the mechanism.
- **Squeaking Noises:** Check all joints and fasteners for tightness. Applying a small amount of silicone lubricant to pivot points may help.



SPECIFICATIONS



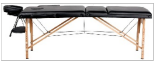

Feature	Specification
Brand	Yaheetech
Material (Table)	Aluminium, Alloy, PVC Leather, High-Density Foam
Material (Stool)	Chrome Iron, PU Leather, High-Density Sponge
Product Dimensions (Table)	84 x 31.5 x 35 inches (L x W x H, fully extended)
Adjustable Table Height	26.6 - 35 inches (8 positions)
Product Dimensions (Stool)	13.6 x 13.6 x 20.9-26.6 inches (Diameter x Height)
Item Weight	29.4 pounds (total for set)
Color	Black

WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the manufacturer's official website or contact the seller directly through your purchase platform. Keep your proof of purchase for any warranty claims.

Related Documents - Aluminium 3 Folding Massage Table with Rolling Stool

	<p><a href="#">Yaheetech Adjustable Massage Bed W/3 Sections (Model # 611035) Assembly Instructions</a></p> <p>Detailed assembly instructions for the Yaheetech Adjustable Massage Bed W/3 Sections, Model # 611035. Learn how to quickly, safely, and correctly assemble your new massage bed.</p>
	<p><a href="#">Yaheetech Portable Massage Table 610920 Installation Guide</a></p> <p>Step-by-step installation instructions for the Yaheetech 610920 portable massage table. Includes parts list, safety warnings, and assembly guidance.</p>

<div><div>Yaheetech</div><div>592032</div><div></div><div><div>1.</div><div>2.</div><div>3.</div><div>4.</div><div>5.</div><div>6.</div><div>7.</div></div><div>1. Please read the instruction manual carefully before assembly. 2. Please read the safety instructions carefully before assembly. 3. Please read the safety instructions carefully before assembly. 4. Please read the safety instructions carefully before assembly. 5. Please read the safety instructions carefully before assembly. 6. Please read the safety instructions carefully before assembly. 7. Please read the safety instructions carefully before assembly.</div><div>1</div></div>	<p><a href="#">Yaheetech 592032 Office Chair Assembly Instructions</a></p> <p>Comprehensive assembly instructions for the Yaheetech 592032 office chair, including parts list, step-by-step assembly guidance, and maintenance tips. Available in English, German, Spanish, French, and Italian.</p>
<div><div>Yaheetech</div><div>614812/614813</div><div></div><div><div>1.</div><div>2.</div><div>3.</div><div>4.</div><div>5.</div><div>6.</div><div>7.</div></div><div>1. Please read the instruction manual carefully before assembly. 2. Please read the safety instructions carefully before assembly. 3. Please read the safety instructions carefully before assembly. 4. Please read the safety instructions carefully before assembly. 5. Please read the safety instructions carefully before assembly. 6. Please read the safety instructions carefully before assembly. 7. Please read the safety instructions carefully before assembly.</div><div>1</div></div>	<p><a href="#">Yaheetech Ergonomic Gaming Chair Installation Guide</a></p> <p>Step-by-step instructions and parts list for assembling your Yaheetech ergonomic gaming chair, featuring massage lumbar support and adjustable footrest. Ensure proper setup for optimal comfort and performance.</p>
<div><div>Yaheetech</div><div>615206</div><div></div><div><div>1.</div><div>2.</div><div>3.</div><div>4.</div><div>5.</div><div>6.</div><div>7.</div></div><div>1. Please read the instruction manual carefully before assembly. 2. Please read the safety instructions carefully before assembly. 3. Please read the safety instructions carefully before assembly. 4. Please read the safety instructions carefully before assembly. 5. Please read the safety instructions carefully before assembly. 6. Please read the safety instructions carefully before assembly. 7. Please read the safety instructions carefully before assembly.</div><div>1</div></div>	<p><a href="#">Yaheetech 615206 Portable Massage Table Assembly and Care Instructions</a></p> <p>Comprehensive assembly guide and care instructions for the Yaheetech 615206 portable massage table. Learn how to set up, use, and maintain your massage bed for optimal performance and longevity.</p>
<div><div>Yaheetech</div><div>592328</div><div></div><div><div>1.</div><div>2.</div><div>3.</div><div>4.</div><div>5.</div><div>6.</div><div>7.</div></div><div>1. Please read the instruction manual carefully before assembly. 2. Please read the safety instructions carefully before assembly. 3. Please read the safety instructions carefully before assembly. 4. Please read the safety instructions carefully before assembly. 5. Please read the safety instructions carefully before assembly. 6. Please read the safety instructions carefully before assembly. 7. Please read the safety instructions carefully before assembly.</div><div>1</div></div>	<p><a href="#">Yaheetech 592328 Ergonomic Office Chair Assembly and Maintenance Guide</a></p> <p>Comprehensive assembly instructions, maintenance tips, and user guide for the Yaheetech 592328 ergonomic office chair. Learn how to assemble, care for, and adjust your chair for optimal comfort and longevity.</p>