Manuals+

Q & A | Deep Search | Upload

manuals.plus /

- Bistos
- » Bistos Bt-410 Head Lamp User Manual

Bistos Bt-410

Bistos Bt-410 Head Lamp User Manual

Model: Bt-410

1. Introduction

Thank you for choosing the Bistos Bt-410 Head Lamp. This ultra-bright LED head lamp is designed for various applications requiring focused and hands-free illumination. Featuring a high-density LED, an easy-to-adjust head strap, and an extensive LED lifetime, the Bt-410 provides reliable performance. This manual provides essential information for the safe and effective use, setup, operation, and maintenance of your device. Please read this manual thoroughly before using the product and retain it for future reference.

2. SAFETY INFORMATION

Always observe the following safety precautions to prevent injury or damage to the device.

- Eye Safety: Do not stare directly into the LED light. The ultra-bright light can cause temporary vision impairment or permanent eye damage.
- Heat: The LED light source may become warm during prolonged use. Avoid direct contact with the light housing when it is hot.
- Battery Safety: Ensure proper battery polarity when inserting. Do not mix old and new batteries, or different types of batteries.

 Dispose of batteries according to local regulations.
- Water Exposure: This device is not waterproof. Avoid exposing the head lamp to water or excessive moisture to prevent electrical shock or damage.
- Cleaning: Use only a soft, dry cloth for cleaning. Do not use abrasive cleaners or solvents.
- Children: Keep the device out of reach of children.

3. SETUP

3.1 Adjusting the Head Strap

The Bt-410 Head Lamp features an easy-to-adjust head strap for a comfortable and secure fit. Place the head lamp on your head and adjust the strap using the buckles until it feels snug but not too tight. Ensure the main unit rests comfortably on your forehead.



The Bistos Bt-410 Head Lamp worn on a head, demonstrating its adjustable strap and the position of the ultra-bright LED light.

3.2 Attaching the Loupe (Optional)

If you have the optional attachable loupe, carefully align it with the designated attachment points on the head lamp's main unit and secure it according to the loupe's specific instructions. Ensure it is firmly attached before use.



A side view of the Bistos Bt-410 Head Lamp, highlighting the main unit, adjustable light head, and the comfortable head strap.

3.3 Charging the Battery

The Bt-410 Head Lamp is designed for easy recharging. Connect the provided charging cable to the charging port on the head lamp and plug the other end into a suitable power source. An indicator light will typically show the charging status. Ensure the device is fully charged before first use to maximize battery life. The device provides more than 4 hours of continuous use on a full charge.

4.1 Powering On/Off

Locate the power button on the main unit of the head lamp. Press the button once to turn the light on. Press it again to turn the light off.

4.2 Adjusting the Light Angle

The LED light head is designed to be adjustable. Gently tilt the light head up or down to direct the beam to your desired angle. Do not force the adjustment beyond its natural range of motion.

5. MAINTENANCE

5.1 Cleaning

To clean the head lamp, wipe the exterior surfaces with a soft, dry, lint-free cloth. Do not use liquid cleaners, aerosols, or solvents. Ensure the device is powered off and disconnected from any charging source before cleaning.

5.2 Battery Care and Storage

To prolong battery life, avoid fully discharging the battery frequently. Recharge the device when the light begins to dim. If storing the head lamp for an extended period, ensure it is partially charged (around 50%) and store it in a cool, dry place away from direct sunlight. Recharge it every few months to prevent deep discharge.

5.3 LED Lifetime

The Bt-410 Head Lamp features an extensive LED lifetime of more than 50,000 hours, ensuring long-term performance under normal operating conditions.

6. TROUBLESHOOTING

If you encounter issues with your Bt-410 Head Lamp, refer to the following common problems and solutions:

. Light does not turn on:

- Ensure the battery is charged. Connect the head lamp to its charger and allow sufficient time for charging.
- Check if the power button is functioning correctly.

. Light is dim:

- The battery may be low. Recharge the head lamp.
- Ensure the LED lens is clean and free from dirt or smudges.

· Head strap is uncomfortable:

- Adjust the strap tension to find a comfortable fit.
- Ensure the main unit is positioned correctly on your forehead.

If the problem persists after attempting these solutions, please contact Bistos customer support for further assistance.


7. SPECIFICATIONS

Feature	Specification
Model	Bt-410
Light Source	Ultra bright LED with high density LED
LED Lifetime	More than 50,000 hours
Continuous Use Time	More than 4 hours (on full charge)
Adjustability	Easy to adjust head strap and light angle
Optional Accessories	Attachable loupe
Manufacturer	Bistos
Package Dimensions	32.8 x 23.6 x 9.2 cm
Package Weight	2.54 kg
ASIN	B091BMBRNM
Date First Available	29 March 2021
Origin	Made in Korea

8. WARRANTY AND SUPPORT

Specific warranty information for the Bistos Bt-410 Head Lamp is not provided within this manual. For details regarding warranty coverage, terms, and conditions, please refer to the documentation included with your purchase or visit the official Bistos website. For technical support, service, or inquiries, please contact Bistos customer service directly.

© 2023 Bistos. All rights reserved.

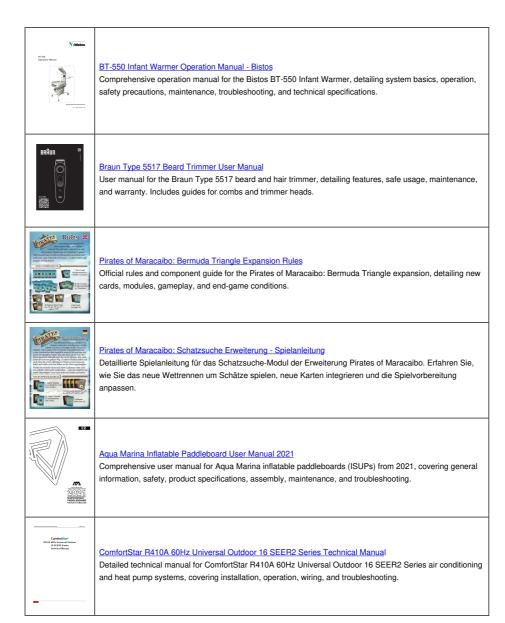

Documents - Bistos - Bt-410



[pdf] User Manual

BT 410 Operation Manual 2 Bistos Hewlett Packard BISTOS 01 MANUAL drive google u 0 uc id 1e9Bh7pUPFTWezAC tZISX0TJKHSTgTVc export |||

BT-410 Head-worn light Operation Manual Keep this manual for future reference P/N: 410-ENG-OPM-EUR-R10 BT-410 Operation descriptions contained in this manual are the property of Bistos Co. Ltd. and may not be copied, reproduced, disseminated or dis... lang:en score:45 filesize: 855.15 K page_count: 23 document date: 2019-07-23



[pdf] User Manual Owner's Manual Instructions

OHJ Bistos BT 410 Owners Manual lib store yahoo net yhst 47024838256514 $\mid\mid\mid$

BT-410 Head-worn light Operator Manual BT-410 P/N: OPMBT-410 EN, Rev.05 BT-410 Operator Manual 1 Proprietary Material Inf the property of Bistos Corporation and may not be copied, reproduced, disseminated or distributed without express wri... lang:en score:44 filesize: 509.96 K page_count: 19 document date: 2013-04-11



[pdf] User Manual

137222 bistos led head lamp bt410f user manual Download Bistos LED Head Lamp BT410F Praxisdienst Shop praxisdienst out media ||| ||| BT-410F and BT-410A reference P/N: 410-ENG-OPM-EUR-R12 BT-410F and BT-410A Operation manual 1 Proprietary Material Information and descriptions contained in this manual an BT-410F and BT-410A Head-worn light Operation Manual Keep this manual for future reference P/N: 410-ENG-OPM-EUR-R12 BT-Material Information and descriptions contained in this manual are the property of Bistos Co. Ltd. and may not be copied, reprod... lang:en score:33 filesize: 1.4 M page_count: 17 document date: 2023-04-26



[pdf] Specifications Warranty Borchure

Bistos BT 410 Surgical Headlight Brochure lib store yahoo net yhst 47024838256514 |||

Medical Head Lamp BT - 410F BT - 410A Medical Head Lamp **BT-410** F Illumination Uniformity Multi-lens structure ensure high leve Light source is located between two eyes is able to cast least shadow on the surface. Also, the luminance angle move... lang:en score:32 filesize: 524.12 K page_count: 3 document date: 2017-09-05



[pdf] User Manual

148160 bistos bt 410a led headlamp operation manual Download Lámpara LED Frontal Bistos BT 410A praxisdienst es out media |||


BT-410F and BT-410A Head-worn light Operation Manual Keep this manual for future reference P/N: 410-ENG-OPM-EUR-R12 BT-Material Information and descriptions contained in this manual are the property of Bistos Co. Ltd. and may not be copied, reprod... lang:en score:32 filesize: 1.29 M page count: 17 document date: 2025-01-08



[pdf] User Manual Instructions

BT 400 2 Bistos RA Users Manual Yahoo this Operator's Only the AC cord supplied with is approved for use Unit Monitoring patient status while receiving photothe stored in radiation area Do touch equipment wet hand bt users manual lib store yahoo net yhst 47024838256514

BT-400 Neonatal Phototherapy unit Operation Manual BT-400 Keep this manual for future reference P/N: ... - guidance Portable mot closer to any part of the **BT-410**, including cables, than the recommended separation distance calculated from the equation appl... lang:en score:30 filesize: 796.8 K page_count: 31 document date: 2019-08-16

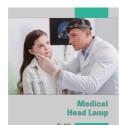

[pdf] Specifications Warranty

Medical Head LampSingle lens loupe is light weight and simple but provides 2x magnification Page 3 BISTOS Technical Specification Model BT 410 Order codeHe AfmBOopeR3hOws9XBB1FN1K oJ1ox3WBJEr6nkJwqqXIVgfW1uCu4GMk |||

Medical Head Lamp BT - 410F BT - 410A Medical Head Lamp BT-410 F Illumination Uniformity Multi-lens structure ensure high leve Light source is located between two eyes is able to cast least shadow on the surface. Also, the luminance angle move... lang:en score:29 filesize: 526.95 K page_count: 3 document date: 2018-05-17



[pdf] User Manual Catalog

Bio Signal Total Solution npcdn net BT 750 Patient Monitor 710 Pulse Oximeter 740 700 Pre warm Baby Manual mode Swivel Head Left right LED Bisto general C 6cabc39636115d55f8c9c1732bd24354 |||

Bio Signal Total Solution Fetal Monitor Fetal Doppler Infant Incubator Infant Warmer Phototherapy Pa ... evention BPA / Latex free R bottle set Others **BT-410** F Head Lamp Power LED light source 6,000 k Color temperature Adjustable Illumination Coaxial ... lang:en score:28 filesize: 1.39 M page_count: 5 document date: 2018-10-19



[pdf] Specifications Warranty

Head lamp BT 410 halomedicals 2020 08 |||

Medical Head Lamp BT - 410F BT - 410A Medical Head Lamp BT-410 F Illumination Uniformity Multi-lens structure ensure high leve Light source is located between two eyes is able to cast least shadow on the surface. Also, the luminance angle move... lang:en score:26 filesize: 526.95 K page_count: 3 document date: 2018-05-17



[pdf] Decleration of Conformity

PowerPoint Dasom Bistos Co Ltd HKTDC Sourcing 9401c889491243e4940365244a9c2260 sourcing media hktdc doc |||

Table of Contents 01 Prologue 02 Company Introduction 03 Core Competence Management 04 Products 2 ... re FM1 Arugus AFM F GT-1500 BT-200V BT-400 BT-700 BT-200 BT-200 BT-300 BT-250 BT-200T BT-500 Babycare KR-100 Pumpables I lang:en score:23 filesize: 1.61 M page_count: 29 document date: 2021-01-11



[pdf]

137222 bistos led kopflampe bt410f gebrauchsanweisung BT 410F und 410A Kopflampe Gebrauchsanweisung praxisdienst de out media |||

BT-410F und BT-410A Kopflampe Gebrauchsanweisung Bewahren Sie dieses Handbuch zum spteren Nachschlagen auf P/N: 410-Gebrauchsanweisung 1 Urheberrechtlich geschtztes Material Die in diesem Handbuch enthaltenen Informationen und Beschreibung lang:de score:21 filesize: 1011.36 K page_count: 18 document date: 2023-04-26



[pdf] User Manual Borchure

Bistos Co Ltd E Katalog LKPP prieš 6 dienas — Operator s manual 1EA Astral Lamp Optional Accessary Loupe 1 5 x 3 Battery Module GEA MEDICAL bistos EYE

f6950576e56b32f9fd9393d75e4ace40b513216c0cdd22672937a1ad8c4eebe73bf9dd019f8ed7261000e6a26e88d5b7bfb239967825064fbd4bb19c1182711afdb213

lang:de score:18 filesize: 3.73 M page count: 2 document date: 2016-05-31



[pdf] Decleration of Conformity

General~2020~indd~Bistos~Co~Ltd~HKTDC~Sourcing~571ed 1a3ef 1245a3963ef 8f 0c7a4243d~sourcing~media~hktdc~doc~|||

BT-380 Fetal Maternal Monitor 1MHz Frequency with Twin Doppler 10.1 Color Touch Screen Fetal Mon ... unction Blue/Red Color M Measurement High Accuracy BT-410A BT-410F BT-410 Head Lamp Power LED Light with 6,000K Color Temperature Adjustable ill lang:en score:13 filesize: 998.42 K page_count: 8 document date: 2020-12-10