



Manuals.plus /

› GROWNSY /

› GROWNSY Nasal Aspirator for Baby | Electric Nose Suction | Model BC024 Instruction Manual

GROWNSY BC024

GROWNSY Nasal Aspirator for Baby | Electric Nose Suction | Model BC024 Instruction Manual

Your guide to safe and effective use

1. INTRODUCTION

Thank you for choosing the GROWNSY Electric Nasal Aspirator. This device is designed to provide gentle and effective relief from nasal congestion for infants and toddlers. Featuring adjustable suction modes, soothing music and lights, and an ergonomic design, it aims to make nasal clearing a comfortable experience for your child. Please read this manual thoroughly before first use and keep it for future reference.

2. PRODUCT OVERVIEW



Figure 1: GROWNSY Electric Nasal Aspirator (Model BC024)

The GROWNSY Nasal Aspirator is a compact, handheld device designed for ease of use and portability. It features a main body with control buttons, a transparent mucus collection chamber, and interchangeable silicone tips. An LED display indicates suction level and battery status.

3. WHAT'S IN THE BOX

Supporting Sustainable Development



ClimatePartner
certified



Global Recycled
Standard

✓ Carbon impact	✓ Safer chemicals
✓ Recycled materials	✓ Worker well-being

This product has sustainability features recognized by trusted certifications.

Figure 2: Package Contents

- Nasal Aspirator Unit × 1
- USB-C Charging Cable × 1
- Storage Box × 1
- Silicone Nasal Tips (different sizes) × 3
- Instruction Manual × 1

4. KEY FEATURES

4.1. Three Suction Modes



Figure 3: Three Adjustable Suction Modes

The aspirator offers three adjustable suction levels (Weak, Moderate, Strong) to effectively clear different consistencies of mucus. Adjust the mode based on your baby's age and congestion level.

4.2. Ergonomic Design

ERGONOMIC BUTTONS & GRIP

Engineered with exclusive ergonomic design, ensuring effortless handling and optimal comfort.



Figure 4: Ergonomic Buttons & Grip

Engineered with an exclusive ergonomic design, the device ensures effortless handling and optimal comfort during use. The scientific grip helps prevent accidental drops.

4.3. Ambient Light & Music Function

Ambient Light & Music Function

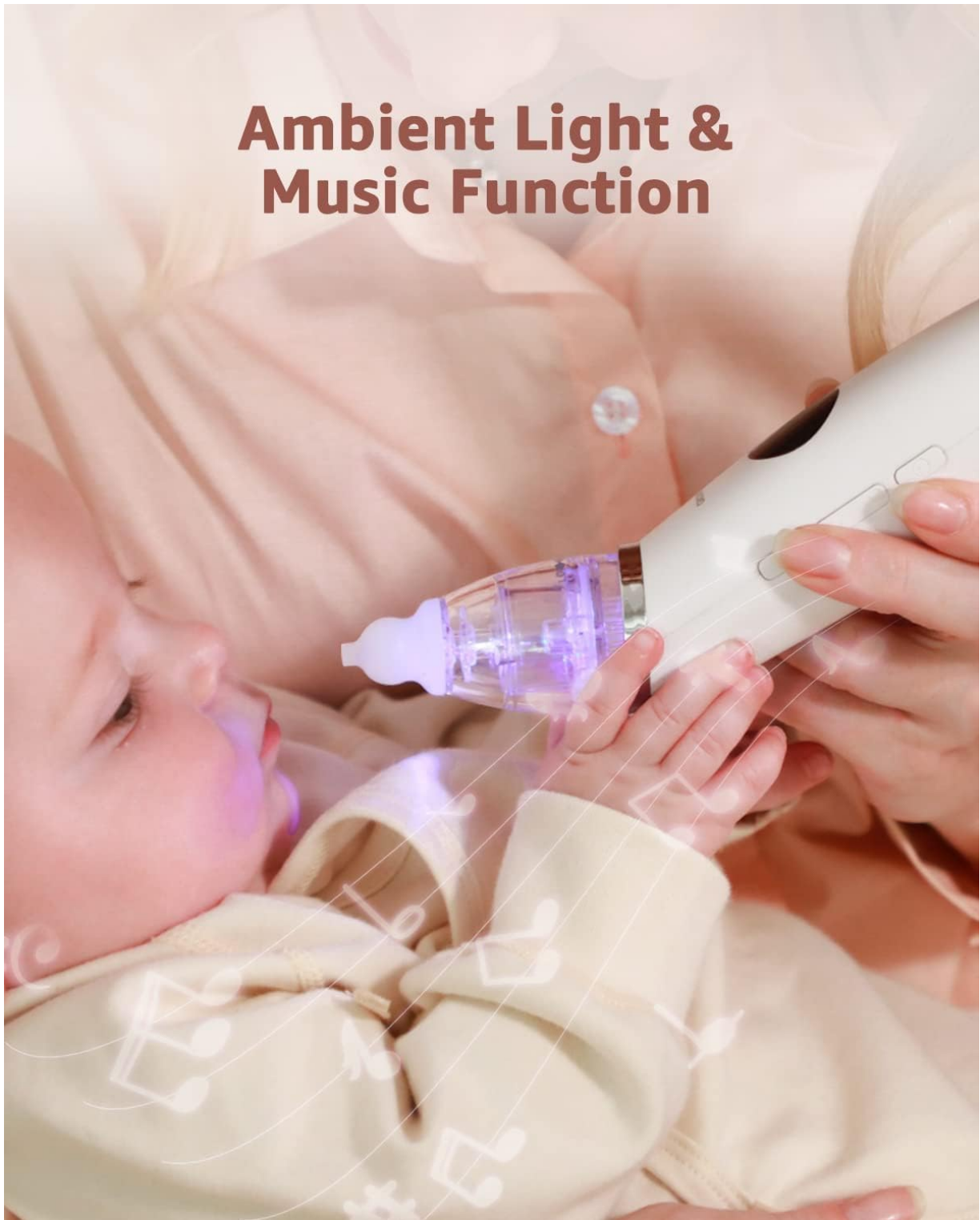


Figure 5: Ambient Light & Music Function

To help soothe and distract your baby during nasal clearing, the aspirator includes an LED display with ambient light and a music function. This can make the experience more pleasant for your child.

4.4. USB Rechargeable

Rechargeable

 FSA or HSA eligible



Figure 6: USB Charging Port

The device is USB rechargeable, offering approximately 120 minutes of use on a full charge. This ensures portability and convenience. Note: The aspirator will not operate while charging.

4.5. Low Noise Operation

LOW NOISE



Figure 7: Low Noise Operation

Designed for quiet operation, the low noise level helps prevent startling your baby, making the suction process calmer.

4.6. Anti-Backflow Design

ANTI-BACKFLOW

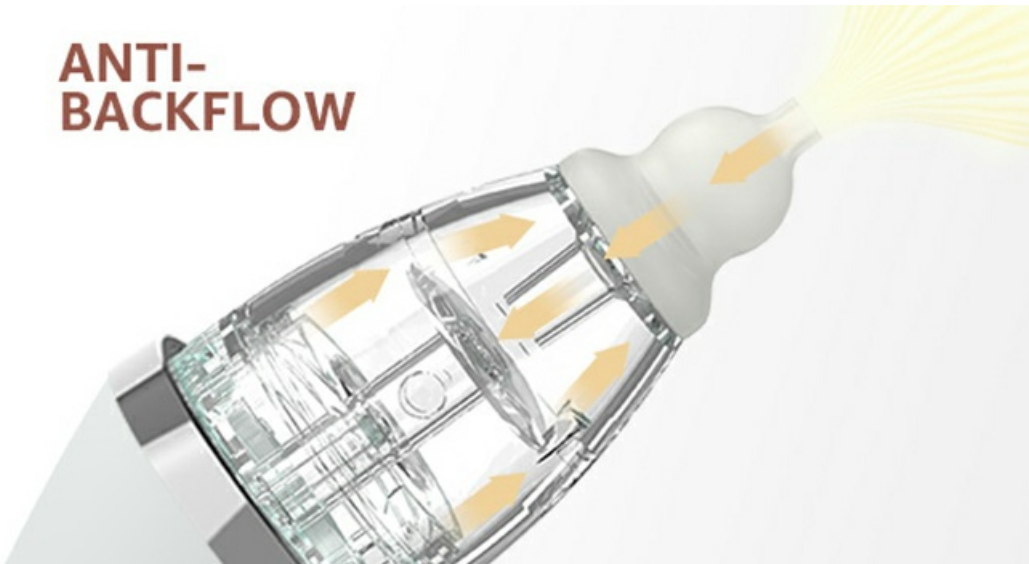


Figure 8: Anti-Backflow Mechanism

A unique anti-backflow design prevents mucus from flowing back into the device or motor, ensuring hygienic operation and prolonging product life.

4.7. BPA-Free Materials

BPA-FREE



Figure 9: BPA-Free Silicone Tips

The nasal tips are made from soft, food-grade silicone and are BPA-free, ensuring ultimate safety and comfort for your baby's delicate nose.

5. SETUP

1. **Charge the Device:** Fully charge the nasal aspirator before its first use.
2. **Select Nasal Tip:** Choose an appropriate soft, food-grade silicone tip based on your baby's age and the size of their nostrils.
3. **Assemble Suction Chamber:** Insert the selected silicone tip onto the suction chamber. Ensure the sealant rubber ring is fitted around the suction chamber base.
4. **Attach Chamber to Main Body:** Twist the assembled suction chamber clockwise onto the main body until it clicks securely into place. Ensure it is tight to prevent air leakage and maintain optimal suction.

6. OPERATING INSTRUCTIONS

1. **Prepare Mucus (if needed):** For stubborn or dry mucus, soften with saline nasal spray or drops. Wait 30

seconds to 1 minute before using the aspirator.

2. **Power On:** Press and hold the power button to turn on the device.
3. **Adjust Suction Mode:** Press the mode button to cycle through the 3 suction levels (Weak, Moderate, Strong) to find the most suitable setting.
4. **Activate Light/Music (Optional):** Press the light button to turn on the ambient light or the music button to play soothing melodies, helping to distract and calm your baby.
5. **Position the Aspirator:** Hold the aspirator in an upright position. Gently place the silicone tip into your baby's nostril, ensuring a snug fit to create a seal.
6. **Begin Suction:** The device will automatically begin suctioning. Move the tip slightly to ensure all mucus is cleared.
7. **Repeat for Other Nostril:** Repeat the process for the other nostril as needed.
8. **Power Off:** Turn off the device immediately after suctioning is complete to prevent mucus backflow into the motor.

7. CLEANING AND MAINTENANCE



Figure 10: Easy-to-Clean Components

1. **Disassemble:** After each use, remove the suction chamber and silicone tip from the main body.
2. **Clean Components:** Wash the reusable nasal suction tips and the collection chamber with warm water and mild soap. Ensure all mucus residue is removed.
3. **Dry Components:** Air dry all cleaned parts completely before reassembling or storing.
4. **Clean Main Body:** Do not wash the main body of the machine. Wipe it with a clean, damp cloth if necessary.
5. **Address Liquid Entry:** If liquid accidentally enters the machine, remove the storage chamber and run the machine for 5-10 minutes to drain any liquid from the motor. Wipe away any stains afterwards. This helps prevent damage and prolong device life.
6. **Battery Maintenance:** To maintain battery health, charge the battery to 100% every two months if the device is not used for an extended period.

8. TROUBLESHOOTING

What's in the box?

1. Nasal aspirator × 1
2. USB-C cable × 1
3. Storage box × 1
4. Nozzles × 3
5. Manual × 1



Figure 11: Troubleshooting Guide

Malfunction	Cause of Issue	Troubleshooting
Weak suction / No suction / No power after period of use	<ol style="list-style-type: none"> 1. The sealing ring has fallen off. 2. The suction chamber is not installed tightly. 3. The suction chamber is full of mucus. 4. The device is not charged. 	<ol style="list-style-type: none"> 1. Re-install the sealing ring. 2. Twist the suction chamber clockwise until it clicks to ensure a tight fit. 3. Clean up mucus in the storage chamber and clean it thoroughly. 4. Charge the device.
Not working / Not turning on after a charge	<ol style="list-style-type: none"> 1. The nasal aspirator can NOT be turned on while charging. 2. The device is not charged. 	<ol style="list-style-type: none"> 1. Disconnect from charging cable before use. 2. Ensure the device is fully charged.
Can't charge	<ol style="list-style-type: none"> 1. Charging voltage over 5V. 2. Charging cable is damaged. 	<ol style="list-style-type: none"> 1. Please ONLY use 5V adapter. Voltage over 5V can NOT charge the device. 2. Try a different charging cable.

9. SPECIFICATIONS

Feature	Detail
Product Dimensions	4.13 x 4.13 x 6.69 inches
Item Weight	11.29 ounces
Model Number	BC024
Brand Name	GROWNSY
Material Type	Silicone (Phthalate Free)
Color	White
Power Source	Battery Powered (USB Rechargeable)
Water Resistance Level	Waterproof (for washable parts)

10. WARRANTY AND SUPPORT

GROWNSY offers a 12-month warranty on this product. For any product issues, please do not hesitate to contact us. We will respond within 24 hours.

Email: service@grownsyofficial.com

SMS/WhatsApp: (833) 524-0498

Thank you again for your custom, and we truly hope you love your new nasal aspirator as much as we do.

Best Regards,

Randell

Founder of GROWNSY