

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [step2bed](#) /

› [step2bed XL Bedside Step Stool with LED Light Guided Rails - User Manual](#)

## step2bed Step2Bed XL

# step2bed XL Bedside Step Stool with LED Light Guided Rails

Model: Step2Bed XL

## INTRODUCTION

The step2bed XL is a robust and adjustable bedside step stool designed to provide assistance for individuals requiring support when getting in and out of bed. It features LED light-guided rails for enhanced visibility and safety, and is constructed from stainless steel to ensure durability and stability. This manual provides essential information for assembly, setup, operation, and maintenance of your step2bed XL.

## SAFETY INFORMATION

- Always ensure all bolts and fasteners are securely tightened before each use.
- Verify that the step stool is stable and does not wobble.
- Do not exceed the maximum weight capacity of 800 lbs.
- Ensure the motion-activated LED light is functioning correctly for nighttime visibility.
- Use the provided self-adhering straps to secure the step stool to the bed frame to prevent movement.
- Keep the step surface clean and dry to prevent slips.
- Consult a healthcare professional if you have specific mobility concerns or conditions.

## PACKAGE CONTENTS

Before assembly, ensure all components are present:

- Main step platform
- Tall grab bar (with foam padding)
- Short grab bar
- Crossbar
- Adjustable legs (4)
- Motion-activated LED light

- Self-adhering Velcro straps (2)
- Assembly hardware (bolts, nuts, washers)
- Instruction manual

## ASSEMBLY INSTRUCTIONS

---

Follow these steps to assemble your step2bed XL. It is recommended to have a second person assist with assembly.

1. Insert the tall and short grab bars into the designated holes on the main step platform.
2. Attach the crossbar between the two grab bars, ensuring proper alignment. Secure with bolts and tighten.
3. Insert the four adjustable legs into the bottom of the main step platform.
4. Turn the step2bed XL on its side to access the leg adjustments.
5. Adjust each leg to the desired height by pressing the pin and sliding the inner tube. Ensure all four legs are set to the same height for stability.
6. Secure the motion-activated LED light to the underside of the step platform using the provided adhesive or brackets, ensuring the sensor faces outwards.
7. Stand the step2bed XL upright and ensure all connections are secure.



Figure 1: Visual guide for step2bed XL assembly. This diagram illustrates the sequence of attaching the grab bars, crossbar, and legs to the main step platform, and positioning the LED light.

## SETUP AND ADJUSTMENT

---

The step2bed XL offers adjustable height settings and can be secured to your bed for added stability.

### Height Adjustment

The step height can be adjusted from 2.5 inches to 7.5 inches. To adjust, press the spring-loaded pin on each leg and slide the inner tube to the desired height. Ensure all four legs are adjusted to the same level to maintain stability and prevent wobbling.

### Securing to Bed

Use the two self-adhering Velcro straps to attach the step2bed XL to your bed frame. Wrap the straps around the bed frame and the step stool's frame, tightening them securely. This prevents the step stool from shifting during use.

# Adjustable Height Step With Non-Slip Tape

Height adjustment settings from 2.5" - 7.5"  
for any size bed. Weight of the step2bed XL is 22 pounds  
and weight capacity is 800 lbs.



Figure 2: Dimensions of the step2bed XL, highlighting the adjustable height range from 2.5 to 7.5 inches and overall width of 32 inches.



Figure 3: Key features of the step2bed XL, including soft foam padding on grab bars, motion-activated LED light, Velcro straps for securing, reinforced steel base, and height adjustment points.

## OPERATING INSTRUCTIONS

The step2bed XL is designed to facilitate safe and independent bed transfers.

### Getting Into Bed

1. Approach the step2bed XL, turning sideways towards the bed.
2. Place both hands on the tall grab bar for support.
3. Step onto the platform with one foot, then bring the other foot up.
4. With hands still on the tall grab bar, sit down on the edge of the bed.
5. Lay down into a resting position, lifting your legs over the short rail and into bed.

### Getting Out of Bed

1. Reach for the closest arm of the tall grab bar.

2. Use the rail to lift your shoulders off the bed while swinging your legs over the short rail and placing your feet on the step platform.
3. Take 15-30 seconds while sitting up to regain your composure.
4. When ready, use one hand on the tall grab bar or two hands (one on the tall rail, one on the short rail) and stand up.
5. While keeping one hand on the tall grab bar, slowly step off the step2bed onto the floor.



Figure 4: Step-by-step visual instructions for safely entering and exiting a bed using the step2bed XL.



Figure 5: The motion-activated LED light provides illumination for the step and surrounding floor area, enhancing safety in low-light conditions.

## MAINTENANCE

- **Cleaning:** Wipe down the stainless steel frame and step platform with a damp cloth and mild detergent. Avoid abrasive cleaners that could damage the finish.
- **Inspection:** Regularly inspect all bolts, nuts, and connections to ensure they are tight. Check the rubber feet for wear and tear.
- **LED Light:** Ensure the motion-activated LED light is clean and unobstructed. Replace batteries as needed.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Step stool wobbles	Uneven leg adjustment; loose bolts; not secured to bed.	Ensure all legs are at the same height. Tighten all assembly bolts. Use Velcro straps to secure to bed frame.
LED light not working	Dead batteries; sensor obstructed; light switch off.	Replace batteries. Clear any obstructions from the sensor. Check if the light has an on/off switch and ensure it is in the 'on' position.
Difficulty getting in/out of bed	Incorrect step height; improper technique.	Adjust step height to a comfortable level. Review operating instructions and practice transfers.

## SPECIFICATIONS

Feature	Detail
Product Dimensions	18"D x 32"W x 32"H (Maximum Height)
Material	Stainless Steel
Color	White
Special Features	Adjustable Height, Durable, Extra Wide, Portable, Motion-Activated LED Light
Adjustable Step Height	2.5" - 7.5"
Weight Limit	800 Pounds
Item Weight	23.7 pounds (10.8 Kilograms)
Assembly Required	Yes
Number of Steps	1
UPC	860250001031

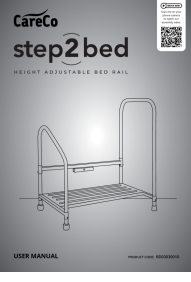

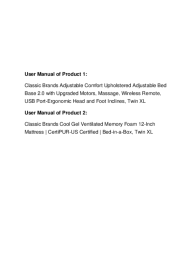


## WARRANTY INFORMATION

The step2bed XL comes with a 30-day warranty on our elderly assistance products. This warranty covers manufacturing defects and ensures the product meets quality standards. For warranty claims or further details, please contact customer support.

## SUPPORT

For any questions, concerns, or assistance with your step2bed XL, please refer to the contact information provided with your purchase or visit the official step2bed website. Our customer service team is available to help with assembly, troubleshooting, and general inquiries.

### Related Documents - Step2Bed XL

	<p><a href="#">CareCo Step2Bed Height Adjustable Bed Rail User Manual</a></p> <p>Comprehensive user manual for the CareCo Step2Bed Height Adjustable Bed Rail, covering setup, usage, maintenance, safety warnings, and warranty information. Learn how to safely install and use this bed rail.</p>
	<p><a href="#">STRUCTUBE LOFT Adjustable Twin XL Bed Base Assembly and Operation Guide</a></p> <p>Comprehensive assembly and operation manual for the STRUCTUBE LOFT Lifestyle Adjustable Twin XL Bed Base. Learn setup, features, remote functions, and troubleshooting. Visit <a href="http://www.structube.com">www.structube.com</a> for the latest version.</p>
	<p><a href="#">Classic Brands Adjustable Comfort Bed Base 2.0 User Manual and Mattress Care Guide</a></p> <p>Comprehensive user manual for the Classic Brands Adjustable Comfort Upholstered Adjustable Bed Base 2.0 with upgraded motors, massage, wireless remote, and USB port. Includes assembly instructions, remote operation, troubleshooting, and warranty information. Also covers care and unwrapping for Classic Brands Cool Gel Ventilated Memory Foam Mattresses.</p>
	<p><a href="#">Motorized Adjustable Bed Base Frame Assembly and Operation Guide</a></p> <p>This guide provides comprehensive instructions for assembling, operating, and troubleshooting the Flash Furniture Motorized Adjustable Bed Base Frame. It includes a parts list, safety precautions, warranty information, and remote control functions.</p>
	<p><a href="#">Renanim Adjustable Bed Owner's Manual - SerenityFlex Model</a></p> <p>Comprehensive owner's manual for the Renanim Adjustable Bed, Model SerenityFlex. Includes safety precautions, parts list, installation guide, remote control functions, Bluetooth app guide, troubleshooting, and warranty information.</p>



### [Ebern Designs BW00354 Twin XL Loft Bed Assembly Instructions](#)

Comprehensive assembly guide for the Ebern Designs BW00354 Twin XL Loft Bed, featuring a slide and ladder. This manual provides detailed parts lists, hardware specifications, and step-by-step instructions with textual descriptions of diagrams.