

## Oticon 200294

# Oticon Charger 1.0 User Manual

Model: 200294 | Brand: Oticon

## INTRODUCTION

The Oticon Charger 1.0 miniRITE R is a desktop charger specifically designed for charging miniRITE R styles of hearing aids. Utilizing inductive technology, it wirelessly charges hearing instruments to full capacity within approximately three hours. This charger is engineered to simplify the daily routine of charging, making it smooth and effortless. It is ideal for users who primarily charge their hearing aids at home. The charger is intended to be permanently connected to a power supply. A green LED indicator on the charger confirms it is ready for use, allowing users to easily place their hearing aids in the charger for overnight charging or quick recharges as needed.

## PRODUCT OVERVIEW



Figure 1: The Oticon Charger 1.0, a compact desktop unit designed for charging Oticon miniRITE R hearing aids.

## SETUP INSTRUCTIONS

Follow these steps to set up your Oticon Charger 1.0:

- 1. Connect to Power:** Plug the charger's power adapter into a standard electrical outlet. The charger is designed to remain connected to the power supply.
- 2. Verify Power Indicator:** A green LED light on the charger will illuminate, indicating that the charger is powered on and ready for use.
- 3. Prepare Hearing Aids:** Ensure your Oticon miniRITE R hearing aids are clean and dry before placing them in the charger.

Your browser does not support the video tag.

Video 1: Demonstrates the setup and initial use of the Oticon Charger 1.0, including connecting it to power and placing hearing aids for charging.

## OPERATING INSTRUCTIONS

To charge your hearing aids:

- Place Hearing Aids:** Gently insert each hearing aid into its designated slot in the charger. Ensure they are seated correctly.
- Charging Indication:** The hearing aids will begin charging automatically. The indicator lights on the hearing aids themselves will show their charging status (refer to your hearing aid manual for specific light patterns).
- Full Charge:** A full charge typically takes approximately 3 hours. Once fully charged, the hearing aids are ready for use.
- Removing Hearing Aids:** Carefully lift the hearing aids from their slots when charging is complete or when you need to use them.

## MAINTENANCE

To ensure the longevity and optimal performance of your Oticon Charger 1.0, follow these maintenance guidelines:

- Cleaning:** Regularly wipe the charger with a soft, dry cloth. Do not use liquids or abrasive cleaners.
- Storage:** Store the charger in a cool, dry place away from direct sunlight and extreme temperatures.
- Avoid Moisture:** Keep the charger away from water and other liquids. If it comes into contact with moisture, disconnect it immediately and allow it to dry completely before re-using.
- Handle with Care:** Avoid dropping the charger or subjecting it to strong impacts.

## TROUBLESHOOTING

If you encounter issues with your Oticon Charger 1.0, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Charger LED is not green.	Not connected to power or power outlet issue.	Ensure the power adapter is securely plugged into both the charger and a working electrical outlet. Try a different outlet.
Hearing aids are not charging.	Hearing aids not properly seated or charger not powered.	Re-seat the hearing aids firmly in their slots. Check if the charger's green LED is illuminated. Ensure hearing aids are compatible with this charger model.
Charging takes longer than 3 hours.	Battery degradation in hearing aids or environmental factors.	Ensure the charging environment is within recommended temperature ranges. If the issue persists, consult your hearing aid professional as the hearing aid batteries may need service.

If the problem persists after trying these solutions, please contact Oticon customer support or your hearing care professional.

## SPECIFICATIONS

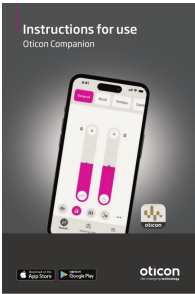
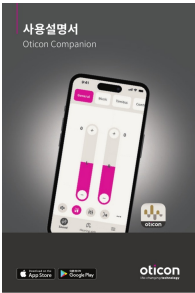
Feature	Detail
Model Number	200294
Product Dimensions	19 x 6 x 14 inches


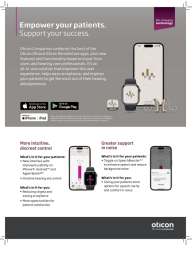


Feature	Detail
Item Weight	10.2 ounces
Brand	Oticon
Connectivity Technology	USB
Connector Type	USB
Compatible Devices	Oticon More, OPN & Ruby Hearing aids, Mini RITE R Hearing Aids, Bernafon, Philips, Sonic Brand Rite Model Hearing Aids
Included Components	Power supply
Special Features	Travel, Magnetic, Light, Fast Charging
Mounting Type	Wall Mount
Amperage	0.5 Amps
Total USB Ports	1

WARRANTY AND SUPPORT

For information regarding warranty coverage, please refer to the documentation provided with your original purchase or visit the official Oticon website. For technical support, troubleshooting assistance, or to inquire about replacement parts, please contact Oticon customer service or your authorized hearing care professional. You can find more information and support resources on the [Oticon official website](#).

Related Documents - 200294

	<p><a href="#">Oticon Companion App: User Guide for Hearing Aid Control and Remote Care</a></p> <p>Explore the Oticon Companion app's features for managing your hearing aids, including volume control, program selection, RemoteCare visits, HearingFitness tracking, and troubleshooting. Compatible with iOS and Android devices.</p>
	<p><a href="#">Oticon Companion</a></p> <p>Oticon Companion . , RemoteCare , Tinnitus SoundSupport Oticon .</p>

	<p><a href="#">Oticon Genie 2 2021.2: Fitting Software Instructions for Hearing Care Professionals</a></p> <p>Comprehensive instructions for using Oticon Genie 2 2021.2 fitting software. Learn about installation, system requirements, fitting steps, tools like In-situ Audiometry and RemoteCare, and important warnings for hearing care professionals.</p>
	<p><a href="#">Oticon Companion App: Empowering Patients and Supporting Hearing Care Professionals</a></p> <p>Discover the Oticon Companion app, designed to enhance the hearing aid experience for patients with intuitive control, noise management, personalization, and remote support. Learn about its benefits for both users and hearing care professionals.</p>
	<p><a href="#">Oticon Charger 1.0: Instructions for Use</a></p> <p>Comprehensive instructions for using and maintaining the Oticon Charger 1.0, compatible with various Oticon hearing aid models including Real, More, Play PX, Zircon, Jet PX, Opn S, Opn Play, Ruby, and CROS PX. Learn about charging, troubleshooting, and safety.</p>
	<p><a href="#">1.0</a></p> <p>1.0 (miniRITE R, miniBTE R )</p>