

[Manuals.plus](#) /

› [Audinell](#) /

› Audinell PerfectDry LUX Hearing Aid Dryer Instruction Manual

## Audinell PerfectDry LUX

# Audinell PerfectDry LUX Hearing Aid Dryer Instruction Manual

Model: PerfectDry LUX

## PRODUCT OVERVIEW

---

The PerfectDry LUX is a compact electronic dehumidifier designed to remove sweat and moisture from hearing aids, wireless earbuds, and ear amplifiers. It features automatic temperature regulation, pulsed hot air dehumidification, controlled ventilation, and a 360° UV-C lamp for sanitization. This device helps prolong the life and improve the performance of your listening devices.



Figure 1: PerfectDry LUX Hearing Aid Dryer

## SETUP GUIDE

---

### Unboxing and Initial Inspection

Carefully remove the PerfectDry LUX from its packaging. Ensure all components, including the device and power cable, are present and undamaged.

### Power Connection

Connect the provided USB cable to the PerfectDry LUX and then to a suitable power adapter or USB port.



Figure 2: PerfectDry LUX with USB cable and power adapter

## Preparing Your Device

Before placing your hearing aids or other devices into the PerfectDry LUX, ensure you remove any batteries. If your device has a removable dome, it is recommended to remove it for thorough drying and cleaning.

## OPERATING INSTRUCTIONS

---

### Placing Devices

Open the lid of the PerfectDry LUX. Place your hearing aids, earbuds, or ear amplifiers inside the drying chamber. The device is designed to accommodate various types and brands.



Figure 3: Open PerfectDry LUX drying chamber

### **Starting a Drying Cycle**

Close the lid. Press the power button located on the front of the device once. A blue light will illuminate, indicating that the drying and sanitization cycle has begun.

Video 1: Visual overview of the PerfectDry LUX features and operation, including UV-C and heat drying.

### **Cycle Duration**

The PerfectDry LUX completes a full drying and sanitization cycle in approximately 45 minutes. The device will automatically shut off at the end of the cycle.

The UV-C rays are generated by a short wavelength lamp whose 360 ° annular arrangement makes it possible to reach all the surfaces of the hearing aid.

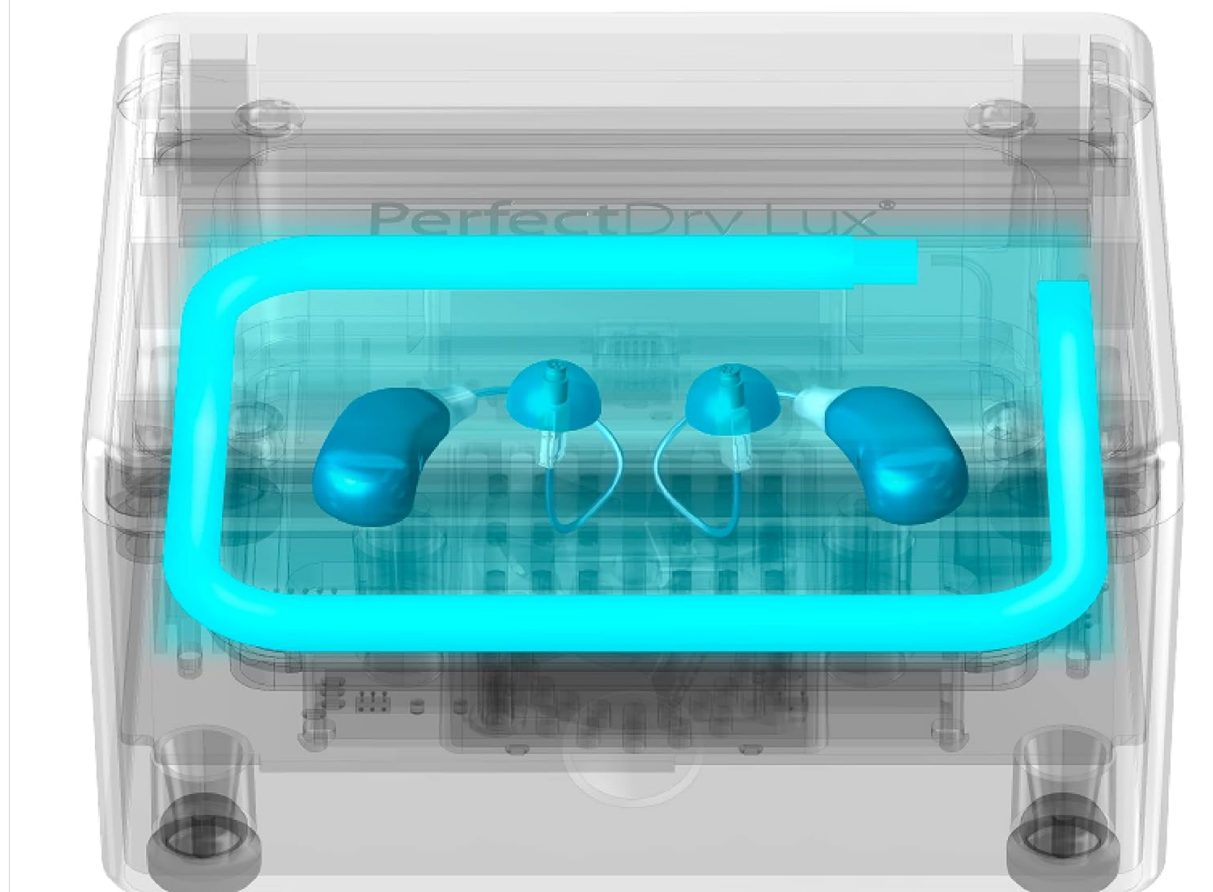


Figure 4: Diagram showing 360° UV-C light coverage inside the PerfectDry LUX

- Pulsed hot air dehumidification system for an optimized, even flow.
- Automatic temperature regulation ( $<35^{\circ}\text{C} \pm 3^{\circ}\text{C}$ ).
- Controlled ventilation system.

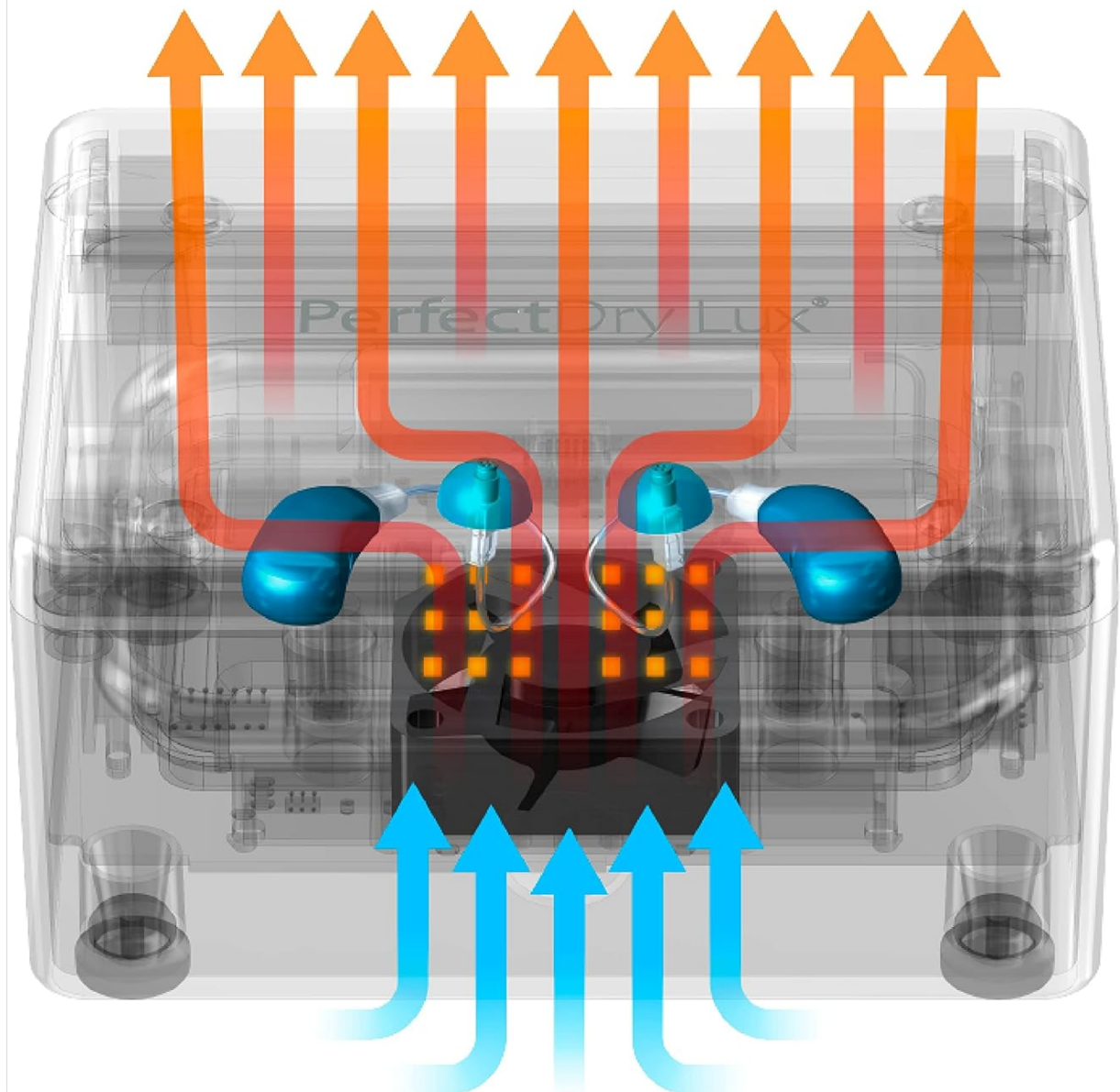


Figure 5: Diagram illustrating pulsed hot air dehumidification within the PerfectDry LUX

## CARE AND MAINTENANCE

### Cleaning the Device

Periodically wipe the exterior of the PerfectDry LUX with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.

### Internal Cleaning

Ensure the drying chamber is free of debris. If necessary, gently wipe the interior with a soft, slightly damp cloth.

## PRODUCT SPECIFICATIONS

- **Brand:** Audinell

- **Model Name:** PerfectDry LUX
- **Product Dimensions:** 3"D x 3.6"W x 1.9"H
- **Item Weight:** 0.52 Pounds (8.3 ounces)
- **Material:** Plastic, Glass
- **Operation Mode:** Automatic
- **Upper Temperature Rating:** 95 Degrees Fahrenheit (35°C)
- **Special Feature:** Portable
- **UPC:** 855347004378
- **Manufacturer:** MG Development



Figure 6: Dimensions of the PerfectDry LUX device

## TROUBLESHOOTING

### Device Not Turning On

- Ensure the USB cable is securely connected to both the PerfectDry LUX and the power source.
- Verify that the power source is active.

### Drying Cycle Not Starting

- Make sure the lid is fully closed before pressing the power button. The device will not operate if the lid is open.

### Ineffective Drying

- Ensure batteries are removed from hearing aids before placing them in the dryer.
- For optimal results, remove any domes or ear tips from hearing aids to allow for better air circulation.

- Avoid overloading the drying chamber.

## WARRANTY AND CUSTOMER SUPPORT

---

The PerfectDry LUX comes with an industry-leading 3-year warranty. For warranty claims or technical assistance, please contact Audinell customer support.

**Manufacturer:** MG Development

**Brand Store:** [Visit the Audinell Store](#)