

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

Manuals.plus /

- › [Braun](#) /
- › Braun Manual Nasal Aspirator Filters (Model BNF040US) - Instruction Manual

Braun BNF040US

Braun Manual Nasal Aspirator Filters (Model BNF040US) - Instruction Manual

Brand: Braun | Model: BNF040US

1. INTRODUCTION

Thank you for choosing Braun Manual Nasal Aspirator Filters. These filters are designed to ensure hygienic and effective use of your Braun Manual Nasal Aspirator (BNA050), sold separately. This manual provides instructions for the proper installation, use, and maintenance of the filters and the aspirator system.

2. PRODUCT OVERVIEW

The Braun Manual Nasal Aspirator Filters are disposable and hygienic, designed to prevent mucus from entering the tube during aspiration, ensuring a clean and safe experience for both user and child. Each package contains 40 BPA-free filters.

BRAUN

Nasal aspirator filters

Use only with Braun manual nasal aspirator

40 Count package

Single use filters

Hygienic & safe



Image: Braun Nasal Aspirator Filters 40 Count Package. This image displays the product packaging for the Braun Nasal Aspirator Filters, highlighting that it is a 40-count package of single-use, hygienic filters for use with the Braun manual nasal aspirator.

3. SETUP AND FILTER INSTALLATION

The filters are an essential component for the hygienic operation of your Braun Manual Nasal Aspirator. Always ensure a new filter is installed before each use.

3.1 Installing a New Filter

1. Remove the nasal aspirator from its packaging.
2. If a filter is not already installed, or if replacing an old one, unscrew the top and bottom chambers of the aspirator.
3. Insert a new filter into the designated holding area within the mucus chamber.
4. Reassemble the aspirator by screwing the chambers back together securely.

Easy-to-change filters

Simply unscrew the top and bottom chambers, insert a new filter, and reassemble.

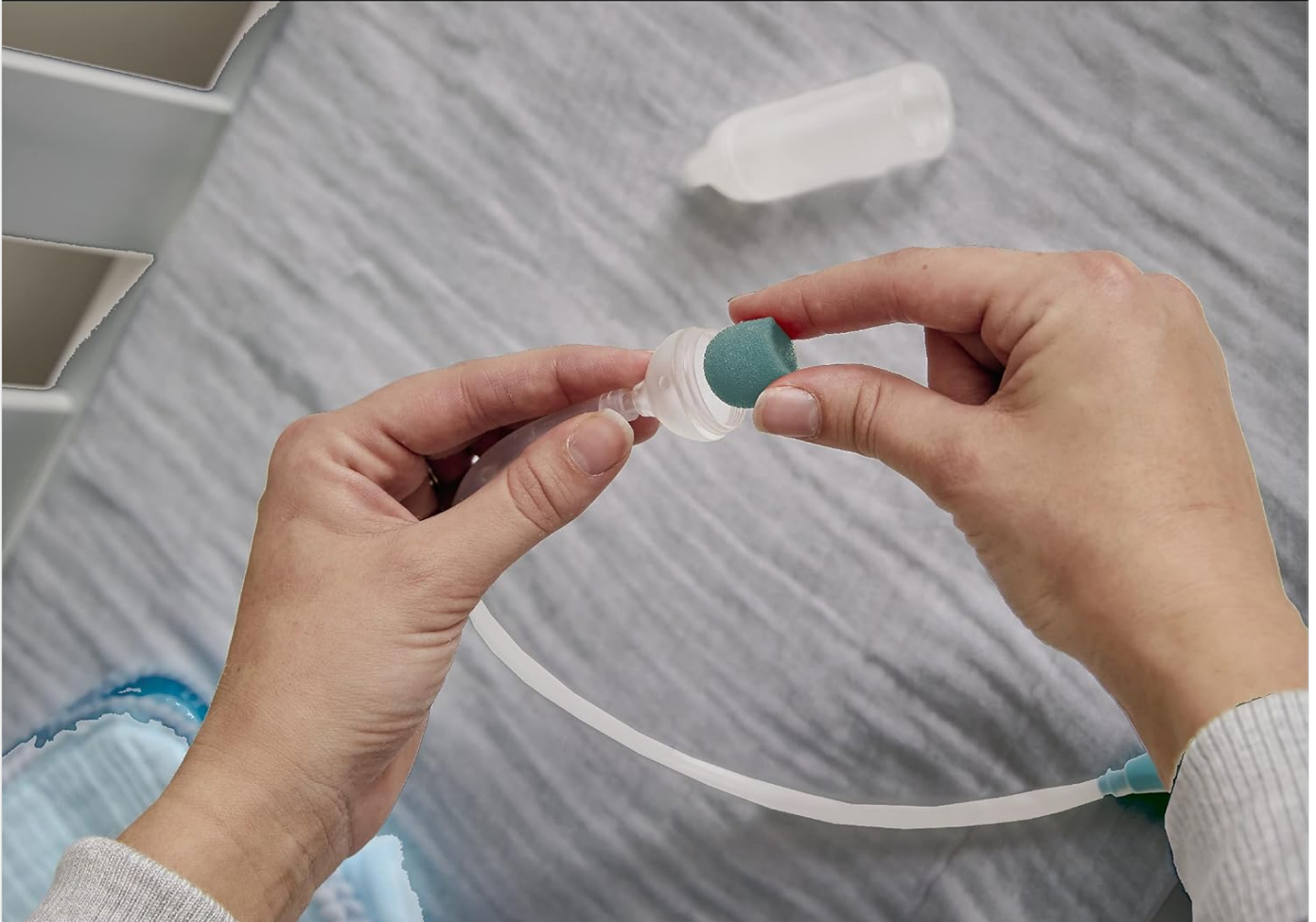


Image: Easy-to-change filters. This image shows hands unscrewing the top and bottom chambers of the nasal aspirator to insert a new filter, illustrating the simple filter replacement process.

3.2 Getting Started with the Aspirator

For a visual guide on assembling and getting started with your Braun Manual Nasal Aspirator, please watch the video below:

Video: Braun BNA050 Manual Nasal Aspirator - Getting Started. This video demonstrates the initial setup and assembly of the Braun Manual Nasal Aspirator, including how to prepare it for use.

4. OPERATING INSTRUCTIONS

The Braun Manual Nasal Aspirator, used with these filters, is designed to gently clear stuffy noses in infants and toddlers.

1. Ensure a new filter is correctly installed in the aspirator.
2. Select the correct nose tip for your child to ensure a proper seal.
3. Place the soft tip of the aspirator gently into one nostril of the child.

4. Place the mouthpiece in your mouth and apply gentle, continuous suction. The mucus will be pulled from the child's nose into the device and captured by the filter in the mucus chamber.
5. Remove the nasal aspirator from the child's nostril.
6. Repeat for the other nostril if necessary.



Image: Mother using Braun manual nasal aspirator on baby. This image shows a mother gently using the Braun manual nasal aspirator to clear her baby's stuffy nose, demonstrating the product in action.



Image: Baby with Braun nasal aspirator, showing filter function. This image illustrates how the filters ensure hygienic and effective use of the Braun nasal aspirator by capturing mucus.

5. MAINTENANCE AND CLEANING

Proper cleaning of your nasal aspirator after each use is crucial for hygiene and longevity. The filters are disposable and should be discarded after every use.

5.1 Disposing of Used Filters

After each use, disassemble the aspirator and gently shake the used filter out of the filter holder. Dispose of the used filter in a waste bin. Do not reuse filters.

5.2 Cleaning the Aspirator Components

The aspirator components (mucus chamber, nose tip, mouthpiece, and tube) can be cleaned by hand washing in hot water or by placing them in the top rack of a dishwasher. Ensure all parts are thoroughly cleaned and dried before reassembly.

Video: Braun BNA050 Manual Nasal Aspirator - How to Clean. This video provides detailed instructions on how to disassemble and clean all components of the Braun Manual Nasal Aspirator after use.

6. SAFETY INFORMATION

- Always use a new, clean filter for each aspiration session to maintain hygiene and prevent cross-contamination.
- Ensure the aspirator is fully assembled and the filter is correctly seated before use.
- Do not force the nose tip into the child's nostril. Apply gentle pressure only.
- Keep all components out of reach of children when not in use.
- This product is intended for clearing nasal congestion only. Do not use for other purposes.
- If your child experiences discomfort or if congestion persists, consult a healthcare professional.

7. SPECIFICATIONS

Feature	Detail
Model Number	BNF040US
Unit Count	40 Count
Compatibility	Braun Manual Nasal Aspirator (BNA050)
Material Type Free	Bisphenol-A (BPA) Free
Manufacturer	Kaz, Inc
Recommended Age	0 - 8 years (for aspirator use)
UPC	328785008005

8. WARRANTY AND SUPPORT

For warranty information, product registration, or further support, please refer to the documentation included with your Braun Manual Nasal Aspirator (BNA050) or visit the official Braun Healthcare website. Keep your purchase receipt for warranty claims.