

Mighty Well PICCPerfect PICC Line Cover

Mighty Well PICCPerfect PICC Line Cover Instruction Manual

Model: PICCPerfect PICC Line Cover

1. PRODUCT OVERVIEW

The Mighty Well PICCPerfect PICC Line Cover is designed to provide comfort, security, and discretion for individuals using a Peripherally Inserted Central Catheter (PICC) line, midline, or other medical dressings on the arm. This cover features a double-layer design and is made from advanced fabric technology.



Image 1.1: The Mighty Well PICCPPerfect PICC Line Cover in Olive, showcasing the brand logo.

2. KEY FEATURES

- **Double-Layer Security:** The cover features a double-layer design to securely cover PICC lines, midlines, casts, wounds, and dressings. This design helps keep tubing away from the skin.
- **Dual Access Points:** Equipped with two access points, allowing easy access to single and double

lumens for various placement heights on the arm.

- **Advanced Fabric Technology:** Made with moisture-wicking and anti-odor materials for enhanced comfort. The fabric is also antimicrobial.
- **Secure Fit:** Bonded anti-slip edges are incorporated to prevent the cover from sliding down the arm.
- **Designed by a PICC Line User:** Developed with direct experience to address common challenges faced by PICC line users.

PICCPerfect vs. Others		
	PICCPerfect	others
Double-layer of security	✓	✗
Designed by a PICC Line Patient	✓	✗
Dual Access Points for Versatility	✓	✗
Antimicrobial, Moisture-Wicking, Anti-Odor	✓	✗
Wide Range of Sizes	✓	✗
MIGHTYWELL®		

Image 2.1: A comparison highlighting the unique features of the PICCPerfect cover, including double-layer security, dual access points, antimicrobial properties, and a wide range of sizes.

Your browser does not support the video tag.

Video 2.1: An overview of the PICCPerfect Antimicrobial PICC Line Cover by Mighty Well, demonstrating its moisture-wicking, antimicrobial, and breathable fabric, along with its unique fold-up design.

3. SIZING GUIDE


To ensure a proper fit, measure the circumference of your bicep where your PICC line is placed. Keep your arm relaxed and straight down your side during measurement. If your measurement falls between two sizes, select the smaller size for a snug fit.

PICCPerfect® size guide

measure around your bicep where your PICC line is placed

see the guide below to select your size.

if you are in between two sizes, select the smaller size.

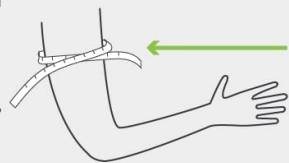

PICCPerfect®	range (in):	range (cm):
XXXS	6" - 7"	15 CM - 17.5 CM
XXS	7 - 8"	17.5 CM - 20 CM
XS	8" - 9"	20 CM - 23 CM
S	9 - 10.5"	23 CM - 26.5 CM
M	10.5" - 12.5"	26.5 CM - 32 CM
L	12.5" - 17.5"	32 CM - 37 CM
XL	17.5" - 20"	37 CM - 46 CM

MIGHTYWELL®

Image 3.1: Detailed sizing chart for the PICCPerfect cover, showing recommended bicep circumference ranges for sizes XXXS to XL in both inches and centimeters.

PICCPerfect Sizing Chart

Size	Range (in)	Range (cm)
XXXS	6" - 7"	15 CM - 17.5 CM
XXS	7" - 8"	17.5 CM - 20 CM
XS	8" - 9"	20 CM - 23 CM
S	9" - 10.5"	23 CM - 26.5 CM
M	10.5" - 12.5"	26.5 CM - 32 CM
L	12.5" - 17.5"	32 CM - 37 CM
XL	17.5" - 20"	37 CM - 46 CM

4. USAGE INSTRUCTIONS

Follow these steps to properly put on and use your PICCPerfect PICC Line Cover:

1. **Step 1: Measure for Size.** Before applying, ensure you have the correct size by measuring your bicep circumference as described in the Sizing Guide.
2. **Step 2: Prepare the Cover.** Open the PICCPerfect cover so the inner layer is visible, with the coated slits positioned on top.
3. **Step 3: Position on Arm.** Gently pull the cover up your arm, aligning the slits with your PICC line insertion site.
4. **Step 4: Thread Lumens.** Thread the lumen(s) through one of the slits. Ensure the external line is completely through and not touching the skin.
5. **Step 5: Secure Tubing (if applicable).** If your PICC has extension tubing, carefully wrap the tubing around the arm, on top of the PICCPerfect's inner layer.
6. **Step 6: Create Pocket.** Fold the outer layer up to create a pocket, which helps secure the PICC and prevents accidental pulling.

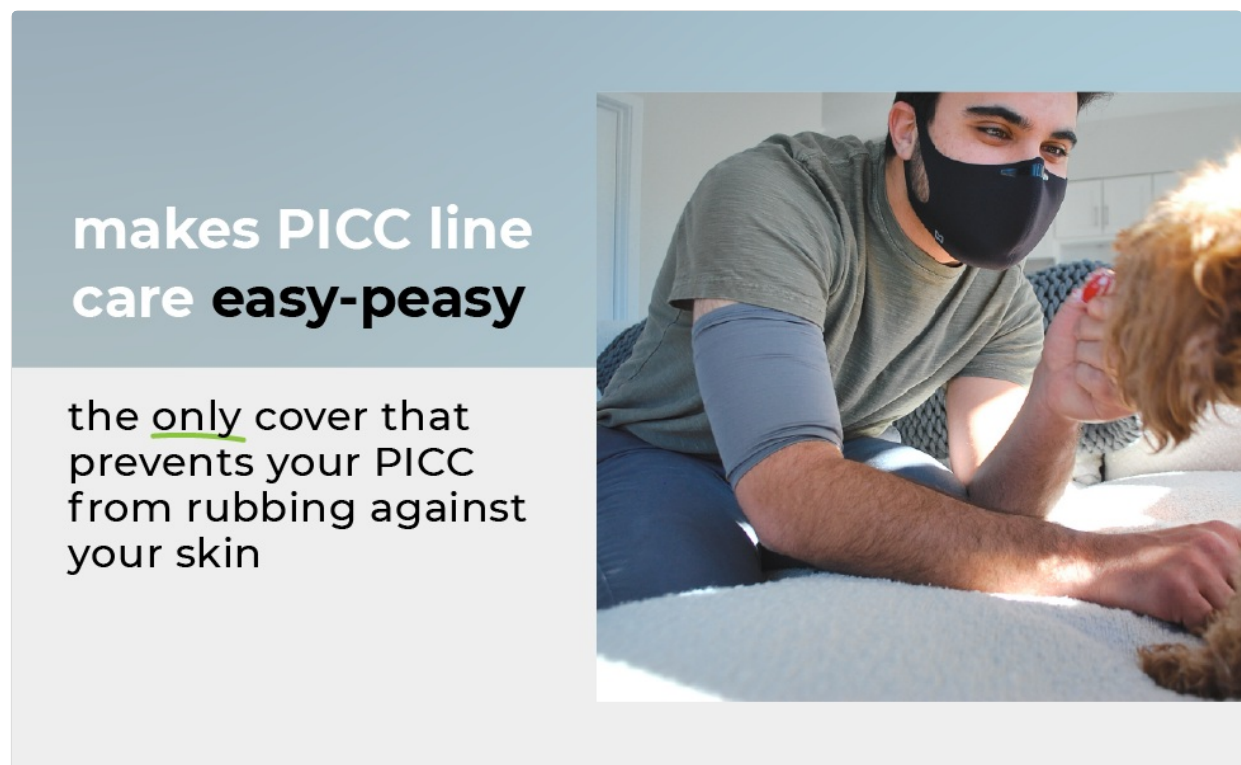


Image 4.1: A visual step-by-step guide demonstrating the application of the PICCPerfect PICC Line Cover, from measuring to securing the PICC line.

Your browser does not support the video tag.

Video 4.1: A demonstration of the PICCPerfect PICC Line Cover's features, including its easy application, dual-layer design, double access points, antimicrobial and moisture-wicking properties, and bonded anti-slip edges.

5. CARE AND MAINTENANCE

To maintain the quality and effectiveness of your PICCPerfect PICC Line Cover, follow these care instructions:

- **Washing:** The cover is machine washable. Use a gentle cycle with cold water.
- **Drying:** Air dry or tumble dry on low heat. Avoid high heat to preserve fabric integrity and elasticity.
- **Avoid Bleach:** Do not use bleach or harsh detergents, as these can damage the fabric's antimicrobial

and moisture-wicking properties.

- **Storage:** Store in a clean, dry place when not in use.

6. SPECIFICATIONS

- **Product Dimensions:** 6.3 x 0.2 x 4.33 inches
- **Item Weight:** 0.63 ounces
- **Material:** Advanced fabric with moisture-wicking, anti-odor, and antimicrobial properties.
- **Design:** Double-layer with dual access points and bonded anti-slip edges.
- **Manufacturer:** Mighty Well

7. TROUBLESHOOTING

This product is a non-electronic garment. Common issues are typically related to fit or wear. If you experience any of the following, consider these solutions:

- **Cover slipping down:** Ensure you have selected the correct size based on the sizing guide. If between sizes, opt for the smaller size. The bonded anti-slip edges are designed to prevent this.
- **Discomfort or irritation:** Verify that the PICC line and any tubing are properly positioned within the cover and not directly against the skin. Ensure the cover is clean and dry.
- **Difficulty accessing lumens:** Utilize the dual access points by gently adjusting the cover to expose the necessary lumen(s).

For persistent issues or concerns not addressed here, please contact Mighty Well customer support.

8. WARRANTY AND SUPPORT

For information regarding product warranty, returns, or customer support, please refer to the official Mighty Well website or contact their customer service directly. Amazon's return policy for this item allows returns within a specified period from the date of purchase. Please check your purchase details for specific return window information.

Manufacturer Contact: Mighty Well

Website: [Visit the Mighty Well Store on Amazon](#)