

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [AURAGLOW](#) /
- › [Auraglow Teeth Whitening Trays Instruction Manual](#)

AURAGLOW Teeth Whitening Trays

Auraglow Teeth Whitening Trays Instruction Manual

Brand: AURAGLOW

PRODUCT OVERVIEW

Auraglow teeth whitening trays offer an effective and comfortable alternative to custom dental whitening trays. Crafted from a soft thermoform material, these trays can be precisely molded to your teeth by dipping them in warm water. Featuring a 20% slimmer design, they provide a comfortable fit without the bulk often associated with other trays. These trays are designed for optimal use with Auraglow teeth whitening gel (sold separately) for convenient at-home whitening.

Auraglow



20%
SLIMMER
DESIGN

Teeth Whitening Trays

Moldable slim fit trays for teeth whitening at home.

4 MOLDABLE TRAYS AND STORAGE CASE

Image: Auraglow Teeth Whitening Trays product packaging.

Custom Moldable Teeth Whitening Trays

- ✓ 4 moldable teeth whitening trays
- ✓ Slim design for a better fit
- ✓ Great for teeth whitening
- ✓ Not bulky or thick fitting
- ✓ Easy to mold to teeth
- ✓ Great alternative to expensive dental trays



Image: Features of Auraglow Custom Moldable Teeth Whitening Trays.

Easy to mold teeth whitening trays

Designed for comfort and a close fit like dental trays.

- ✓ Reusable trays
- ✓ EVA material
- ✓ BPA and latex free
- ✓ No-slip tray design
- ✓ Upper and lower tray
- ✓ Use with Auraglow gel
(not included)



*Final molded
tray example*

Image: Molded trays and key material features.

WHAT'S INCLUDED


Your Auraglow Teeth Whitening Trays package includes everything you need to get started with your at-home whitening routine:

- **4 Thermoform Trays:** Made from soft, moldable material for a custom fit.
- **1 Storage Case:** For hygienic storage of your molded trays.
- **Detailed Instructions:** A guide for molding and using your trays effectively.

What's included



1 Retainer Case

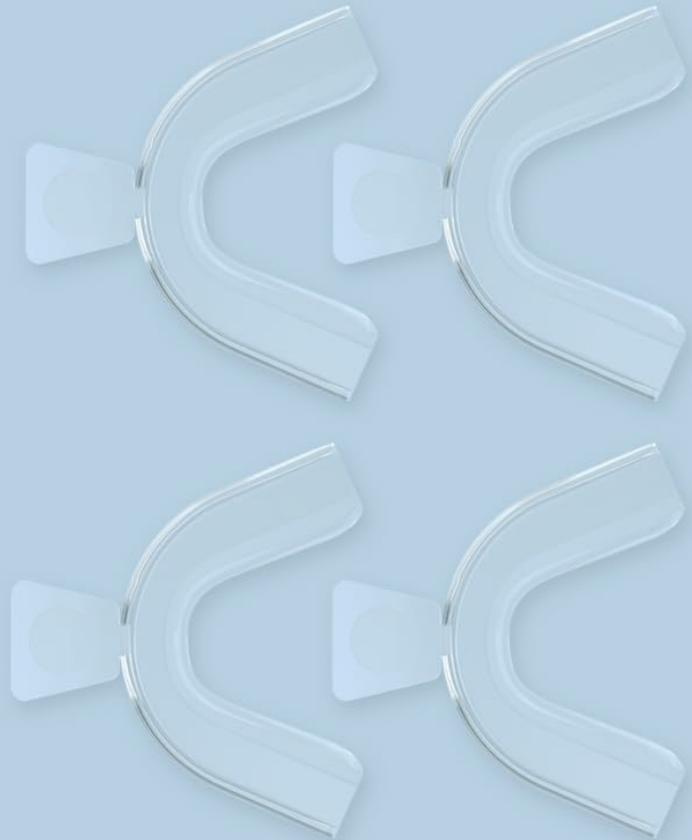

4 Slim Moldable Trays

Image: Contents of the package.

IMPORTANT SAFETY INFORMATION

Please read and adhere to the following safety guidelines before using the Auraglow Teeth Whitening Trays:

- Always use these trays with Auraglow teeth whitening gel (not included).
- Do not use for children. This product is intended for adult use only.
- Keep out of reach of children.
- Trays are BPA free, latex free, and silicone free.

SETUP: MOLDING YOUR AURAGLOW TEETH WHITENING TRAYS

Achieving a custom fit is crucial for effective whitening. Follow these steps carefully to mold your trays:

Your browser does not support the video tag.


Video: Official Auraglow Teeth Whitening Trays Molding Instructions. This video demonstrates the step-by-step process of molding the trays for a custom fit.

Step-by-Step Molding Instructions:

Instructions for molding



Step 1

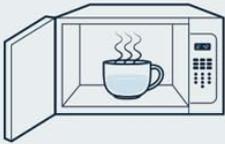


Bring a cup of water to a low boil, deep enough to fully submerge the trays.

Step 2



Hold the first tray by the tab and dip it into the hot water for approximately 3-4 seconds.

Step 3

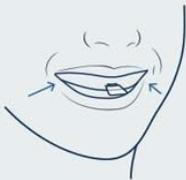


Remove the tray from the water, place it in your mouth and mold with your teeth. Repeat step 2 for the other tray.

Step 4



Try both trays on together to see how they fit. If necessary, dip in the water to mold again or use the extra set of trays provided.

Step 5

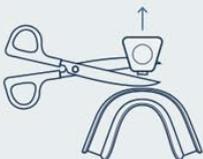

Trim off the front tabs. You can also trim the tray area if there is excess material for a better fit.

Step 6



Fill the tray with Auraglow whitening gel and wear based on the instructions.

Image: Visual guide for molding instructions.

- 1. Prepare Hot Water:** Bring a cup of water to a low boil. Ensure the cup is deep enough to fully submerge the trays.
- 2. Soften the Tray:** Hold the first tray by the tab and dip it into the hot water for approximately 3-4 seconds. The tray will begin to wilt and become soft. *Tip: The edges of the trays should not fold over onto each other.*
- 3. Mold to Your Teeth:** Remove the softened tray from the water. Carefully place the tray over your upper teeth. Bite gently onto the tray and suck the air out by puckering your lips tightly together to create a vacuum and mold it to your teeth.
- 4. Repeat for Other Tray & Test Fit:** Repeat steps 2 and 3 for the other tray (for your lower teeth). Once both are molded, try both trays together to see how they fit. If the fit isn't quite right, you can re-dip them in hot water and repeat the fitting process. You can also use the extra set of trays provided if needed.

5. **Trim for Comfort:** When you are satisfied with the fitment, trim off the front tabs of the trays with scissors. You can also trim any excess material from the tray area for a better, more comfortable fit.

OPERATING: USING YOUR AURAGLOW TEETH WHITENING TRAYS

Once your trays are molded and trimmed, they are ready for use with your Auraglow teeth whitening gel (sold separately).

1. **Apply Whitening Gel:** Apply a small, thin line of Auraglow teeth whitening gel into the front surface of the molded tray, ensuring it covers the area that will contact your teeth. Do not overfill.
2. **Insert Trays:** Carefully place the trays onto your upper and lower teeth.
3. **Wear Time:** Refer to the instructions provided with your Auraglow teeth whitening gel for recommended wear time.
4. **After Use:** After each use, remove the trays and rinse them thoroughly with cool water.

Use with Auraglow whitening gel (Sold separately)



Image: Auraglow whitening gel (sold separately).

CARE AND MAINTENANCE

Proper care will extend the life of your Auraglow Teeth Whitening Trays:

- **Cleaning:** After each use, rinse your trays with cool water to remove any residual gel. You can gently brush them with a soft toothbrush and water if needed. Avoid hot water as it may alter the tray's shape.
- **Storage:** Store your clean, dry trays in the provided storage case to protect them from dust and damage.
- **Avoid Harsh Chemicals:** Do not use abrasive cleaners or harsh chemicals on your trays, as this can damage the material.

TROUBLESHOOTING

Here are some common issues you might encounter during the molding process and how to resolve them:

Issue	Solution
Tray feels hard when molding.	The tray was not held in the hot water long enough. Try holding it in the hot water for a slightly longer duration (e.g., 5-6 seconds) until it becomes soft and pliable.
Tray shrunk or folded onto itself.	The tray was held in the hot water for too long. This can cause the material to over-soften and deform. Use one of the extra trays and try again, reducing the time in hot water.
The fit isn't quite right after molding.	You can re-dip the tray in hot water for a few seconds to soften it again and repeat the molding process. Alternatively, use one of the two extra trays provided to start fresh.

PRODUCT SPECIFICATIONS

- **Brand:** AURAGLOW
- **Product Benefits:** Whitening
- **Item Form:** Gel (Note: Trays are for use with gel, not gel themselves)
- **Number of Items:** 1 (package containing 4 trays and case)
- **Material Type Free:** BPA Free, Latex Free, Silicone Free
- **Package Dimensions:** 3.5 x 3.35 x 2.24 inches
- **Item Weight:** 2.4 ounces
- **Manufacturer:** Auraglow
- **ASIN:** B09KDYZNFV

WARRANTY INFORMATION

No specific warranty information is provided with this product. For details regarding potential warranties or guarantees, please refer to the manufacturer's official website or contact their customer support directly.

CUSTOMER SUPPORT

For further assistance, product inquiries, or to explore other Auraglow products, please visit the official AURAGLOW store:

Visit the AURAGLOW
Store

