



[Manuals.plus](#) /

- › [Besst Survivor](#) /
- › Besst Survivor Comprehensive First Aid Kit (324 Piece Set Max) User Manual

## Besst Survivor B09KMTLS9X

# Besst Survivor Comprehensive First Aid Kit (324 Piece Set Max) User Manual

Model: B09KMTLS9X

## INTRODUCTION

---

This manual provides essential information for the proper use and maintenance of your Besst Survivor Comprehensive First Aid Kit (324 Piece Set Max). This kit is designed to assist in managing minor injuries and providing initial response to more serious emergencies in various settings, including home, car, camping, and outdoor activities. Familiarizing yourself with its contents and instructions is crucial for effective emergency preparedness.



- **Additional Space:** Allows for customization with personal items.



## Organized & Labeled:

Our color-coded and labeled design lets you find what you need in seconds, not minutes.

## First Aid Handbook:

Know that professional-level care is at your fingertips.

**Image Description:** The Besst Survivor First Aid Kit opened, revealing its highly organized interior with color-coded and labeled mesh compartments. Various medical supplies like bandages, shears, and gloves are neatly stored, demonstrating the quick access design.

### Kit Contents (324 Pieces):

The kit is organized into sections for easy identification. Below is a detailed list of included items:

#### Severe Bleeding

- 1 x Israeli Bandage 4"
- 1 x Compressed Gauze 4.1 yards
- 2 x Trauma Pads

#### Sprain & Fracture

- 1 x First Aid Splint 18"
- 1 x Ice Pack

- 1 x Elastic Wrap Bandage
- 1 x Triangular Bandage

### **Burn Care**

- 3 x Burn Cream Packets 3.5g
- 6 x Non-Adherent Pads 3"x4"
- 6 x Non-Adherent Pads 2"x3"

### **Eye Care**

- 6 x Eye Pads

### **CPR Kit**

- 1 x CPR Face Shield
- 1 x Gloves (pair)
- 2 x Towelettes
- 1 x Mini CPR Keychain

### **Light & Communication**

- 1 x Emergency Bright Stick
- 1 x Tactical Flashlight
- 1 x Survival Whistle

### **First Aid Instruments**

- 1 x Scissor
- 1 x Tweezer
- 2 x Gloves (pair)
- 10 x Safety Pins
- 30 x Cotton Tipped Applicators
- 2 x Plastic Bags (2 sizes)
- 1 x First Aid Handbook
- 1 x Pencil

### **Tapes, Wraps & Gauze**

- 8 x Sterile Gauze Pads 3"x3"
- 8 x Sterile Gauze Pads 2"x2"
- 2 x Gauze Rolls 3"
- 2 x Gauze Rolls 2"
- 1 x Adhesive Tape

### **Adhesive Bandages**

- 50 x Standard
- 50 x Mini
- 10 x Knuckle
- 10 x Fingertip
- 8 x Knee/Elbow

### **Outdoor Skin Protection**

- 8 x Sting Relief Pads

### **Personal Space**

- 1 x Organizer Box for Pills
- 1 x Labels

### **Clean and Prep**

- 8 x Towelettes
- 20 x Cleaning Wipes

### **Shelter and Warmth**

- 1 x Emergency Sleeping Bag
- 1 x Disposable Silver Rescue Blanket

### **Others**

- 1 x First Aid Bag
- 1 x Adjustable Shoulder Strap



## 324 Pieces - Comprehensive First Aid Kit

### Severe Bleeding:

- 1 Israeli Bandage 4"
- 1 Compressed Gauze 4.1yards
- 2 Trauma Pads

### Sprains & Fractures:

- 1 First Aid Splint 18"
- 1 Ice Pack
- 1 Elastic Wrap Bandage
- 1 Triangular Bandage

### Burn Care:

- 3 Burn Cream Packets
- 6 Non-Adherent Pads 3\*4"
- 6 Non-Adherent Pads 2\*3"

### Eye Care:

- 6 Eye Pads

### CPR Kit:

- 1 CPR Face Shield
- 1 Gloves (pair)
- 2 Towelettes
- 1 Mini CPR Keychain

### Personal Space:

- 1 Organizer Box for Pill
- 1 Labels

### Hygiene:

- 8 Towelettes
- 20 Cleaning Wipes

### Outdoor Skin Protection:

- 8 Sting Relief Pads

### Tapes,Wraps & Gauze:

- 8 Sterile Gauze Pads 3\*3"
- 8 Sterile Gauze Pads 2\*2"
- 2 Gauze Roll Bandages3"
- 1 Gauze Roll Bandage 2"
- 1 Adhesive Tape

### Adhesive Bandages:

- 50 Standard
- 50 Mini
- 10 Knuckle
- 10 Fingertip
- 8 Knee/Elbow

### Shelter & Warmth:

- 1 Emergency Sleeping Bag
- 1 Disposable Silver Rescue Blanket

### Light & Communication:

- 1 Emergency Bright Stick
- 1 Tactical Flashlight
- 1 Survival Whistle

### First Aid Tools:

- 1 Scissor
- 1 Tweezer
- 2 Gloves (pair)
- 20 Safty Pins
- 70 Cotton Tipped Applicators
- 2 Plastic Bags (2 sizes)

### First Aid Instruments:

- 1 First Aid Handbook
- 1 Pencli

**Image Description:** A comprehensive display of all 324 items included in the first aid kit, neatly arranged and categorized. This image provides a visual inventory of the kit's extensive contents, from bandages and gauze to specialized trauma items and tools.

## SETUP AND ORGANIZATION

Upon receiving your Besst Survivor First Aid Kit, take time to familiarize yourself with its layout and contents. The kit is designed with color-coded and labeled compartments to facilitate quick access during an emergency.

1. **Inspect Contents:** Open the kit and verify all 324 pieces are present according to the list in the "Product Overview" section.
2. **Personalize:** Utilize the "Personal Space" organizer box and labels to add any specific medications or items you may require.
3. **Attachment Options:**
  - **Shoulder Strap:** Attach the adjustable shoulder strap for easy carrying.
  - **MOLLE System:** The removable MOLLE-compatible panel allows attachment to tactical backpacks, car seats, or other gear. Secure the D-rings to anchor points for stability.
  - **Handles:** Use the top or side handles for quick grab-and-go situations.

4. **Placement:** Store the kit in an easily accessible location at home, in your vehicle, or within your camping/hiking gear. Ensure it is protected from extreme temperatures and moisture.

## Portable, Easy Access:

It to be attached to various platforms like cars, bikes, or backpacks, enhancing its portability and accessibility.



Removable Shoulder Strap



MOLLE Compatible



Removable Velcro Back Panel

**Image Description:** This image illustrates the versatility of the Best Survivor First Aid Kit's carrying and attachment options. It shows the kit being worn with a shoulder strap, secured to a backpack using its MOLLE system, and attached to the back of a car seat via its removable velcro panel.

## OPERATING INSTRUCTIONS: GENERAL FIRST AID GUIDANCE

This section provides general guidance on using the items in your first aid kit. For detailed instructions on specific medical procedures, refer to the included First Aid Handbook and seek professional medical training. Always prioritize safety and call for professional medical help when necessary.

### 1. Severe Bleeding Management

- **Israeli Bandage:** Apply directly over the wound. Wrap tightly and secure the pressure bar to maintain direct pressure.
- **Compressed Gauze:** Pack into deep wounds to control bleeding. Apply direct pressure.
- **Trauma Pads:** Use for large wounds to absorb fluid and apply pressure.
- *Always apply direct pressure to the wound and elevate if possible. Seek immediate professional medical*

attention for severe bleeding.

## 2. Sprains & Fractures

- **First Aid Splint:** Mold the splint to immobilize injured limbs. Secure with elastic wrap or tape.
- **Ice Pack:** Activate and apply to sprains or strains to reduce swelling and pain. Do not apply directly to skin; wrap in cloth.
- **Elastic Wrap Bandage:** Use to provide compression and support for sprains or to secure splints.
- **Triangular Bandage:** Can be used as a sling for arm injuries or to secure dressings.
- *Immobilize the injured area and seek medical evaluation for suspected fractures.*

## 3. Burn Care

- **Burn Cream Packets:** Apply a thin layer to minor burns after cooling the area with water.
- **Non-Adherent Pads:** Cover burns gently with these pads to prevent sticking and protect the wound.
- *For severe burns, cool with water, cover loosely, and seek immediate medical attention. Do not apply ice.*

## 4. Eye Care

- **Eye Pads:** Gently cover an injured eye to protect it from further irritation. Secure with tape.
- *Do not attempt to remove embedded objects from the eye. Seek immediate medical attention for eye injuries.*

## 5. CPR Kit Usage

- **CPR Face Shield:** Place over the victim's mouth and nose during rescue breaths to create a barrier.
- **Gloves:** Wear gloves to protect yourself from bodily fluids.
- **Towelettes:** Use for cleaning hands or the victim's face if necessary.
- *Only perform CPR if you are trained. Call emergency services immediately.*

## 6. Wound Cleaning & General Care

- **Gloves:** Always wear gloves when treating wounds to prevent infection.
- **Cleaning Wipes/Towelettes:** Clean the area around the wound gently.
- **Tweezers:** Carefully remove splinters or small debris from wounds.
- **Scissors:** Cut bandages, tape, or clothing as needed.
- **Adhesive Bandages:** Apply to minor cuts and scrapes after cleaning.
- **Gauze Pads/Rolls:** Cover larger wounds after cleaning. Secure with adhesive tape.

## 7. Light & Communication

- **Tactical Flashlight:** For illumination in dark environments. Check batteries regularly.
- **Emergency Bright Stick:** Activate for temporary light or signaling.
- **Survival Whistle:** Use to signal for help, especially in outdoor settings.

## 8. Outdoor Skin Protection

- **Sting Relief Pads:** Apply to insect bites or stings to alleviate discomfort.

## 9. Shelter & Warmth

- **Emergency Sleeping Bag:** Provides insulation and protection from the elements in survival situations.
- **Disposable Silver Rescue Blanket:** Helps retain body heat in cold conditions or prevent overheating in hot conditions.

# Comprehensive First Aid


A comprehensive set of professional-grade medical supplies covering various emergency needs, ensuring you're well-prepared.



Severe Bleeding:  
Israeli Bandage



Sprains & Fractures:  
Splint

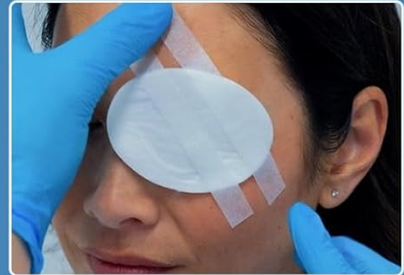

Eye Care



Severe Bleeding:  
Compressed Gauze



CPR Kit:  
Instructions Attached To Kit



Burn Care



Travel Pills Organizer



Emergency Sleeping Bag



First Aid Guide

**Image Description:** A visual guide showcasing the application of various first aid items for different emergencies. It includes scenes of applying an Israeli bandage for severe bleeding, using a splint for fractures, eye care, the CPR kit, burn care, organizing pills, deploying an emergency sleeping bag, and consulting the first aid guide.

## MAINTENANCE

Regular maintenance ensures your first aid kit remains effective and ready for use.

- **Regular Inspection:** At least once every 6-12 months, open your kit and inspect all contents.
- **Check Expiration Dates:** Dispose of and replace any expired medications or sterile items.
- **Restock Used Items:** After any use, immediately replace items to ensure the kit is complete for the next emergency.
- **Clean the Bag:** Wipe down the exterior of the bag with a damp cloth as needed. Ensure it is completely dry before zipping it closed.
- **Battery Check:** For the tactical flashlight, check and replace batteries annually or as needed.
- **Storage:** Store the kit in a cool, dry place away from direct sunlight and extreme temperatures to preserve the integrity of the supplies.

## TROUBLESHOOTING

---

While a first aid kit is straightforward, here are common concerns and their resolutions:

- **Missing Items:** If you find items missing upon initial inspection, contact Besst Survivor customer support. For items used during an emergency, restock them promptly.
- **Damaged Packaging:** If sterile items have damaged packaging, they are no longer considered sterile. Replace them immediately.
- **Difficulty Finding Items:** Re-familiarize yourself with the kit's labeled compartments. The color-coded system is designed for quick identification.
- **Flashlight Not Working:** Check the batteries and replace them if necessary. Ensure the battery cap is securely tightened.

## SPECIFICATIONS

---

<b>Product Dimensions</b>	10.9 x 7.5 x 4.6 inches
<b>Weight</b>	3.5 Pounds
<b>Material</b>	600D Rip-resistant Polyester
<b>Number of Pieces</b>	324
<b>Model Number (ASIN)</b>	B09KMTLS9X
<b>Manufacturer</b>	Besst Survivor
<b>Special Features</b>	Compact, Mountable, Portable, Reusable, Reflective Piping, Labeled Compartments, MOLLE System



600D Ripstop Polyester



MOLLE Compatible Removable Velcro Panel



D- Ring



Adjustable Strap



Reflective Design



Detachable Shoulder Strap



**Weight:**  
3.5 lbs

**Image Description:** A visual representation of the Besst Survivor First Aid Kit's physical specifications, including its dimensions (10.9 x 7.5 x 4.6 inches) and weight (3.5 lbs). Key features like 600D Ripstop Polyester, MOLLE compatibility, D-rings, adjustable strap, reflective design, and detachable shoulder strap are highlighted.

## WARRANTY AND SUPPORT

Specific warranty details are not provided in the product information. For any inquiries regarding your Besst Survivor First Aid Kit, including missing parts, defects, or general support, please contact the manufacturer directly.

**Manufacturer:** Besst Survivor

Please refer to the retailer's website or your purchase documentation for contact information or return policies.

## Important Information

**Legal Disclaimer:** Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This first aid kit is for emergency use and is not a substitute for professional medical care. Always seek qualified medical assistance when needed.

© 2025 Besst Survivor. All rights reserved.

