

[manuals.plus](#) /

- › [Audien Hearing](#) /
- › [Audien ATOM Rechargeable Hearing Amplifier to Aid and Assist Hearing, Premium Comfort Design and Nearly Invisible](#)

Audien Hearing Atom

Audien ATOM Hearing Amplifier User Manual

Model: Atom | Brand: Audien Hearing

INTRODUCTION

The Audien ATOM Rechargeable Hearing Amplifier is designed to aid and assist hearing with a premium comfort design and nearly invisible fit. This manual provides essential information for setting up, operating, and maintaining your device to ensure optimal performance and longevity.

✓ **500,000+ HAPPY CUSTOMERS**

✓ **NEARLY INVISIBLE IN THE EAR**

✓ **SUPERIOR SOUND**

✓ **100% RECHARGEABLE**

✓ **AFFORDABLE**

✓ **COMFORTABLE**



Image: Audien ATOM Hearing Amplifiers highlighting key features like nearly invisible, superior sound, 100% rechargeable, affordable, and comfortable.

WHAT'S IN THE BOX

Your Audien ATOM package includes the following components:

- 2 hearing amplifiers
- 1 charging dock
- 1 Charging Cable and Wall Plug
- 1 carrying case
- 4 wax guards
- 1 ear tips set (3 different sizes: Small, Medium, Large)
- User Manual & Getting Started Guide



Image: Contents of the Audien ATOM package including two hearing aids, charging dock, cables, and accessories.

SETUP

Initial Charging

Before first use, fully charge your Audien ATOM hearing amplifiers. Place both devices into the sleek magnetic charging dock. Connect the charging cable to the dock and plug it into a wall outlet using the provided wall plug. A full charge takes approximately 4-6 hours and provides up to 20 hours of battery life.



**NO BATTERIES
TO REPLACE... EVER**

Image: Audien ATOM charging dock with hearing aids, showing 4-6 hour charging time and 20 hour battery life.

Fitting the Ear Tips

Select the ear tip size (Small, Medium, or Large) that provides the most comfortable and secure fit within your ear canal. The correct fit is crucial for optimal sound quality and to minimize whistling or feedback. Gently push the chosen ear tip onto the amplifier until it is firmly attached.



Image: Close-up view of Audien ATOM hearing amplifiers and charging dock, illustrating their compact size and design.

OPERATING INSTRUCTIONS

Power On/Off

To power on the device, locate the small switch on the amplifier and slide it to the "On" position. To power off, slide the switch to the "Off" position. The devices are designed to be left on when placed in the charging dock.

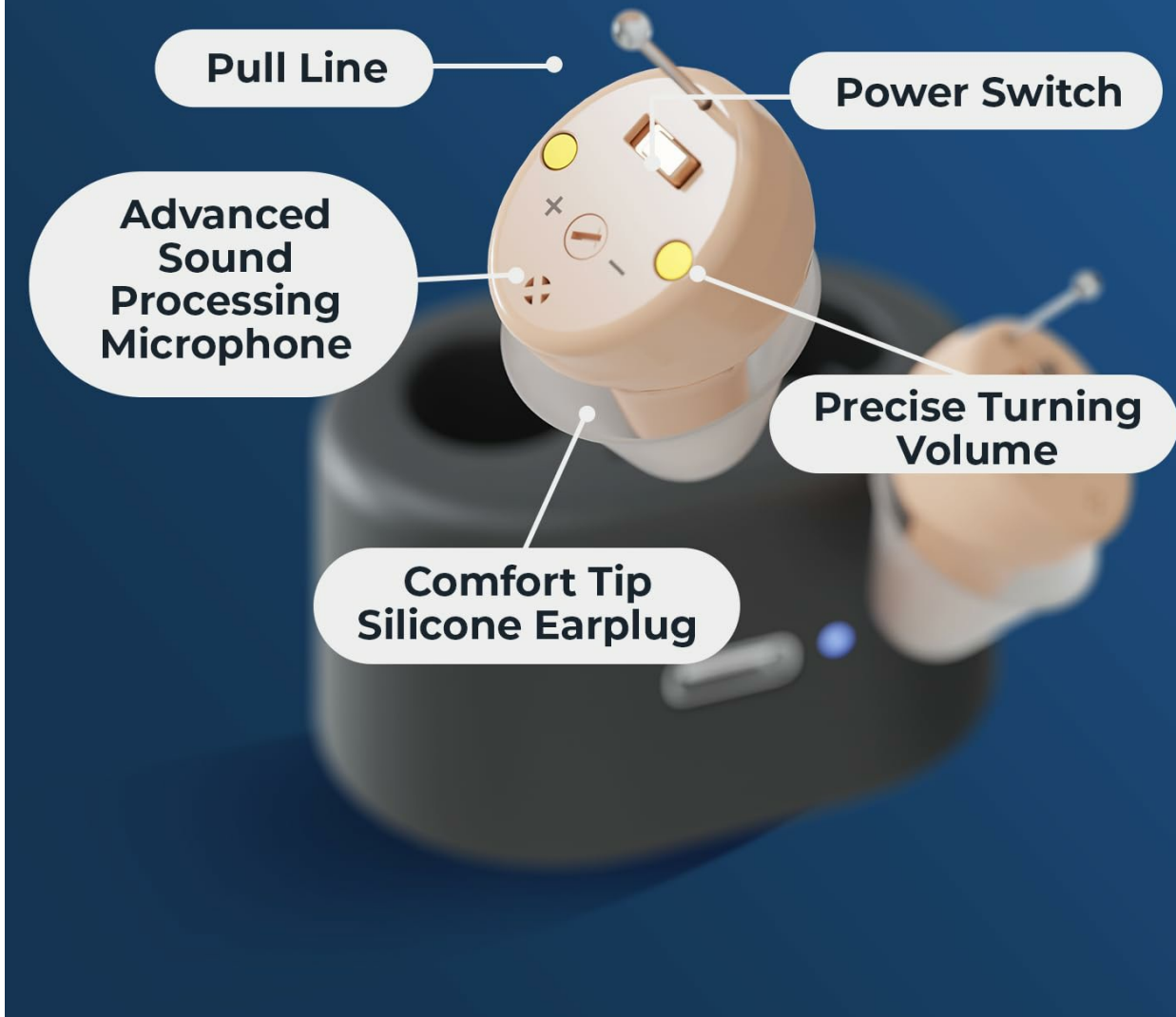


Image: Diagram of Audien ATOM hearing amplifier components, including pull line, power switch, microphone, and volume control.

Volume Adjustment

The Audien ATOM features a precise turning volume control. Use the small screwdriver tool provided to adjust the volume to your individual needs. Start with a lower setting and gradually increase until the sound is comfortable and clear. The devices amplify speech while reducing background noise, making it easier to hear in various settings.

For a visual guide on operation, please refer to the official product video below:

Your browser does not support the video tag.

Video: Official Audien Hearing Atom product overview, demonstrating features and usage.

MAINTENANCE

Cleaning Your Amplifiers

Regular cleaning is essential to maintain the performance and longevity of your Audien ATOM hearing amplifiers. Use the

provided cleaning brush to gently remove any earwax or debris from the ear tips and the device's surface. Ensure the built-in, replaceable wax guards are clear to prevent sound blockage.

Replacing Wax Guards

The devices include built-in, replaceable wax guards. If you notice a decrease in sound quality or blockage, replace the wax guards using the spares provided in your kit. Refer to the Getting Started Guide for detailed instructions on wax guard replacement.

TROUBLESHOOTING

If you experience issues with your Audien ATOM hearing amplifiers, consider the following common troubleshooting steps:

- **No Sound:** Ensure the device is powered on and fully charged. Check that the ear tip is not blocked by earwax.
- **Weak Sound:** Adjust the volume using the provided tool. Clean the ear tip and wax guard. Ensure the ear tip is properly seated in your ear canal.
- **Whistling/Feedback:** This often indicates an improper fit. Try a different size ear tip from the included set to ensure a snug seal in your ear canal.
- **Intermittent Sound:** Ensure the device is fully charged. Clean the charging contacts on both the device and the dock.

For further assistance, please refer to the comprehensive User Manual or contact Audien Hearing customer support.

SPECIFICATIONS

Feature	Detail
Model Name	Atom
Power Source	Battery Powered
Batteries	2 Lithium Ion batteries required (included)
Battery Life	Up to 20 hours
Charging Time	4-6 hours
Material	Silicone
Color	Black
Package Dimensions	5.98 x 3.78 x 1.61 inches; 7.05 ounces





IMPORTANT SAFETY INFORMATION



Hearing amplifiers are not a replacement for hearing aids. Misuse can cause hearing damage. Consult a healthcare professional before use.
This product is from a small business brand. Support small.

WARRANTY AND SUPPORT


Audien Hearing is committed to ensuring your experience is exceptional. For any questions, troubleshooting, or support needs, please contact Audien Hearing customer support. The brand offers USA-based support to assist you. Please refer to the contact information provided in your Getting Started Guide for direct support channels.

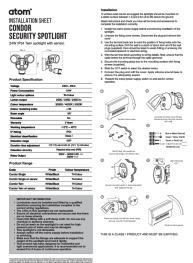Related Documents - Atom

	<p>Atom BILI-THERAPY Spot Type LED Phototherapy Unit - Neonatal Jaundice Treatment</p> <p>Learn about the Atom BILI-THERAPY Spot Type, an advanced LED phototherapy unit designed for effective neonatal jaundice treatment. This brochure details its features, space-saving design, switchable examination lamp function, product lineup, and technical specifications. Provided by Parker Healthcare.</p>
	<p>Atom Massage Gun User Manual - Model 648459</p> <p>Comprehensive user manual for the Atom massage gun (Model 648459), detailing its features, included heads, usage instructions, charging, safety guidelines, and medical disclaimers. Available in multiple languages.</p>
	<p>Atom Condor Security Spotlight Installation and User Guide</p> <p>Comprehensive guide for installing and operating the Atom Condor Security Spotlight, featuring a 24W IP54 twin spotlight with a PIR sensor. Includes product specifications, installation steps, PIR control adjustments, troubleshooting, and warranty information.</p>
	<p>Atom AT9090/TRI LED Downlight Installation Sheet & Specifications</p> <p>Comprehensive installation guide and product specifications for the Atom AT9090/TRI 8W low profile LED downlight. Includes safety information, warranty details, and technical data.</p>

 <p>624160</p> <p>atom Binocular · Kiiikarit · Kikare Binokkel · Binoklis Binoklis · Биноколь</p>	<p>Atom Binoculars 624160 User Manual - English, Finnish, Swedish, Estonian, Latvian, Lithuanian, Russian</p> <p>Instruction manual for Atom binoculars (Model 624160), providing detailed guidance on setting interpupillary distance, focusing, diopter adjustment, use with eyeglasses, zoom function, and essential cleaning and maintenance procedures. Includes important safety warnings regarding sun observation.</p>
	<p>Atom Eclipse 240mm Ducted Exhaust Fan - White Product Specifications</p> <p>Explore the Atom Eclipse 240mm ducted exhaust fan, featuring a powerful ball bearing motor, quiet operation, energy efficiency, and a 3-year warranty. Includes detailed specifications, dimensions, and accessories.</p>


Documents - Audien Hearing – Atom

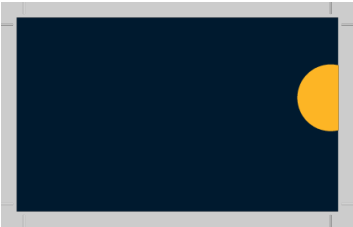

[Audien Atom 2 User Manual: Your Guide to Enhanced Hearing](#)

Comprehensive user manual for the Audien Atom 2 hearing aids, covering setup, usage, charging, cleaning, troubleshooting, and technical specifications. Learn how to fit your device, adjust volume and listening modes, and maintain your hearing aids for optimal performance.


lang:en score:28 filesize: 2.03 M page_count: 24 document date: 2024-05-12

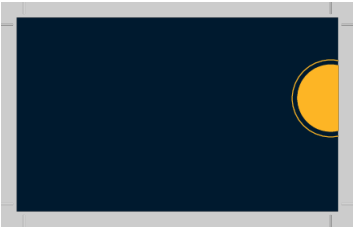

[Audien Atom Pro 2 User Manual: Your Guide to Enhanced Hearing](#)

Comprehensive user manual for the Audien Atom Pro 2 hearing aid, covering setup, usage, cleaning, troubleshooting, and technical specifications. Learn how to fit, adjust, and maintain your device for optimal hearing.

lang:en score:28 filesize: 2.07 M page_count: 24 document date: 2024-05-12

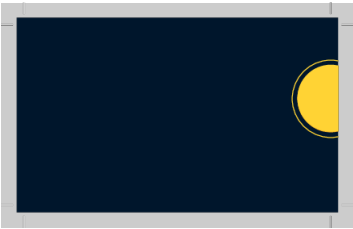

[Audien Atom X User Manual: Comprehensive Guide to Your Hearing Aids](#)

Explore the Audien Atom X user manual for setup, operation, troubleshooting, and maintenance of your OTC hearing aids. Learn about Bluetooth connectivity, charging, cleaning, and more.

lang:en score:22 filesize: 3.49 M page_count: 36 document date: 2025-07-03

Audien Hearing 1-Year Limited Warranty

Warranty: Audien Hearing warrants that our hearing aids shall be free from material defects in material and workmanship, conform to applicable specifications, and perform in accordance with the product documentation under normal usage for a period of 1 Year from the date of purchase or delivery to the Customer.

Remedy: In the event of a material defect or non-conforming product, Audien Hearing shall, at its sole option and expense, promptly replace the defective product or refund the purchase price. Such remedy shall be the Customer's sole and exclusive remedy under this Warranty.

Exclusions: The Warranty does not cover

- (a) Consumable parts such as Wax Guards, Cleaning Brushes, and Ear Domes;
- (b) Damage by accident, abuse, misuse, liquid contact, or other external cause;
- (c) Damage caused by operating the product outside of the product's published guidelines;
- (d) Defects caused by normal wear and tear or otherwise due to the normal aging of the Product.

Procedure: To file a Warranty claim, the Customer must contact Audien Hearing's customer support center and specify the nature of the defect. Audien will arrange for the return of the defective product at Audien's cost for inspection and subsequent Customer Remedy.

[Audien Hearing 1-Year Limited Warranty](#)

Details the 1-year limited warranty for Audien Hearing hearing aids, covering terms, remedies, exclusions, and claim procedures.

lang:en score:17 filesize: 64.5 K page_count: 1 document date: 0000-00-00

