

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [KMINA](#) /

› [KMINA PRO Shower Chair with Wheels User Manual](#)

KMINA 0712

KMINA PRO Shower Chair with Wheels User Manual

MODEL: 0712

1. Introduction

The KMINA PRO Shower Chair with Wheels is designed to provide comfort, safety, and independence for individuals with limited mobility. This versatile chair can be used as a shower chair, a commode chair, or for assisted transport, simplifying daily hygiene and movement tasks for both users and caregivers. Its robust aluminum construction ensures durability, while thoughtful features like adjustable components and an improved braking system enhance user experience and safety.

2. Safety Information

- Always ensure all components are securely fastened before use.
- Do not exceed the maximum user weight of 265 lb (120 kg).
- Engage all wheel brakes when the chair is stationary, especially during transfers or when using it as a commode.
- Ensure the floor surface is clear of obstacles and dry to prevent slips during movement.
- Seek assistance for transfers if you have significant mobility limitations.
- Regularly inspect the chair for any signs of wear, damage, or loose parts. Discontinue use if any damage is found and contact customer support.
- Keep hands and fingers clear of moving parts during adjustments.

3. Components and Package Contents

Upon opening the box, you will find the following main components:

- Main frame with seat and backrest
- Two armrests
- Horizontal tube for armrest connection
- Two footrests
- Commode bucket with lid
- Assembly hardware (knobs, screws)

For a visual guide on unboxing and assembly, please refer to the assembly video in the next section.



Figure 3.1: KMINA PRO Shower Chair with Wheels, fully assembled.

4. Setup and Assembly

The KMINA PRO Shower Chair requires minimal assembly. Follow these steps:

1. **Attach the Armrests:** Insert each armrest into the designated vertical tubes on the main frame. Adjust the height to one of the 6 available positions for optimal comfort and secure it with the provided knobs.
2. **Insert the Horizontal Tube:** Connect the two armrests using the horizontal tube. This tube can be adjusted in two angles (0°/45°) and two positions (68/75 cm) to fit various toilet sizes and user preferences. Secure it firmly.
3. **Attach the Footrests:** Slide the footrests into their respective slots on the front of the chair. They are removable and can be folded up when not in use.
4. **Place the Commode Bucket:** Position the commode bucket securely under the seat opening. Ensure it is properly seated before use.

Video 4.1: KMINA PRO - ASSEMBLY | Rolling Shower Chair. This video demonstrates the step-by-step assembly process of the shower chair.



Figure 4.1: Illustration of removable footrests and folding armrests for user convenience.

5. Operating Instructions

5.1. Multifunctional Use

The KMINA PRO Shower Chair is designed for various uses:

- **Shower Chair:** Use in the shower for personal hygiene. The aluminum frame is resistant to water.
- **Commode Chair:** Utilize the incorporated commode bucket for toileting needs.
- **Over-the-Toilet Use:** Position the chair directly over a standard WC. Ensure the toilet's width (measurement A) is less than or equal to 14.75 inches and height (measurement B) is less than or equal to 17 inches for proper fit.
- **Transport Chair:** The directional wheels allow for easy movement and transfer of the user.

5.2. Braking System

The chair features an improved braking system on all four wheels for maximum safety and grip. To engage the brake, press down on the red lever located on each wheel. To release, lift the lever up. Always ensure brakes are engaged when the chair is stationary or during transfers.

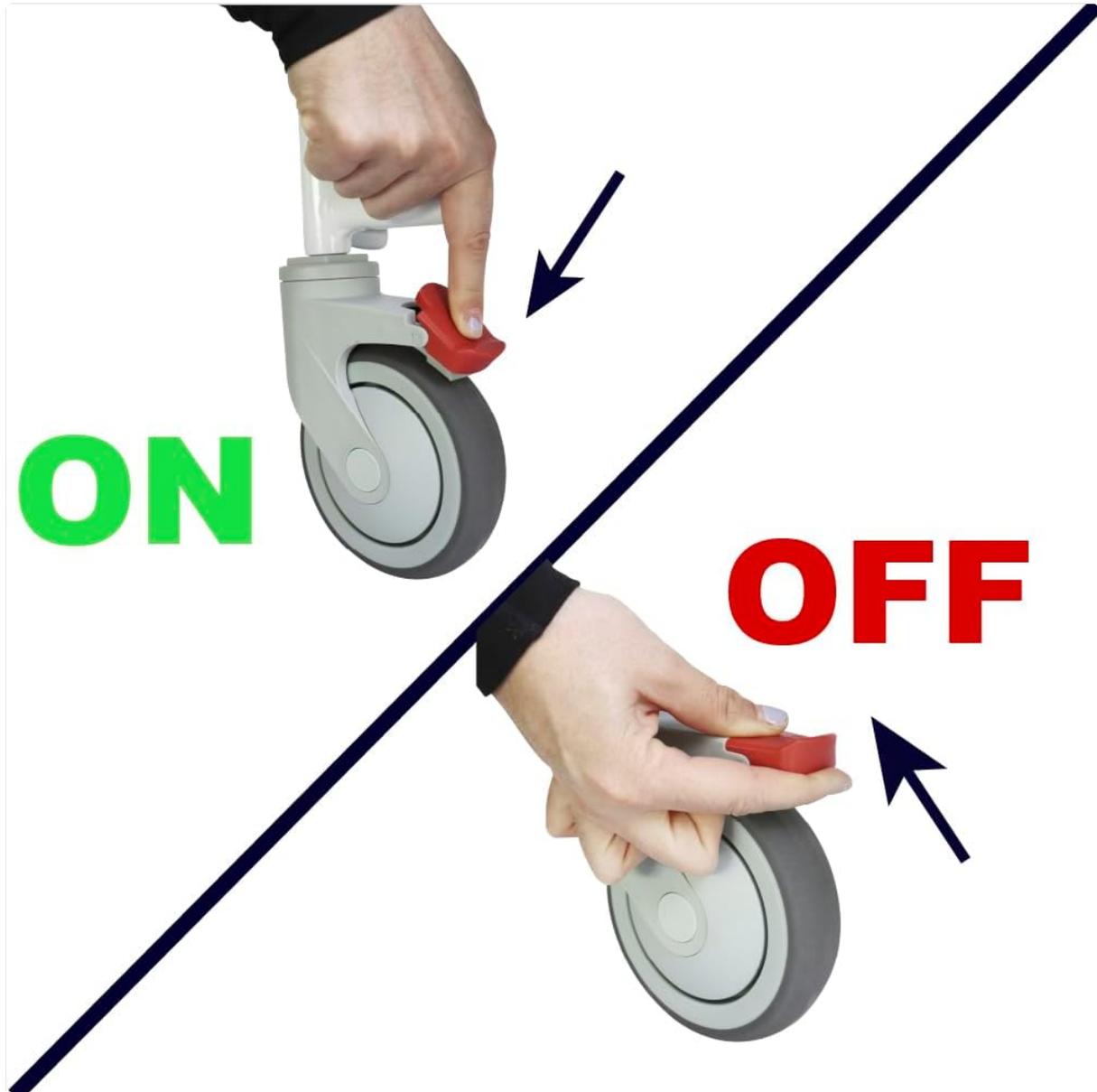


Figure 5.1: How to engage (ON) and disengage (OFF) the wheel brakes.

5.3. Adjustments

- **Armrest Height:** Adjust the armrests to one of the 6 available height settings by loosening the knob, sliding the armrest to the desired height, and re-tightening the knob.
- **Footrests:** The footrests can be removed or folded up to facilitate transfers or when not needed.

Video 5.2: KMINA PRO - USE | Rolling Shower Chair. This video demonstrates the various uses and features of the shower chair.



Figure 5.3: The KMINA PRO Shower Chair positioned over a toilet.



+1.000.000 More than a million people trust KMINA.

+30 Countries Available in +30 countries.

PATENTED +20 patented products.

MADE IN EUROPE 40% manufactured in Europe.

Figure 5.4: Toilet width (A) must be less than or equal to 14.75 inches.



Figure 5.5: Toilet height (B) must be less than or equal to 17 inches.

6. Maintenance

- **Cleaning:** Clean the chair regularly with mild soap and water. Rinse thoroughly and dry completely to prevent water spots and corrosion. Do not use abrasive cleaners or solvents.
- **Inspection:** Periodically check all screws, bolts, and knobs to ensure they are tight. Inspect the wheels and brakes for smooth operation and wear.

- **Storage:** Store the chair in a dry, cool place away from direct sunlight when not in use for extended periods.

7. Troubleshooting

- **Difficulty Engaging Brakes:** Ensure no debris is obstructing the brake mechanism. Check if the brake lever is fully depressed. If issues persist, inspect the wheel and brake for damage.
- **Chair Does Not Roll Smoothly:** Check wheels for hair, lint, or other obstructions. Ensure all wheels are properly aligned and not damaged.
- **Instability:** Verify that all assembly components (armrests, horizontal tube) are securely tightened. Ensure the chair is on a flat, even surface.

If you encounter problems not listed here or if troubleshooting steps do not resolve the issue, please contact KMINA customer support for assistance.

8. Specifications

Product Dimensions	21.85"D x 21.85"W x 20.87"H
Item Model Number	0712
Item Weight	11.5 Kilograms (25.35 Pounds)
Maximum User Weight	265 lb (120 kg)
Material	Aluminum frame, PU padded seat
Color	Black

9. Warranty and Support

For detailed warranty information and customer support, please refer to the official KMINA website or the documentation included with your product. You may also find additional resources and a comprehensive user manual in PDF format at the following link:

[Download Full User Manual \(PDF\)](#)

For direct assistance, please contact KMINA customer service through their official channels.

