



[Manuals.plus](#) /

> [COMHOMA](#) /

> COMHOMA Power Lift Recliner H7135-BLK User Manual

COMHOMA H7135-BLK

COMHOMA Power Lift Recliner User Manual

Model: H7135-BLK

1. INTRODUCTION

Thank you for choosing the COMHOMA Power Lift Recliner. This manual provides essential information for the safe and efficient use of your new recliner. Please read it thoroughly before assembly and operation, and retain it for future reference.



Image 1.1: COMHOMA Power Lift Recliner in black, featuring highlighted areas indicating the massage and heating functions on the backrest, lumbar, and leg rest.

2. SAFETY INFORMATION

Your safety is paramount. Please adhere to the following guidelines:

- Always use the hand controller to operate the chair's lift and recline functions.
- Do not attempt to stand up or sit down while the chair is in motion.
- Ensure the chair is on a stable, level surface before use.
- Do not exceed the maximum weight capacity of **300 lbs.**
- Keep the chair's power cord away from heat sources and liquids to prevent fire or electrical hazards.
- Always follow the assembly instructions carefully to ensure the chair is properly assembled and stable.
- Keep children and pets away from the chair's moving parts.
- Unplug the chair from the power outlet before cleaning or maintenance.

3. ASSEMBLY INSTRUCTIONS

The COMHOMA Power Lift Recliner is designed for straightforward assembly. All necessary tools are included with the chair parts.

3.1 Unpacking

1. Carefully remove all components from the packaging.
2. Verify that all parts listed in the packing list are present. Contact customer support if any parts are missing or damaged.
3. Place the main chair base on a clean, flat, and stable surface.

3.2 Component Assembly

1. Slide the backrest into the designated slots on the chair base until it clicks securely into place.
2. Connect any electrical cables from the backrest to the chair base, ensuring a snug fit.
3. Attach the armrests by aligning them with the corresponding brackets and pushing down firmly until they lock.
4. Connect the power adapter to the chair's power input and then plug it into a wall outlet.

Refer to the included pictorial guide for detailed visual instructions during assembly.

4. OPERATING INSTRUCTIONS

Your recliner is equipped with a user-friendly remote control to manage all functions.

4.1 Power Lift Function

LIFT RECLINER CHAIR



Image 4.1: A sequence showing a person using the power lift function, transitioning from a seated position to a standing position with the chair gently elevating.

- To **lift** the chair: Press and hold the 'UP' button on the remote control. The chair will slowly tilt forward and elevate, assisting you in standing up.
- To **lower** the chair: Press and hold the 'DOWN' button on the remote control. The chair will gently lower back to a seated position.

4.2 Recline Function



Image 4.2: An image demonstrating the adjustable backrest, showing three positions: upright for chatting, reclined for watching TV, and fully flat for sleeping.

- To **recline**: Use the 'RECLINE' button (or similar directional buttons) on the remote. The backrest will move backward and the footrest will extend.
- To **return to upright**: Use the 'UPRIGHT' button (or similar directional buttons) on the remote. The footrest will retract, and the backrest will return to its original position.
- The chair offers an **infinite number of reclining positions**, allowing you to stop at any desired angle.

4.3 Massage and Heating Functions



Image 4.3: This graphic highlights the key features of the recliner's massage and heating system, including 5 massage modes, 4 massage parts, lumbar heating, time function, and vibration strength.

- **Massage Modes:** Select from 5 different massage modes using the 'MODE' button on the remote.
- **Massage Parts:** Target 4 specific massage areas (back, lumbar, thighs, calves) by pressing the corresponding buttons.
- **Lumbar Heating:** Activate the lumbar heating function for soothing warmth by pressing the 'HEAT' button.
- **Time Function:** Set a massage duration using the 'TIMER' button.
- **Vibration Strength:** Adjust the intensity of the massage using the 'INTENSITY' buttons.

5. MAINTENANCE

Proper care will extend the life and appearance of your recliner.

- **Cleaning:** For the PU leather surface, wipe with a soft, damp cloth. Avoid harsh chemicals or abrasive cleaners. For spills, blot immediately with a clean, dry cloth.
- **Electrical Components:** Regularly check the power cord for any damage. Ensure all connections are secure. Do not attempt to repair electrical components yourself.

- **Moving Parts:** Periodically inspect the lift and recline mechanisms for any obstructions or unusual noises. Lubricate moving parts if necessary, using a silicone-based lubricant.
- **General Care:** Avoid placing the chair in direct sunlight or near heat sources to prevent material degradation.

6. TROUBLESHOOTING

If you encounter issues with your recliner, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Chair does not move/respond to remote.	Power cord unplugged; loose connection; remote battery low/dead; power outage.	Check power cord connection to wall and chair. Ensure wall outlet has power. Replace remote batteries if applicable.
Lift or recline mechanism is stuck.	Obstruction in mechanism; motor issue.	Inspect for any objects blocking the mechanism. Do not force movement. If issue persists, contact customer support.
Massage/Heat function not working.	Function not activated; loose internal connection.	Ensure the specific function is activated via the remote. Check power connection.
Unusual noises during operation.	Lack of lubrication; loose components.	Apply silicone lubricant to moving parts. Check for loose screws or bolts and tighten if necessary.

If these steps do not resolve the issue, please contact COMHOMA customer support.

7. SPECIFICATIONS

Detailed specifications for your COMHOMA Power Lift Recliner:



Image 7.1: A technical diagram illustrating the various dimensions of the recliner, including height, width, and depth in both upright and reclined positions.

- **Brand:** COMHOMA
- **Model Name:** Recliner (H7135-BLK)
- **Color:** Black
- **Product Dimensions (Upright):** 28.7"D x 33.4"W x 39.3"H
- **Overall Size:** 32.5 inches x 28.7 inches x 39.9 inches
- **Material:** PU Leather (Furniture Finish), Sponge (Fill Material)
- **Maximum Weight Recommendation:** 300 Pounds
- **Item Weight:** 88.4 pounds
- **Reclining Position Count:** Infinite
- **Special Features:** Power Lift, Massage (5 modes, 4 parts), Lumbar Heating, Ergonomic Design, Side Pockets

- **Included Components:** Arm Pad
- **UPC:** 810154191178

8. WARRANTY AND SUPPORT

COMHOMA is committed to providing quality products. For warranty information or technical support, please refer to your purchase documentation or contact COMHOMA customer service directly.

You may also visit the official COMHOMA store for additional product information and support resources: [COMHOMA Store](#)

8.1 Additional Comfort Features




Spacious Seat



Expandable Footrest



Side Pockets

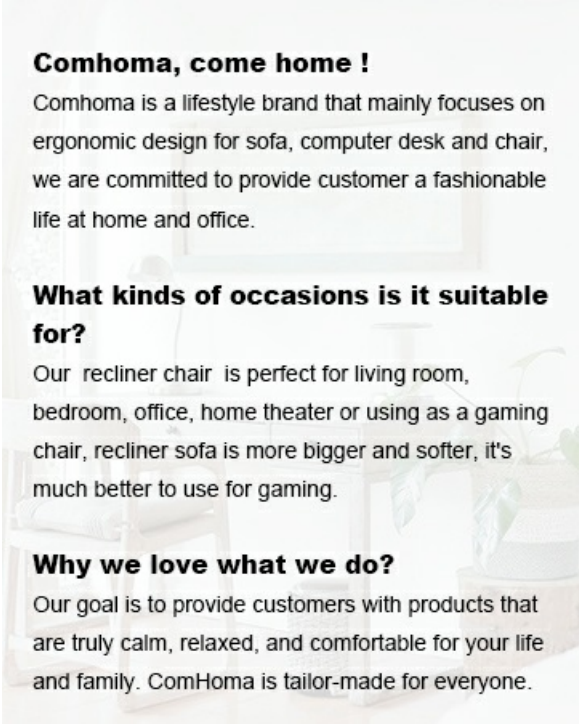

Comhoma, come home !

Comhoma is a lifestyle brand that mainly focuses on ergonomic design for sofa, computer desk and chair, we are committed to provide customer a fashionable life at home and office.

What kinds of occasions is it suitable for?

Our recliner chair is perfect for living room, bedroom, office, home theater or using as a gaming chair, recliner sofa is more bigger and softer, it's much better to use for gaming.

Why we love what we do?

Our goal is to provide customers with products that are truly calm, relaxed, and comfortable for your life and family. ComHoma is tailor-made for everyone.

Padded Armrests