

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [ELENKER](#) /

› [ELENKER Adjustable Steerable Knee Scooter User Manual](#)

## ELENKER B0CWQFDCNP

# ELENKER Adjustable Steerable Knee Scooter User Manual

Model: B0CWQFDCNP

## INTRODUCTION

---

The ELENKER Adjustable Steerable Knee Scooter is designed to provide a stable and comfortable mobility solution for individuals recovering from foot injuries, ankle surgery, or other lower leg conditions. It serves as an effective alternative to traditional crutches, offering enhanced stability and ease of use for daily activities.



Image: The ELENKER Knee Scooter in light blue and white, featuring a comfortable knee pad, steerable handlebars, and four wheels, with a black storage bag attached.

## SAFETY INFORMATION

---

Please read all safety instructions carefully before using the ELENKER Knee Scooter to ensure safe operation and prevent injury.

- **Speed Control:** Do not slide at high speed. The recommended safe speed is 1 meter per second (1m/s).
- **Braking System:** Always ensure the dual adjustable locking handbrakes and dual rear on-wheel brakes are functioning correctly before each use. Familiarize yourself with their operation.
- **Surface Use:** Use the knee scooter on smooth, dry, and level surfaces. Avoid uneven terrain, stairs, or slippery areas.
- **Stability:** Maintain proper balance while operating the scooter. Do not lean excessively or attempt sharp turns at speed.
- **Adjustments:** Ensure all adjustable components (handlebars, knee rest) are securely locked in place after adjustment and

before use.

- **Weight Capacity:** Adhere to the specified weight capacity of the device to ensure stability and prevent damage.

**Legal Disclaimer:** Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is a medical device and should be used as directed by a healthcare professional.

## COMPONENTS

---

The ELENKER Knee Scooter consists of the following main components:

- **Handlebars:** For steering and control, equipped with handgrips and brake levers.
- **Knee Rest Pad:** Padded support for the injured leg, adjustable for height.
- **Frame:** Sturdy iron frame providing the main structure.
- **Wheels:** 8-inch PVC wheels (front and rear) for smooth movement.
- **Braking System:** Dual handbrakes and rear on-wheel brakes for stopping and parking.
- **Storage Bag:** Detachable cloth bag for personal belongings.



Image: Detachable Cloth Bag. This convenient bag allows you to carry personal items while using the scooter.

# Knee Pad

comfortable and soft



Image: Knee Pad. The comfortable and soft knee pad provides support and reduces pressure during use.

## SETUP AND ASSEMBLY

Follow these steps to assemble your ELENKER Knee Scooter:

1. **Unpack:** Carefully remove all components from the packaging.
2. **Assemble Frame:** Connect the main frame components as per the included quick-start guide (if applicable).
3. **Attach Handlebars:** Insert the handlebar post into the main frame and secure it with the locking mechanism.
4. **Attach Knee Rest:** Mount the knee rest pad onto its support post and secure it.
5. **Adjust Height:** Adjust the height of the handlebars and knee rest to your comfortable and safe level. Ensure all quick-release levers or knobs are tightened securely.
6. **Check Brakes:** Verify that both handbrakes engage properly and effectively lock the wheels.
7. **Attach Bag:** Secure the detachable cloth bag to the handlebar post.

For storage or transport, the knee scooter can be folded:

# Save Space

Compact design for easy transport or storage




01



Fold down steering column and keep it in alignment.

02

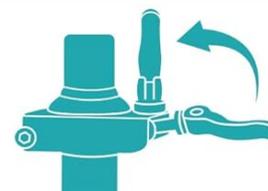

Fold the lock lever up

03



Turn the handle to left or right to fasten the part.

Image: Save Space. The compact design allows for easy transport or storage. Steps include folding down the steering column, folding up the lock lever, and turning the handle to fasten.

## OPERATING INSTRUCTIONS

Using your ELENKER Knee Scooter effectively and safely:

- **Positioning:** Place the injured leg's knee comfortably on the knee rest pad. Use your uninjured leg to propel the scooter.
- **Steering:** Use the handlebars to steer the scooter in the desired direction. The steerable design allows for easy maneuverability.
- **Braking:** To slow down or stop, gently squeeze the handbrake levers. For parking, engage the locking mechanism on the handbrakes.

# Dual Braking System



Simply squeeze the brake handles to engage the lock.



Please press the unlocked handle according to the arrow position in the picture.



Image: Dual Braking System. Simply squeeze the brake handles to engage the lock. To release, press the unlocked handle according to the arrow position.

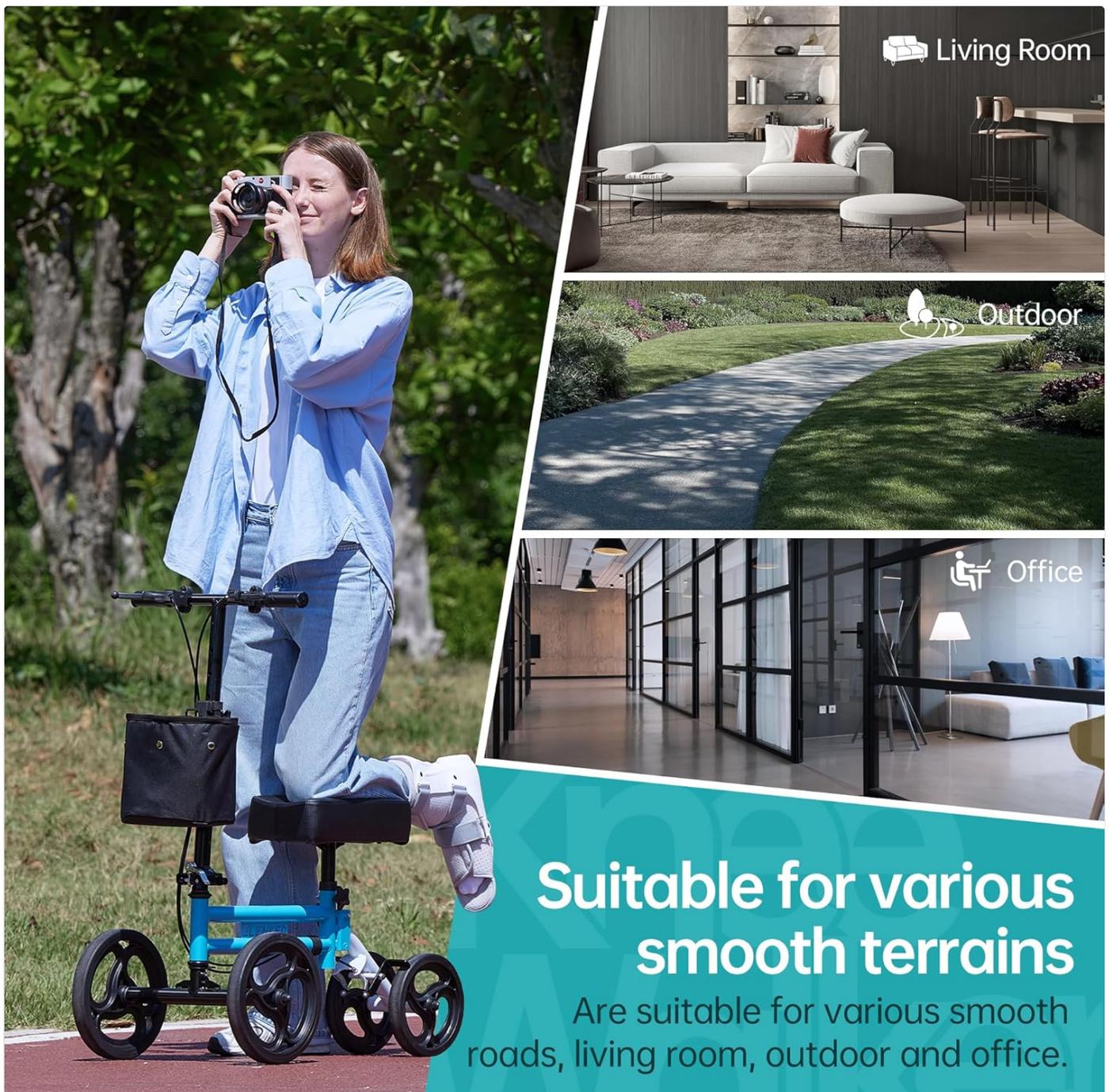


Image: Suitable for various smooth terrains. The knee scooter is designed for use on smooth roads, in living rooms, outdoors, and in office environments.

## MAINTENANCE

Regular maintenance will ensure the longevity and safe operation of your knee scooter:

- **Cleaning:** Wipe down the frame and knee pad with a damp cloth and mild detergent. Avoid abrasive cleaners.
- **Wheel Inspection:** Periodically check the wheels for wear and tear. Ensure they rotate freely and are free of debris.
- **Brake Check:** Regularly test the brakes to ensure they are responsive and hold securely. Adjust if necessary (refer to a qualified technician if unsure).
- **Fasteners:** Check all screws, bolts, and quick-release mechanisms to ensure they are tight and secure.
- **Storage:** Store the knee scooter in a dry, cool place away from direct sunlight when not in use. Its foldable design makes storage convenient.

## TROUBLESHOOTING

---

Common issues and their potential solutions:

| Problem                     | Possible Cause  | Solution  |
|-----------------------------|---|---|
| Brakes not engaging fully   | Cable tension too loose; worn brake pads.                     | Adjust brake cable tension; inspect and replace brake pads if worn.   |
| Scooter feels unstable      | Loose fasteners; uneven surface; incorrect height adjustment. | Check and tighten all fasteners; use on smooth, level surfaces; ensure handlebars and knee rest are properly adjusted and locked. |
| Wheels not rolling smoothly | Debris in wheels; worn bearings.                              | Remove any debris from wheels; if bearings are worn, contact customer support for replacement.                                    |

## SPECIFICATIONS

---

| Feature                 | Detail   |
|-------------------------|--|
| Brand                   | ELENKER  |
| Model                   | B0CWQFDCNP   |
| Color                   | Light Blue+white   |
| Material                | Iron   |
| Wheel Size              | 8" PVC wheels  |
| Braking System          | Dual adjustable locking handbrakes and dual rear on-wheel brakes |
| Foldable                | Yes, for easy storage and transportation                         |
| Package Dimensions      | 20 x 18.5 x 10.4 inches  |
| Item Weight             | 24.75 Pounds   |
| Age Range (Description) | Adult  |

# Adjustable Design

Adjustable knee pad and handlebars to fit you well


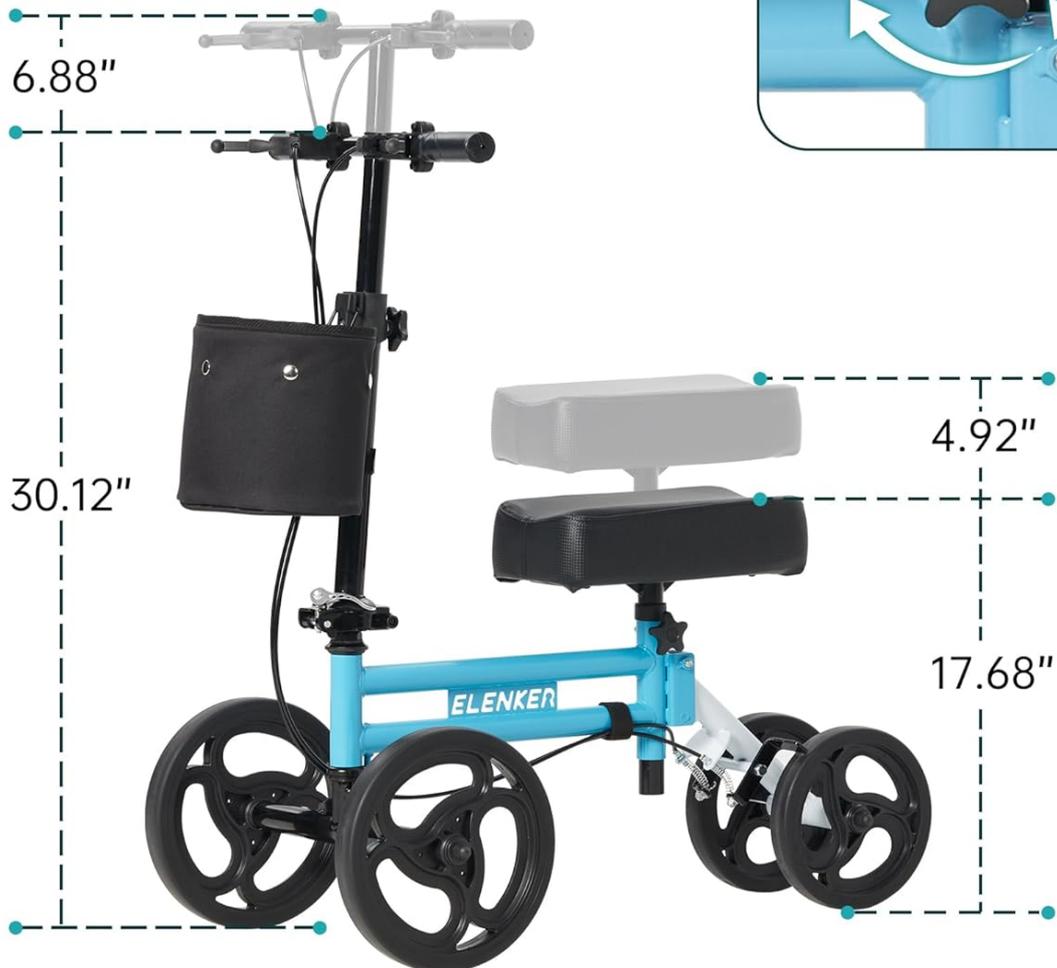


Image: Adjustable Design. Adjustable knee pad and handlebars to fit you well, with key dimensions indicated.

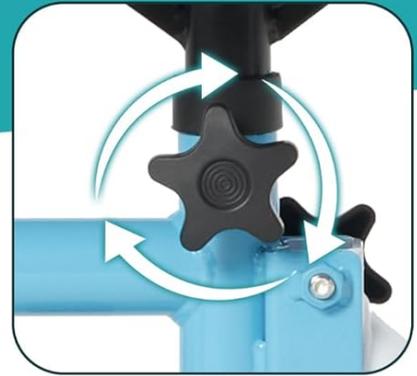

## 10" Front Wheels

Wheels with a wider contact surface, it is more non-slip, durable and provide greater stability for those with poor balance.

Image: 10" Front Wheels. Wheels with a wider contact surface provide non-slip durability and greater stability, especially for those with poor balance.

### CUSTOMER SUPPORT

---

For any questions, concerns, or assistance with your ELENKER Knee Scooter, please contact our customer support team.

- **Email Support:** Available 24/7 for your convenience.
- **Hotline:** 8-hour hotline for direct assistance.

ELENKER is committed to providing high-quality products and excellent customer satisfaction. We continuously update our designs and features to improve user experience.



Image: ELENKER's Commitment. Highlighting customer service, trusted brand status, high-quality products, and high customer satisfaction.



© 2024 ELENKER. All rights reserved.

For more information, visit [ELENKER Store](#)