

Rollz Rollz Motion 2.1

Rollz Motion 2.1 Rollator Walker Wheelchair Combo User Manual

Model: Dark Purple Regular

INTRODUCTION

The Rollz Motion 2.1 is an innovative 2-in-1 mobility aid designed to function as both a rollator walker and a transport wheelchair. This device offers flexibility, allowing users to transition between supported walking and being pushed in a wheelchair, enhancing independence and comfort. It features a sleek, foldable design and is crafted from high-quality aluminum for durability and ease of use.

Key features include:

- **2-in-1 Functionality:** Seamlessly converts from a rollator walker to a transport wheelchair.
- **Adjustable Handles:** Accommodates users of various heights for comfortable upright walking.
- **Foldable Design:** Easily folds into a compact unit for transport and storage.
- **Safe Braking System:** Ensures secure stopping in diverse conditions.
- **Curb Assistant:** Facilitates navigating small obstacles.

SETUP

Unfolding the Rollz Motion 2.1

To unfold the Rollz Motion 2.1, ensure the brakes are engaged. Gently pull the handles outwards until the frame locks into place. Confirm all locking mechanisms are secure before use.

Your browser does not support the video tag.

Video: Demonstration of folding and unfolding the Rollz Motion 2.1 rollator.

Adjusting Handle Height

The handles can be adjusted to suit your height for optimal comfort and posture. Locate the height adjustment buttons or levers on the handle shafts. Press and hold to extend or retract the handles to the desired position. Release the buttons/levers to lock them in place. Ensure both handles are set to the same height.

Your browser does not support the video tag.

Video: Guide on how to set up the handle height for the Rollz Motion 2.1.

Available in 2 sizes

Adjustable handles

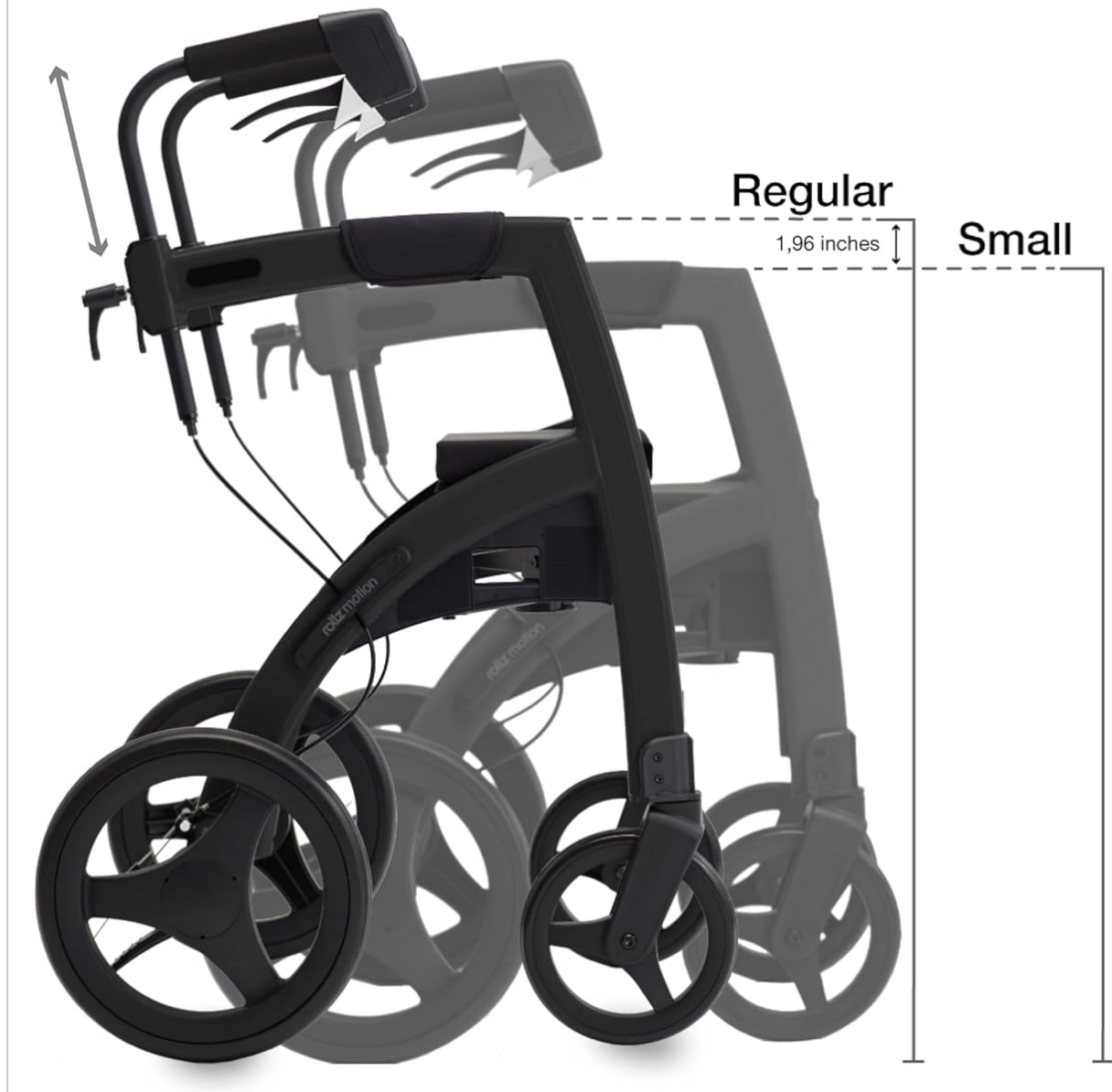


Image: The Rollz Motion 2.1 with adjustable handles, illustrating the height range for both Regular and Small sizes.

Converting to Wheelchair Mode

To convert the Rollz Motion 2.1 from a rollator to a wheelchair:

1. Engage the parking brakes.
2. Remove the backrest from its rollator position.
3. Flip the seat cushion forward to reveal the wheelchair seat.
4. Attach the footrests to the designated slots at the front of the frame.
5. Re-attach the backrest to the wheelchair position, ensuring it is securely locked.

Your browser does not support the video tag.

Video: Step-by-step demonstration of converting the Rollz Motion 2.1 from rollator to wheelchair and back.

Effortlessly transforms from rollator walker to wheelchair and back



Image: Visual representation of the Rollz Motion 2.1's transformation capability between rollator and wheelchair modes.

OPERATING THE ROLLZ MOTION 2.1

Using as a Rollator

When using the Rollz Motion 2.1 as a rollator, stand upright between the handles. Push the rollator forward, maintaining a comfortable distance. Use the handbrakes to control speed and stop. For temporary seating, engage the parking brakes and sit on the padded seat.



Always the right support with you

Upright rollator
walker for the most
convenient walking
journey

Image: A user demonstrating upright walking with the Rollz Motion 2.1 rollator.

Using as a Wheelchair

In wheelchair mode, the user sits on the seat with their feet on the footrests. The caregiver pushes the wheelchair using the handles. Ensure the user is seated comfortably and securely, and the footrests are properly positioned. Always engage the parking brakes when stopping or transferring the user.



Image: The Rollz Motion 2.1 in wheelchair mode, with a user being assisted.

Braking System

The Rollz Motion 2.1 is equipped with a reliable braking system. Squeeze the handbrakes to slow down or stop. Push the brake levers down to engage the parking brakes, which will lock the wheels for stability when sitting or transferring.

Curb Assistant

The integrated curb assistant allows for easier navigation over small curbs or obstacles. Place your foot on the curb assistant lever located near the rear wheels and gently tilt the rollator/wheelchair backward to lift the front wheels over the obstacle.

rollz motion²



Image: An annotated image of the Rollz Motion 2.1, pointing out the curb assistant and other key features.

MAINTENANCE

Regular maintenance ensures the longevity and safe operation of your Rollz Motion 2.1.

- **Cleaning:** Wipe down the frame with a damp cloth and mild detergent. Avoid abrasive cleaners.
- **Brakes:** Periodically check the brake cables and levers for proper function. Ensure they engage smoothly and hold securely. Adjust if necessary, or consult a professional.
- **Wheels and Tires:** Inspect wheels for wear and tear. The soft tires are designed to absorb vibrations; ensure they are free from damage.
- **Folding Mechanism:** Keep the folding mechanism clean and free of debris to ensure smooth operation.
- **Storage:** When not in use, store the Rollz Motion 2.1 in a dry, clean place, preferably folded to save space.

TROUBLESHOOTING

If you encounter any issues with your Rollz Motion 2.1, refer to the following common solutions:

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Brakes not engaging properly	Loose brake cables or worn brake pads.	Check cable tension and adjust. Inspect brake pads for wear and replace if necessary. Consult a professional if unsure.
Difficulty folding/unfolding	Debris in the mechanism or stiff joints.	Clean the folding mechanism. Apply a silicone-based lubricant to moving parts if needed.
Wheels not rolling smoothly	Dirt in wheel bearings or damaged wheels.	Clean wheels and axles. Check for any obstructions. If damage is present, contact customer support.
Unstable when walking	Handles not adjusted correctly or uneven ground.	Ensure handles are at the correct height and locked. Use on even surfaces.

If the issue persists after attempting these solutions, please contact Rollz customer support for further assistance.

SPECIFICATIONS

The following specifications apply to the Rollz Motion 2.1 Dark Purple Regular model:

- **Material Type:** Aluminum
- **Walker Type:** Rollator Walker
- **Is Foldable:** Yes
- **Maximum Height Recommendation:** 76.75 Inches
- **Color:** Dark Purple
- **Weight Capacity Maximum:** 275 Pounds
- **Item Weight:** 23.7 Pounds
- **Model Name:** Rollz Motion 2.1
- **Manufacturer Part Number:** 8718347120634

WARRANTY AND SUPPORT

The Rollz Motion 2.1 comes with a 30-day easy return policy. For specific warranty details, product support, or to purchase replacement parts, please refer to the official Rollz website or contact your authorized dealer. Rollz International is committed to providing quality mobility solutions and customer satisfaction.

For additional information and resources, visit the [Rollz Store on Amazon](#).