

Medline G2-105RX1

Medline Premium Series Shower Chair with Backrest (Model G2-105RX1) - User Manual

For safe and comfortable use in shower and bathtub.

1. INTRODUCTION

The Medline Premium Series Shower Chair with Backrest is designed to provide enhanced safety and comfort during showering or bathing. This plastic shower stool features slip-resistant feet and adjustable height settings, offering stability and confidence for individuals requiring assistance, including seniors, the elderly, and those with disabilities.

2. SAFETY INFORMATION

Please read all instructions carefully before assembly and use. Failure to follow these instructions may result in injury.

- **Weight Capacity:** This shower chair is designed to support a maximum user weight of **350 lbs (158 kg)**. Do not exceed this limit.
- **Stability:** Ensure all four legs are evenly adjusted and firmly on the floor before use. The slip-resistant rubber feet are crucial for preventing movement.
- **Assembly:** Verify that all components are securely fastened after assembly and before each use.
- **Inspection:** Regularly inspect the chair for any signs of wear, damage, or loose parts. Discontinue use if any damage is found and contact customer support.
- **Surface:** Use the chair on a flat, stable surface. Avoid uneven or slippery surfaces.
- **Adjustments:** Do not attempt to adjust the chair height while seated on it.

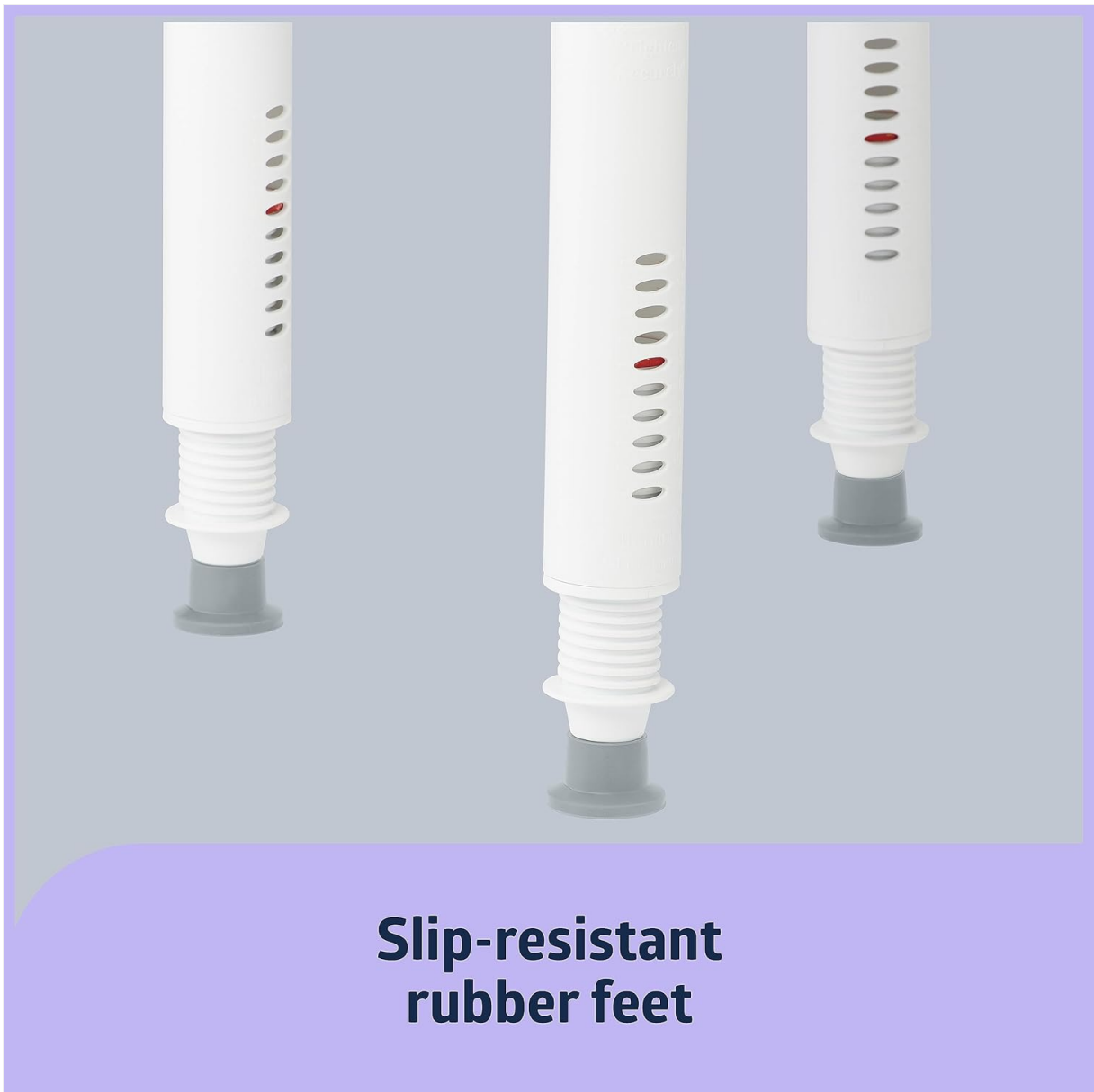


Image: Close-up view of the shower chair's legs with slip-resistant rubber feet, highlighting their design for stability.

3. PACKAGE CONTENTS

Before assembly, ensure all parts are present and undamaged. The package should contain the following components:

- 1x Seat Panel
- 1x Backrest
- 2x Armrests
- 4x Leg Extensions with Rubber Feet



Image: All individual components of the shower chair, including the seat, backrest, armrests, and four leg extensions, laid out on a flat surface, indicating "Easy to assemble - no tools required."

4. SETUP AND ASSEMBLY

The Medline Premium Series Shower Chair is designed for easy, tool-free assembly. Follow these steps:

1. **Attach Legs to Seat:** Insert each of the four leg extensions into the designated openings on the underside of the seat panel. Ensure they click securely into place.
2. **Attach Backrest:** Slide the backrest into the slots at the rear of the seat panel until it locks firmly.
3. **Attach Armrests:** Insert the armrests into the side slots of the seat panel, ensuring they are fully engaged.
4. **Adjust Height:** To adjust the chair's height, press the push-pins on each leg and slide the inner leg extension up or down to the desired height. Ensure the red lines on the legs are aligned and all four legs are set to the same height for stability. The seat height can be adjusted from 16 inches to 20.5 inches.
5. **Final Check:** Before first use, gently shake the chair to confirm all connections are secure and the chair is stable.

5. OPERATING INSTRUCTIONS (USAGE)

Using your Medline Shower Chair safely and effectively:

- **Placement:** Position the shower chair securely in your shower or bathtub on a flat, non-slippery surface. Ensure the rubber feet have good contact with the floor.
- **Sitting:** Approach the chair carefully. Use the armrests for support when sitting down or standing up. Ensure you are fully seated before beginning to shower.
- **Backrest Support:** The backrest provides additional support and comfort. Lean back gently as needed.
- **Stability:** Avoid sudden movements or leaning excessively to one side, which could compromise stability.

**Easy to assemble—
no tools required**



Image: The Medline shower chair with backrest and armrests positioned inside a shower stall, demonstrating its intended use for reliable safety and stability.

Healthcare's **#1**
choice in bath safety aids



Comfortable and sturdy backrest for extra support

Image: A full view of the Medline shower chair, highlighting the comfortable and sturdy backrest and armrests designed for extra support.

6. MAINTENANCE AND CLEANING

Proper maintenance ensures the longevity and hygiene of your shower chair.

- **Cleaning:** Clean the chair regularly with mild soap and water. Rinse thoroughly and wipe dry to prevent water spots and mildew. Do not use abrasive cleaners or harsh chemicals.
- **Drying:** Allow the chair to air dry completely after each use to prevent mold and mildew growth.
- **Storage:** Store the chair in a dry place when not in use for extended periods.
- **Inspection:** Periodically check all connections, especially the leg extensions and rubber feet, for wear or damage. Replace worn parts immediately.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Chair feels unstable or wobbly.	Legs are not adjusted to the same height or not fully extended/locked. Connections are loose.	Ensure all four legs are adjusted to the same height and the push-pins are fully engaged. Check all assembly connections for tightness. Ensure the chair is on a flat surface.
Leg adjustment mechanism is stiff or difficult to operate.	Accumulation of soap residue or mineral deposits.	Clean the leg extensions and push-pin areas with warm, soapy water. Ensure no debris is obstructing the mechanism.
Rubber feet are slipping.	Feet are worn, dirty, or the shower floor is excessively slippery.	Clean the rubber feet and the shower floor. If feet are worn, they may need replacement. Ensure the chair is placed on a non-slip surface.

8. SPECIFICATIONS

Feature	Detail
Model Number	G2-105RX1
Brand	Medline
Material	Plastic, Rubber
Color	White
Maximum Weight Capacity	350 Pounds (158 kg)
Adjustable Seat Height	16 to 20.5 Inches
Seat Depth	15 Inches
UPC	193489038715

Reliable safety and stability



Image: A diagram illustrating the key dimensions of the shower chair, including adjustable height (16-20.5 inches), seat width (22.05 inches), and base width (19.95 inches).

9. WARRANTY AND SUPPORT

For warranty information, replacement parts, or technical support, please contact Medline customer service. Keep your purchase receipt and model number (G2-105RX1) handy when contacting support. You can visit the official Medline website for more information: www.medline.com