

CONTEC NE-M01VET

CONTEC NE-M01VET Nebulizer Machine User Manual

Model: NE-M01VET | Brand: CONTEC



Image: The CONTEC NE-M01VET Nebulizer Machine, a compact white device with a clear mask, positioned next to a cat, highlighting its veterinary application.

1. INTRODUCTION

The CONTEC NE-M01VET Nebulizer Machine is a portable and efficient device designed to deliver medication in the form of a fine mist for veterinary applications. Its compact size and quiet operation make it suitable for use at home, during travel, or in outdoor settings, providing a comfortable experience for animals requiring respiratory support. This manual provides essential information for the safe and effective operation, maintenance, and troubleshooting of your nebulizer.



Image: A banner showcasing the CONTEC NE-M01VET Pet Nebulizer with icons representing its key features: Low Noise, Softer Mist, Easy to Use, AC/Battery Charging, and User-friendly Material.

2. SAFETY INFORMATION

Please read all safety instructions carefully before using the nebulizer. Failure to follow these instructions may result in injury or damage to the device.

- Always consult a veterinarian before administering any medication to your pet using this device.
- Use only solutions prescribed or recommended by a veterinarian.
- Keep the device out of reach of children and unsupervised pets.
- Do not immerse the main unit in water or any other liquid.
- Ensure proper ventilation during use.
- Clean the device thoroughly after each use to prevent contamination and ensure optimal performance.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

3. PRODUCT OVERVIEW AND COMPONENTS

The CONTEC NE-M01VET Nebulizer is designed for ease of use and effective mist delivery. Familiarize

yourself with its components:

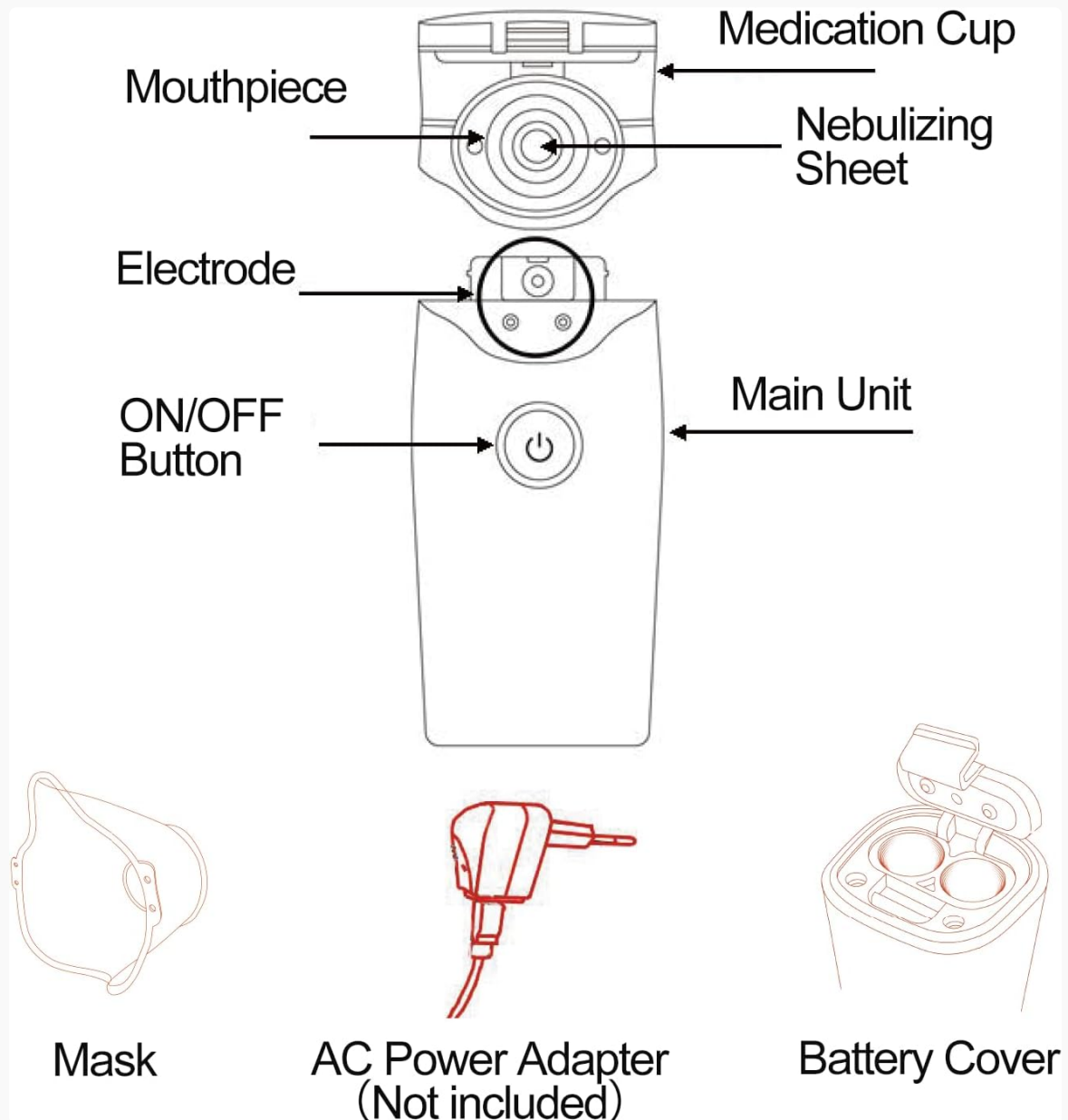


Image: An exploded diagram illustrating the main components of the nebulizer, including the main unit, medication cup, nebulizing sheet, electrode, ON/OFF button, mask, and battery cover.

Key Features:

- **Medical Nebulization:** Designed for effective delivery of prescribed solutions.
- **Fine Particle Size:** Produces a tiny spray with particles around 5 μm , easy to absorb for faster relief. Over 60% of mist particles are below 5 μm .
- **Low-Noise Design:** Operates virtually noiselessly, with a noise level as low as 25 dB, reducing anxiety for pets.
- **Automatic Switch-off:** The nebulizer automatically turns off when the medicine/liquid is almost used up.
- **Compact and Convenient:** Small in size and light in weight, making it highly portable for travel or home use.
- **One-button Operation:** Simple and convenient to use.

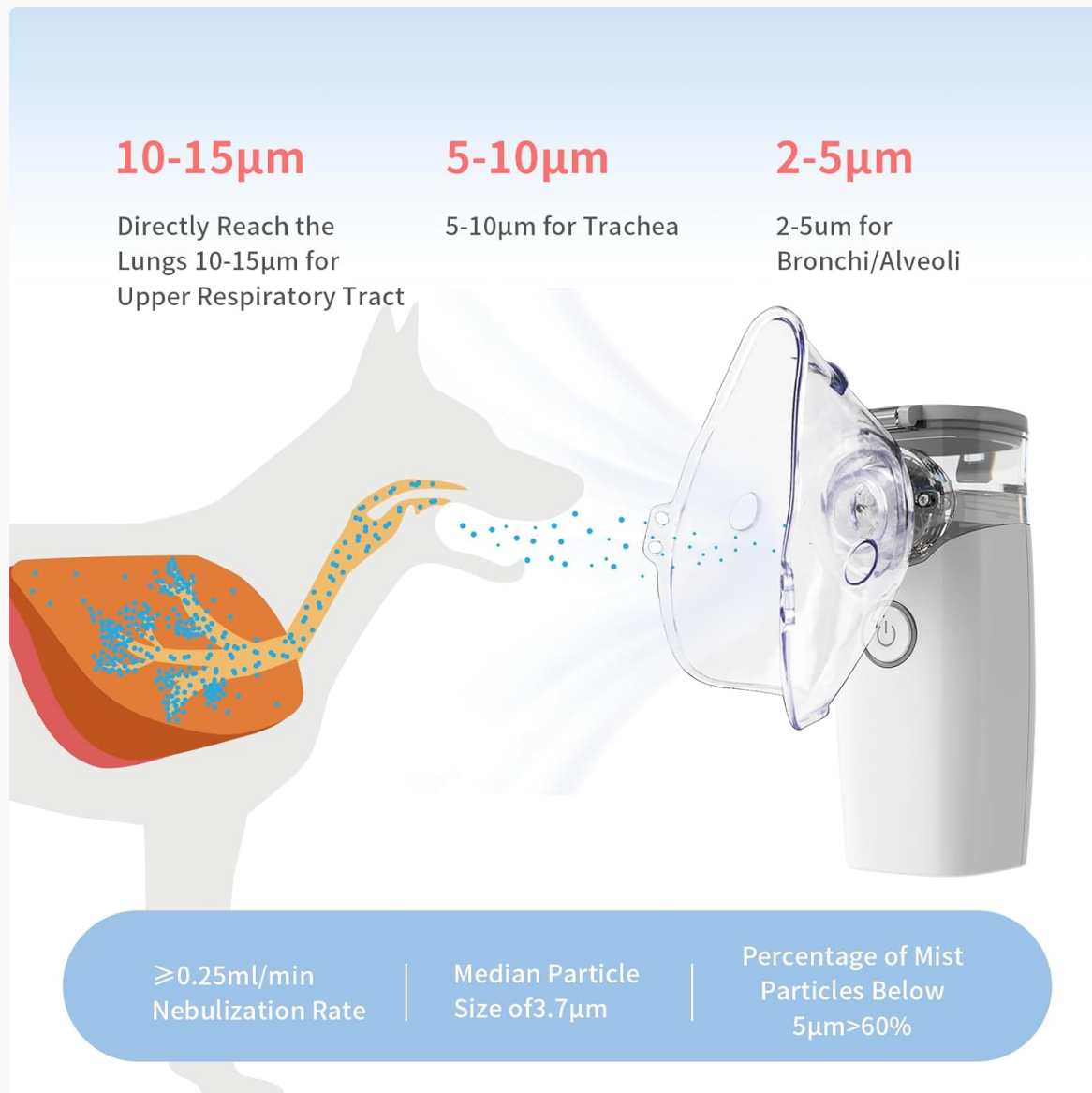


Image: A visual representation of the nebulizer's mist particle sizes and their target areas within the respiratory system, along with key performance metrics like nebulization rate and median particle size.

4. SETUP

Prepare your nebulizer for first use by following these steps:

Power Supply:

The NE-M01VET can be powered in two ways:

- **USB Cable:** Connect the provided USB cable to the device and a compatible power source (e.g., USB adapter, power bank).
- **Batteries:** Insert two AA batteries into the battery compartment located at the bottom of the device. Ensure correct polarity.

Two Power Supply Mode



powered by plug-in adapter



powered by 2AA batteries
(not included)



Image: A visual guide demonstrating the two power supply options for the nebulizer: via a plug-in adapter or using two AA batteries.

Assembling the Device:

1. Open the medication cup by sliding the top cover.
2. Pour the prescribed liquid medication into the medication cup. The maximum capacity is 10ml.
3. Close the medication cup securely.
4. Attach the appropriate mask (large or small) or mouthpiece to the nebulizer unit. The masks are designed to fit securely.

Installation and Use



Image: A series of six images demonstrating the installation and use steps, including battery insertion, liquid filling, and mask attachment.

5. OPERATING INSTRUCTIONS

Once assembled and powered, operating the nebulizer is straightforward:

1. Press the ON/OFF button to start the nebulization process. The device will begin producing a fine mist.
2. Position the mask or mouthpiece appropriately for your pet to inhale the mist.
3. The nebulizer operates quietly, minimizing stress for the animal during treatment.
4. The device features an automatic switch-off function, turning off when the medication is nearly depleted.

Video: An official CONTEC video demonstrating the nebulizer's mist generation and quiet operation, highlighting its efficiency and user-friendly design.

Video: An official CONTEC video showcasing the nebulizer from various angles and demonstrating its simple assembly process.

6. MAINTENANCE AND CLEANING

Proper maintenance and cleaning are crucial for the longevity and hygienic operation of your nebulizer.

- **After Each Use:** Disassemble the medication cup, mask/mouthpiece, and clean them thoroughly. It is recommended to use rubbing alcohol and distilled water for cleaning.
- **Main Unit:** Wipe the main unit with a soft, damp cloth. Do not submerge it in water.
- **Drying:** Ensure all components are completely dry before reassembling or storing the device.
- **Storage:** Store the nebulizer in a clean, dry place away from direct sunlight and extreme temperatures.

7. TROUBLESHOOTING

If you encounter issues with your CONTEC NE-M01VET Nebulizer, please refer to the following common solutions:

- **No Mist Production:**
 - Ensure the device is properly charged or that fresh batteries are installed with correct polarity.
 - Check if the medication cup is securely attached and contains sufficient liquid.
 - Verify that the nebulizing sheet is clean and free from blockages.
 - Some users have reported that aligning the cup correctly by following specific button pressing procedures (often detailed in a small slip of paper in the box) can resolve issues.
- **Weak Mist:**
 - Ensure the medication cup is not overfilled or underfilled.
 - Clean the nebulizing sheet thoroughly.
- **Device Not Turning On:**
 - Check battery installation or USB connection.
 - Ensure the ON/OFF button is pressed firmly.

For persistent issues not resolved by these steps, please contact CONTEC customer support.

8. SPECIFICATIONS

Specification	Detail
Product Dimensions	5.31 x 4.13 x 3.54 inches
Item Weight	13.12 ounces
Model Number	NE-M01VET
Manufacturer	CONTEC MEDICAL SYSTEMS CO., LTD
First Available	April 26, 2022

Image: A diagram illustrating the nebulizer's advanced absorption technology, featuring SUS alloy mesh, 1700 holes, a low nebulization rate, and quiet operation.

9. WARRANTY AND SUPPORT

For warranty information, technical support, or any inquiries regarding your CONTEC NE-M01VET Nebulizer

Machine, please contact CONTEC Medical Systems Co., Ltd. directly. Refer to the product packaging or the official CONTEC website for the most up-to-date contact details.

© 2024 CONTEC Medical Systems Co., Ltd. All rights reserved.