Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Easy@Home](#) /
- › [Easy@Home Portable Nebulizer User Manual](#)

Easy@Home YS35

Easy@Home Portable Nebulizer User Manual

Model: YS35 | Brand: Easy@Home

INTRODUCTION

The Easy@Home Portable Nebulizer is designed to provide efficient and quiet nebulization for both adults and children. Utilizing advanced micro-mesh technology, it converts liquid medication into fine mist particles, allowing for effective delivery to the respiratory system. Its compact and lightweight design makes it ideal for use at home or while traveling. This manual provides detailed instructions for the proper setup, operation, maintenance, and troubleshooting of your nebulizer.

Image: The Easy@Home Portable Nebulizer, shown with its included adult mask, child mask, and mouthpiece.

PRODUCT OVERVIEW & WHAT'S IN THE BOX

Your Easy@Home Portable Nebulizer package includes all necessary components for immediate use. Please verify that all items listed below are present and in good condition upon unboxing.

- Main Nebulizer Unit
- Medicine Cup
- USB Charging Cable
- Mouthpiece
- Child Mask
- Adult Mask

- Carry Pouch
- User Manual (this document)

Image: A complete set of Easy@Home Portable Nebulizer components, including the main unit, masks, mouthpiece, USB cable, and carry pouch.

SETUP

1. **Initial Charge:** Before first use, fully charge the nebulizer using the provided USB cable. Connect the USB cable to the nebulizer and a compatible USB power source (DC 3-5V, 1-2A), such as a laptop, power bank, or wall adapter (adapter not included). The charging indicator light will show the charging status.

USB Rechargeable

Output: DC3~5V,1~2A
(Adaptor not included)

Image: The nebulizer being charged via USB, illustrating various compatible power sources like laptops, power banks, and wall sockets.

2. **Prepare Medicine Cup:** Open the medicine cup by gently twisting or sliding it off the main unit, depending on the model. Ensure it is clean and dry before adding medication.
3. **Add Medication:** Pour the prescribed liquid medication into the medicine cup. Do not overfill; adhere to the maximum fill line indicated on the cup.
4. **Attach Accessory:** Securely attach either the mouthpiece, adult mask, or child mask to the nebulizer unit, depending on the user's preference and age. Ensure a snug fit to prevent leakage during operation.

OPERATING INSTRUCTIONS

1. **Power On:** Press and hold the power button for a few seconds until the indicator light illuminates and the nebulizer begins to produce mist.
2. **Nebulization:** Place the mouthpiece in your mouth or position the mask securely over your nose and mouth.

Breathe normally and deeply to inhale the medication mist. The nebulizer operates with a super silent design, ensuring minimal disturbance, even for sleeping children.

Super Quiet

Quiet operation,
baby sleep is not disturbed.

Image: The nebulizer's super quiet operation ensures that even a sleeping baby is not disturbed during use.

3. **Micro-Mesh Technology:** The device utilizes advanced micro-mesh technology to produce fine particles (2.5-5 microns), which are optimal for reaching the lower respiratory tract for effective treatment.

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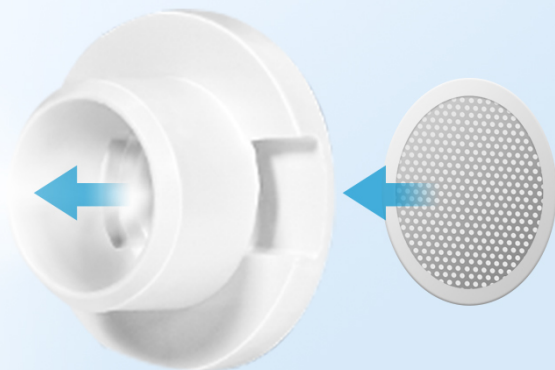
Image: A visual representation of the mesh vibration technology, showing how liquid medication is atomized into fine particles for efficient delivery.

- 4. Automatic Shut-off:** The nebulizer is equipped with an automatic shut-off feature that activates after 20 minutes of continuous use. This conserves energy and protects the device's performance.

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Image: The nebulizer automatically turns off after 20 minutes of use to save energy and protect the device.

- 5. Power Off:** Once the medication is fully nebulized or your treatment is complete, press and hold the power button again to turn off the device.

Image: A child demonstrating the proper use of the nebulizer with the child-sized mask.

Image: An adult demonstrating the proper use of the nebulizer with the adult-sized mask.

Image: An adult demonstrating the proper use of the nebulizer with the mouthpiece.

MAINTENANCE

Proper cleaning and maintenance are crucial for the longevity and hygienic operation of your nebulizer. Always clean the device after each use.

1. **Disassemble:** After each use, detach the medicine cup and the mask/mouthpiece from the main unit.
2. **Rinse Components:** Rinse the medicine cup, mask, and mouthpiece thoroughly with clean, warm water. Shake off excess water and allow them to air dry completely.
3. **Automatic Cleaning Mode:** For thorough cleaning of the internal mesh, pour pure water into the medicine cup. With the nebulizer turned off, press and hold the power button for five seconds. The nebulizer will enter automatic cleaning mode, indicated by the illumination of three indicator lights. This mode helps prevent residue buildup.

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Image: Activating the one-button automatic cleaning mode by pressing and holding the power button.

4. **Sterilization:** Periodically, you may sterilize the medicine cup and accessories using a solution of alcohol or vinegar, followed by a thorough rinse with pure water.
5. **Clean Main Unit:** Wipe the main unit with a clean, damp cloth. Do not immerse the main unit in water.

6. **Storage:** Store the nebulizer and its accessories in a clean, dry place, preferably in the provided carry pouch, when not in use.

TROUBLESHOOTING

If you encounter issues with your Easy@Home Portable Nebulizer, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
No mist produced	Low battery; Medicine cup empty; Mesh clogged; Improper assembly.	Charge the device; Refill medicine cup; Perform automatic cleaning; Reassemble correctly.
Weak mist production	Low battery; Insufficient medication; Mesh partially clogged.	Charge the device; Ensure adequate medication; Perform automatic cleaning.
Device not turning on	Battery depleted; Power button not pressed correctly.	Charge the device fully; Press and hold the power button firmly for a few seconds.
Leakage from medicine cup	Medicine cup not properly sealed; Overfilled.	Ensure medicine cup is securely attached; Do not exceed the maximum fill line.

If the problem persists after attempting these solutions, please contact Easy@Home customer support for further assistance.

SPECIFICATIONS

Detailed technical specifications for the Easy@Home Portable Nebulizer (Model YS35):

Feature	Specification
Model Number	YS35
Brand	Easy@Home
Units	1.00 unit
Country of Origin	China
Battery Type	1 Lithium-ion (included)
Package Dimensions	14.7 x 14.2 x 6.8 cm; 300 grams
Particle Size	2.5 - 5 microns (Micro-Mesh Technology)
Charging Input	DC 3-5V, 1-2A (USB)

Feature	Specification
Automatic Shut-off	After 20 minutes

WARRANTY AND SUPPORT

Easy@Home is committed to providing high-quality products and customer satisfaction. For any questions, concerns, or warranty claims regarding your Easy@Home Portable Nebulizer, please contact our customer support team. Refer to the contact information provided on the product packaging or the official Easy@Home website for the most up-to-date support details.

Please retain your purchase receipt as proof of purchase for warranty purposes.