

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [YUNCHI](#) /

› [YUNCHI Waterproof Extra Wide Leg Cast Cover for Shower - Adult Instruction Manual](#)

YUNCHI Extra-wide leg cast cover for adults

YUNCHI Waterproof Extra Wide Leg Cast Cover for Shower - Adult Instruction Manual

Model: Extra-wide leg cast cover for adults

1. INTRODUCTION

The YUNCHI Waterproof Extra Wide Leg Cast Cover is designed to provide a watertight seal, protecting casts, bandages, wounds, and splints on the leg and foot during showering. This manual provides essential information for safe and effective use.



Figure 1: YUNCHI Waterproof Extra Wide Leg Cast Cover in use.

2. SAFETY INFORMATION

- Due to the potential risk of choking and suffocation, children should use this product under parental supervision.
- Made for therapeutic showering and bathing only. Not for sustained or vigorous water activities like diving.
- Keep away from fire and forbid drying from the dryer.
- Please air-dry in the shade between uses to maintain elasticity.
- Please use cautiously when walking on wet surfaces in case of slipperiness.

- Do not wear the cover for longer than 20 minutes to avoid cutting off blood circulation with extended wear.

3. PRODUCT FEATURES

- **Extra Large Size:** Designed for overweight/plus-size adults, thicker large casts, splints, or medical fracture boots.
- **Reusable & Durable Material:** Made of medical-grade PVC bag, thick, waterproof, and durable. The silicone hole is upgraded with an anti-tear strip for repeated use.
- **Comfortable & Safe Watertight Seal:** Non-latex material that is skin-friendly and forms a watertight seal below the knee to prevent wound infections and keep casts dry.
- **Easy to Put On and Off:** Designed to stretch easily over large casts without touching the wounded part. Hard ring design allows for easy one-handed application and removal.

4. PRODUCT DIMENSIONS & SIZING

This XL-Extra Wide model features a widened hard ring diameter of 9.5 inches and a center hole diameter of 3.7 inches. It fits calf circumferences from 13 to 28 inches. The overall length of the cover is 23 inches.



Figure 2: Product dimensions for the XL-Extra Wide Leg Cast Cover.

Note: This product is NOT SUITABLE FOR KNEE AND THIGH WOUNDS. Please measure your calf and cast size before purchasing to ensure proper fit.

EXTRA WIDE OVER 2.5" DESIGNED FOR LARGER CASTS/LEGS

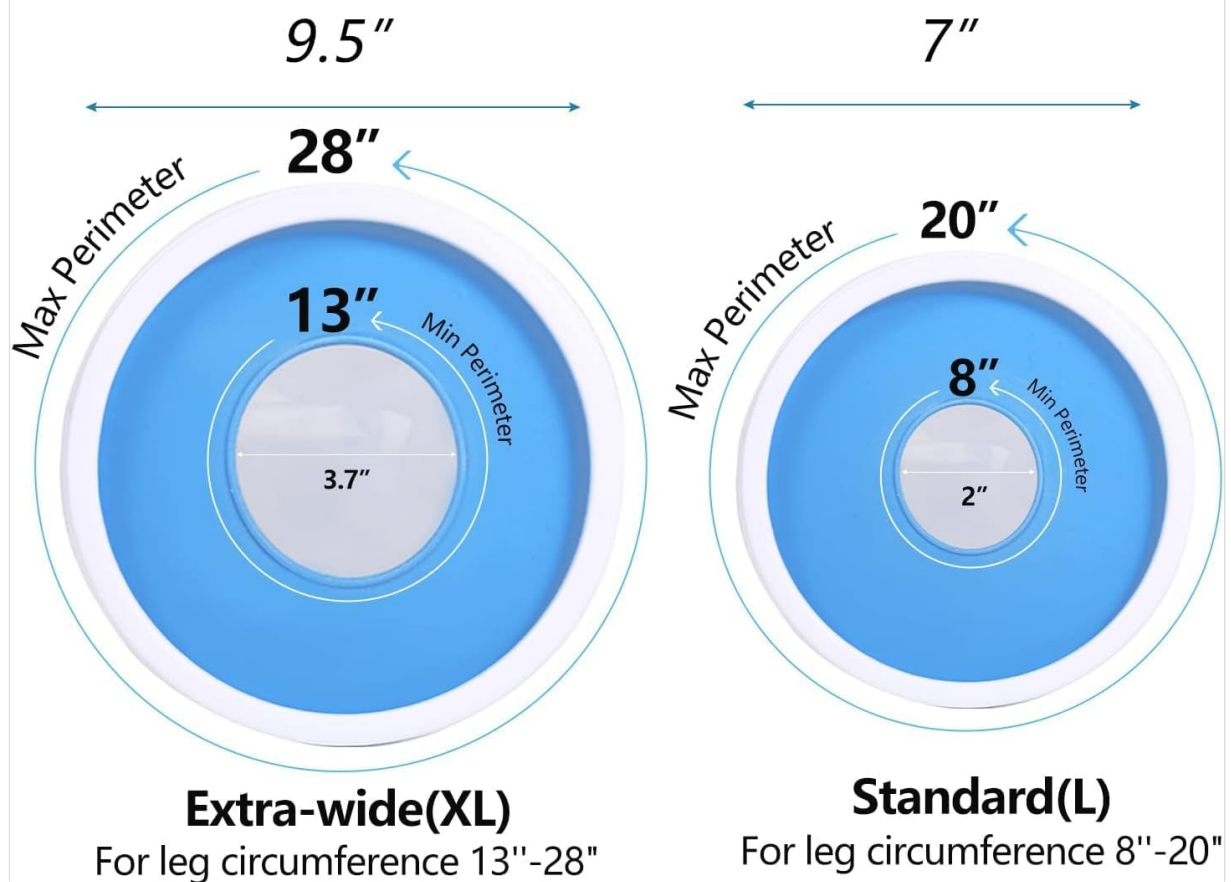


Figure 3: Size comparison between Extra-wide (XL) and Standard (L) models.

5. USAGE INSTRUCTIONS

5.1. Preparation

1. Before initial use, perform a waterproof test after wearing the cast cover to ensure it does not leak before bathing or showering.

5.2. Applying the Cast Cover

2. **Expand the Silicone Hole:** Gently stretch the silicone opening to create enough space for your leg.



Figure 4: Expanding the silicone opening.

3. **Slowly Extend Your Leg:** Carefully insert your leg into the cover, guiding your foot and cast through the expanded silicone hole. Ensure the cast or bandage does not snag on the silicone.

XL / EXTRA WIDE LEG CAST COVER

Wider than regular size, design for large legs/casts and boots

EXTRA WIDTH 2.8"
6" - 7"
EXTRA WIDTH 2.8"
9.8"

EXTRA-WIDE SIZE

VS

REGULAR SIZE

X

Figure 5: Slowly extending the leg into the cover.

4. **Adjust the Sealing Ring:** Once your leg is fully inside, adjust the white sealing ring to fit snugly against your skin, just below the knee. This creates the watertight seal.



Figure 6: Adjusting the sealing ring for a watertight fit.

5.3. During Use

5. **Safely Take a Bath or Shower:** With the cast cover securely in place, you can now safely shower or bathe without wetting your cast or wound.

Video 1: Demonstration of the YUNCHI Waterproof Extra Wide Cast Cover in a shower, showing its effectiveness in keeping the leg dry.

Video 2: A user demonstrating the YUNCHI cast cover in the shower, highlighting its ability to keep the cast dry.

5.4. Removal and Drying

6. **Remove the Cover:** Gently stretch the silicone opening again and slide the cover off your leg.
7. **Drying:** Use the provided S-hook to hang the cast cover in the shower or a well-ventilated area to air-dry. Ensure it is completely dry before storing.



Figure 7: Hanging the cast cover to air-dry using the provided S-hook.

6. CARE & MAINTENANCE

- After each use, rinse the cast cover with clean water.
- Air-dry the cover in the shade. Avoid direct sunlight or using a dryer, as this can damage the material and reduce its elasticity.
- Store in a cool, dry place when not in use.

7. IMPORTANT CONSIDERATIONS

- Always ensure the cast cover forms a complete watertight seal before entering the shower or bath.
- Do not submerge the entire leg and cover in water for extended periods.
- If you experience any discomfort, numbness, or changes in skin color while wearing the cover, remove it immediately.

8. SPECIFICATIONS

Feature	Detail
Brand	YUNCHI
Model Number	Extra-wide leg cast cover for adults
Size	XL-Extra wide
Calf Circumference Fit	13 to 28 inches
Hard Ring Diameter	9.5 inches
Center Hole Diameter	3.7 inches
Product Dimensions	11.81 x 10.24 x 1.18 inches; 11.68 ounces
Material	Medical-grade PVC bag, non-latex silicone
Specific Uses	Fracture, wounds, abrasions, rashes, burns, cuts
Use For	Ankle, Foot, Shins, Toes

9. WARRANTY AND SUPPORT

YUNCHI offers a 30-day refund policy for this product. For any inquiries or assistance, 24-hour online quick response after-sale service is available.