

HALIDODO KN95 Face Mask

HALIDODO 5-Layer KN95 Face Mask User Manual

Model: KN95 Face Mask

Introduction

Safety Information

Product
Overview

Troubleshooting

Specifications

Setup &
Wearing
Support

Care & Maintenance

1. INTRODUCTION

Thank you for choosing HALIDODO 5-Layer KN95 Face Masks. This manual provides essential information for the proper use, care, and maintenance of your face masks to ensure optimal protection and comfort. Please read this manual thoroughly before using the product.



A box containing 50 HALIDODO KN95 face masks, showcasing various colors like sky blue, pink, grey, wine red, and royal blue. The packaging highlights the 50 PCS quantity and 'Foldable Protective Mask' description.

2. IMPORTANT SAFETY INFORMATION

Disclaimer: This KN95 mask is NOT an N95 mask. For more information about KN95 masks, please refer to official guidelines.

- Keep the face mask clean before use.
- Do not touch the rear side of the face mask during use.
- Do not wash the face mask with water. Water will destroy the filtering ability and protective nature of the mask.
- Do not steam the face mask or put it in the microwave.
- Discard the face mask when it becomes visibly dirty or breathing resistance increases remarkably.
- Ensure a proper seal around the nose and mouth for effective filtration.

3. PRODUCT OVERVIEW AND FEATURES

The HALIDODO KN95 Face Mask is designed for enhanced protection and comfort, featuring a 5-layer

filtration system.

3.1 5-Ply Protective Structure

Each mask is constructed with five distinct layers:

1. **Outer Non-Woven Layer:** Provides initial filtration and splash resistance.
2. **First Melt-Blown Layer:** Blocks more than 95% of particles.
3. **Middle Non-Woven Fabric:** Absorbs moisture and provides structural support.
4. **Second Melt-Blown Layer:** Blocks more than 95% of particles. (Note: These two melt-blown layers may stick together during production; they can be separated by rubbing.)
5. **Inner Non-Woven Layer:** Soft and breathable for skin comfort.

This 5-ply design offers superior protection compared to standard 3-ply disposable masks.



An exploded view diagram showing the five distinct layers of the KN95 face mask: an outer non-woven fabric, two melt-blown fabric layers, a middle non-woven fabric, and an inner non-woven fabric, all contributing to a filter efficiency of $\geq 95\%$.

3.2 High Filtration Efficiency

The masks provide a minimum filter efficiency of 95%, effectively protecting against dust, PM 2.5 particles, haze, smoke, and automobile exhaust.

3.3 Comfortable & Breathable Design

- **High-Elastic Ear Loops:** Designed for a comfortable fit without straining the ears.
- **Adjustable Nose Clip:** Ensures a secure and customized fit, helping to prevent glasses from fogging.
- **Lightweight and Foldable:** Easy to carry and store.
- **Full Coverage:** Effectively covers the mouth, nose, and chin.


4. SETUP AND WEARING INSTRUCTIONS

Follow these steps to properly wear your HALIDODO KN95 Face Mask:

1. **Unfold the Mask:** Gently unfold the mask, ensuring the inner surface remains untouched.
2. **Position the Mask:** Hold the mask with the nose clip facing upwards. Place the mask over your nose and mouth.
3. **Secure Ear Loops:** Place an ear loop around each ear. Adjust them for a snug but comfortable fit. For a more sealed wearing, you can cross the earloops before placing them over your ears.
4. **Adjust Nose Clip:** Press down on the adjustable nose clip with both index fingers to mold it firmly to the shape of your nose bridge. This creates a tight seal and helps prevent glasses from fogging.
5. **Check for Seal:** Inhale sharply. If air leaks around the edges, readjust the mask and nose clip until a proper seal is achieved.

One Size Fits All

Cross earloop for more sealed wearing

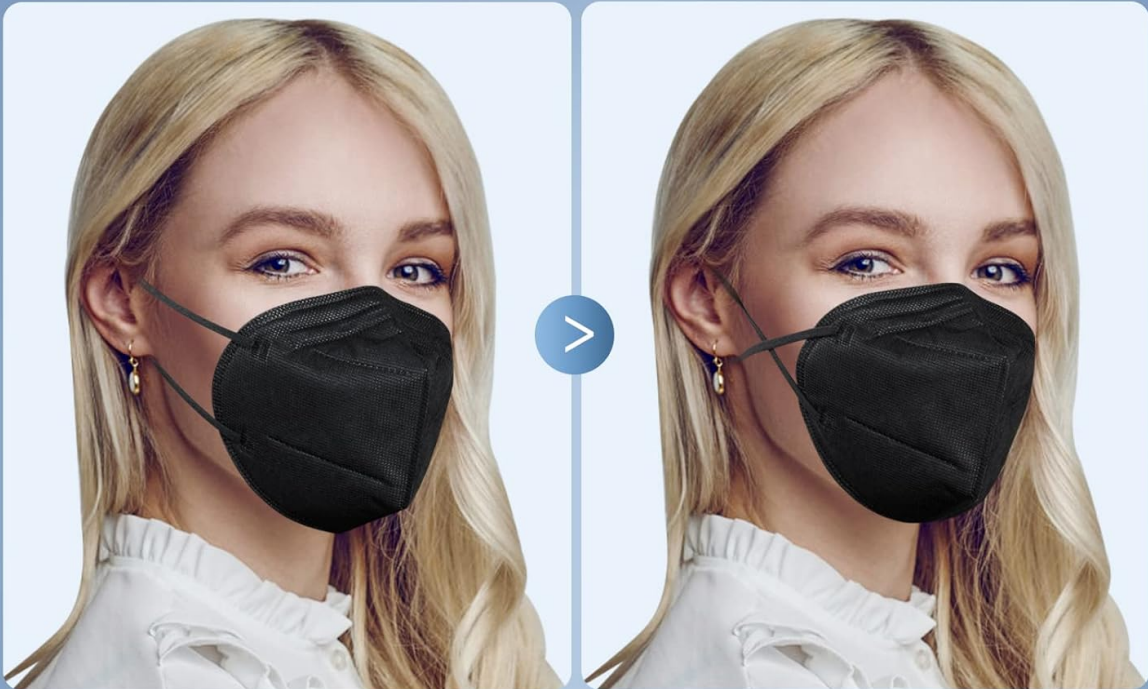

Normal

Tightened

SPECIFICATIONS

Length : 6.5inch

Width : 5.2inch

Range of earloop length : From 3.54 to 8.46inch

A comparison of a mask worn normally and tightened with cross earloops, illustrating how 'One Size Fits All' with a sealed wearing. Specifications are listed: Length 6.5 inches, Width 5.2 inches, and earloop length range from 3.54 to 8.46 inches.

5. OPERATING (USAGE) GUIDELINES

Once properly fitted, the HALIDODO KN95 Face Mask can be used in various environments where respiratory protection is desired. These masks are suitable for daily use by families, office staff, and warehouse workers.

- Wear the mask continuously in environments where protection is needed.
- Avoid touching the mask surface while wearing it. If you must touch it, wash or sanitize your hands immediately afterward.
- If the mask becomes damp, damaged, or breathing becomes difficult, replace it with a new one.

6. CARE AND MAINTENANCE

Proper care ensures the effectiveness and hygiene of your mask:

- **Storage:** Store unused masks in a clean, dry place, away from direct sunlight and contaminants.

- **Disposal:** Dispose of used masks responsibly in a closed waste bin. Do not reuse disposable masks.
- **Cleaning:** Do not wash, steam, or microwave the mask. These methods will compromise the filter material and reduce its protective capabilities.
- **Replacement:** Replace the mask if it becomes visibly dirty, damaged, or if you experience increased breathing resistance.

7. TROUBLESHOOTING COMMON ISSUES

7.1 Glasses Fogging

If your glasses fog up while wearing the mask, it indicates air is escaping from the top of the mask and reaching your lenses.

- **Solution:** Re-adjust the nose clip to ensure it is tightly molded to the bridge of your nose. Ensure the top edge of the mask is snug against your face.



Two side-by-side images demonstrating the anti-fogging feature. On the left, a person wearing a HALIDODO KN95 mask has clear glasses. On the right, a person wearing a different KN95 mask has fogged glasses, indicating better breathability and seal with the HALIDODO mask.

7.2 Discomfort Around Ears

If the ear loops cause discomfort, try adjusting their position or using a mask extender (not included) to relieve pressure.

7.3 Difficulty Breathing

If you experience significant difficulty breathing, ensure the mask is not wet or visibly dirty. If it is, replace it immediately. If the issue persists with a new, clean mask, consult a healthcare professional.

8. SPECIFICATIONS

Feature	Detail
Product Type	KN95 Face Mask
Layers	5-Ply (2 non-woven outer, 2 melt-blown inner, 1 non-woven middle)
Filter Efficiency	≥95%
Dimensions (Length)	Approximately 6.5 inches
Dimensions (Width)	Approximately 5.2 inches
Earloop Length Range	3.54 to 8.46 inches
Origin	Imported
Model ASIN	B0FKMQPWWG

9. WARRANTY INFORMATION

HALIDODO disposable face masks are intended for single use and do not come with a specific warranty. Please inspect the product upon receipt for any manufacturing defects. If you encounter any issues, please contact customer support immediately.

10. CUSTOMER SUPPORT

For further assistance, questions, or concerns regarding your HALIDODO KN95 Face Masks, please visit the official HALIDODO store on Amazon or contact their customer service directly.

HALIDODO Store: [Visit the HALIDODO Store on Amazon](#)